
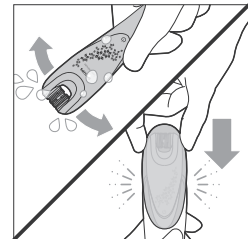
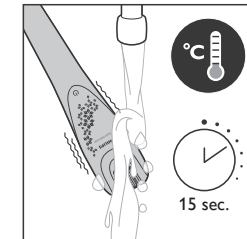
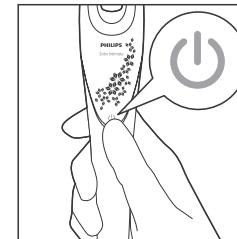
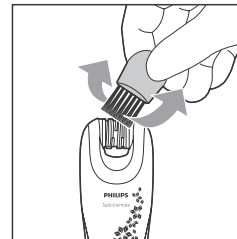
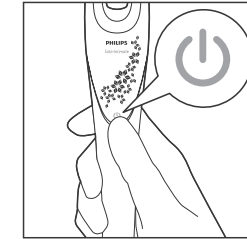
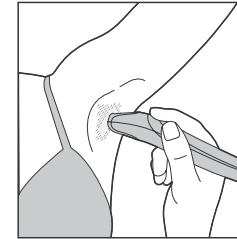
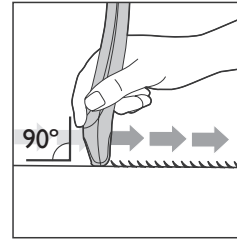
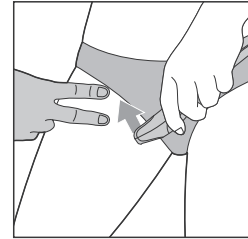
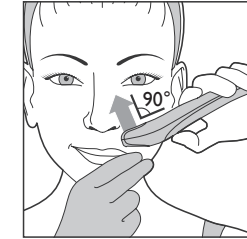
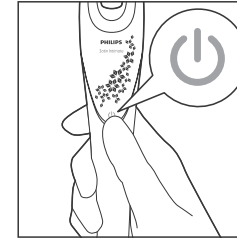
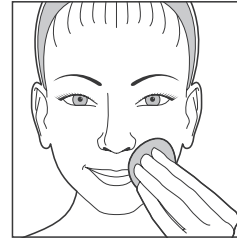
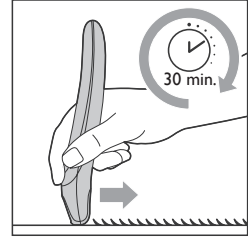
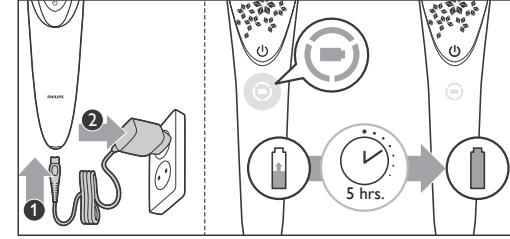
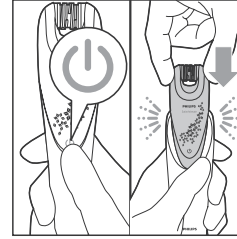
PHILIPS

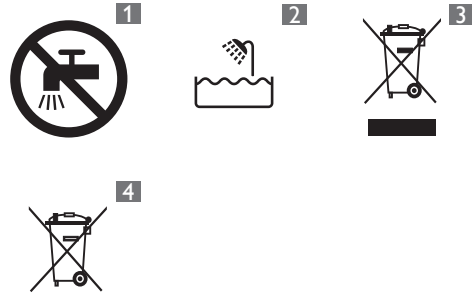Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Important

Read this important information carefully before you use the delicate area epilator and save it for future reference.

Danger

- Keep the adapter dry (Fig. 1).
- If you charge the epilator in the bathroom, do not use an extension cord.

Warning

- Do not use the appliance if it is damaged.
- The adapter contains a transformer. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- This appliance can be used by children aged from 8 years and above persons with reduced physical, sensory or mental capabilities or lack experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Only use the epilator in combination with the adapter supplied.

Caution

- The epilator is intended for epilating delicate areas such as the bikini area, bikini line, underarms and lower face.
- For hygienic reasons, the epilator should only be used by one person.
- To prevent damage and injuries, keep the epilator away from scalp hair, eyebrows, eyelashes, clothes, threads, cords, brushes etc.
- Always check the appliance before you use it. Do not use an appliance, adapter or any other part if it is damaged, as this may cause injury. Always replace a damaged adapter or part with one of the original type.
- Do not use the epilator on irritated skin or skin with varicose veins, rashes, spots, moles (with hairs) or wounds without consulting your doctor first. People with a reduced immune response or people who suffer from diabetes mellitus, haemophilia or immunodeficiency should also consult their doctor first.

- Your skin may become a little red and irritated the first few times you use the epilator. This phenomenon is absolutely normal and disappears. As you epilate more often, your skin gets used to epilation, skin irritation decreases and regrowth becomes thinner and softer. If the irritation has not disappeared within three days, we advise you to consult a doctor.
- To prevent from burning, do not clean the epilator with water that is hotter than shower temperature (max. 40°C).
- Never use compressed air, scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliances.

Electromagnetic fields (EMF)

This Philips epilator complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

General

- This symbol means that the delicate area epilator is suitable for use in a bath or shower (Fig. 2).
- Maximum noise level of the epilator:
Lc = 77dB(A)

Preparing for use

Charging

Charging the delicate area epilator takes approx. 5 hours. When the epilator is fully charged, it has a cordless operating time of up to 30 minutes. If the epilator is not fully charged, the operating time is shorter.

Note: No disposable batteries needed.

- 1 Switch off the appliance.
- 2 Place the appliance in the charging case.

Note: Make sure the charging pins of the epilator are in full contact with the charging case. You hear a click when full contact has been established.

- 3 Insert the small plug in the bottom of the charging case (1) and insert the adapter into the wall socket (2).
- ▶ The charging light lights up green continuously to indicate that the appliance is charging.

Note: When the battery is fully charged, the light goes out.

Using the delicate area epilator

The epilator is primarily intended to remove hairs by the root in the delicate areas such as bikini area, bikini line, underarms and lower face. You can also use it to epilate hard-to-reach areas such as knees and ankles. The appliance can be used wet or dry.

Epilating tips

- **Wet:** We advise you to start epilating on wet skin, as water relaxes the skin and makes epilation less painful and gentler. Use your common shower gel or foam for a smoother epilation.
- **Dry:** If you epilate dry, epilation is easier just after taking a bath or shower. Make sure your skin is completely dry when you start epilating.
- **General:** Your skin should be entirely clean and free of oil. We advise you to epilate in the evening, as skin irritation usually disappears overnight.
- **General:** Epilation is easier and more comfortable when the hairs have a length of 3-4 mm. If the hairs are longer, we recommend that you shave first and epilate the shorter regrowing hairs after 1 or 2 weeks. Alternatively, you can trim the hairs to 3-4 mm before you start epilating.

Epilating your delicate areas

To get used to epilation in delicate areas, for instance your face, first epilate a small area. Make sure your skin is stretched. It is recommended to use the delicate area epilator in the evening so that any irritation can disappear overnight. To protect your skin, the protective bar on the back of the epilating head helps you to keep the epilator at the right angle.

Caution: Do not use this appliance to epilate eyebrows or eyelashes.

Caution: Do not use this appliance to remove hairs that grow out of moles.

- 1 Thoroughly clean the area to be epilated. Remove any residues of deodorant, make-up, cream etc.

Note: Check if the appliance is cleaned thoroughly before every use.

- 2 Press the on/off button to switch on the appliance.
- 3 To epilate your upper lip, hold the appliance with the on/off button pointing upwards. Place the appliance on the skin above your upper lip at an angle of 90°. Stretch the skin

with your other hand and move the appliance upwards.

- 4 To epilate the bikini line, stretch the skin with one hand to make the hairs stand upright and reduce skin pulling. Place the appliance onto the skin at an angle of 90° and slowly move it against the direction of hair growth.
- 5 To epilate the underarms, lift your arm to stretch the skin. Place the appliance onto the skin at an angle of 90° and move the appliance in different directions to catch all hairs.
- 6 When you have finished epilating, switch off the appliance and clean the epilating head according to the instructions in chapter 'Cleaning and maintenance'.

Tip: To relax the skin, we advise you to apply a moisturising cream or a calming gel, immediately or a few hours after epilation.

Cleaning and maintenance

Cleaning the delicate area epilator

- 1 Remove loose hairs with the cleaning brush.
- 2 Switch on the appliance.
- 3 Rinse the appliance while it is switched on, with cold water for 15 seconds after every use.
- 4 Shake off excess water and let the appliance dry.

Storage

- For protection during storage or travel, store the appliance in the charging case.

Note: Make sure that all parts are dry when you store them.

Recycling

- This symbol on a product means that the product is covered by European Directive 2012/19/EU (Fig. 3).
- This symbol means that the product contains a built-in rechargeable battery covered by European Directive 2006/66/EC which cannot be disposed of with normal household waste. We strongly advise you to take your product to an official collection point or a Philips service centre to have a professional remove the rechargeable battery (Fig. 4).
- Inform yourself about the local separate collection system for electrical and electronic products and rechargeable batteries. Follow local rules and never dispose of the product

and rechargeable batteries with normal household waste. Correct disposal of old products and rechargeable batteries helps prevent negative consequences for the environment and human health.

Removing the rechargeable batteries

Only remove the batteries when they are completely empty to avoid a hazard. Never let the batteries come in contact with metal.

- 1 Make sure the batteries are completely empty.
- 2 Put the device in the freezer for at least 24 hours to make the material brittle.
- 3 Wrap the device in a towel.
- 4 Stamp on it with your foot or use a hammer to break the device.

Caution: Only do this while wearing solid shoes. Attempting this with bare feet could cause injuries.

- 5 Take out the batteries and remove the printed circuit boards from the rechargeable batteries with a screwdriver.

Caution: Be careful when using a hammer to avoid a hazard.

Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/ support for a list of frequently asked questions or contact the Consumer Care Centre in your country.

Problem	Possible cause	Solution
The delicate area epilator does not work.	The appliance is not charged.	Charge the appliance. See chapter 'Preparing for use' section 'Charging'.

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