

# PHILIPS

## sonicare

Sonic electric  
toothbrush with app

DiamondClean  
Smart

Personalized coaching

Smart brush head sensors

Smart brush head recognition

5 modes, 3 intensities



HX9985/18



## Complete system for a healthier mouth

Our best toothbrush with smart sensor technology

Get the ultimate clean with DiamondClean Smart. Eliminate guesswork with real time feedback and personalized coaching for superior teeth, gum, whitening & fresh breath. It's more than just a toothbrush - it's a complete oral care solution.

### Proven to improve oral health

- Enjoy long-lasting fresh breath
- Up to 7x healthier gums\*\*\*\* in just 2 weeks
- All-in-One brush head for complete care

### Personalized in-app coaching to perfect brushing

- Takes you back to any missed spots for 100% coverage
- Alerts you if you're pressing too hard
- Start taking better care of your trouble areas
- Guides you to brush and highlights missed areas
- Guides you to reduce your scrubbing

### Optimize your brushing

- Brush heads automatically select the optimal settings
- Choose from 5 modes and 3 intensity settings
- Make sure you're getting the most from your brush head
- Premium travel case and charging glass


# Highlights

## Premium All-in-One brush head



Clean without compromise with the A3 Premium All-in-One brush head. Our all-time best has angled bristles to help remove up to 20x more plaque even in hard-to-reach areas\*. Triangular tips work to remove up to 100% more stains in less than two days\*\*\*. And longer bristles clean deep for gums that are up to 15x healthier in just two weeks\*\*. All that without ever switching brush heads. Did you know: brush heads become less effective after 3 months of use? Our BrushSync™ feature can remind you when it's almost time for a replacement.

## Care for your gums

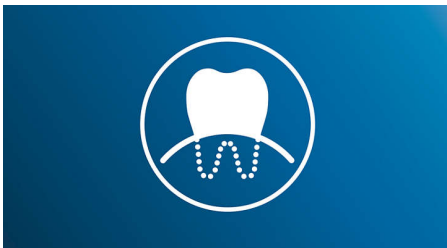

Click on the G3 Premium Gum Care brush head to improve your gum health. Its smaller size and targeted gum line bristles deliver gentle yet effective cleaning along the gum line, where gum disease starts. It's clinically proven to deliver up to 100% less gum inflammation\*\*\*\* and up to 7x healthier gums in just two weeks.\*\*\*\*


## TongueCare+ tongue brush



Click on the TongueCare+ tongue brush to gently remove odor causing bacteria from the pores of your tongue. Its 240 specially


designed MicroBristles get between all your tongue's ridges and grooves to remove bacteria and debris that lead to bad breath. Team with our antibacterial BreathRx tongue spray for superior cleaning and super fresh breath.

## Focus areas and goal-setting

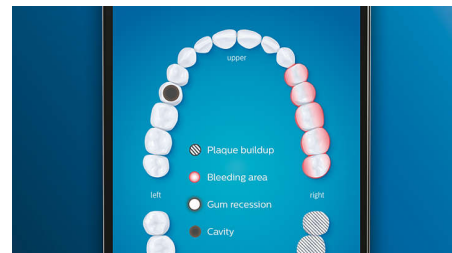


Any trouble areas your dentist has pointed out? We'll highlight them on your in-app 3D mouth map, and remind you to pay extra attention to these areas.

## Never miss a spot

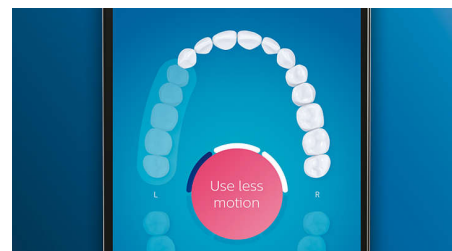

The Philips Sonicare app lets you know when you've achieved a thorough clean, and coaches you to be a better, more mindful brusher.

## Too much pressure?

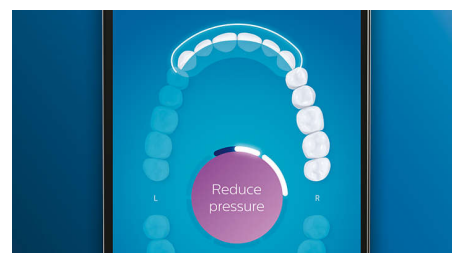


You might not notice when you brush too hard, but your DiamondClean Smart will. If you use too much pressure, the light ring on the end of your handle will flash. This is a gentle reminder to ease off the pressure, and let your brush head do the work. 7 out of 10 people found this feature helped them become a better brusher.

## Easy does it




With an electric toothbrush, you let the brush do the work. We've built a scrubbing sensor into your handle as a gentle reminder to reduce scrubbing. This way, you can improve your technique and get a gentler, more effective clean.

## TouchUp feature

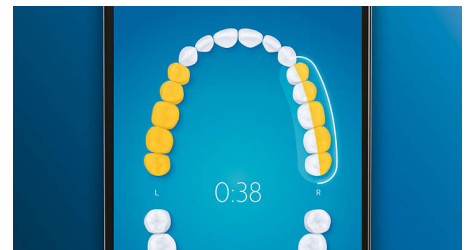

If you happen to miss a spot while brushing your teeth, your app's TouchUp feature will show you. You can then go back for a second pass, and be confident you are getting a complete clean, every time.

## BrushSync mode pairing

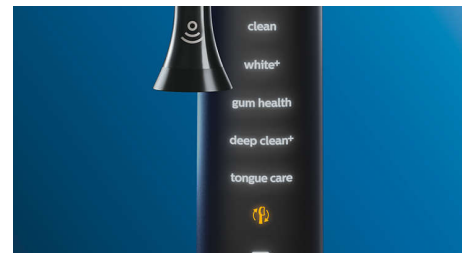

Wondering which mode and intensity you should use? Wonder no longer. A microchip in your brush head tells your DiamondClean Smart which one you're using. So if you click on a gum-care brush head, for example, your toothbrush knows to select the optimal mode and intensity to provide a gentle yet effective care for your gums. All you need to do is to press the power button.

# Specifications

## Connectivity

**Bluetooth® wireless technology:** Connected brushing app

## Power

**Voltage:** 110–220 V

## Technical specifications

**Operating time (full to empty):** 14 days\*\*\*\*\*

**Battery:** Rechargeable

**Battery type:** Lithium ION

## Design and finishing

**Color:** Black

## Service

**Warranty:** 2-year limited warranty

## Compatibility

**Android compatibility:** Android phones, Bluetooth 4.0 enabled tablets

**iOS compatibility:** iPad 3rd Gen or higher, iPhone 4S or higher, with iOS7 or higher

## Ease of use

**Handle:** Slim ergonomic design

**Battery indicator:** Illuminated icon indicates battery life

**Handle compatibility:** Easy click-on brush heads

**Display:** Illuminated display

**Brush head recognition:** Syncs with optimal mode

**Replacement reminder:** To always ensure best results, reminder icon lights up

## Items included

**Brush heads:** 2 A3 Premium All-in-One, 1 G3 Premium Gum Care, 1 TongueCare+ tongue brush

**Travel case:** Premium travel case

**Handle:** 1 DiamondClean Smart

**Glass charger:** 1

## Cleaning performance

**Performance:** Removes up to 20x more plaque\*

**Health benefits:** Up to 15x healthier gums\*\*

**Speed:** Up to 62000 brush movement/min

**Whitening benefits:** Up to 100% less stains\*\*\*

**Timer:** BrushPacer and SmarTimer

**Pressure feedback:** Ring lights up in purple, Vibrates handle to alert user

**TouchUp:** Ensures 100% coverage

## Modes

**3 intensities:** Low, Medium, High

**Clean:** For exceptional everyday clean

**Deep Clean+:** For an invigorating deep clean

**Gum Health:** Special attention to molars

**TongueCare+:** For long lasting clean breath

**White+:** To remove surface stains

## BrushSync Mode Pairing

**C3 Premium Plaque Control:** Pairs with Clean mode

**G3 Premium Gum Care brush head:** Pairs with Gum Health mode

**W3 Premium White brush head:** Pairs with White+ mode

**TongueCare+ tongue brush:** Pairs with TongueCare mode

## Smart sensor technology

**Pressure sensor:** Alerts when brushing too hard

**Location sensor:** Tracks & improves coverage

**Scrubbing sensor:** Guides to reduce scrubbing

**3D mouth map:** Highlights trouble areas

**TouchUp:** To address missed spots



\* vs. a manual toothbrush

\*\* vs. a manual toothbrush in 6 weeks

\*\*\* in lab test vs. a manual toothbrush

\*\*\*\* in Gum Health mode vs a manual toothbrush at 2 weeks

\*\*\*\*\* based on two periods of two-minute brushings per day on standard mode

\*\*\*\*\* based on two periods of two-minute brushings per day on standard mode

\*\*\*\*\* Recommended with interdental cleaning