

PHILIPS

sonicare

*DiamondClean
Essential*
7000 Series/系列

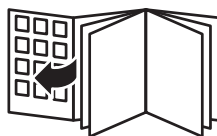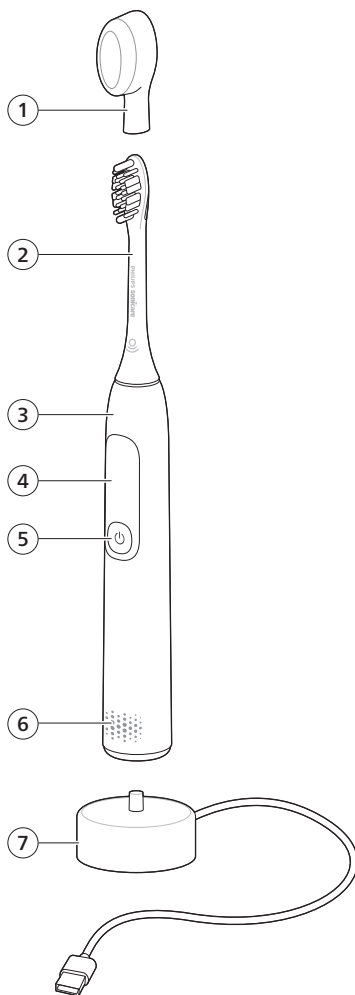English 6

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Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **www.philips.com/welcome**.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.

Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support').
- If the USB cable or USB power supply unit is damaged, stop using it.
- Keep the charger away from water.
- After cleaning, make sure the charger is completely dry before you connect it.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not use adapters and chargers outdoors or near heated surfaces.

- This appliance has only been designed for cleaning teeth, gums and tongue.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- The maximum altitude of use is 4500 meters.

Medical warnings

- This appliance is a personal care device and is not intended to be shared by multiple users.
- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.

- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5 V, ≥ 2 W. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country (see the international warranty leaflet for contact details). Using a non-certified supply unit may cause hazards or serious injuries.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).

Storage and transportation conditions

Temperature -20°C to +60°C (-4 °F to 140 °F).

Electromagnetic fields (EMF)


This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

LIST OF SYMBOLS USED

The following symbols may appear on the product:
This symbol means read the instructions for use, before using the appliance.

This symbol means to use the indicated power supply.

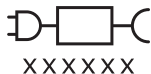Intended use

Sonicare power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Your Philips Sonicare toothbrush

- 1 Brush head cap(s)
- 2 Smart brush head(s)
- 3 Handle
- 4 Smart display
- 5 Power on/off and mode button
- 6 Pressure sensor indicator
- 7 USB charger

Note: Wall adapter not included. You can use any compatible IPX4 splash-proof adapter (identifiable by the "IPX4" text or symbols on the adapter).

Getting started

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.


Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using

a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.


To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gradually increases the power over the first 14 brushing sessions to help you adapt to brushing with your Sonicare toothbrush (see 'EasyStart').

Brushing instructions

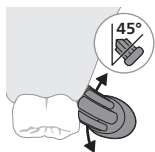
- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Press the power on/off button to turn on the Philips Sonicare. Then select the brushing mode.
- 3 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

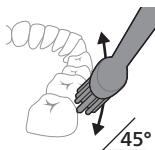
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

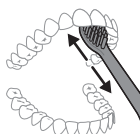
- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the QuadPacer feature (see chapter 'Features').

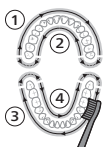- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

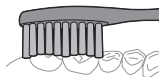
- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare toothbrush should be safe to use on braced (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Note: When the Philips Sonicare toothbrush is used in clinical studies, it must be used in Clean mode with the handle fully charged and the EasyStart feature deactivated.

Charging and battery status

Charging your Philips Sonicare

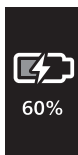- 1 Attach the USB charger to a wall adapter and plug the wall adapter into an electrical outlet.
- 2 Place the handle on the charger.
 - The battery icon with lightning symbol is shown on the display and you hear two beeps. The light ring lights up white for 5 seconds.

Note: It can take up to 16 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.


Note: This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day).

Battery status


The battery level is shown on the display:

- When the battery level is below 20%, the battery icon flashes in red color and you hear two short beeps.

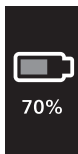- At the end of every brushing session.
- When you remove the toothbrush from the charger.

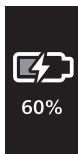- When you press the power on/off button while the toothbrush is charging.

Features

- Brushing modes
- EasyStart
- Brush head replacement reminder
- Pressure sensor
- QuadPacer
- SmarTimer
- Smart display
- Brushing report

Brushing modes

Based on the smart brush head you use; the optimal mode will be selected automatically by the BrushSync technology. To change modes, press the mode button while the handle is turned on. Your power toothbrush is equipped with five different modes:

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Clean

For exceptional everyday cleaning (2 minutes).

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Deep Clean

For invigorating and thorough cleaning (3 minutes).

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Gum Pro

For exceptional cleaning and gentle gum massaging (3 minutes).

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Sensitive

For gentle cleaning of gums (2 minutes).

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White

For removing stains from teeth (2 minutes 30 seconds). After brushing your whole mouth (2 minutes), spend the remaining time brushing the front teeth.

Note: When you use the toothbrush with the smart brush head for the first time, the default mode is 'Gum Pro'.

To toggle between the brushing modes:

- Press the power on/off button once to turn on the toothbrush.
- Press a second time **within 5 seconds** to change the brushing mode.

- After the toothbrush has been on 5 seconds, if you want to change modes, you'll need to turn the toothbrush off, then on and press the power on/off button **within 5 seconds**.

EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. To deactivate EasyStart see 'Activating or deactivating features'.

Brush head replacement reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When a new smart brush head is attached for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and the display shows the brush head symbol.

- 2 Over time, based on the pressure you apply and how long the brush head has been used, the handle will track brush head wear in order to determine the optimal time to change your brush head. When the life of the brush head ends, the display shows the brush head replacement reminder symbol and indicates that you should replace the brush head as soon as possible.

Pressure Sensor

Your Philips Sonicare is equipped with an advanced sensor that measures the pressure you apply while brushing. If you apply too much pressure, the toothbrush provides immediate feedback to indicate that you need to adjust the pressure. When you apply too much pressure, the light ring at the bottom of the handle lights up orange.

QuadPacer

The QuadPacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. This Philips Sonicare toothbrush comes with the QuadPacer feature activated.

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush will pause. After a pause of 15 seconds, the SmarTimer resets.

Brushing reports

Your Philips Sonicare regularly provides a report on your brushing behavior.

After brushing the display shows a report with the brushing time.

- Brushing time of 2 minutes or more

- Brushing time of less than 2 minutes

- During brushing the display shows the current brushing mode and brushing time

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder

To activate or deactivate these features, follow the instructions below:

EasyStart

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- 1 Put the handle on the plugged-in charger.
 - 2 Keep the power on/off button pressed until you hear a single short beep (after 3 seconds).
 - 3 Release the power on/off button.
- Triple tone of low-medium-high means the EasyStart feature has been activated. The icon on the display shows that the EasyStart has been activated.

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- Triple tone of high-medium-low means the EasyStart feature has been deactivated. The icon on the display shows that the EasyStart has been deactivated.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

Brush head replacement reminder

Note: You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

- 1 Put the handle on the plugged-in charger.
- 2 Press and hold the power on/off button while the handle remains on the charger.
- 3 Keep the power on/off button pressed until you hear two short beeps.
- 4 Release the power on/off button.
 - Triple tone of low-medium-high means the brush head replacement reminder feature has been activated. The icon on the display shows that the brush head replacement reminder has been activated.
 - Triple tone of high-medium-low means the brush head replacement reminder feature has been deactivated. The icon on the display shows that the brush head replacement reminder has been deactivated.

Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean the product or accessories with sharp objects or in the dishwasher, microwave, with chemicals or in boiling hot water

Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

Brush head

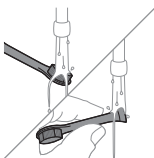- 1 Rinse the brush head and bristles after each use and allow to air dry.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Storage

If you are not going to use the product for a long time, unplug it, and clean it (see chapter 'Cleaning'). Store it in a cool and dry place away from direct sunlight. Charge the toothbrush regularly to prevent the battery from discharging completely. A completely discharged battery may fail to charge again.

Brush head replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm). When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX379X).

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

**CONTAINS LITHIUM-ION RECHARGEABLE CELL.
MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.


To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come

into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.


Before removing the battery, make sure that the battery is completely empty.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Remove the film at the bottom of handle, then release the screw.

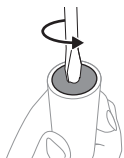
- 3 Remove the end cap from the toothbrush handle.


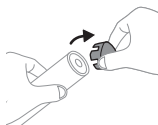
- 4 Holding the handle upside down, press the shaft down on a hard surface.

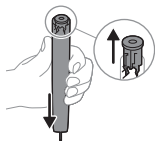- 5 Wedge the screwdriver between the battery and the white frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the white frame.

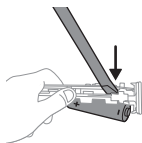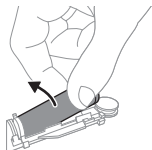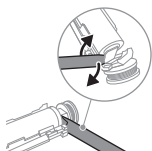- 6 Insert the screwdriver between the bottom of the battery and the white frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

- 7 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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商標	型號	額定電壓/電流	消耗電功率	輸出負載		
				商標	牙刷握柄型號	牙刷握柄額定輸入
PHILIPS	HX6110 AA3 HX6110 ABA3 HX6110 ACA3	5 V $\overline{\text{---}}$ 0.3 A	1.5 W	PHILIPS	HX379B 此牙刷握柄 為專屬專用	3.6 V $\overline{\text{---}}$ 2.8 W
設備名稱：飛利浦 USB 無線充電器 型號（型式）：HX6110 AA3, HX6110 ABA3, HX6110 ACA3 Equipment name Type designation (Type)						
單元Unit	限用物質及其化學符號 Restricted substances and its chemical symbols					
	鉛 Lead (Pb)	汞 Mercury (Hg)	鎘 Cadmium (Cd)	六價鉻 Hexavalent (Cr ^{VI})	多溴聯苯 Polybrominated biphenyls (PBB)	多溴二苯醚 Polybrominated diphenyl ethers (PBDE)
電路板 PCB's	-	O	O	O	O	O
塑膠部件 Plastics parts	O	O	O	O	O	O
USB 電線 USB Cable Unit	O	O	O	O	O	O
金屬部件 Internal Metal parts	O	O	O	O	O	O
備考1. "超出0.1 wt %"及"超出0.01 wt %"係指限用物質之百分比含量超出百分比含量基準值。 Note: 1. "Exceeding 0.1 wt %" and "exceeding 0.01 wt %" indicate that the percentage content of the restricted substance exceeds the reference percentage value of presence condition. 備考2. "O"係指該項限用物質之百分比含量未超出百分比含量基準值。 Note 2: "O" indicates that the percentage content of the restricted substance does not exceed the percentage of reference value of presence. 備考3. "-"係指該項限用物質為排除項目。 Note 3: The "-" indicates that the restricted substance corresponds to the exemption.						

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IPX7

www.philips.com/Sonicare

©2025 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.

©2025 Koninklijke Philips N.V. (KPNV)。保留所有权利。飞利浦及飞利浦盾牌标志是 KPNV 的商标。Sonicare 和 Sonicare 徽标是 Philips Oral Healthcare, LLC 和/或 KPNV 的商标。

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Lifestyle B.V.