PHILIPS

sonicare


9000

DiamondClean

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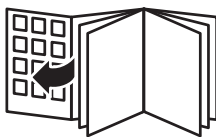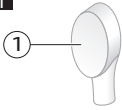
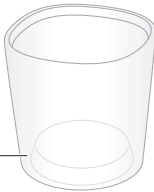
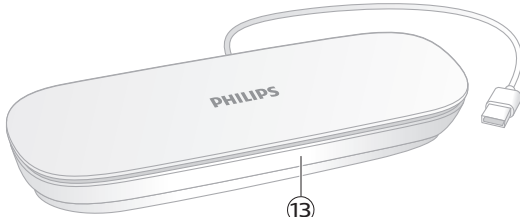
English 6
Français 27
Türkçe 50
العربية 71

Introduction

Congratulations on your new Philips Sonicare power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: www.philips.com/welcome.

Welcome to the Philips Sonicare family!

Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Any serious incident that has occurred in relation to this device should be reported to Philips (www.philips.com/support) and to the competent authority of the member state in which the user and / or patient is established.

Note: The accessories and features of different models vary. For example, some models come with a sanitizer or are equipped with Bluetooth or BrushSync.

Warnings

- Keep the USB wall adapter and charging base away from water.
- Do not immerse the USB wall adapter or charging base in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to the wall socket.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use brush heads recommended by Philips Sonicare. Only use the charger provided with the product. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support'). If the cord is damaged, replace the charger and/or sanitizer with a new one.
- Do not use the charger outdoors or near heated surfaces.

- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Use by children should be supervised by an adult.
- To charge the battery, only use the USB wall adapter, charging base and charging travel case provided with the product: TCAXXXX, CBAXXXX. Do not use an extension cord.
- Charge and use the product at a temperature between 5 °C and 40 °C (41 °F and 104 °F).
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes equipped with a radio frequency interface (Bluetooth® or NFC) are in compliance with Directive 2014/53/EU.

- The Bluetooth radio frequency interface in applicable products operates at 2.4 GHz.
- The maximum radio frequency power of products equipped with Bluetooth is 3 dBm.
- The NFC radio frequency interface in applicable products operates at 13.56 MHz.
- The maximum radio frequency power of products equipped with NFC is 30.16 dBm.

The full text of the EU declaration of conformity is available at the following internet address: **www.philips.com/support**

Note: The features of different models may vary. Some models may not be equipped with Bluetooth or NFC.

Intended Use

DiamondClean series power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. DiamondClean series power toothbrushes are intended for consumer home use.

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush heads*
- 3 Handle (HX991W, HX991B)
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Modes
- 8 Brush head replacement reminder light
- 9 Battery indication
- 10 Charging glass
- 11 Charging stand
- 12 Charging base (100-240 V~, 50/60 Hz, 1.3W or 0.05 A)
- 13 (Charging) travel case (4.75-5.25 Vdc, 1.5 W)*
- 14 USB wall adapter*(100-240V~ 50/60Hz, 3.5W)

* Note: The content of the box may vary based on the model purchased.

The following symbols may appear on the product or package.

Read Operator's Manual.

Detachable Power Supply Part: The "xxxxxx" indicates the model number of the approved chargers to be used.

Consult instructions for use

Sonicare app – Getting Started

The Sonicare app pairs with your toothbrush to provide you a connected experience. By connecting your toothbrush to your Sonicare app account, you will be able to:

- Easily review your progress report for ways to improve your brushing routine
- Receive personalized tips and recommendations for managing your oral health

To start using the Sonicare app:

- 1 Download the Sonicare app to your phone
- 2 Open the app and follow the guided steps
- 3 Pair your toothbrush with the app
- 4 Create your account
- 5 Complete firmware update, if prompted, to access the latest improvements and features.
- 6 Brush regularly

You are ready to start your connected experience. Be sure to brush regularly. You can receive updates to help improve your oral healthcare.

Note: Make sure your phone's Bluetooth is turned on when using the app. Your toothbrush uses the Bluetooth connection to transfer your brushing data to the app. If you have questions about why your brushing data is collected, be sure to review the Sonicare Privacy Statement, available throughout the app setup process.

Sonicare app – Track and improve your oral health care

Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session
- Applied brushing pressure
- Brush head replacement reminder (based on actual use)

The Sonicare app collects the data provided by the smart sensors and tracks your brushing habits over time.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

Each of the features available in the app are designed to help you understand your brushing habits and provide you with simple recommendations for improving and maintaining your oral health. You will receive personalized brushing information in the Sonicare app.

Brush heads

To see what type of brush head you have, check the label at the bottom of the brush head. The symbol at the bottom of the brush head indicates that it is a smart brush head, equipped with the BrushSync technology.

Note: The content of the box may vary based on the model purchased.


Brush head name	Premium Plaque Defense	Premium White	Premium Gum Care
Brush head label	C3	W3	G3
Benefit	Plaque removal	Plaque and stain removal	Healthier gums

To explore our full range of brush heads, go to www.philips.com/toothbrush-heads for more information.

Using your Philips Sonicare

Brushing instructions

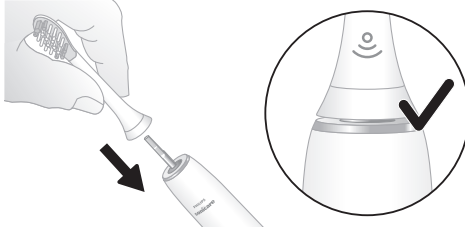- 1 Push the brush head firmly onto the handle. Wet the bristles and apply toothpaste.


Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

12 English


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), apply light pressure to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush head in contact with the teeth at all times.

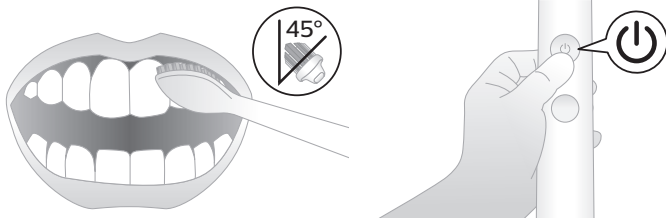- 3 Press the power on/off button to turn on the Philips Sonicare.

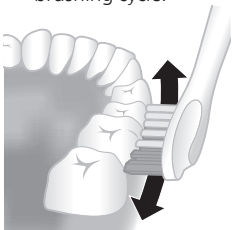
- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.

Note: Different modes are available (see 'Brushing modes').

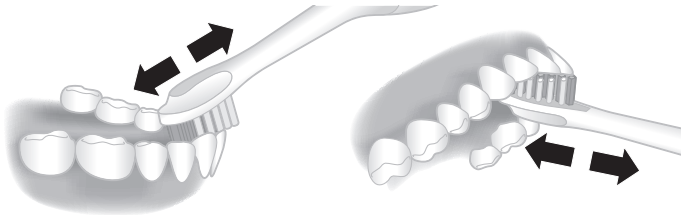- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

**Note:**

- A change in vibration from the handle alerts you when you apply too much pressure (see 'Pressure sensor').
- Your Philips Sonicare toothbrush is safe to use on braces, dental modifications and dental restorations (fillings, crowns, veneers)
- Do not use your power toothbrush in the shower.

BrushSync technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. The symbol at the bottom of the brush head indicates that it is equipped with the BrushSync technology. This enables the following features:

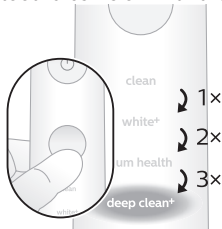
- Brush head mode pairing
- Brush head replacement reminder

Brushing modes

Your power toothbrush is equipped with 4 different modes: clean, white+, gum health and deep clean+.

Based on the smart brush head you use, the optimal mode and intensity will be selected automatically by the BrushSync technology.

To manually change the mode, press the mode/intensity button while the toothbrush is off with the brush head attached.

* Note:

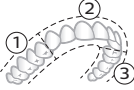
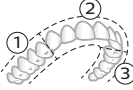
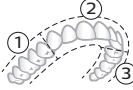
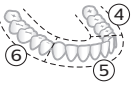
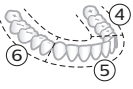
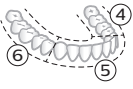
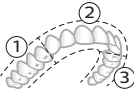
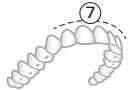
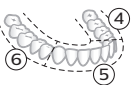
- All brush heads work with each mode.
- When the Philips Sonicare toothbrush is used in clinical studies, it must be used at high intensity, with the EasyStart feature turned off.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates to the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the brush head will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended mode/intensity in the following table.

	Clean and Deep Clean+	White+	Gum health
Benefit	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage
Recommended brush head	Premium Plaque Defense (C3)	Premium White (W3)	Premium Gum Care (G3)
Recommended intensity level	3	3	3
Total brushing time	Clean: 2 minutes Deep Clean+: 3 minutes*	2 minutes and 40 seconds	3 minutes and 20 seconds

	Clean and Deep Clean+	White+	Gum health
How to brush	Clean: Brush each segment for 20 seconds	Brush each segment for 20 seconds,	Brush each segment for 20 seconds,
			
			
	Deep Clean+: Brush each segment for 30 seconds*	then brush segments 7 and 8 for 20 seconds each.	then brush segments 7, 8, 9 and 10 for 20 seconds each.
			
			

*Note: When brushing with the Philips Home+ app in 'Deep Clean+' mode, the brushing time will be 3 minutes.

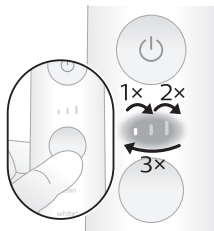
Intensity settings

Your power toothbrush comes with 3 different intensity settings:

- low intensity (one bar)
- medium intensity (two bars)
- high intensity (three bars)

When attached to the handle, the smart brush head will automatically select the recommended intensity. To manually change the intensity, turn on the handle and press the mode/intensity button.

While toothbrush is ON, press this button to change the intensity setting

Features

- BrushPacer
- Pressure sensor
- Brush head replacement reminder
- EasyStart

BrushPacer

The BrushPacer features indicates when you have to move to the next section of your teeth with a brief change in vibration.

The toothbrush turns off automatically at the end of the brushing session.

Pressure sensor

Your Philips Sonicare measures the pressure you apply while brushing to protect your gums and teeth from damage.

If you apply too much pressure, the handle will change its vibration until you reduce the pressure.

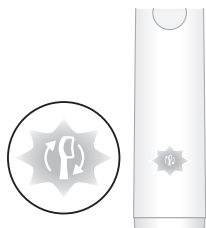
Note: The Pressure Sensor comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

Brush head replacement reminder

Over time, the handle will track the wear of your smart brush head by measuring:

- The overall pressure you apply while brushing
- The total time you have brushed with your brush head

Once your brush head is no longer effective, the brush head replacement reminder will blink amber and the handle will make a series of beeps and tones.

Note: The brush head replacement reminder comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

EasyStart

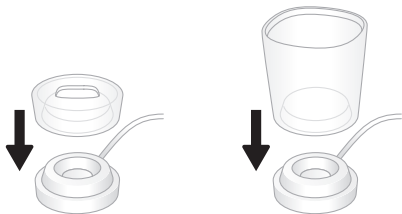
The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power of the Philips Sonicare. The EasyStart feature works with all modes.

Note: The EasyStart feature comes deactivated with your product. To activate EasyStart, see 'Activating or deactivating features'.


Charging and battery status

Charging on charger base

- 1 Plug the cord of the charging base into an electrical outlet.
- 2 Place the charging stand or the charging glass on the charging base.

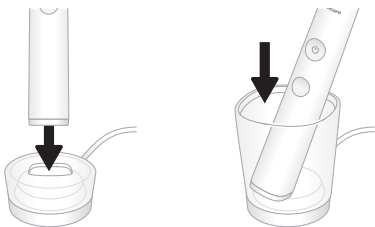- 3 Place the handle on the charging stand or in the charging glass.

- 4** The charger will charge your toothbrush until it is fully charged. The handle is charging when the battery indication is blinking white. When the handle is fully charged, the battery indication illuminates white for 30 seconds and then turns off.

Note: Your handle comes pre-charged for first use. After first use, charge for at least 24 hours.

Charging with the travel case*

*Note: The charging travel case may not be included depending on the model purchased.

- 1 Take the USB cable from the bottom of the charging travel case.
- 2 Connect the charging travel case to 5V DC input power either from:
 - The USB adapter provided
 - a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
 - a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
 - a UL-listed or certified 1310 Class 2 power supply.

Important: The wall adapter must be a safe extra low voltage source, compatible with the electrical ratings marked on the bottom of the device.

- 3** The travel case charges your toothbrush until your toothbrush is fully charged.

Battery Status (handle on plugged in charger)

When the handle is placed on the charger or in the travel case, the battery indication will communicate the battery level.

Battery status	Battery indication
Charging	Blinks in white
Full	Illuminates white for 30 seconds and then turns off.

Battery Status (when handle is not placed on charger)

When you turn on the handle, the battery light at the bottom of the handle will indicate the status of the battery.

Battery status	Battery indication	Sound
Full	Illuminated in white	-
Partially full	Blinks white	-
Low	Blinks amber	The handle emits 3 beeps
Empty	Blinks amber	The handle emits 2 sets of 5 beeps

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

Step 1: Place the handle on the charging stand or in the charging glass.

Step 2: Press and hold power button for:

Easy start

Brush head
replacement
reminder

Pressure sensor

3 seconds

Up to 5 seconds

Up to 7 seconds

Step 3: Let go power button when you hear

1 beep

1 beep and then
2 beeps

1 beep, 2 beeps and
then 3 beeps

If you see the battery indication blink white 2 times and hear 3 tones from low to high, then the feature has been activated.

OR

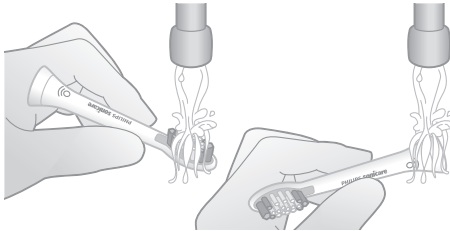
If you see battery indication blink white once and hear 3 tones high to low, then the feature has been deactivated.

Cleaning


- All parts can be cleaned with a damp cloth or rinsed with warm water.


- Remove the brush head from the handle and rinse it thoroughly.

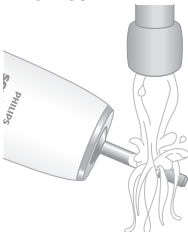- Rinse the entire handle, especially the brush head connection, at least once a week.

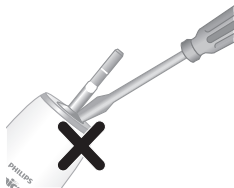Caution: Do not push on the rubber seal at the top of the handle. Gently clean around the rubber seal.

**Cautions:**

- Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.
- Do not use isopropyl alcohol, vinegar, bleach, or any other house hold cleaning products, to clean the handle as this may cause discoloration.
- Unplug the charger before you clean it.
- Do not use essential oils to clean the brush heads as this can cause damage.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Frequently asked questions

Question	Answer
Can I leave my Sonicare on the charger between brushings?	Yes, you can always keep your Philips Sonicare toothbrush on a plugged-in charger. This does not affect the battery lifetime.
Will other brush heads fit my Philips Sonicare toothbrush?	All Philips Sonicare electric toothbrush brush heads that snap on will fit your toothbrush. However, not all are equipped with BrushSync technology. You can identify the brush heads with BrushSync technology by the icon at the base of the brush head.

Question**Answer**

How often do I need to replace my brush head?

We advise you to replace your Philips Sonicare brush head after 3 months of use (based on brushing twice a day for 2 minutes each session). The brush head replacement reminder will also notify you when it is time to replace your brush head.

Is my Philips Sonicare toothbrush waterproof?

Yes, your Philips Sonicare toothbrush is waterproof. However, never put your handle in a dishwasher to clean it. And do not use it in the shower.

What phones are compatible with the Philips Sonicare app?

The Philips Sonicare app support Android OS, and Apple iOS. For the most up to date version of the Philips Sonicare app check out the Android Play Store or Apple App Store. Mobile devices must have Bluetooth 4.0 or newer.

What other Philips Sonicare toothbrushes can I use with the Philips Sonicare app?

You can only use Philips Sonicare Flexcare Platinum Connected, DiamondClean Smart, DiamondClean 9000 and ExpertClean toothbrushes with the Sonicare app. Look for the Bluetooth symbol on the front panel of the package.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removing the rechargeable battery

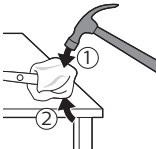
Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. When you handle batteries, make sure that your hands, the product and the batteries are dry.

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.


- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

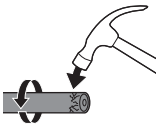
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

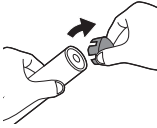
Note: You may have to hit on the end several times to break the internal snap connections.

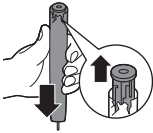
- 4** Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

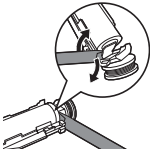- 5** Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.
- 6** Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

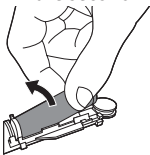
- 7** Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.
- 8** Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

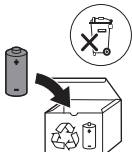- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings).

Note: Do not wrap batteries in aluminium foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

Note: Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

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© 2021 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.
Philips Oral Healthcare, LLC.