PHILIPS

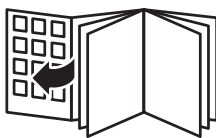
sonicare

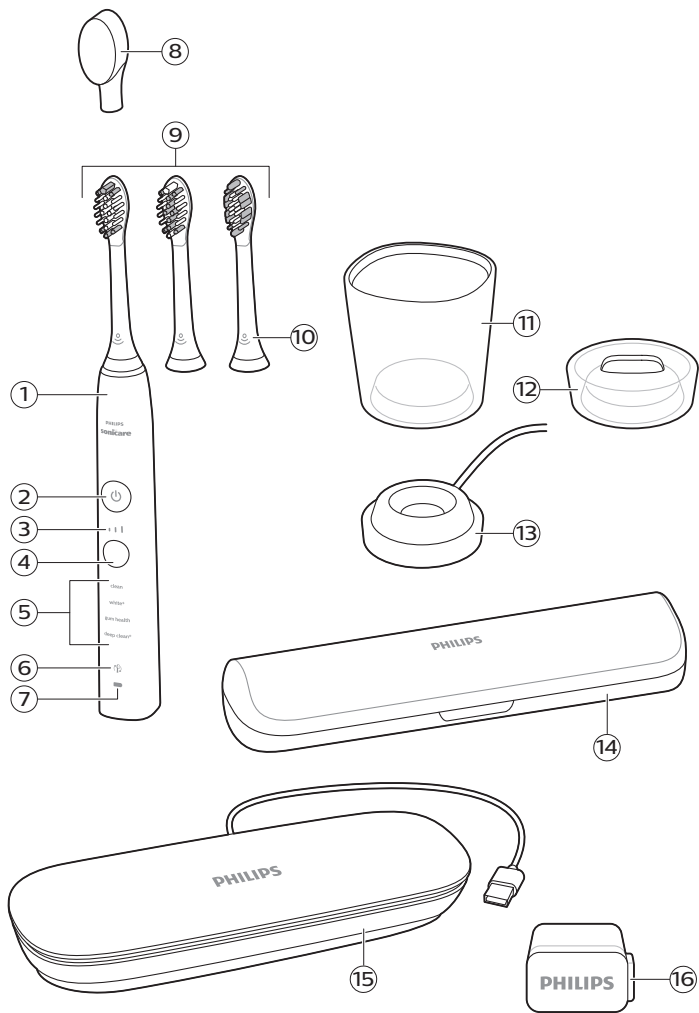
9000

DiamondClean

9000 Series

English 6
Español 35
Français (Canada) 67

Introduction

Welcome to the Philips Sonicare family!

This toothbrush enables you to obtain superior plaque removal, whiter teeth and healthier gums. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time.

Further support and product registration are available to you at: www.philips.com/support

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USE

Read this user manual carefully before you use the appliance and save it for future reference.

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

WARNINGS

- Keep the charger and travel-charging-case away from water.
- After cleaning, make sure the charger and/or sanitizer is completely dry before you connect USB / wallsocket Charger.

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use charger provided with the product. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country. Devices with cord only: If the cord is damaged, replace the charger and/or with sanitizer with a new one.
- Do not use the charger outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

MEDICAL WARNINGS

- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

8 English

- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Battery safety instructions

- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables. Only use detachable supply unit TCAxxxx and/or CBAxxxx. Only use the rechargeable battery provided with the product.
- Charge and use the product at a temperature between 41°F and 104°F (5 °C and 40 °C)
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eye. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.

- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap batteries in aluminium foil. Tape battery terminals or put batteries in a plastic bag before you discard them."

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC Compliance

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).

- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.

ISED Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1 This device may not cause interference.
- 2 This device must accept any interference, including interference that may cause undesired operation of the device.

SAVE THESE INSTRUCTIONS

Intended use

Sonicare DiamondClean 9000 power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare DiamondClean 9000 power toothbrushes are intended for personal use. Use by children should be with adult supervision.

Your Philips Sonicare (Fig 1.)

- 1 Handle
- 2 Power on/off button
- 3 Intensity indicator
- 4 Mode/intensity button
- 5 Mode indicator
- 6 Brush head replacement reminder light
- 7 Battery indicator

Accessories:

- 8 Brush head cap
- 9 Smart brush heads
- 10 BrushSync symbol
- 11 Charging glass
- 12 Charging stand
- 13 Charging base
- 14 Travel case
- 15 Charging travel case
- 16 USB(-A) wall adapter

Note: Accessories included may vary based on the model purchased.

Sonicare app - Getting started

The Sonicare app pairs with your toothbrush to provide you a connected experience. The app is compatible with a wide range of iPhone and Android™ smartphones.

To start using the Sonicare app:

- 1 Download the Sonicare app to your phone.
- 2 Ensure your phone's Bluetooth is turned on.
- 3 Pick up your toothbrush to ensure it is active (lights on).
- 4 Open the app and follow the steps.
- 5 Pair your toothbrush with the app.
- 6 Create your account via the app. Complete firmware update, if prompted, to access the latest improvements and features.

7 Brush regularly with your toothbrush connected to the app. When you regularly sync the toothbrush with the app, you will receive updates to help improve your oral healthcare.

8 Sync regularly.

9 Keep your phone close to your toothbrush for optimal data transfer.

To sync manually: Pair/connect your toothbrush with the app every two weeks to benefit from the app features/capabilities.

To sync automatically: Allow location permission while setting up the app. By allowing location permission, your phone knows when it is in the connection range of your toothbrush and can refresh your brushing data in the app to provide the latest insights and recommendations.

Note: Make sure your phone's Bluetooth is switched on when using the app so that your toothbrush can transfer and update your brushing data in the app. If you have questions about why your brushing data is collected, be sure to review the Sonicare Privacy Statement, available throughout the app setup process.

App features

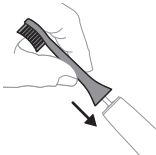
Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session

- Applied brushing pressure
- Brush head replacement reminder (based on actual use)

The Sonicare app collects the data provided by the smart sensors and tracks your brushing habits over time.

Getting started

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

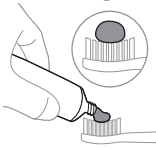
Note: It is normal to see a slight gap between the brush head and the handle.

Using your Philips Sonicare toothbrush


If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, use a lower setting for the first 1-2 weeks then transition to a higher setting.

Brushing instructions

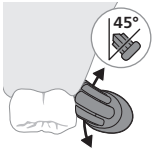
- 1 Wet the bristles and apply a small amount of toothpaste.


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

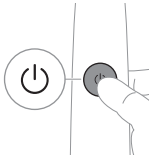Note: Keep the center of the brush in contact with the teeth at all times.

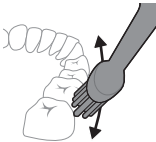
- 3 Press the power on/off button to switch on the Philips Sonicare.

- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

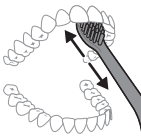
Note: If you apply too much pressure, you will be alerted by a change in vibration of the handle and the brush head replacement reminder flashing amber.

- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

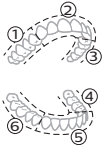BrushPacer

To make sure you brush evenly throughout the mouth, divide your mouth into 6 sections using the BrushPacer feature. The BrushPacer divides the brushing time into 6 equal segments and indicates when you should move to the next area. The segments are indicated with a brief pause in vibration. The toothbrush automatically stops at the end of the brushing session.

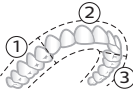
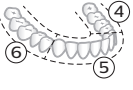
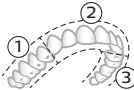
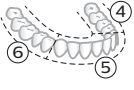
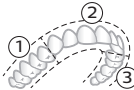
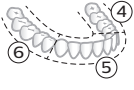
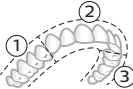
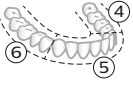
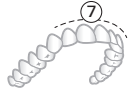
Note: Additional segments may be added when you use the White+ or Gum Health Modes.

Modes

Your power toothbrush is equipped with 4 different modes: Clean, White+, Gum Health and


Deep Clean+. Based on the smart brush head you use, the optimal mode and intensity will be selected automatically by the BrushSync technology. To change modes, press the mode/intensity button while the handle is turned off.

	Clean and Deep Clean+	White+	Gum health
Benefit	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage
Recommended intensity level	High	High	High
Total brushing time	Clean: 2 minutes Deep Clean+: 3 minutes	2 minutes and 40 seconds	3 minutes and 20 seconds

	Clean and Deep Clean+	White+	Gum health
How to brush	Clean: Brush each segment for 20 seconds	Brush each segment for 20 seconds,	Brush each segment for 20 seconds,
	 	 	 
	Deep Clean+: Brush each segment for 30 seconds*	then brush segments 7 and 8 for 20 seconds each	then brush segments 7, 8, 9 and 10 for 20 seconds each
	 	 	 

Note: When brushing with the Philips Sonicare app in 'Deep Clean+' mode, the brushing time will be 2 minutes when connected to the app and 3 minutes when not connected to the app.

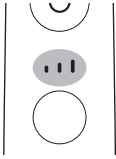Intensities

Your power toothbrush comes with 3 different intensity settings:

- low intensity (one bar)
- medium intensity (two bars)
- high intensity (three bars)

When attached to the handle, the smart brush head will automatically select the recommended intensity. If you wish to change to a different intensity, you can do so by pressing the mode/intensity button while brushing.

The intensity cannot be changed as long as the handle is powered off or paused.

BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The ☺ symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- BrushSync mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to

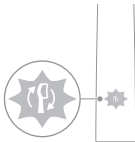
www.philips.com/toothbrush-heads for more information.

Features

- Pressure Sensor Feedback
- BrushPacer
- SmarTimer
- EasyStart
- Brush Head Replacement Reminder
- Brush Head Mode Pairing

Pressure Sensor Feedback

Your Sonicare toothbrush measures the pressure you apply while brushing, to protect your gums and teeth from damage.

If you apply excess pressure, the handle changes its vibration and the brush head replacement reminder icon flashes amber until you reduce the pressure.

Note: Your product comes with the Pressure Sensor activated. To deactivate this feature (see 'Activating or deactivating features').

BrushPacer

The BrushPacer is an interval timer that helps you to brush evenly through the mouth by dividing your mouth into 6 sections.

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. After a pause of 30 seconds, the SmarTimer resets.

EasyStart

Note: To activate EasyStart. (see 'Activating or deactivating features')

This Philips Sonicare model comes with the EasyStart feature deactivated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare.

The EasyStart feature works with all modes.

Disclaimer

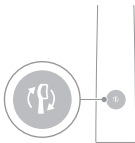
When the Philips Sonicare toothbrush is used in clinical studies, it must be used at high intensity with the handle fully charged and the EasyStart feature deactivated.

Brush Head Replacement Reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When you attach a new smart brush head for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and begins tracking brush head wear.
- 2 Over time, based on the pressure you apply and the total time the brush head has been used, the handle tracks brush head wear in order to determine the optimal time to replace your brush head. This feature guarantees the best cleaning and care of your teeth.
- 3 When the Brush Head Replacement Reminder indicator lights up amber, you should replace your brush head.

Note: All smart brush heads come with the brush head replacement reminder feature activated. To deactivate the brush head replacement reminder see 'Activating or deactivating features'.

Brush Head Mode Pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates with the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the brush head will remember the last used setting.

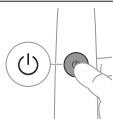
Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush Head Replacement Reminder
- Pressure Sensor Feedback

Note: You can only activate or deactivate the Brush Head Replacement Reminder when a smart brush head is on the handle.

Step 1: Place the handle on the charging stand or in the charging glass.

Step 2: Press and hold power on/off button for:

Easy start	Brush Head Replacement Reminder	Pressure Sensor Feedback
3 seconds	Up to 5 seconds	Up to 7 seconds


Step 3: Let go of the power on/off button when you hear:

1 beep	1 beep and then 2 beeps	1 beep, 2 beeps and then 3 beeps
--------	-------------------------	----------------------------------

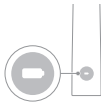
Together with:

If you see the battery indicator flash 2 times and hear 3 tones from low to high, then the feature has been activated.

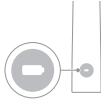OR

If you see battery indicator flash amber 1 time and hear 3 tones high to low, then the feature has been deactivated.

Charging and battery status

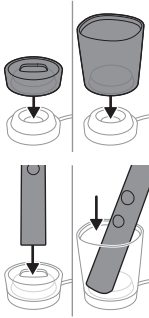
This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

Charging

Charging on charging base

- 1 Plug the cord of the charging base into an electrical outlet.
- 2 Place the charging stand or the charging glass on the charging base.

3 Place the handle on the charging stand or in the charging glass.

- The battery indicator lights up and the handle beeps twice. This indicates the toothbrush is charging.

4 The charger will charge your toothbrush until it is fully charged.

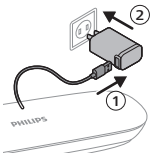
The handle is charging when the battery indicator is flashing. When the handle is fully charged, the battery indicator lights up for 30 seconds and then switches off.

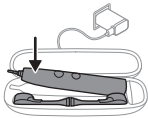
Charging with the travel case

Note: Accessories included may vary based on the model purchased.

- 1** Take the USB cable from the bottom of the charging travel case.
- 2** Connect the charging travel case via the USB cable to 5V DC input power either from:

- The USB adapter provided
- a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
- a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
- a UL-listed or certified 1310 Class 2 power supply.

- 3 Place the handle in the travel case. The travel case charges your toothbrush until your toothbrush is fully charged.

Battery status

Battery status (handle on plugged in charger or in charging travel case)

When the handle is placed on the charger or in the charging travel case, the battery indicator shows the battery level.

Battery status	Battery indicator
Charging	Flashes
Full	Lights up for 30 seconds and then switches off

Battery status (when handle is not placed on charger or in charging travel case)

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	Battery indicator	Sound
Full	Lights up	-
Partially full	Flashes	-
Low	Flashes amber	The handle emits 3 beeps after brushing

Battery status	Battery indicator	Sound
Empty	Flashes amber	The handle emits 2 sets of 5 beeps after brushing

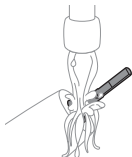
Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush on the charger when not in use.

Cleaning

Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

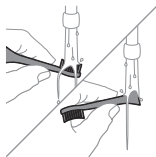Brush head

- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the brush head cap as often as needed.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Travel case

Rinse with warm water and use a damp cloth to clean inside of travel case.

Charging travel case

- 1 Unplug the charging travel case before cleaning.
- 2 Use damp cloth to wipe the surface of the travel case.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit

www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This is normal. A small gap is necessary to allow movement in order to deliver the vibrations needed to clean your teeth.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the handle.	Remove the brush head from the handle and re-set it on the metal shaft leaving small gap. (1-2mm)
	The battery is low.	Charge the toothbrush handle on the charger.
I have charged my toothbrush, but it only works for a short time.	The toothbrush handle is not placed vertically on the charger, resulting in no charge.	Place handle vertically on charger base and make sure the charger is plugged in.

Problem	Possible cause	Solution
The toothbrush stops working.	The battery is empty.	Charge the toothbrush handle on the charger.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX99XX).

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided.

Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 (North America), outside North America contact your local Philips Consumer Care Center. Internet information:

www.sonicare.com (North America) or
www.philips.com/support (outside North America).

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.

- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Disposal

CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

- Contact your local town or city officials for battery disposal information. You can also call 1-800-822-8837 or visit www.call2recycle.org for battery drop-off locations.
- For assistance, visit our website www.philips.com/support or call 1-800-682-7664 toll free.
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit www.recycle.philips.com.
- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

Removing the rechargeable battery

The built-in rechargeable battery **must only be removed by a qualified professional** when the

appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

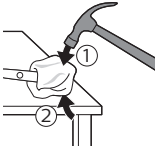
To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.


- 2** Remove and discard the brush head. Cover the entire handle with a towel or cloth.

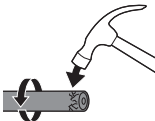
- 3** Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

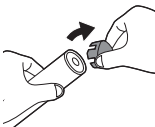
Note: You may have to hit on the end several times to break the internal snap connections.

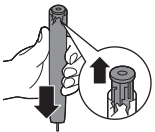
- 4** Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

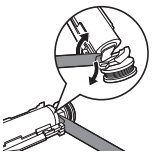- 5** Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

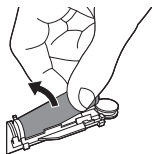- 6** Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.


- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.
- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

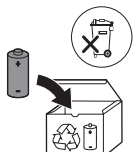The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

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©2022 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.
Philips Oral Healthcare, LLC.