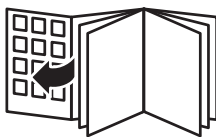
PHILIPS


Colgate

SonicPro 70

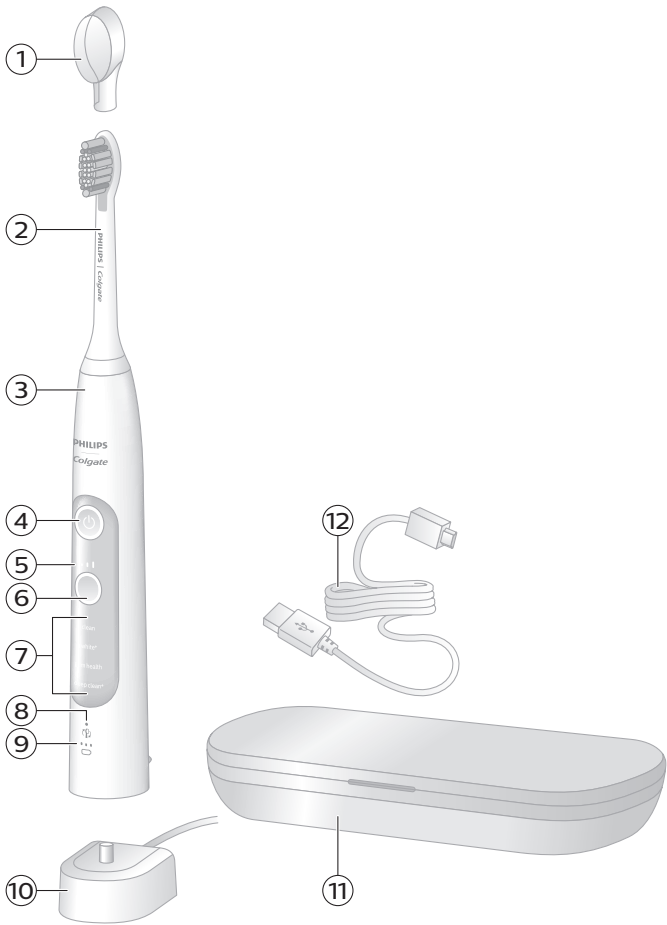1

— |

| —

— |

| —

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6 English

Introduction

Congratulations on your new Philips Colgate SonicPro power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using a combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: www.philips.com/welcome. Welcome to the Philips Colgate SonicPro family!

Intended Use

Philips Colgate SonicPro power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. They are intended for consumer home use. Use by children should be with adult supervision.

Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Danger

- Keep the charger away from water. Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to the wall socket.

Warning

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- If the appliance is damaged in any way (brush head, toothbrush handle and charger), stop using it.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.

- Always replace the charger with one of the original type in order to avoid a hazard.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see 'Warranty and support').
- Keep the cord away from heated surfaces.
- Do not use the charger outdoors or near heated surfaces.

Caution

- Do not clean the brush head, handle, or charger in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.
- This Philips Colgate SonicPro toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not use other brush heads than the ones recommended by Philips Colgate SonicPro.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

Battery safety instructions

- Only use original Philips accessories or consumables. Only use detachable supply unit HX6100.
- Charge, use and store the product at a temperature between 0 °C and 40 °C.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.

8 English

- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio frequency interface (Bluetooth) are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: **www.philips.com/support**

The radio frequency interface in this product operates at 13.56 MHz.

The maximum RF power transmitted by the appliance is 30.16 dBm.

Your SonicPro (Fig. 1)

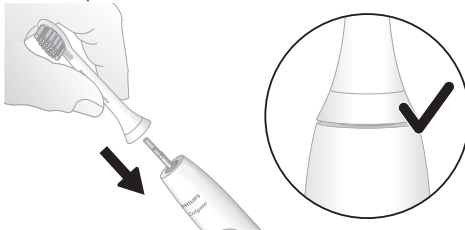
- 1 Brush head cap
- 2 Smart brush head(s)
- 3 Handle
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Modes
- 8 Brush head replacement reminder light
- 9 Battery indicator lights
- 10 Charger
- 11 Travel case *
- 12 Micro USB cable *

* Note: The content of the box may vary based on the model purchased.

Using your SonicPro

Brushing instructions

- 1 Push the brush head firmly onto the handle. Wet bristles and apply toothpaste.


Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.


Note: Keep the center of the brush in contact with the teeth at all times.

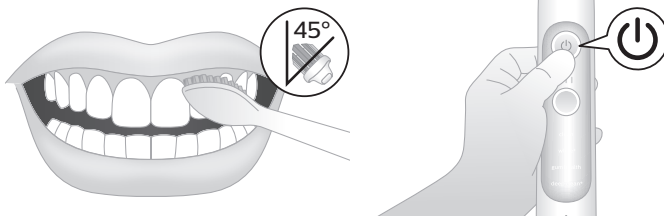- 3 Press the power on/off button to turn on the toothbrush.


- 4 Apply light pressure to maximize SonicPro's effectiveness and let the toothbrush do the brushing for you. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.

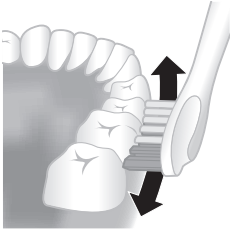10 English

Note: To learn how to brush with specific modes, see 'Modes'.

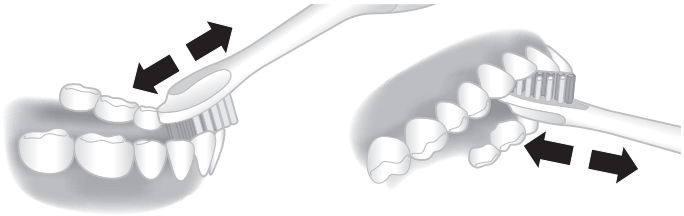- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note:

- A change in vibration from the handle alerts you when you apply too much pressure (see Pressure sensor).
- Your toothbrush is safe to use on braces, dental moderations and dental restorations (fillings, crowns, veneers).
- Do not use your power toothbrush in the shower.

Brush head

Your SonicPro comes with the 2-in-1 brush head for healthier gums and a deep clean. This brush head is equipped with the BrushSync technology.

Note: Your toothbrush is compatible with all other Philips Colgate SonicPro brush heads.

BrushSync Technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. This enables the following features:

- Brush head mode pairing
- Brush head replacement reminder

Brushing modes

Your power toothbrush is equipped with up to 4 different modes*: Clean, Deep Clean+, White+* and Gum health*.

Based on the smart brush head you use, the mode will be selected automatically.

To manually change the mode, press the mode/intensity button while the toothbrush is off.

* Note:

- The type and number of modes vary based on the model purchased.
- All brush heads work with each mode.

12 English

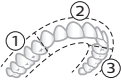
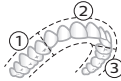
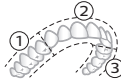
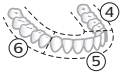
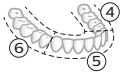
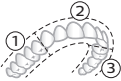
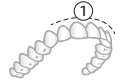
- When the Philips Colgate SonicPro toothbrush is used in clinical studies, it must be used at high intensity, with the EasyStart feature turned off.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates to the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the brush head will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended mode/intensity in 'Brushing Modes'.

	Clean and Deep Clean+	White+	Gum Health
Benefit	Plaque removal	Plaque and surface stain removal	Plaque removal and gum massage
Recommended intensity level	3	3	3
Total brushing time	Clean: 2 minutes Deep Clean: 3 minutes*	2 minutes and 40 seconds	3 minutes and 20 seconds

	Clean and Deep Clean+	White+	Gum Health
How to brush	Clean: Brush each segment for 20 seconds	Brush each segment for 20 seconds,	Brush each segment for 20 seconds,
			
			
	Deep Clean+: Brush each segment for 30 seconds	then brush segments 1 and 2 for 20 seconds each.	then brush segments 1, 2, 3 and 4 for 20 seconds each.
			
			

***Note:**

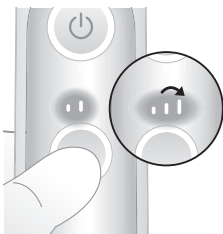
- The modes and brush head type vary based on the model purchased.

- When brushing with the Philips Colgate SonicPro App in 'Deep Clean' mode, the brushing time is 2 minutes.

Intensity settings

Your power toothbrush comes with 3 different intensity settings: low, medium and high. When attached to the handle, the smart brush head will automatically select the recommended intensity. To manually change the intensity, turn on the handle and press the mode/intensity button. While toothbrush is ON, press this button to change intensity setting.

14 English

Features

- EasyStart
- Connected Brush experience
- Brush head replacement reminder
- Pressure sensor
- BrushPacer

EasyStart

The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power.

The EasyStart feature works with all modes.

Note: Your toothbrush comes with the EasyStart feature deactivated. To activate EasyStart, see 'Activating or deactivating features'.

Philips Colgate SonicPro app - Connected brushing experience

The Philips Colgate SonicPro app pairs with your toothbrush to provide you with interactive brushing features. By brushing with your Philips Colgate SonicPro app account, you will be able to:

- Receive brushing suggestions and personalized feedback
- Track your brushing progress over time

To start using the connected brushing features:

- 1 Download the Philips Colgate SonicPro app to your phone.
- 2 Follow the steps in the app setup.
- 3 Pair your toothbrush with the app.

- 4 Create your account.
- 5 Personalize your app by answering questions about your brushing habits and interests.
- 6 Start your connected brushing experience.

Note: Make sure your phone's Bluetooth is turned on when using the connected brushing features.

Note: If you have any questions while setting up your account, be sure to review the Philips Colgate SonicPro Privacy Statement, available throughout the app setup process.

How to brush with the app

- 1 Go to the brushing screen in the app and select 'Start Brushing'.
- 2 First brushing starts on the upper right molar. You can customize your starting location in 'Settings'.
- 3 Follow the brushing guidance.

Tip: To ensure the best performance of the brushing sensors, keep your head facing forward.

- 4 Gently begin brushing. The app will highlight the correct segment. The timer will help you keep track.
- 5 Move your brush head when the highlighted segment changes.
- 6 Repeat, until you finish brushing all segments.
- 7 When you have finished brushing, review and touch up any missed areas. Follow the in-app instructions on how to use the touch-up feature.

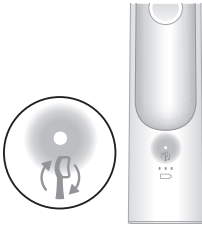
Brush head replacement reminder

Over time, the handle tracks the wear of your smart brush head by measuring:

- The overall pressure you apply while brushing.
- The total time you have brushed with your brush head.

Once your brush head is no longer effective, the brush head replacement reminder will flash amber. This indicates that it is recommended to replace your brush head.

16 English

Note: Your toothbrush comes with the brush head replacement reminder activated. To deactivate this feature, see 'Activating or deactivating features'.

Pressure sensor

Your SonicPro measures the pressure you apply while brushing to protect your gums and teeth from damage.

If you apply too much pressure, the handle changes its vibration until you reduce the pressure.

Note: Your toothbrush comes with pressure sensor activated. To deactivate this feature, see 'Activating or deactivating features'.

BrushPacer

The BrushPacer features indicates when you are supposed to move onto the next section of your teeth via a brief change in vibration.

The toothbrush turns off automatically at the end of the brushing session.

Charging and battery status

Charging on the charger

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
 - When fully charged, all battery indicator lights light up green.

Note: Your handle comes pre-charged for first use. After first use, charge for at least 24 hours.

Charging with the travel case

Note: Only the following SonicPro handles can be charged in the charging travel case: HX972G, HX972K, HX972V, HX972P, HX962C. The HX number can be found on the bottom of your handle.

Note: The charging travel case requires a USB adapter to charge the handle. The USB adapter is not included.

- 1 Connect the charging travel case to 5V DC input power either from:
 - a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
 - a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
 - a UL-listed or certified 1310 Class 2 power supply.

The USB adapter must be a safe extra low voltage source, compatible with the electrical ratings marked on the bottom of the device.

- 2 The travel case charges your toothbrush until it is fully charged.
 - When fully charged, all battery lights light up green.

Note: Charging the handle with the charging travel case may take longer than on the charger.

Battery status (handle on plugged in charger)

When the handle is placed on the charger or travel case, the battery indicator lights flash green based on the current battery level.

Number of flashing lights	Battery charge level
1	Low
2	Medium
3	High

Note: When the handle is fully charged, the three battery indicator lights light up green. After 30 seconds the battery indicator lights switch off.

Battery status (when handle is not placed on charger)

When removing the handle from the charger or travel case, the battery indicator lights indicate the status of the battery.

LED color	Number of LEDs	Battery charge level
Solid green	3	High
Solid green	2	Medium
Solid green	1	Low
Flashing amber	1 (with beeps)	Empty or near empty

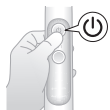
Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

Step 1: Place the handle on the charger

Step 2: Press and hold power button for:

EasyStart

Brush head
replacement
reminder

Pressure sensor

Up to 3 seconds

Up to 5 seconds

Up to 7 seconds

Step 3: Let go of the power button when you hear

1 beep

1 beep, 2 beeps
and then 3 beeps

1 beep and then 2
beeps

If you see battery indicator lights flashing green 2 times and hear low to high triple tones, then the feature has been activated.

OR

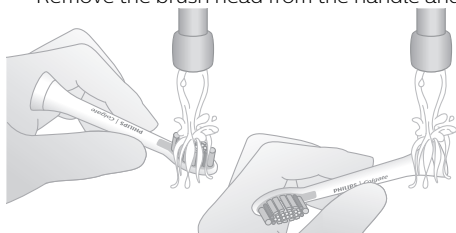
If you see battery indicator lights flashing amber once and hear high to low triple tones, then the feature has been deactivated.

Cleaning


- All parts can be cleaned with a damp cloth or rinsed with warm water.


- Remove the brush head from the handle and rinse it thoroughly.

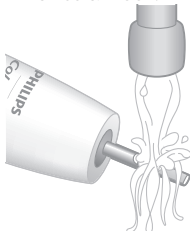- Rinse the entire handle, especially the brush head connection, at least once a week.

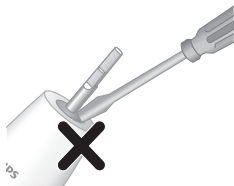Caution: Do not push on the rubber seal at the top of the handle. Gently clean around the rubber seal.

Cautions:

- Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.
- Do not use isopropyl alcohol, vinegar, bleach, or any other household cleaning products, to clean the handle as this may cause discoloration.
- Unplug the charger before you clean it.
- Do not use essential oils to clean the brush heads as this can cause damage.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Frequently asked questions

Question	Answer
Can I leave my toothbrush on the charger between brushings?	Yes, you can always keep your toothbrush on a plugged-in charger. This will not affect the battery lifetime.
Will other brush heads fit my Philips Colgate SonicPro toothbrush?	All Philips Colgate SonicPro electric toothbrush brush heads that snap on will fit your toothbrush.

Question	Answer
How often do I need to replace my brush head?	We recommend you to replace your Philips Colgate SonicPro brush head after 3 months of use (based on brushing twice a day for 2 minutes each session).The brush head replacement reminder will also notify you when it is time to change your brush head.
Is my toothbrush waterproof?	Yes, your toothbrush is waterproof. However, never put your handle in a dishwasher to clean it. And do not use it in the shower.
What phones are compatible with the Philips Colgate SonicPro app?	The Philips Colgate SonicPro app supports Android OS 5.0 and newer, and iPhone 9.0 and newer. Mobile devices must also have Bluetooth 4.0 and newer.
Can I use the Philips Colgate SonicPro app without Wi-Fi?	Wi-Fi is not needed to use the Philips Colgate SonicPro app.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

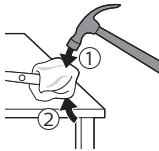
To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver.

Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.


- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

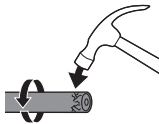Note: You may have to hit on the end several times to break the internal snap connections.

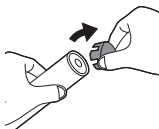
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.


24 English

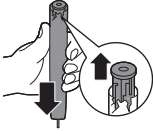- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

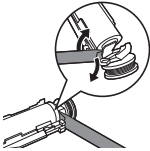
- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.


- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.
- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

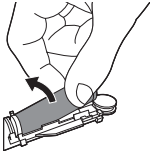Caution: Be aware of the sharp edges of the battery tabs to avoid injury to your fingers.

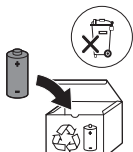- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

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© 2021 Koninklijke Philips N.V. (KPNV) All rights reserved. Philips and the Philips shield are trademarks of KPNV.
Colgate and the Colgate logo are the registered trademarks of Colgate-Palmolive Company.