PHILIPS

sonicare

For Kids

Design a Pet Edition

English 6
Azərbaycanca 21
Čeština 38
Dansk 54
Deutsch 70
Eesti 88
Español 104
Français 121
Hrvatski 138
Italiano 153
Latviešu 170
Lietuviškai 185
Magyar 202
Nederlands 219
Norsk 236
O'zbek 252
Polski 269
Português 286
Română 303
Shqip 320
Slovenščina 337
Slovensky 353
Srpski 370
Suomi 386
Svenska 402
Türkmençe 418
Ελληνικά 435
Български 453
Кыргызча 471
Македонски 489
Русский 507
Тоҷикӣ 525
Українська 543

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Introduction

This toothbrush enables you to obtain superior plaque removal, whiter teeth and healthier gums. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time.

Further support and product registration are available to you at: **www.philips.com/support**
Welcome to the Philips Sonicare family!

Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Warnings

- Keep the charger away from water.
- After cleaning, make sure the charger is completely dry before you connect.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use the charger provided with the product. If supplied, only use the original USB cable and USB power supply unit.
- Do not use the charger outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.

- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables. Only use detachable supply unit HX6110
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.

- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on recycling.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings).


Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Storage and transportation conditions

Temperature -10°C to 60°C (14 °F to 140 °F).

The following symbols may appear on the product:
Read Operator's Manual.

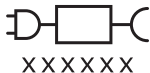Detachable Power Supply Part: The "xxxxxx" indicates the model number of the approved chargers to be used.

Intended use

Sonicare for Kids - Design a Pet Edition power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare for Kids - Design a Pet Edition power toothbrushes are intended

for personal use. Use by children should be with adult supervision.

Your Philips Sonicare (Fig 1.)

- 1 Handle
- 2 Power on/off button
- 3 Battery indicator

Accessories

- 4 Brush head cap
- 5 Brush head(s)
- 6 USB charger (wall adapter not included)
- 7 Stickers
- 8 Travel case

Note: Accessories included may vary based on the model purchased.

Note: Wall adapter not included. Only use 5V IPX4 adapter (USB Type A).

Sonicare app - Getting started

Download the Sonicare for Kids App to support healthy habit building. The App is compatible with wide range of iOS and Android (™) smart devices.

Note: If your device is not Bluetooth capable you can still use the Sonicare for Kids app on iOS and Android (TM) smart devices to enhance your brushing experience.

Getting started

Applying your favorite sticker


The Sonicare for Kids - Design a Pet Edition comes with several stickers to choose from. The stickers can be moved around to create the design you choose.

- 1 Apply your favorite sticker(s) to the handle

Note: Make sure the handle is clean and dry, and that stickers are fully pressed onto the handle with no air gaps to ensure full adhesion.

Note: Kids should not be left unattended when applying the stickers.

Attaching the brush head

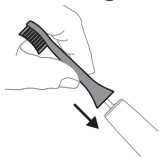- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

Note: It is normal to see a slight gap between the brush head and the handle.

Using your Philips Sonicare toothbrush


If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience. To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gently increases the power over the first 14 brushings to help you adapt to brushing with a Philips Sonicare toothbrush.

Brushing instructions

Note: For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.


- 1 Wet the bristles and apply a small amount of toothpaste.

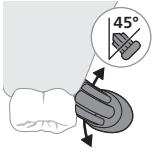
- Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

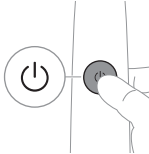
- Press the power on/off button to switch on the Philips Sonicare.

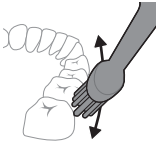- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

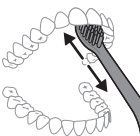- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

KidPacer

To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the KidPacer feature. The KidPacer divides the brushing time into 4 equal segments and indicates when you should move to the next area. The segments are indicated with a short tune. The toothbrush automatically stops at the end of a brushing session with a congratulatory tune.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush switches off.

Intensities

Your power toothbrush is equipped with 2 different intensities: gentle and extra gentle.

Note: When you use the toothbrush for the first time, the default setting is extra gentle.

To change intensities:

- Press the power on/off button once to turn on the toothbrush
- Press a second time within 2 seconds to change the intensity

Press a third time within 2 seconds to pause. After 2 seconds of brushing, pressing the power on/off button again will pause the toothbrush

	Extra gentle	Gentle
Benefit	Reduced power intensity used for training and introducing younger children to the sonic experience.	Full-strength intensity designed for older/independent brushers.
Recommended brush heads	Sonicare for Kids Compact or Standard brush head	Sonicare for Kids Compact or Standard brush head
Total brushing time	2 minutes	2 minutes

How to brush

Features

- KidPacer
- SmarTimer
- EasyStart

KidPacer

The KidPacer is an interval timer that helps your child to brush evenly through the mouth by dividing your mouth into 4 sections.

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. After a pause of 30 seconds, the SmarTimer resets.

EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare.


To deactivate EasyStart (see 'Activating or deactivating features').

Activating or deactivating features


You can activate or deactivate the following features of your toothbrush:

- EasyStart

Step 1: Place the handle on the charger.

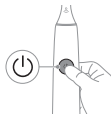
Step 2: Press and hold power on/off button for up to 3 seconds.

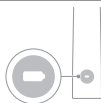If you see the battery indicator flash green 2 times and hear 3 tones from low to high, then the feature has been activated

OR

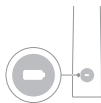If you see the battery indicator flash amber 1 time and hear 3 tones high to low, then the feature has been deactivated.


Note: Using the EasyStart feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.

Charging and battery status

This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

Charging

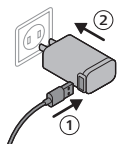USB charging

Note: Accessories included may vary based on the model purchased.

- 1 Plug the USB cord of the charging stand into a wall adapter that is connected to a electrical outlet.
- 2 Place the handle on the charger.
 - The battery indicator lights up and the handle beeps twice. This indicates that the toothbrush is charging.

Battery status

When the handle is placed on the charger, the battery indicator shows the battery level.

Battery status	Battery indicator
Charging	Flashes green
Full	Solid green for 30 seconds then switches off

Battery status (when handle is not placed on charger)

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	Battery indicator	Sound
Full	Solid green	-
Partially full	Flashes green	-
Low	Flashes amber	The handle emits 3 beeps after brushing
Empty	Flashes amber	The handle emits 2 sets of 5 beeps after brushing


Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush on the charger when not in use.

Cleaning

Toothbrush handle

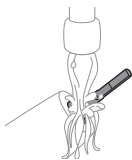
- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

Brush head

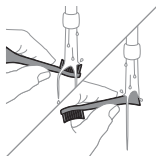- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Travel case

Rinse with warm water and use a damp cloth to clean the travel case.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Only use Sonicare for Kids Compact or Standard replacement brush heads.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX360XX).

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty exclusions

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removing the rechargeable battery

CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to



protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.


To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.

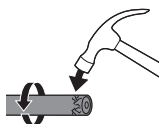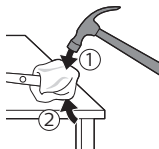
- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.


Note: You may have to hit on the end several times to break the internal snap connections.

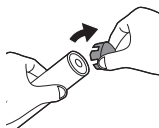- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

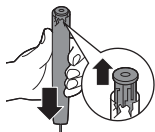
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

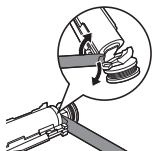- 6 Wedge the screwdriver between the battery and the frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the frame.

- 7 Insert the screwdriver between the bottom of the battery and the frame to break the metal tab connecting the battery to the printed circuit board. This will release the bottom end of the battery from the frame.

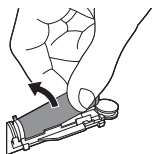
- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

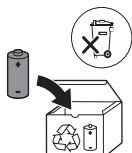Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

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©2023 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.