
PHILIPS

One

HY1200

English	3	Română	39
Azərbaycanca	4	Shqip	41
Čeština	6	Slovenščina	43
Dansk	8	Slovensky	45
Deutsch	10	Srpski	47
Eesti	13	Suomi	48
Español	15	Svenska	50
Français	17	Türkmençe	52
Hrvatski	19	Ελληνικά	54
Italiano	21	Български	56
Latviešu	23	Кыргызча	58
Lietuviškai	25	Македонски	61
Magyar	27	Русский	63
Nederlands	29	Тоҷикӣ	65
Norsk	31	Українська	67
O'zbek	33	Қазақша	69
Polski	35	Հայերեն	71
Português	37	ქართული	73

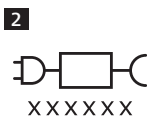Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Warnings

- Keep the charging cable away from water.
- After cleaning, make sure the charging cable is completely dry before you connect it to a USB power supply unit.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use the charging cable provided with the product with 5V DC input. If supplied, only use the original USB power supply unit.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country.
- Do not use near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.

Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.

- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on recycling.

Storage and transportation conditions

Temperature -10°C to 60°C (14 °F to 140 °F).

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Explanation of symbols

The following symbols may appear on the product:

- Read Operator's Manual (Fig. 1).
- Detachable Power Supply Part (Fig. 2): The "xxxxxx" indicates the model number of the approved chargers to be used.

Support

For all product support such as frequently asked questions, please visit www.philips.com/support.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 3).
- Follow your country's rules for the separate collection of electrical products and batteries.

CONTAINS NICKEL METAL HYDRIDE. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

Removal of built-in rechargeable battery

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Instructions for removal of built-in rechargeable batteries can be found on www.philips.com/support.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

If batteries are damaged or leaking, avoid contact with the skin or eye. If this occurs, immediately rinse well with water and seek medical care.

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All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.