

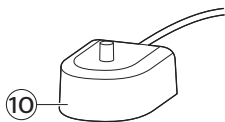
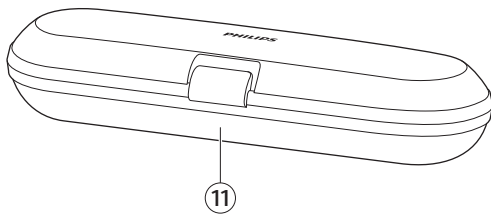
PHILIPS
sonicare

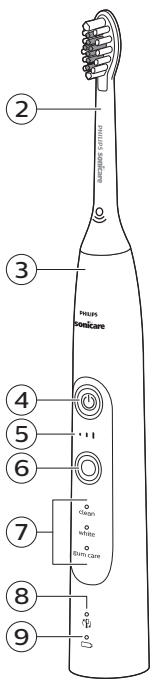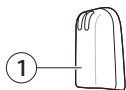English 6

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Introduction

Congratulations on your new Philips Sonicare power toothbrush! Superior plaque removal, whiter teeth and healthier gums are at your fingertips. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time. Further support and product registration are available to you at: www.philips.com/welcome. Welcome to the Philips Sonicare family!

Intended use

ExpertResults series power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. ExpertResults series power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Important safety information

Read this important information carefully before you use the appliance and save it for future reference.

Note: The accessories and features of different models vary. For example, some models come with a sanitizer or are equipped with Bluetooth.

Danger

- Keep the charger and/or sanitizer away from water. Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not immerse the charger and/or sanitizer in water or any other liquid.
- After cleaning, make sure the charger and/or sanitizer is completely dry before you connect it to the wall socket.

Warning

- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Keep the sanitizer out of the reach of children at all times.
- Stop using the sanitizer if the UV light bulb remains on when the door is open. UV light can be harmful to the human eye and skin.
- If the appliance is damaged in any way (brush head, toothbrush handle, charger and/or sanitizer), stop using it.
- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger and/or sanitizer.
- Always replace the charger and/or sanitizer with one of the original type in order to avoid a hazard.
- This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see 'Warranty and support').
- Keep the cord away from heated surfaces.
- Do not use the charger and/or sanitizer outdoors or near heated surfaces.

Caution

- Do not clean the brush head, handle, charger, charger cover or sanitizer in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- Consult your dentist if excessive bleeding occurs after using this appliance or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use this appliance.

- This Philips appliance complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your doctor or the appliance manufacturer of the implanted appliance prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Do not use other brush heads than the ones recommended by Philips Sonicare.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The UV light bulb of the sanitizer is hot during and immediately after the sanitizing cycle. Do not touch the UV light bulb when it is hot.
- Do not operate the sanitizer without the protective screen in place to avoid contact with a hot UV light bulb.
- If you have used the sanitizer for 3 UV clean cycles in a row, switch it off for at least 30 minutes before you start another UV clean cycle.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio frequency interface (Bluetooth) are in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address:

www.philips.com/support

The radio frequency interface in this product operates at 13.56 MHz.

The maximum RF power transmitted by the appliance is 30.16 dBm.

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush head(s)*
- 3 Handle
- 4 Power on/off button
- 5 Intensity light
- 6 Mode/intensity button
- 7 Modes
- 8 Brush head replacement reminder light
- 9 Battery indication
- 10 Charger
- 11 Travel case

*Note: The content of the box may vary based on the model purchased.

Using your Philips Sonicare

Brushing instructions

- 1 Push the brush head firmly onto the handle.

Note: It is normal to see a slight gap between the brush head and the handle.

- 2 Wet the bristles and apply toothpaste.
- 3 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

- 4 Press the power on/off button to turn on the Philips Sonicare.
- 5 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.


Note: To learn how to brush with specific modes, see 'Modes'.

Note: The bristles should slightly flare. Do not scrub.

Note: A change in vibration from the handle alerts you when you apply too much pressure while brushing.

Note: Your handle will briefly pause its vibration when its time to brush the next section, see 'QuadPacer'.

- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
 - Braces (brush heads wear out sooner when used on braces)
 - Dental restorations (fillings, crowns, veneers)

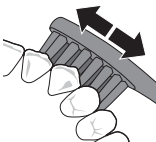Brush heads

Your brush heads

To see what type of brush head you have, check the brush head label at the bottom.

Brush head name	DiamondClean	Premium Plaque Control
Brush head label	W	C3
Benefit	Stain removal	Plaque removal


To explore our full range of brush heads, go to www.philips.com/toothbrush-heads for more information.

BrushSync Technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. The symbol at the bottom of the brush head indicates that it is equipped with the BrushSync technology.



This enables the following features:

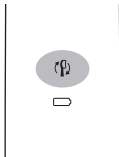- Brush head mode pairing
- Brush head replacement reminder

Brushing modes

Your power toothbrush is equipped with 3 different modes: Clean, White, and Gum Care. Based on the smart brush head you use, the mode will be selected automatically. To manually change the mode, press the mode/intensity button while the toothbrush is off.

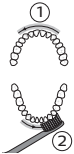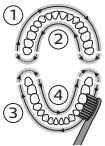	Clean	White	Gum Care
Benefit	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage
Re-recommended brush head	Premium Plaque Control (C3)	Diamond Clean (W)	Optimal Gum Care (G2)
Re-recommended intensity level	3	3	3
Total brushing time	2 minutes	2 minutes and 30 seconds	3 minutes
How to brush	<p>Brush each quadrant for 30 seconds</p> 	<p>Brush each quadrant for 30 seconds,</p>  <p>then brush the front top and bottom teeth for an additional 30 seconds.</p> 	<p>Brush each quadrant for 30 seconds,</p>  <p>then brush along the gumline for an additional minute.</p> 

Note: After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can

also brush your tongue, with the toothbrush turned on or off, as you prefer.

Note: All brush heads work with each mode.

Note: When the Philips Sonicare is used in clinical studies, it must be used in the Clean mode (for plaque and gum claims) or in White mode (for whitening claims), at high intensity, with the EasyStart feature turned off.

Intensity settings

Your power toothbrush comes with 3 different intensity levels: low, medium and high. When attached to the handle, the smart brush head will automatically select the recommended intensity. To manually change the intensity, turn on the handle and press the mode/intensity button.

Features

- EasyStart
- Brush head replacement reminder
- Pressure sensor
- Brush head mode pairing
- QuadPacer

EasyStart

The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power of the Philips Sonicare.


The EasyStart feature works with all modes.

Note: The EasyStart feature comes deactivated with your product. To activate EasyStart, see 'Activating or deactivating features'.

Brush head replacement reminder

The handle will track the wear of your brush head. It measures:

- The pressure you apply when you brush
- The time you have brushed with the brush head

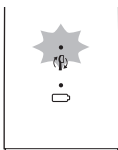Once your brush head is no longer effective, the brush head replacement reminder will blink amber. This indicates that it is recommended to change your brush head.

Note: The brush head replacement reminder comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

Pressure sensor

Your Philips Sonicare measures the pressure you apply while brushing to protect your gums and teeth from damage.

If you apply too much pressure, the handle will change its vibration until you reduce the pressure.

Note: The pressure sensor comes activated with your product. To deactivate this feature, see 'Activating or deactivating features'.

Brush head mode pairing

Brush heads with BrushSync technology are equipped with an RFID chip. The chip communicates to the handle and automatically pairs it with the recommended mode and intensity.

Note: If you change modes/intensity, the handle will remember the last used setting. If you wish to brush with the recommended setting, you can find the recommended modes/intensity in 'Brushing Modes'.

QuadPacer

The QuadPacer feature indicates via a brief change in vibration when you are supposed to move onto the next brushing section.

The toothbrush will automatically turn off at the end of the brushing session.

Charging and battery status

Charging

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
- 3 The charging indicator will flash green until your toothbrush is fully charged.

Battery status (handle on plugged-in charger)

LED color	Battery charging status
Flashing green	Charging
Solid green	Fully charged
Off	Fully charged and in sleep mode

When placing the handle on the charger:

- The handle will emit 2 short beeps to confirm the handle is placed properly and the handle starts charging.
- The battery light will flash green to indicate it is charging.
- The handle will show solid green when it is fully charged. In order to preserve the charge, the green light will switch off 30 seconds after the charging process is complete.

Battery status (when handle is not placed on charger)

When removing the handle from the charger, the battery light at the bottom of the handle will indicate the status of the battery.

LED color	Number of beeps	Status
Solid green	n/a	Full
Flashing green	n/a	Medium
Flashing amber	3	Low
Flashing amber	5	Battery empty

Note: Your toothbrush comes pre-charged for first use. After first use, charge for at least 24 hours.

Note: It is OK to leave the handle on the plugged in charger in between brushings. This will not impact the battery life negatively.

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
 - Brush head replacement reminder
 - Pressure sensor
- 1 Place the handle on the charger.
 - 2 Press and hold the power button for:

Feature	Press the power button for
EasyStart	3 seconds
Brush head replacement reminder	Up to 5 seconds
Pressure sensor	Up to 7 seconds


3 Let go of the power button when you hear:

Feature	You will hear
EasyStart	1 beep
Brush head replacement reminder	1 beep then 2 beeps
Pressure sensor	1 beep, 2 beeps then 3 beeps

If you see the brush head replacement reminder and the battery indicator lights **flashing green** and hear **low to high triple tones**, then the feature had been **activated**.
OR

If you see the brush head replacement reminder and the battery indicator lights **flashing amber** and hear **high to low triple tones**, then the feature had been **deactivated**.

Cleaning

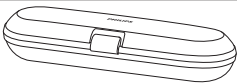			
	✓	✓	✗
	✓	✓	✗
	✓	✓	✗
	✓	✓	✗
	✓	✓	✗

Caution: Unplug the charger before you clean it.

Caution: Do not use essential oils to clean the brush heads as this can cause damage.

Caution: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean the handle as this may cause discoloration.

Caution: Do not push on the rubber seal at the top of the handle.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX751x).

Frequently asked questions

Question	Answer
Is my Philips Sonicare toothbrush water resistant?	Yes. However, do not fully submerge the handle in liquid for an extended period of time.
How often should I clean the connection between the brush head and the handle?	It is recommended to clean this after each brushing, or at least once a week.
Can I leave my Philips Sonicare toothbrush on the charger between brushings?	Yes, you can always keep your Philips Sonicare toothbrush on a plugged-in charger. This will not affect the battery lifetime.
How often can I use my Philips Sonicare toothbrush per charge?	Your Philips Sonicare toothbrush should last 2-3 weeks with standard brushings (2x/day) with a full battery. You do not have to charge your Philips Sonicare toothbrush on a daily basis.

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- This product contains a built-in rechargeable battery which shall not be disposed of with normal household waste. Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- The UV light bulb of the sanitizer contains mercury. Do not throw the UV light bulb away with the normal household waste at the end of its life. Hand it in at an official collection point for recycling.


- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery at the end of lifetime

Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. Any attempt to remove the battery will void the warranty.


To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 Remove and discard the brush head.
- 2 Deplete the battery until the handle no longer turns on.
- 3 Cover the entire handle with a towel or cloth.

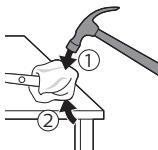
- 4 Hold the top of the handle with one hand.
- 5 With your other hand, use your hammer to strike the handle 0.5 inches above the bottom end on all sides to eject the end cap.


Note: You may have to strike it several times to break the internal cap connections.

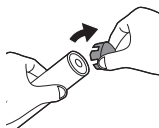
- 6 Remove the end cap. If the end cap does not release easily, repeat step 3 until the end cap is released.

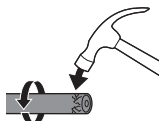
- 7 Holding the handle upside down with a cloth or towel, push the handle down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

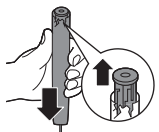
- 8 Remove the rubber battery cover with a screwdriver.

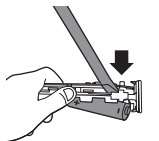- 9 Wedge the screwdriver between the battery and the internal black frame towards the bottom to separate them.

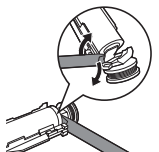
- 10 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

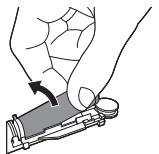- 11 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

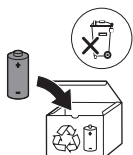- 12** Cover the battery contacts with tape to prevent any electrical shock from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- UV light bulb.

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©2023 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.