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©2024 AquaShield

All rights reserved

Philips and the Philips Shield Emblem are registered trademarks of Koninklijke Philips N.V. and are used under license.

This product has been manufactured by and is sold under the responsibility of Hongkong AquaShield Health Technology Company Limited and Hongkong AquaShield Health Technology Company Limited is the warrantor in relation to this product.

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## 1 Readme

Symbol Description

People Restricted or Prohibited from Use

Precautions for Use

Precautions for Use Places

Precautions for Maintenance

---

## 2 Product Features

---

## 3 Product Overview

---

## 4 How to Use

Operation

Charging

---

## 5 Maintenance

Cleaning

Storage

---

## 6 Technical Parameters

---

## 7 Troubleshooting

---

## 8 Service Guide for Product Warranty

Warranty and Service

Warranty Regulations

Warranty Card

---

## 9 Packing List

---

# 1 Readme

---

Please read this manual carefully before using the product and keep it properly. The manual is available at [www.philips.com](http://www.philips.com)

## Symbol Description

To avoid personal injury or property damage, be sure to observe the contents below.

- The following symbols are used to illustrate what you should observe.

---

 Indicates prohibition

 Indicates the actions that need to be subject to instructions

---

## People Restricted or Prohibited from Use

 The product should be used with caution by the following groups/ under the following circumstances

People who are being treated in medical institutions, or have the following conditions, please be sure to consult a doctor or physiotherapist before use.

- People who are implanted with pacemakers or other medical electronic devices in their bodies
- People with skin trauma
- People with fever or body temperature exceeding 38°C
- People who are easy to suffer from numbness
- People who have difficulty with respirators
- People with hypertension

- People with acute symptoms such as visceral diseases (gastritis, hepatitis, gastroenteritis)
- People who have lost skin temperature sensation
- People who have experienced skin burns
- People who have abnormal body sensations

 The product is prohibited from use by the following groups/ under the following circumstances

People who are prohibited from exercising or stretching activities as advised by doctor

For example, people with thromboembolism, severe aneurysm, acute venous aneurysm, various skin inflammations, or various skin infections are prohibited because they can easily cause deterioration.

- People who cannot express their will or people who cannot operate devices by themselves.
- Children under the age of 6 are not allowed to use. Adolescents aged 6-14 should use it under the guidance of an adult guardian.
- Please stop immediately in case of abnormal reactions such as back pain, dizziness, vomiting, palpitation during use.
- The device is set to automatically shut down after 10 minutes. Do not use it for more than 10 minutes each time. Keep the device rest for 20 minutes after each use to ensure the extension of its service life.
- Do not use the product under the following circumstances:
  - When you are severely tired
  - Within one hour before and after meals
  - When you are drunk

## Precautions for Use

- ❗ The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- ❗ This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- ❗ Children should be supervised to ensure that they do not play with the appliance.
- ❗ This appliance contains batteries that are only replaceable by skilled persons.
- ❗ If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- ❗ Must be only supplied at safety extra low voltage corresponding to the marking on the appliance.
- ❗ Ensure that the device is in good condition before each use. If any damage is found, regardless of size, stop using it immediately.
- ❗ Please remove the dust on the charger regularly. The accumulated dust on the plug can easily cause poor insulation due to humidity, and eventually lead to fire.
- ❗ Before use, be sure to confirm whether the switch is working properly.
- ❗ Please stop using the product immediately when the device cannot be used or you feel abnormal.
- ❗ Unplug the charger plug in case of power outage.
- ❗ Be sure to turn off the power when it is not in use.
- ⊘ Do not load the charger to cause damage.
- ⊘ Do not stand on the device or use the device as a toy so as to avoid accidents.

- ⊘ Do not disassemble, modify or repair the product by yourself, nor unzip the product by yourself.
- ⊘ Do not use the device when your body is wet, otherwise it is easy to get an electric shock.
- ⊘ Do not remove the cover when using the device, for it is easy to cause an accident.
- ⊘ Do not pull the device with force, nor stab it with a knife or a sharp object.

---

## Precautions for Use Places

- ⊘ Do not use the product near a stove, inside an electric heating quilt, or on heating appliances such as an electric blanket, which is easy to cause fire or failure.
- ⊘ Do not cover the product with a towel or cloth during use.
- ⊘ Do not use the device near a bathtub, shower, washbasin or other high-humidity vessel filled with water.
- ⊘ Do not use the device in places with room temperature above 35°C or in places with direct sunlight.
- ⊘ Do not use the device in bathrooms or places with high humidity.

---

## Precautions for Maintenance


- ❗ Be sure to turn off the power and unplug the charger during the maintenance.
- ❗ Be sure to turn off the power and unplug the charger if the device is not in use for a long time.
- ⊘ Do not drop main unit, nor strike it strongly.
- ⊘ Do not exert great force on the device on purpose.
- ⊘ Do not pour water on the product.

- ❗ Be sure to confirm whether the switch is normal in advance before it is used again if it has not been used for a long time.
- ⊘ Do not use volatile oil, thinner, alcohol, strong alkaline detergent, bleach, etc. for cleaning and maintenance.
- ❗ The product is equipped with a simple shockproof facility, and is suitable for air, railway, road and ship transportation. Keep away from rain and snow splashing, inversion and collision during the transportation.
- ❗ Keep the product out of reach of children.
- ❗ If the product is not in use for a long time, store it properly after fully charging and place it in a safe, dry, ventilated and cool place.
- ❗ If the product is not in use for a long time, it is recommended to fully charge it once every 1 month.

## 2 Product Features

This product combines a 2-level hot compress function with 6-point kneading to massage and relax the shoulders, neck and back. A variety of massage modes and intensity levels are available, suitable for different groups of people.

## 3 Product Overview
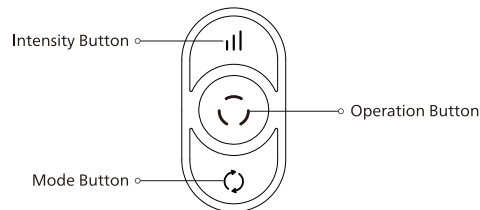
## 4 How to Use

### Preparations before Use

1. Check whether the device and its accessories are complete.
2. Remove the pendants and accessories during use so as to avoid affecting the massage parts.

### Operation

#### Operation button

Operation button - press and hold (for 2 seconds)

- a. In the shutdown state: Press and hold the operation button for 2 seconds to turn on the device (voice prompt: massage on). By default, the device enters the kneading mode, with rst-level strength and rst-level hot compress temperature (button indicator remains pink).
- b. In the working state: Press and hold the operation button for 2 seconds to turn off the device (voice prompt: massage off), and the indicator turns off.

#### Operation button - press

- a. In the shutdown state: The button operation is invalid and the device will not respond.
- b. In the working state: Press the operation button to switch temperature levels  
Cycle through the temperature levels in sequence: second-level hot compress temperature (voice prompt: 45°C; button indicator remains red) → heating off (voice prompt: turn off heating; button indicator remains white) → rst-level hot compress temperature (voice prompt: 40°C; button indicator remains pink)

#### Mode button

- a. In the shutdown state: The button operation is invalid and the device will not respond.
- b. In the working state: Press the mode button to switch the modes circularly.  
Cycle through modes in sequence: Relaxing mode (voice prompt: relaxing mode) → Neck-care mode (voice prompt: neck-care mode) → Kneading mode (voice prompt: kneading mode).

#### Intensity button

- a. In the shutdown state: The button operation is invalid and the device will not respond.
- b. In the working state: Press the intensity button to switch the levels circularly.  
Cycle through the levels in sequence: second-level strength (voice prompt: Intensity level 2) — rst-level strength (voice prompt: Intensity level 1).

### Charging

1. When the battery is low, the button indicator light ashes red (voice prompt: The battery power is low, please charge it in time).
2. The product cannot work when charging, and the button indicator ashes red slowly. After fully charging, the button indicator remains white.

#### Notes:

1. Please use the charging cable provided with the product.
2. When the battery is low, the button indicator light ashes quickly. After 2 seconds, the indicator light goes out and the product will automatically shut down.
3. The device is set to automatically shut down after 10 minutes. If you need to stop midway, press and hold the operation button to turn it off.

## 5 Maintenance

### Cleaning

\* The product cover cannot be removed.

- (1) Wipe off the dirt on the surface with a clean cloth;
- (2) For the deep dirt, dampen a cloth with the neutral detergent and then gently wipe it;
- (3) Prevent the water from entering the product;
- (4) Let it air dry.

❗ Do not put it in places with direct sunlight or in places with high humidity after use. If it will not be used for a long time, after removing the dirt, cover it with a cloth and keep it away from dust.

❗ Be sure not to leave the product in a charging or working state during maintenance. Do not use wet hands to unplug the power supply to prevent electric shock.

❗ If the product will not be used for a long time, it is recommended to fully charge it once every 1 month.

### Disposal

The lithium-ion battery used in this product is a valuable resource that can be recycled. Take out the lithium-ion battery when it is used up and discarded, and then hand it over to a professional institution that can recycle the rechargeable batteries for disposal.

Precautions for battery:

1. Do not replace the battery on your own during the normal use.

2. Be sure to remove the battery from the device before disposal, with the steps as follows:

- Remove the product cover;
- Remove relevant internal wires and then remove the battery;
- The product must be powered off before the battery is removed;
- The battery must be disposed of in a safe manner.

## 6 Technical Parameters

Product Name	Philips Massage Shawl
Product Model	PPM3522
Operating Voltage	7.4V ===
Operating Power	16W
Charging Voltage	5V ===
Charging Current	2A
Battery Capacity	2000mAh
Product Dimension	282x361x115mm
Net Weight	Approx. 1.28kg
Date of Manufacture	See product compliance certificate

## 7 Troubleshooting

This chapter lists the most common problems to the product. If the problems you encounter are not included in the following tables, please call our service hotline or contact the authorized service center in your region.

Fault Description	Possible Cause	Solution
Product does not work	Low battery capacity when it is not in use for long time	After recharging, press and hold the operation button for 2 seconds again to turn on and use
	Did not press the operation button to turned on	Press and hold the operation button for 2 seconds to turn on and use
	The product has just completed a 10-minute massage and automatically switches to the power-off state	If you intend to massage again, press and hold the operation button for 2 seconds again
Suddenly stop working halfway	Low battery capacity	After recharging, press and hold the operation button for 2 seconds again to turn on and use
A sound comes from the product (vibration of massage heads)	The sound is caused by the mechanical structure of product	This is a normal phenomenon
Product smoking	Short circuit in the internal wires of product	Stop using immediately and contact the designated service center for the product

If the product has been left unused for a long time, be sure to confirm that its functions are normal before using it again

Regular check	Check for the following faults	In the event of the above faults, stop using the device immediately and cut off the power supply to prevent accidents. Please send it to our maintenance center.
	<ul style="list-style-type: none"> <li>Whether the charging cable is hot</li> <li>Whether the power cord is deeply scratched or deformed, or has a burnt smell</li> <li>Whether there is abnormal sound or vibration during work</li> <li>Whether the product can work properly after power on</li> <li>Whether there are other abnormalities or faults</li> </ul>	

Note: If the fault cannot be eliminated or the product still cannot be operated satisfactorily, contact the after-sales service personnel, and do not disassemble or attempt to repair it by yourself.

## 8 Service Guideline for Warranty

### Warranty and Service

The product is a household appliance with exquisite design and high-quality component. Within two years from the date of purchase, you will receive free warranty service for any damage caused by the manufacturing process or components under normal operation confirmed by our maintenance service.

The warranty service does not include consumable components, transportation fee and door-to-door service. Please show the proof of purchase to the service personnel during maintenance.

You cannot receive free warranty service under the following conditions:

1. Damage by improper operation, such as improper components connection, improper installation, deprived of manual use, mis-operation.
2. Negligence and damage because of transportation or other accident.
3. Maintenance and modification without approval.
4. Other damage caused by force majeure (such as natural disaster, abnormal voltage, etc.).
5. Damage caused by general use outside home (e.g. industrial and commercial use).
6. Aging, wearing and tearing caused by normal operation of the product without affecting normal operation of the product.

This statement serves as a publicity for free maintenance within the specified time limit and conditions, and does not limit the legal rights of consumers.

Please contact our local maintenance service center for any information or difficulties or you can consult through the hotline.

HongKong

Hotline: 852-3188-4987

Service Email: philips.massage.hk@centurytel.com

Taiwan

Hotline: 0800-661-600

Service Email: Philips.massager@xander.com.tw

### Warranty Card

Product Model		Customer Name	
Seller		Contact Number	
Date of Purchase		Serial Number	
Service Center		Customer Address	
Maintenance Content		Date	service personnel

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