English	3	Português	52
Čeština	6	Română	55
Dansk	9	Shqip	58
Deutsch	12	Slovenščina	62
Eesti	17	Slovensky	65
Español	20	Srpski	68
Français	23	Suomi	71
Hrvatski	26	Svenska	74
Italiano	29	Ελληνικά	77
Latviešu	33	Български	80
Lietuviškai	36	Македонски	84
Magyar	39	Русский	88
Nederlands	41	Українська	91
Norsk	45	Қазақша	95
Polski	48	עברית	98


1

2

3
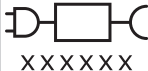
## Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

### Warnings

- Make sure the handle is completely dry before you connect the USB-C cable.
- Keep adapters and chargers away from water.
- After cleaning, make sure the adapters, chargers, and travel cases are completely dry before you connect these.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Only use original Philips accessories or consumables. Only use the USB cable provided with the product.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country. If the USB cable is damaged, stop using it.
- Do not use adapters, chargers and travel cases outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).
- Temperature (Transport and Storage) -10 °C to +60 °C (14°F and 140 °F).
- Make sure brush head and toothbrush are dry before storing in travel case.

## **MEDICAL WARNINGS**

- Stop using the appliance and consult your dentist/ doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

## **BATTERY SAFETY INSTRUCTIONS**

- Charge only with Philips WAA1001, WAA2001, safety approved adapters with input ratings (100 - 240V~; 50/60Hz; 3.5W) and output ratings (5Vdc; 2.5W).
- This appliance contains batteries that are only removable by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on recycling.

## Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Support

For all product support, power consumption details and the EU Declaration of Conformity, please visit [www.philips.com/support](http://www.philips.com/support).

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 1).
- Follow your country's rules for the separate collection of electrical products and batteries.

## Removal of built-in rechargeable battery

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Instructions for removal of built-in rechargeable batteries can be found on [www.philips.com/support](http://www.philips.com/support).

- This appliance contains NiMH rechargeable cell. **MUST BE RECYCLED OR DISPOSED OF PROPERLY.**
- Information on REACH substances, contained in Philips products, can be found on [www.philips.com/REACH](http://www.philips.com/REACH)

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**If batteries are damaged or leaking, avoid contact with the skin or eye. If this occurs, immediately rinse well with water and seek medical care.**

## LIST OF SYMBOLS USED

The following symbols may appear on the product:

- Read Operator's Manual (Fig. 2).
- Detachable Power Supply Part (Fig. 3): The „xxxxxx“ indicates the model number of the approved chargers to be used.

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All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.