

PHILIPS


Lumea

IPL Hair Removal


جهاز إزالة الشعر بتقنية IPL

7000 / 8000 / 9000

BRI930, BRI931, BRI932,
BRI933, BRI937, BRI938,
BRI983, BRI984

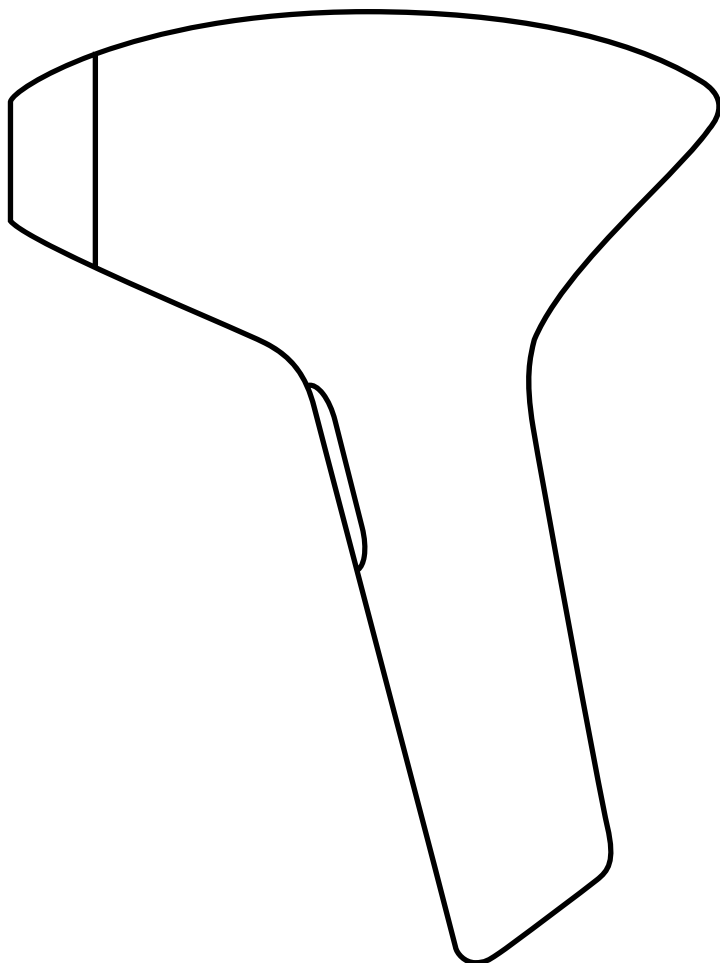

Skin tones - Suitable

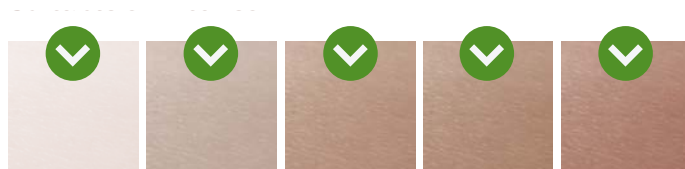

Fair


Light

Beige

Light
brown

Medium
brown

Unsuitable

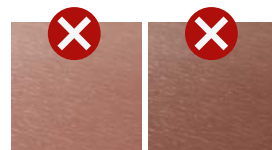

Dark
brown

Dark

Body hair colors - Suitable



Dark
blond

Brown

Dark
brown

Black

Unsuitable



Light
blond

Red

White/
Grey

English	5
العربية	70

Chapters

Welcome to Lumea	6
Intended purpose	6
Important safety instructions	7
Contraindications	7
Warnings	8
Cautions	10
Risks and benefits of using Lumea	11
Potential side effects of IPL treatment	12
What is IPL?	15
Lumea overview	16
<i>SenseIQ</i> 2.0 technology	17
Attachments	18
Treatment modes	19
Your treatment guide	20
Prepare for treatment	21
Test your skin	22
Perform a treatment	24
Clean and store	26
After treatment	27
Plan your next treatments	28
Troubleshooting	29
Additional information	31

Welcome to Lumea

Achieve gentle, fast and effective hair reduction with Lumea.
You can conveniently use Lumea in the comfort of your home.
Are you ready to get started?



Scan with your phone camera, QR app, or visit www.philips.com/myIPL

Note: In this user manual, Lumea refers to the Philips Lumea IPL 7000 Series, 8000 Series and 9000 Series.

Intended purpose

Philips Lumea is a home use device intended for the removal of unwanted hair for a single user. Lumea can be used by lay persons in the age range of 18-65 years with fair up to medium brown skin tones, and naturally dark blond, brown, dark brown and black hair. Lumea can be used on the body and on the female face* (below the cheekbone). Lumea is a handheld device intended for home use.

* Only if your Lumea includes the face attachment (BRI932, BRI933, BRI937, BRI938, BRI983, BRI984).


Is Lumea suitable for you?

Philips Lumea is not suitable for everyone. Check the skin tone and hair color tables at the start of this user manual and read the important safety instructions carefully to see if Lumea is suitable for you.

Unsuitable body hair colors: Treatment with Lumea is not effective if you have light blond, grey, red or white hairs as hairs with these colors do not absorb enough light.

Unsuitable skin tones: You cannot use Lumea if you have dark brown to dark skin tones as these skin tones absorb light more easily which could result in pain, skin redness or burns.

Important safety instructions

 **ALWAYS** read and follow the instructions carefully to fully benefit from the use of your Lumea and to minimize the risk of injury, skin reactions and side effects. Make sure to save these instructions for future reference.

Contraindications

- **DO NOT** use Lumea if you have had chemotherapy or radiation therapy in the last 3 months.
- **DO NOT** use Lumea if your sensitivity to pain is reduced after using medication (e.g. painkiller, numbing cream) or due to a medical condition.
- **DO NOT** use Lumea if you have a dark brown or darker skin tone.
- **DO NOT** use Lumea if you are pregnant, lactating or breastfeeding.
- **DO NOT** use Lumea if you are under the age of 18 years or above 65 years.

Light sensitivity:

- **DO NOT** use Lumea if your skin is extra sensitive to light, for example if you:
 - have a sun allergy;
 - are taking certain medication or herbal remedies. Check the instruction leaflet inside your medication package to see if what you are taking can cause increased sensitivity to sunlight.

Medical conditions:

- **DO NOT** use Lumea – or consult your doctor first – if you have:
 - epilepsy with flashlight sensitivity;
 - an active skin cancer or a history of skin cancer;
 - a skin pigmentation disorder such as vitiligo;
 - a history of poor wound healing;
 - an immunocompromised state caused by a disease such as AIDS or by immunosuppressive drugs;
 - a history of collagen disorders, including a history of keloid scar formation.

Warnings

Thermal safety - To avoid pain, skin redness or burns

- **DO NOT** use Lumea on lips, nostrils, ears, scalp, nipples, labia minora, perineum, around the anus, penis, or scrotum. These areas are highly sensitive and have thinner and darker skin.
- **DO NOT** touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.
- **DO NOT** use Lumea on an active skin condition, as this may worsen the condition or cause additional pain, skin redness or burns. Skin conditions may include:
 - Active bruises, rashes, sores, open wounds, scars or sunburns. Wait until the injury is completely healed before resuming treatment.
 - Active skin irritations, inflammations or infections. Wait until the skin is completely recovered before resuming treatment.
 - Damaged skin after using depilatory (hair removal) creams. Wait at least 2 days before resuming treatment.
 - Damaged skin after using skin peeling, (micro)dermabrasion or skin rejuvenation methods. Wait at least 1 week before resuming treatment.

Eye safety - To avoid eye damage

- **DO NOT** use Lumea close to or on the eye, eyelid or eyebrow. Use Lumea only below the cheekbones.

Electrical safety - To avoid electrical shock

- **DO NOT** use Lumea in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).
- **DO NOT** let water run into Lumea.
- **DO NOT** expose Lumea to large temperature changes, especially in a high humidity environment, without giving Lumea time to acclimatize for 30 minutes to the ambient temperature of the room in which it will be used (see 'Technical specifications').
- **DO NOT** try to open, repair, override or modify Lumea or the power adapter.
- **DO NOT** touch any inner part of a damaged Lumea or the power adapter.
- **DO NOT** use any other attachments or power adapter than supplied with Lumea.
- **ALWAYS** check if the voltage indicated on the power adapter corresponds to the local mains voltage before you connect the equipment.

General safety - To avoid serious injury

- **ALWAYS** check Lumea before you use it. Do not use Lumea or the power adapter if it is damaged. Always replace a damaged part with one of the original type.
- **ALWAYS** unplug Lumea, including the power adapter, after use and in case of a power failure to prevent damage to the equipment.
- **DO NOT** use Lumea directly over or close to any active implants such as a:
 - pacemaker;
 - neurostimulator;
 - insulin pump.Check the active implants' manual for the safe distance to other electronic devices.
- **DO NOT** wear the power adapter cord around your neck to avoid strangulation.
- **DO NOT** use Lumea if you have reduced physical, sensory or mental capabilities as this could result in incorrect usage.
- **DO NOT** use Lumea, including the power adapter cable, within 30 cm of WiFi and/or Bluetooth devices (portable RF communication equipment) or smartphones as it may cause Lumea to malfunction.
- **DO NOT** use Lumea next to or stacked with other equipment. This could result in improper operation.
- **DO NOT** leave Lumea unattended when it is turned on to avoid the risk of fire.
- **KEEP** Lumea out of the reach of children. Children should be supervised to ensure that they do not play with the equipment.
- **ALWAYS** check that the power adapter is connected to a properly installed and working wall socket to avoid the risk of fire.

Cautions

Thermal safety - To avoid pain, skin redness or burns

- **DO NOT** start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills). Do a skin test before starting/resuming your treatment, see section 'Test your skin'.
- **DO NOT** use Lumea with the body or precision attachment on the upper lip or other areas of the face. Only use the face attachment as it has an extra integrated filter for treatment on face below the cheekbones.
- **DO NOT** continue the treatment if you experience pain or a skin reaction that lasts longer than 24 hours. Reduce the light intensity setting or stop using Lumea.
- **DO NOT** expose your skin to sunlight unprotected when using Lumea. During the entire IPL treatment period, and especially in the first 48 hours after treatment, protect your skin from sunlight with clothing or apply sun-blocking products of not less than SPF50 on the non-covered treatment areas.
- **DO NOT** flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase performance.
- **DO NOT** use Lumea on skin that is darker than the surrounding treatment area. These highly pigmented skin areas can lead to excessive heating of the skin. Darker areas may include:
 - dark spots such as moles, warts, large freckles, birthmarks, age spots or dark veins;
 - marks added to the skin, such as tattoos or (permanent) makeup.
- **DO NOT** use Lumea on male neck and face. The hairs in these areas have different characteristics.
- **DO NOT** mark your skin/treatment area with any type of pencil or pen.
- **DO NOT** use Lumea if there is any lotion, gel or visible hair on your skin or on the device. Only treat skin that is clean, dry and hair-free. Make sure Lumea is clean and dry as well.
- **DO NOT** use Lumea if the light exit window or window frame is damaged. Check Lumea for damage before use.
- **DO NOT** subject Lumea to heavy shocks to avoid risk of damage to the lamp.
- **DO NOT** treat more often than once every 2 weeks.

Eye safety - To avoid eye discomfort

- **ALWAYS** use Lumea in a well-lit room. Your eyes are more sensitive to stray light in darker conditions.

General safety - To avoid injury

- **DO NOT** use Lumea over vascular skin disorders (e.g. varicose veins).
- **ENSURE** that the power adapter can be easily disconnected from the electrical outlet in case of a malfunction. The power adapter is considered a disconnection device.
- **DO NOT** share Lumea with other persons for hygienic reasons. Lumea is intended for single person use only.

Notice

Any serious incident that has occurred in relation to this device should be reported to Philips via www.philips.com/support and to the competent authority National Centre for Medical Device Reporting (NCMDR) website <https://ncmdr.sfda.gov.sa/default.aspx>.

Risks and benefits of using Lumea

The safety and effectiveness of Lumea was tested in a clinical study involving 94 women (aged 18-64 years). None of the subjects left the study early due to problems associated with the use of Lumea.

Out of the 94 subjects, 31 experienced one or more mild and transient side effects: pain (26), itchiness (8), redness (5), folliculitis (2), burning sensation (1), superficial crusting (1) and swelling (1).

Results after 4 bi-weekly treatments:

- 75 women received treatments on their upper lip and achieved up to 61.1 % hair reduction.
- 76 women received treatments on their armpit area and achieved up to 60.6 % hair reduction.
- 77 women received treatments on their bikini area and achieved up to 61.9 % hair reduction.
- 76 women received treatments on their legs and achieved up to 93.1 % hair reduction.

35 women completed 8 additional touch-up treatments (once every 4 weeks) and achieved up to 67.0 % hair reduction over all treated areas.

Potential side effects of IPL treatment

The table below lists known side effects that may occur following treatment with IPL devices.

Common side effects	What is observed?	How long does it last and how can it be helped?
Pain/ uncomfortable sensation	Local pain or discomfort on the treated area during or after treatment	If your treatment feels painful/ uncomfortable, reduce the light intensity setting or stop using Lumea. Pain or discomfort after treatment can last up to 1 week.
Itching sensation	Perceived itching in the treated area	These reactions typically resolve within hours but can last up to 1 week. You can cool the area with a cold pack.
Warm/burning sensation	Local warm or burning feeling on the treated area	
Skin redness (erythema)	Redness without a clear edge in the treated area	This reaction typically resolves within hours but can last up to 1 day. You can cool the area with a cold pack. Skin redness lasting longer than 1 day is an uncommon side effect of IPL and a sign of a skin burn. Wait with the next treatment until the skin has healed completely. Make sure you use a lower light intensity for next treatments.

Uncommon side effects	What is observed?	How long does it last and how can it be helped?
Swelling (edema)	Swelling around hair follicles appears as bumps in the treated area, sometimes accompanied by redness.	This reaction typically resolves within hours but can last up to 1 week. You can cool the area with a cold pack.
Skin dryness	Skin appears rough and flaky in the treatment area.	This reaction typically lasts between 1 day and 1 month. You can apply a non-scented moisturizer to the treated area.
Burns and blisters	Extensive skin redness, swelling and possibly blisters, accompanied by pain	These reactions can last between 1 day and 1 month and may require medical treatment or medication. We advise you to consult a doctor.
Superficial crust	A sharply defined brownish and crusty area that can occur on darker skin tones	Wait with the next treatment until the skin has healed completely. Make sure you use a lower light intensity for next treatments.
Folliculitis	Swelling and redness around hair follicles, accompanied by itching and/or pustule formation	Folliculitis may require medical treatment or medication. We advise you to consult a doctor. Folliculitis may last several days up to several weeks. Wait with the next treatment until the skin has healed completely. Make sure you use a clean razorblade when shaving before treatment.
Pigmentary changes/skin discoloration	Skin of the treated area appears darker or lighter than the surrounding skin.	This reaction typically lasts several weeks and, in rare cases, up to a year. Do not treat areas with pigmentary changes/skin discolorations until your skin has regained its normal skin tone. Make sure you use a lower light intensity for next treatments.
Skin infection	Appears as inflamed skin. This is a secondary effect which may occur after a skin burn or other injuries.	Skin infection requires medical treatment or medication. We advise you to consult your doctor. Do not treat over skin infections.

Uncommon side effects	What is observed?	How long does it last and how can it be helped?
Scarring	A scar may occur as a secondary effect after a burn or skin infection.	It is important to consult your doctor if you have a burn or folliculitis to prevent formation of a scar. Do not treat over scar tissue.
Herpes outbreak	Small blisters or open sores, typically near the lip or genitals, experienced by users infected with Herpes Simplex virus	An outbreak typically lasts for several days to weeks and can be shortened by taking anti-viral medications prescribed by your doctor.
Headache and afterimages	Sensitive persons may experience headache or afterimages when the light flash reaches the eye.	These reactions typically last between minutes to several hours. To prevent them, make sure you use Lumea in a well-lit room and press the device firmly onto your skin before flashing.

Hair removal by Intense Pulsed Light can cause increased hair growth in some individuals. Body areas at risk include women's jawline, neck and areas above the cheekbone and men's upper back and shoulders. Based upon current data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage.

What is IPL?



Scan the QR code to watch a short video.



1. IPL stands for Intense Pulsed Light.

Gentle pulses of light are absorbed by the roots of your hair.



2. These pulses of light cause the hair follicles to go into a resting phase.

The hair sheds naturally, and growth is reduced for a certain number of months.




3. Hair does not grow continuously: it grows in cycles.

Not all the hairs in one area are growing at the same time. Since IPL is only effective on hair that is actively growing, it is important to repeat treatments as indicated.

Lumea overview



Scan the QR code to watch a short video.

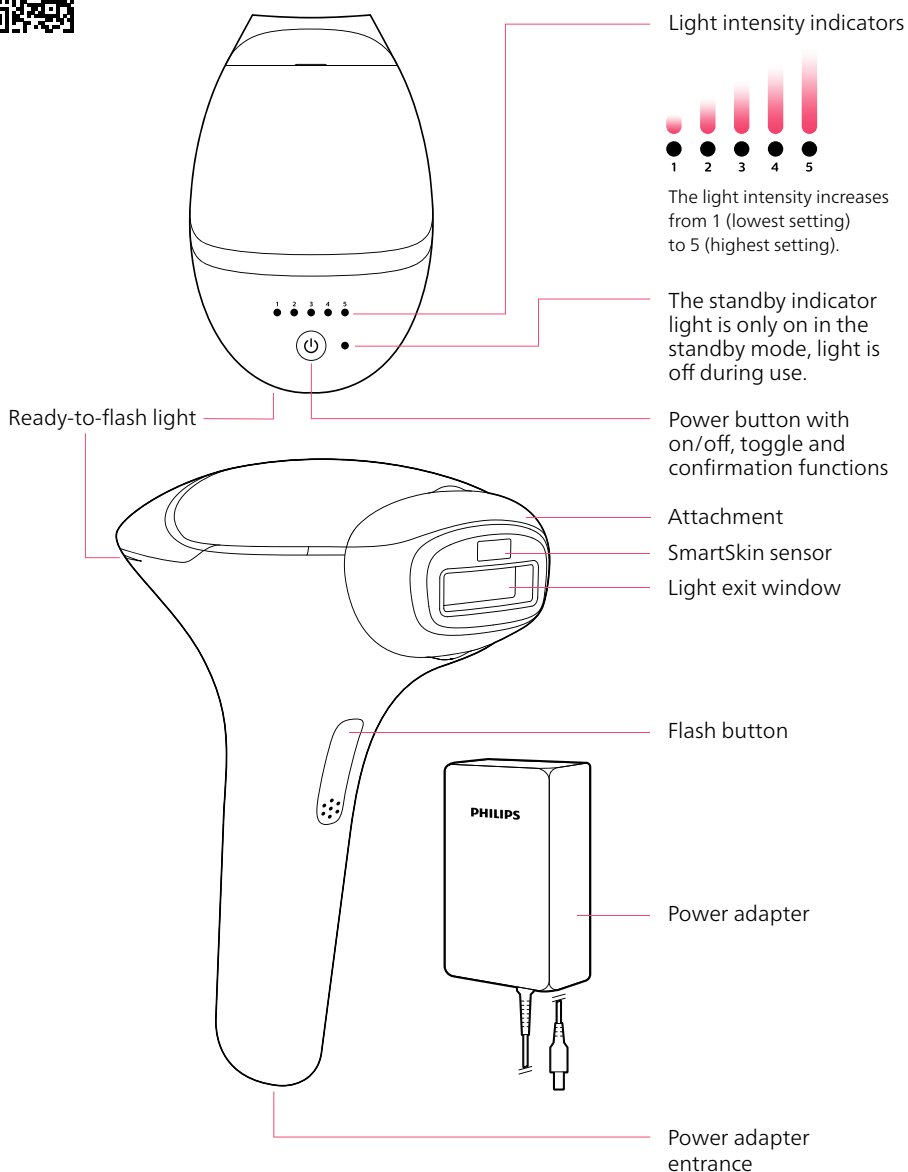

Not shown:
- Cleaning cloth
- Beauty bag or premium beauty pouch, or luxury beauty box

SmartSkin sensor

SmartSkin sensor determines your skin tone and will help you select a comfortable setting by lighting up one of the light intensity indicators on top of your device.

It also automatically prevents you from treating areas of your skin that are too dark for safe use.



Ready-to-flash light

The white ready-to-flash light will let you know when you have good skin contact and your Lumea is ready to flash.

If the ready-to-flash light blinks orange, your skin in this body area is too dark to treat. Treat a lighter colored body area.



SenseIQ 2.0 technology



SenseIQ 2.0
technology

SenseIQ 2.0 technology for a personalized hair removal treatment, our signature combination of:






- Improved smart sensors that measure your skin tone and support you in selecting a setting for a comfortable treatment.
- Unique attachments designed for each curve of your body.
- Philips Lumea IPL app to support your treatment routine.

Attachments



Scan the QR code to watch a short video.

Your Lumea comes with different attachments you can use on different body areas. Select the attachment that is intended for the area you want to treat.


	Attachment type	BR1930	BR1931	BR1932	BR1933	BR1937	BR1938	BR1983	BR1984
	Body attachment <ul style="list-style-type: none"> - largest window size; - a curved-in design to effectively cover and treat areas below the neckline. 	✓	✓	✓	✓	✓	✓	✓	✓
	Face attachment <ul style="list-style-type: none"> - smallest window size for precise treatment below the cheekbones; - small, flat design; - extra integrated filter for treatment on the face below the cheekbones. <p>⚠ DO NOT use Lumea on male neck and face. The hairs in these areas have different characteristics.</p> <p>Tip: Place your tongue between your upper lip and your teeth or puff out your cheeks to make the treatment easier.</p>								
	Precision attachment <ul style="list-style-type: none"> - medium window size; - curved-out design for precise and effective use on bikini and armpit areas. 				✓	✓			
	Precision attachment for bikini <ul style="list-style-type: none"> - medium window size; - curved-out design for precise and effective use on the bikini area. 							✓	✓
	Precision attachment for armpits <ul style="list-style-type: none"> - medium window size; - curved-out design for precise and effective use on the armpit area. 							✓	✓

Treatment modes




Scan the QR code to watch a short video.

Your Lumea has two treatment modes for more convenient use on different body areas.

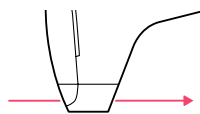

Slide & flash

For larger areas (legs), keep the flash button continuously pressed as you slide Lumea over your skin. Lumea keeps releasing flashes until you release the button or if you lose good skin contact.

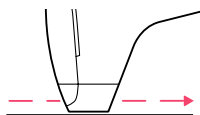

Stamp & flash

For small body areas (bikini, armpit), place Lumea on your skin and press the flash button. Lift Lumea from your skin and place it beside the spot previously treated. Press the flash button again to release another flash.

Your treatment guide

Please read the important safety instructions in the user manual to check if a Lumea treatment is suitable for you, so you can fully benefit from your device.

Before your treatment



1. Prepare for treatment
Remove hair, ensure that your skin is clean and dry



2. Test your skin
Find your comfortable setting



3. Wait 24h for skin reactions

Treatment



4. Perform a treatment



5. Clean and store



6. Plan your next treatments




Scan the QR code and watch our how-to videos for step-by-step guidance on performing a Lumea treatment, along with general videos about your Lumea device.

Prepare for treatment



Scan the QR code to watch a short video.

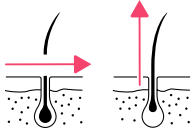

1. Remove hair

Wax, shave or epilate before using Lumea.


Do a treatment before new hair is visible on your skin.

If you choose to wax, wait 24 hours before you use Lumea.



2. Clean and dry skin

No water, oil, cream or body milk should be left on your skin.

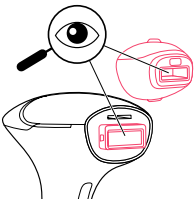

3. Choose the right attachment and check your Lumea for damage and dirt



ALWAYS check Lumea before you use it. Do not use Lumea or the power adapter if it is damaged. Always replace a damaged part with one of the original type.



DO NOT use Lumea if there is any lotion or gel on your skin or on Lumea. Only treat skin that is clean and dry. Make sure Lumea is clean and dry as well.

Test your skin




Scan the QR code to watch a short video.

⚠ DO NOT start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills).

Before your first treatment on each body area or after exposing your skin to sunlight unprotected, test your skin's reaction and find your suitable setting for treatment. It is recommended to use Lumea out of direct sunlight. The sunlight can influence the SmartSkin sensor.


1. Follow the steps in section 'Prepare for treatment'

2. Check your Lumea for damage and dirt

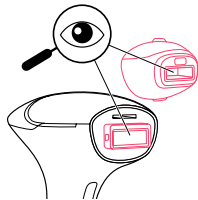


3. Attach the right attachment

Click on the attachment for the area you want to treat.

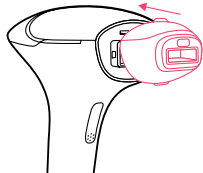

4. Plug in

Plug in Lumea to the wall socket.

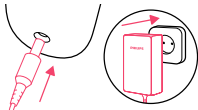


5. Press the power button

The 5 light intensity indicators above the power button will light up continuously from left to right and back, until you place Lumea on your skin.



6. Place Lumea on your skin

It must sit flush against your skin.

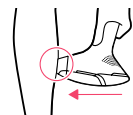



7. Get your light intensity setting

One of light intensity indicators will turn white to indicate the proposed setting for you. The setting is selected automatically after 5 seconds.

Note: If the ready-to-flash light blinks orange your skin may be too dark to treat, please consult chapter 'Troubleshooting'.



8. Find your comfortable setting


Place Lumea back on your skin.

The ready-to-flash light will light up white.



Press the flash button 1 time. If the flash is comfortable, you can increase the setting. If you experience pain, lower the setting. Toggle through the settings by pressing the power button.

Change the setting, move Lumea to the next spot and flash again. Repeat until you find the highest comfortable setting.

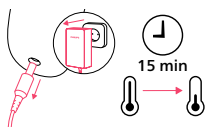


9. Turn off, unplug



DO NOT touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.

Hold the power button for 2 seconds. Unplug your Lumea and the power adapter from the wall socket. Let your Lumea cool down.

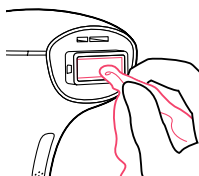

10. Clean all parts of your Lumea after use

It is important to clean all parts of your Lumea. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned. Follow the steps in section 'Clean and store'.



11. Wait 24 hours and check your skin for any reactions

If there are no skin reactions after 24 hours, start the treatment with the highest comfortable setting.

If the highest comfortable setting resulted in a skin reaction, choose a lower setting that did not result in any skin reactions.

If you did not test a lower setting, repeat the skin test.

Perform a treatment




Scan the QR code to watch a short video.

Test your skin's reaction before treating a new body area or after exposing your skin to sunlight unprotected.


1. Follow the steps in section 'Prepare for treatment'

2. Check your Lumea for damage and dirt

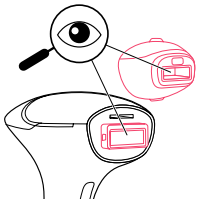

3. Attach the right attachment

Click on the attachment for the area you want to treat.

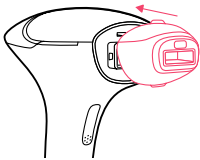

4. Be comfortable

Make yourself comfortable before you start the treatment.

It is recommended to use Lumea out of direct sunlight. The sunlight can influence the SmartSkin sensor.



5. Plug in. Press the power button

The 5 light intensity indicators above the power button will light up continuously from left to right and back, until you place Lumea on your skin.



6. Select setting


Press the power button to toggle through the settings and select your highest comfortable setting that did not result in any skin reactions.



7. Begin treatment

Place Lumea on your skin.

The ready-to-flash light will light up white.

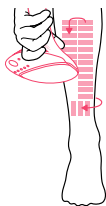


8. Flash. Then move

⚠ DO NOT continue the treatment if you experience pain. Reduce the light intensity setting or stop using Lumea.

DO NOT flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase performance.


Press the flash button. Move to an adjacent spot and flash again.

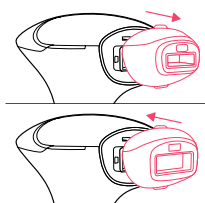

9. Treat another body area

⚠ DO NOT start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills). Do a skin test before starting/resuming your treatment, see section 'Test your skin'.


If you want to treat another body area, make sure you select the right attachment and start again with step 1. If you are treating a body area for the first time or after exposure to sunlight unprotected, make sure you perform a skin test. Wait 24 hours and check your skin for any reactions.

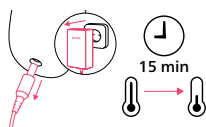

10. Finished. Turn off and unplug

⚠ DO NOT touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.

Hold the power button for 2 seconds. Unplug your Lumea and the power adapter from the wall socket. Let your Lumea cool down.

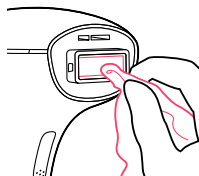

11. Clean all parts of your Lumea after use

It is important to clean all parts of your Lumea. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned. Follow the steps in section 'Clean and store'.




12. Schedule your next treatment

Follow the steps in section 'Plan your next treatments'.

Clean and store




Scan the QR code to watch a short video.

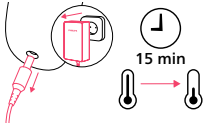

1. Turn off, unplug and cool your Lumea

⚠ DO NOT touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.

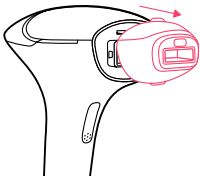

2. Click off the attachment


Gently pull it off.



3. Lightly moisten your Lumea cleaning cloth


⚠ DO NOT use Lumea in wet surroundings.
DO NOT let water run into Lumea.

Dampen the soft Lumea cleaning cloth with a few drops of water.

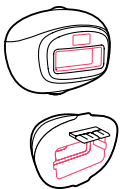


4. Clean the attachments

Clean each part of the attachment you used, including the metal reflectors.

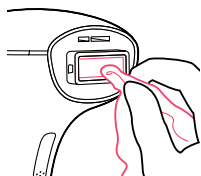


5. Clean the light exit window

Make sure there is no dirt, residue or moisture left on the surface of the light exit window.

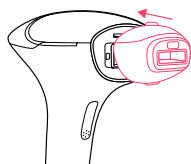

6. Click on an attachment

Make sure the light exit window of your Lumea is protected from damage by clicking on one of the attachments.

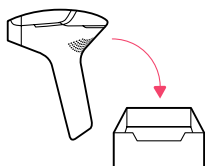

7. Store

Store all parts in the original packaging or in the provided bag, pouch or box in a dust-free and dry place.

After treatment

After use, if you want to apply lotions, creams, deodorants, moisturizers or cosmetics to the treated areas, always check the instruction leaflet of the product you want to apply to make sure it is suitable for the condition of your skin at that moment.

Plan your next treatments



Scan the QR code to watch a short video.


Follow the indicated treatment schedule.



24 hours

Test phase

Before your first treatment on each body area or after exposing your skin to sunlight unprotected, do a skin test by following the steps in section 'Test your skin'. Wait 24 hours after the skin test before you start your treatment phase.

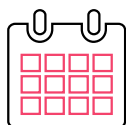

Treatment phase: 4 bi-weekly treatments

During the treatment phase which lasts 4 treatments, perform a treatment every 2 weeks.

Perform a treatment within 3 days before or after the planned treatment date during the treatment phase.



Touch-up phase: once every 4 weeks

After the treatment phase, use your Lumea once every 4 weeks to maintain results.

Perform a treatment within 4 days before or after the planned treatment date during the touch-up phase.

Note: If you miss a treatment, consult chapter 'Troubleshooting' for guidance. If you deviate from the treatment schedule, you may not achieve the intended results.

Troubleshooting

This chapter summarizes the most common problems you could encounter with your Lumea or if you deviate from your treatment schedule. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country as indicated in the included international warranty leaflet.

Common problems while using your Lumea:

Problem

My Lumea is not flashing / the ready-to-flash light turns orange and starts to blink.

My Lumea produces a strange smell while flashing.

The light intensity indicators on top of my Lumea are blinking and the ready-to-flash light is orange and blinking.

Solution

1. Your skin may be too dark. Check the skin tone/hair color tables inside the cover of this booklet.
 2. Make sure your attachment is (properly) connected to the front of your Lumea.
 3. Make sure you have full skin contact for at least 5 seconds. In some areas (sensitive and bony area) it might be difficult to make good skin contact. No light should be able to enter between Lumea and your skin. Try repositioning your Lumea.
 4. Make sure you press the flash button on the inside handle to treat an area of your skin.
-

1. Make sure you have removed all hairs in the area you want to treat, before using your Lumea. Burned hairs can cause smells.
 2. Make sure your skin is clean and dry before starting a treatment.
 3. Make sure your Lumea and its attachments are clean and dry before starting a treatment. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned.
-

1. Lumea is in error mode. Press the power button for more than 2 seconds, Lumea should go back to standby. Or unplug Lumea from the mains, it should turn off.
2. Wait 5 minutes. Lumea will automatically go back to standby after 5 minutes without help.

What to do if you deviate from your treatment schedule:

Phase	Problem	Solution
Treatment phase	Missed treatment. <i>I did my treatment more than 3 days after the originally planned date.</i>	Try to stay closer to your planned treatment date. You can continue with the original planned schedule. You may add one additional bi-weekly treatment to your treatment phase. If you deviate from the treatment schedule, you may not achieve the intended results.
	Last bi-weekly treatment. <i>I finished my treatment phase, but I feel I did not reach the claimed results.</i>	You may choose to add an additional bi-weekly treatment to your treatment phase, but claimed results are based on the indicated treatment schedule and individual results vary.
Touch-up phase	Out of treatment window. <i>I did not treat within 4 days of my scheduled treatment date.</i>	Please follow the correct schedule. If you deviate from the treatment schedule, you may not achieve the intended results. You can continue with the original planned schedule or you may choose to restart the treatment phase, but claimed results are based on the indicated treatment schedule and individual results vary.

Additional information

Disposal

Follow your country’s rules for the separate collection of electrical and electronic products. Correct disposal helps prevent negative consequences for the environment and human health.

Electromagnetic compatibility

This product has been tested and found to comply with IEC 60601-1-2 for electromagnetic compatibility (EMC) as Class B according to CISPR 11. Tests levels are listed in the accompanying tables. Only the power adapter that is supplied with Lumea should be used with your device.

Guidance and manufacturer’s declaration on electromagnetic emissions and immunity

Lumea is intended for use in the electromagnetic environment listed here. The user of the device should assure that it is used in such an environment to prevent improper operation.

Emissions	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	Lumea does not use RF energy for its intended function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	Lumea is suitable for use in locations in residential environments and in establishments directly connected to a low voltage power supply network which supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

Immunity test	Immunity test level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	Floors can be a source of static electricity. When the air is very dry (<20 % RH) the user is advised to use Lumea in a room with concrete or wooden floor (instead of using it in a room with a synthetic floor).
Radiated RF EM fields IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	±8 kV contact ±15 kV air
Proximity fields from RF wireless Communications Equipment IEC 61000-4-3	385 (27), 450 (28), 710 (9), 745 (9), 780 (9), 810 (28), 870 (28), 930 (28), 1720 (28), 1845 (28), 1970 (28), 2450 (28), 5240 (9), 5500 (9), 5785 (9) MHz (V/m)	Lumea is immune against any RF wireless communication device such as smart phone, Wi-Fi or Bluetooth device at a distance of 30 cm or more. Avoid a closer distance as this may result in improper operation.
RATED power frequency magnetic fields IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	
Proximity magnetic fields IEC 61000-4-39	30 kHz, CW, 8 A/m 134.2 kHz, Pulse modulation 2.1 kHz, 65 A/m 13.56 MHz, Pulse modulation 50 kHz, 7.5 A/m	
Electrical fast transient/burst IEC 61000-4-4	±2 kV 100 kHz repetition frequency	Mains power quality expected of a typical commercial or residential environment.
Surges IEC 61000-4-5	±0.5 kV, ±1 kV Line-to-line	Mains power quality expected of a typical commercial or residential environment.
Conducted disturbances induced by RF fields IEC 61000-4-6	3 V _{rms} , before modulation is applied 0.15 MHz – 80 MHz 6 V _{rms} , before modulation is applied in ISM and amateur radio bands between 0.15 MHz and 80 MHz 80 % AM at 1 kHz	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase at 0° 0 % UT; 250/300 cycle	Mains power quality expected of a typical commercial or residential environment. At power interruptions Lumea will switch off.

Note: UT is the AC mains voltage prior to application of the test level.


If in the unlikely event Lumea is locked inoperative in a function by disturbances from other devices beyond the levels stated in the table above, Lumea will shut-off completely and must be restarted by the power button. In that case make sure the distance to the other device causing this disturbance is enlarged.

Electromagnetic fields (EMF)

This device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Technical specifications

Electrical specification	The combination of the power adapter and Lumea is considered: Class II ME EQUIPMENT
Attachment	Type BF applied part The applied part (attachment) can become hot during use <60 °C. If this occurs, let the device cool down and limit the contact duration to less than 1 minute.
Mode of operation	Continuous operation
Intended operator	The user is the intended operator
IP rating power adapter	IP22
Light output	
Risk group classification	Exempt
Wavelength range of emitted optical radiation	Body and precision attachment: 520 nm – 1200 nm Face attachment: 590 nm – 1200 nm
Spot size	Body attachment: 4.1 cm ² (32.7 mm x 12.6 mm) Precision attachment: 3 cm ² (30.0 mm x 10.0 mm) Face attachment: 2 cm ² (25.0 mm x 8.0 mm)
Optical exposure	Setting 1: $2.5 \pm 0.80 \text{ J/cm}^2$ Setting 2: $3.1 \pm 0.85 \text{ J/cm}^2$ Setting 3: $3.8 \pm 0.91 \text{ J/cm}^2$ Setting 4: $4.4 \pm 0.96 \text{ J/cm}^2$ Setting 5: $5.0 \pm 1.00 \text{ J/cm}^2$
Maximum optical energy	Body attachment: 23.0 J Precision attachment: 16.8 J Face attachment: 11.2 J
Optical homogeneity	Max. +/-30 % deviation from average optical exposure in treatment area
Pulse sequence	Single pulse, for all settings
Pulse interval	1.0 – 2.4 s, for all settings
Pulse width	1.25 +/-0.4 ms (FWHM)
Transport and storage conditions	Lumea remains operational in NORMAL USE within its specification after transport or storage in the following environmental range: - a temperature range of -25 °C to 70 °C; and - a relative humidity up to 95 % RH, noncondensing; and - an atmospheric pressure range of 700 hPa to 1060 hPa. If the device has been stored in a hot or cold environment, place it in an environment with a temperature of 20 °C for 30 minutes to let it reach a temperature within the usage conditions (5 °C to 40 °C) before you use it.
Operating conditions	Lumea complies with its specifications when operated in NORMAL USE under the following environmental operating conditions: - a temperature range of 5 °C to 40 °C; - a relative humidity range of 15 % to 90 % RH, noncondensing; and - an atmospheric pressure range of 700 hPa to 1060 hPa.

Power adapter	
Power adapter input	100-240 V ~ 50/60 Hz 1.5 A
Power adapter output	19.5 V  3.33 A 65.0 W
Average active efficiency	≥88.0 %
Efficiency at low load (10 %)	≥79.0 %
No-load power consumption	<0.21 W
Device weight	Power adapter ~ 480 g Handpiece ~ 430 g
Device size	Power adapter ~ 10.8 x 6.5 x 6.8 cm Handpiece ~ 20.0 x 14.5 x 9.5 cm
Service life	5 years, based on Lumea treatment schedule with 14 full body treatments per year.
Materials	
Housing	PC (Polycarbonate), ABS (Acrylonitrile-butadiene-styrene), PA (Polyamide)
Cable	TPE (Thermoplastic elastomer)

Explanation of symbols

The safety signs and symbols are essential to ensure that you use this product safely and correctly and to protect you and others from injury. Below you find the meaning of the safety signs and symbols on the label and in the user manual.



IEC 60417-5333 Type BF applied part
This symbol identifies the type BF applied part complying with IEC 60601-1.



ISO 7000-2620 Humidity limitation
This symbol indicates the acceptable upper and lower limits of relative humidity for transport and storage.



ISO 7000-2621 Atmospheric pressure limitation
This symbol indicates the acceptable upper and lower limits of atmospheric pressure for transport and storage.



ISO 7000-0632 Temperature limit
This symbol indicates the maximum and minimum temperature limits at which the item shall be stored, transported or used.



ISO 7000-0626 Keep dry
This symbol indicates that the device shall be kept dry.

IP22

This symbol indicates the protection against ingress of solid foreign objects ≥ 12.5 mm and against harmful effects due to the ingress of dripping water when tilted at 15 degrees.



IEC 60417-5172 Class II equipment
This symbol indicates 'Class II Equipment'. The power adapter is double-insulated.



ISO 7000-3082 Manufacturer
This symbol indicates the manufacturer of the product.



IEC 60417-6049 Country of manufacture
This symbol indicates the country of manufacture and the manufacturing date. 'CC' is the two-letter country code.



ISO 7000-2492 Batch number
This symbol indicates the manufacturer's batch number.



ISO 7000-2493 Catalogue number
This symbol indicates the manufacturer's catalogue number.



IEC 60417-5032 Alternating current
This symbol indicates alternating current.



IEC 60417-5031 Direct current
This symbol indicates direct current.



ISO 7000-1641 Refer to consult instructions for use
This symbol indicates that operating instructions should be considered when operating the device.



ISO 7010-W017 Warning; Hot surface
This symbol warns about a hot surface. The background color is yellow, the color of the symbol is black.



IEC 61558-1 Rated maximum ambient temperature
This symbol indicates the maximum ambient temperature at which the transformer may be operated continuously under normal conditions of use.



IEC 60417-5957 For indoor use only
This symbol means that the device is designed for indoor use.



IEC 60417-5947 Safety isolating transformer, short-circuit proof
This symbol means that the transformer is designed to be short-circuit proof.



IEC 60417-6190 Power supply unit, switch mode
This symbol means that the transformer is considered a Switch Mode Power Supply (SMPS) unit.



IEC 60417-6343 Maximum altitude
This symbol indicates that the appliance is intended to be usable up to the maximum altitude 3000 m.



ISO 7000-3706 Single patient multiple use
This symbol indicates that the device may be used multiple times on a single user.



SenseIQ 2.0 technology
This symbol indicates that the device is equipped with SenseIQ 2.0 technology.



ISO/IEC 16022 Information technology – Automatic identification and data capture technologies – Data Matrix barcode symbology specification. The data matrix (QR-code) is used to capture information in a machine-readable form.



ISO 15223-1 Unique Device Identifier
Indicates the Unique Device Identifier that is used to identify a specific product.

**ISO 7010-W001 General warning sign**

This symbol means that the user should consult the instructions for use for important safety information such as contraindications and warnings. On the device, the background color is yellow, the color of the symbol is black.

**ISO 7000-0434A Caution**

This symbol means that the user should consult the instructions for use for important safety information on cautions.



A recycling symbol to identify the material of which an item is made, to facilitate recycling or other reprocessing. The symbol includes a number to identify the material code.

**IEC 60417-5009 Stand-by**

This symbol identifies the switch or switch position by means of which part of the equipment is switched on in order to bring it into the stand-by condition. The symbol also identifies the control to shift to or to indicate the state of low power consumption.



This symbol indicates Forest Stewardship Council. The FSC trademarks enable consumers to choose products that support forest conservation, offer social benefits, and enable the market to provide an incentive for better forest management.

**ISO 15223-1 Medical device**

This symbol indicates that this is a medical device.



If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.



ISO 15223-1 معرف الجهاز الفريد

يشير إلى معرف الجهاز الفريد المستخدم لتحديد منتج معين.



ISO 7010-W001 علامة تحذير عامة

يعني هذا الرمز أنه يجب على المستخدم الرجوع إلى تعليمات الاستخدام للحصول على معلومات أمان مهمة مثل موانع الاستخدام والتحذيرات. في الجهاز، تظهر الخلفية باللون الأصفر والرمز باللون الأسود.



ISO 7000-0434A تحذير

يعني هذا الرمز أنه يجب على المستخدم الرجوع إلى تعليمات الاستخدام للحصول على معلومات أمان مهمة حول التنبيهات.



رمز إعادة التدوير لتحديد المادة التي صنع منها الجهاز، لتسهيل إعادة التدوير أو إعادة المعالجة الأخرى. يحتوي هذا الرمز على رقم يشير إلى رمز المادة.



IEC 60417-5009 حالة الاستعداد

يحدد هذا الرمز موضع المفتاح أو المفتاح من خلال أي جزء من الجهاز يتم تشغيله لإدخاله في حالة الاستعداد. كما يحدد هذا الرمز عنصر التحكم للانتقال إليه أو للإشارة إلى حالة انخفاض استهلاك الطاقة.



يشير هذا الرمز إلى مجلس رعاية الغابات. تمكن العلامات التجارية لمجلس رعاية الغابات المستهلكين من اختيار المنتجات التي تدعم الحفاظ على الغابات، وتقدم المنافع الاجتماعية، وتمكن السوق من توفير حافز لإدارة الغابات بشكل أفضل.



ISO 15223-1 الجهاز الطبي

يشير هذا الرمز إلى أنه جهاز طبي.



إذا كنت بحاجة إلى معلومات أو دعم، فيرجى زيارة www.philips.com/support أو قراءة نشرة الضمان الدولي.

تفسير الرموز

تعد العلامات والرموز المتعلقة بالسلامة ضرورية لضمان استخدام هذا المنتج بأمان وبشكل صحيح ولحمايتك أنت والآخرين من الإصابة. تجددين أدناه معنى العلامات والرموز الخاصة بالسلامة على الملصق وفي دليل المستخدم.

IEC 60417-5031 التيار المباشر	
يشير هذا الرمز إلى التيار المباشر.	
ISO 7000-1641 راجعي تعليمات الاستخدام.	
يشير هذا الرمز إلى ضرورة الاطلاع على تعليمات التشغيل في أثناء تشغيل الجهاز.	
تحذير ISO 7010-W017 سطح ساخن	
يحذر هذا الرمز من الأسطح الساخنة، يكون لون الخلفية أصفر ولون الرمز أسود.	
IEC 61558-1 درجة الحرارة المحيطة القصوى المقدرة	ta
يشير هذا الرمز إلى أقصى درجة حرارة محيطة يمكن عندها تشغيل المهايئ بشكل مستمر في ظل ظروف الاستخدام العادية.	
IEC 60417-5957 للاستخدام الداخلي فقط	
يعني هذا الرمز أن الجهاز مصمم للاستخدام الداخلي.	
محول العزل الآمن IEC 60417-5947، مقاوم لحدوث لدائرة كهربائية قصيرة.	
هذا الرمز يعني أن المحول مصمم ليكون مضاد لحدوث قصر الدائرة.	
IEC 60417-6190 وحدة إمداد الطاقة، وضع التبديل	
يعني هذا الرمز أن المحول يعتبر وحدة تبديل وضع امدادات الطاقة (SMPS).	
IEC 60417-6343 أقصى ارتفاع	
يشير هذا الرمز إلى أن الجهاز مصمم للاستخدام حتى أقصى ارتفاع يبلغ 3000 م.	≤ 3000m
ISO 7000-3706 مريض واحد متعدد الاستخدامات	
يشير هذا الرمز إلى أنه يمكن استخدام الجهاز عدة مرات لمستخدم واحد.	
تقنية SenseIQ 2.0	
يشير هذا الرمز إلى أن الجهاز مزود بتقنية SenseIQ 2.0.	
ISO/IEC 16022 تقنية المعلومات - تقنيات التعرف التلقائي للبيانات وجمعها - مواصفات ترميز الباركود	
لمصفوفة البيانات، تستخدم مصفوفة البيانات (رمز الاستجابة السريعة) لجمع المعلومات في شكل يمكن قراءته آلياً.	

IEC 60417-5333 النوع BF الجزء المطبق	
يشير هذا الرمز إلى نوع الجزء المنطبق من النوع BF المتوافق مع IEC 60601-1.	
ISO 7000-2620 تحديد الرطوبة	
يحدد الرمز الحد الأدنى الأقصى والأدنى المقبولين للرطوبة النسبية للنقل والتخزين.	
ISO 7000-2621 تحديد الضغط الجوي	
يحدد الرمز الحد الأدنى الأقصى والأدنى المقبولين للضغط الجوي للنقل والتخزين.	
ISO 7000-632 حد درجة الحرارة	
يشير هذا الرمز إلى الحدود القصوى والدنيا لدرجة الحرارة التي يجب أن يتم تخزين المنتج، أو نقله، أو استخدامه عندها.	
ISO 7000-0626 الحفاظ على الجفاف	
يشير هذا الرمز إلى أن الجهاز يجب أن يظل جافاً.	
IP22	
يشير هذا الرمز إلى الحماية من دخول أجسام غريبة صلبة <12.5 مم وضد الآثار الضارة الناجمة عن دخول الماء المتقطر عند إمالاته إلى 15 درجة.	
IEC 60417-5172 معدات من الفئة الثانية	
يشير هذا الرمز إلى „معدات من الفئة II“، محول طاقة بعزل مزدوج.	
ISO 7000-3082 الشركة المصنعة	
يشير هذا الرمز إلى الشركة المصنعة للمنتج.	
IEC 60417-6049 بلد الصنع	
يشير هذا الرمز إلى بلد الصنع وتاريخ التصنيع. يشير حرفاً „CC“ إلى رمز البلد.	
ISO 7000-2492 رقم الدفعة	LOT
يشير هذا الرمز إلى رقم الدفعة للشركة المصنعة.	
ISO 7000-2493 رقم الكatalog	REF
يشير هذا الرمز إلى رقم كatalog الشركة المصنعة.	
IEC 60417-5032 التيار المتردد	
يشير هذا الرمز إلى التيار المتردد.	

مهايئ الطاقة	
مدخل محول الطاقة	من 100 إلى 240 فولت ~ 50/60 هرتز 1.5 أمبير
مخرج محول الطاقة	19.5 فولت  3.33 أمبير 65.0 واط
متوسط الكفاءة النشطة	≤ 88.0 %
الكفاءة عند التحميل المنخفض (10%)	≤ 79.0 %
استهلاك الطاقة بدون حمل	> 0.21 واط
وزن الجهاز	محول الطاقة ~ 480 جم المقبض ~ 430 جم
حجم الجهاز	محول الطاقة ~ 10.8 × 6.5 × 6.8 سم المقبض ~ 20.0 × 14.5 × 9.5 سم
مدة الخدمة	5 أعوام، وفقًا لجدول الجلسات الخاص بجهاز Lumea الذي يتضمن 14 جلسة إزالة الشعر من الجسم كاملاً في العام.
المواد	
الغلاف	PC (بولي كربونات)، ABS (أكريلونيتريل بوتادين ستايرين)، PA (بولي أميد)
الكابل	TPE (المطاط الصناعي الحراري)

المواصفات الفنية

المواصفات الكهربائية	يكون الجمع بين مجول الطاقة وجهاز Lumea على النحو التالي: معدات من الدرجة الثانية ME
المرفق	جزء مستخدم من النوع BF
طريقة العمل	قد يصبح الجزء المستخدم (الملحق) سائخا في أثناء الاستخدام بدرجة حرارة >60 درجة مئوية. إذا حدث ذلك، فاترك الجهاز ليبرد وقلل مدة الاتصال لأقل من دقيقة واحدة.
المشغل المقصود	العملية المستمرة
محول طاقة بتصنيف IP	المستخدم هو المشغل المقصود
ناتج الضوء	IP22
تصنيف مجموعة المخاطر	معفى
نطاق الطول الموجي للإشعاع الضوئي المنبعث	مرفق الجسم الدقة: 520 نانو متر - 1200 نانو متر ملحق الوجه: 590 نانو متر - 1200 نانو متر
حجم المنطقة	ملحق الجسم: 4.1 سم ² (32.7 × 12.6 ملم) ملحق الدقة: 3 سم ² (30.0 × 10.0 ملم) ملحق الوجه: 2 سم ² (25.0 × 8.0 ملم)
التعرض البصري	الإعداد 1: 2.5 ± 0.80 جول/سم ² الإعداد 2: 3.1 ± 0.85 جول/سم ² الإعداد 3: 3.8 ± 0.91 جول/سم ² الإعداد 4: 4.4 ± 0.96 جول/سم ² الإعداد 5: 5.0 ± 1.00 جول/سم ²
الطاقة الضوئية القصوى	ملحق الجسم: 23.0 جول ملحق الدقة: 16.8 جول ملحق الوجه: 11.2 جول
التجانس البصري	الحد الأقصى +/- 30% انحراف عن متوسط التعرض البصري في منطقة المعالجة
تتابع النبض	نبضة واحدة لجميع الإعدادات
المسافة بين نبضتين	1.0 - 2.4 ثانية، لجميع الإعدادات
عرض النبضة	1.25 +/- 0.4 ملتي ثانية (العرض التام عند منتصف الارتفاع)
شروط النقل والتخزين	يظل جهاز Lumea قيد التشغيل في الاستخدام العادي ضمن مواصفاته بعد النقل أو التخزين في النطاق البيئي التالي: - نطاق درجة حرارة من -25 درجة مئوية إلى 70 درجة مئوية؛ و - رطوبة نسبية تصل إلى 95 % من الرطوبة النسبية، غير المكثفة؛ و - ضغط جوي يتراوح من 700 هيكوباسكال إلى 1060 هيكوباسكال. إذا تم تخزين الجهاز في بيئة سائخة أو باردة، فأنقله إلى بيئة بدرجة حرارة 20 درجة مئوية لمدة 30 دقيقة حتى يصل إلى نطاق درجة حرارة الاستخدام (من 5 درجات مئوية إلى 40 درجة مئوية) قبل أن تستخدمه.
ظروف التشغيل	يتوافق جهاز Lumea مع المواصفات الخاصة به عند تشغيله في الاستخدام العادي في ظل ظروف التشغيل البيئية التالية: - درجة حرارة تتراوح من 5 درجة مئوية إلى 40 درجة مئوية؛ - رطوبة نسبية بين 15 % و 90 % من الرطوبة النسبية، غير المكثفة؛ و - ضغط جوي يتراوح من 700 هيكوباسكال إلى 1060 هيكوباسكال.

اختبار المناعة	مستوى اختبار المناعة	البيئة الكهرومغناطيسية - التوجيه
التفريغ الكهروستاتيكي (ESD) IEC 61000-4-2	± 8 كيلو فولت اتصال ± 15 كيلو فولت هواء	يمكن أن تكون الأرضيات مصدراً للكهرباء الساكنة، عندما يكون الجو جافاً جداً (>20% من الرطوبة النسبية)، يُنصح المستخدم باستخدام Lumea في غرفة ذات أرضية خرسانية أو خشبية (بدلاً من استخدامه في غرفة ذات أرضية صناعية).
حقول الترددات الراديوية المشعة IEC 61000-4-3	10 فولت/م 80 ميغا هرتز - 2.7 جيجا هرتز 80% AM عند 1 كيلو هرتز	± 8 كيلو فولت اتصال ± 15 كيلو فولت هواء
مجالات القرب من معدات الاتصالات اللاسلكية RF IEC 61000-4-3	385 (27)، 450 (28)، 710 (9)، 745 (9)، 780 (9)، 810 (28)، 870 (28)، 930 (28)، 1720 (28)، 1845 (28)، 1970 (28)، 2450 (28)، 5240 (9)، 5500 (9)، 5785 (9) ميغاهرتز (فولت/م)	جهاز Lumea محصن ضد أي جهاز اتصال لاسلكي RF مثل الهاتف الذكي أو Wi-Fi أو جهاز Bluetooth على مسافة 30 سم أو أكثر. تجنبي قرب المسافة لأن هذا قد يؤدي إلى تشغيل غير سليم.
المجالات المغناطيسية ذات تردد الطاقة المصنفة IEC 61000-4-8	30 أمبير/م 50 هرتز أو 60 هرتز	
المجالات المغناطيسية القريبة IEC 61000-4-39	30 كيلو هرتز، الموجة المستمرة، 8 أمبير/م 134.2 كيلو هرتز، تعديل نبضي 2.1 كيلوهرتز، 65 أمبير/م 13.56 ميغا هرتز، تعديل نبضي 50 كيلوهرتز، 7.5 أمبير/م	
عابر/انفجار كهربائي سريع IEC 61000-4-4	تردد تكرر ±2 كيلو فولت 100 كيلو هرتز	جودة الطاقة الرئيسية المتوقعة من بيئة تجارية أو سكنية نموذجية.
الطفرات IEC 61000-4-5	±0.5 كيلو فولت، ±1 كيلو فولت بين الخطيين	جودة الطاقة الرئيسية المتوقعة من بيئة تجارية أو سكنية نموذجية.
الاضطرابات الناجمة عن مجالات التردد الراديوي IEC 61000-4-6	3 فولت جذر متوسط التربيع، قبل تطبيق التشكيل 0.15 ميغا هرتز - 80 ميغا هرتز 6 فولت لجذر متوسط التربيع، ينطبق قبل التعديل على النطاق الترددي (ISM) ونطاقات راديو الهواة بين 0.15 ميغاهرتز و80 ميغاهرتز، و80% للنطاق الترددي (AM) عند 1 كيلو هرتز	
انخفاضات الجهد، الانقطاعات القصيرة وتغيرات الجهد على خطوط إدخال إمداد الطاقة IEC 61000-4-11	0% UT؛ 0.5 دورة عند 0 درجة مئوية، 45 درجة مئوية، 90 درجة مئوية، 135 درجة مئوية، 180 درجة مئوية، 225 درجة مئوية، 270 درجة مئوية و 315 درجة مئوية 0% UT؛ 1 دورة 90% UT؛ 25/30 دورات مرحلة واحدة عند 0 درجة مئوية 0% UT؛ 250/300 دورة	جودة الطاقة الرئيسية المتوقعة من بيئة تجارية أو سكنية نموذجية، عند انقطاع التيار الكهربائي، يتوقف جهاز Lumea.

ملاحظة: UT هو جهد التيار الكهربائي المتردد قبل تطبيق مستوى الاختبار. في بعض الحالات المستبعد حدوثها، إذا تعطلت إحدى وظائف جهاز Lumea بسبب تشوشات ناتجة عن أجهزة أخرى وبمستويات تزيد عما تم تحديده في الجدول أدناه، فسيغلق Lumea تمامًا وتجب إعادة تشغيله بالضغط على زر التشغيل. في هذه الحالة، تأكد من إبعاد الجهاز الآخر الذي يسبب هذا الاضطراب.

المجالات الكهرومغناطيسية (EMF)

يتوافق هذا الجهاز مع جميع المعايير واللوائح المعمول بها فيما يتعلق بالتعرض للحقول الكهرومغناطيسية.

المعلومات الإضافية

التخلص من الجهاز

اتبع القواعد المحلية المعمول بها في بلدك والمعنية بتجميع المنتجات الكهربائية والإلكترونية. يساعد التخلص الصحيح من المنتجات القديمة في منع حدوث عواقب سلبية على البيئة وصحة الإنسان.

التوافق الكهرومغناطيسي

تم اختبار هذا المنتج ووجد أنه يتوافق مع معيار IEC 60601-1-2 للتوافق الكهرومغناطيسي (EMC) ضمن الفئة B وفقاً لـ CISPR 11. وستجدين مستويات الاختبارات مدرجة في الجداول المصاحبة. يجب استخدام محول الطاقة المرفق مع Lumea فقط مع جهازك.

الإرشادات وإعلان الشركة المصنعة حول الانبعاثات الكهرومغناطيسية والحصانة منها

تم تصميم Lumea للاستخدام في البيئة الكهرومغناطيسية المدرجة هنا. يجب أن يتأكد مستخدم الجهاز من أنه يتم استخدامه في مثل هذه البيئة لمنع التشغيل غير السليم.

الانبعاثات	التوافق	البيئة الكهرومغناطيسية - التوجيه
انبعاثات التردد اللاسلكي CISPR 11	المجموعة 1	لا يستخدم جهاز Lumea طاقة التردد اللاسلكي للقيام بالوظيفة المخصص لها. لذلك، فإن انبعاثات التردد اللاسلكي الخاصة بها منخفضة جداً ومن غير المحتمل أن تسبب أي تدخل مع المعدات الإلكترونية القريبة.
انبعاثات التردد اللاسلكي CISPR 11	الفئة (ب)	Lumea مناسب للاستخدام في المواقع في البيئات السكنية وفي المؤسسات المتصلة مباشرة بشبكة إمداد الطاقة ذات الجهد المنخفض والتي تزود المباني المستخدمة للأغراض المنزلية.
انبعاثات توافقية IEC 61000-3-2	الفئة "أ"	
تقلبات الجهد/انبعاثات الوميض IEC 61000-3-3	يتوافق	

ماذا تفعلين إذا لم تلتزمي بجدول جلساتك؟

المرحلة	المشكلة	الحل
مرحلة الجلسات.	جلسة فائتة.	حاولي إجراء الجلسة في موعد قريب من موعد الجلسة المخطط لديك. يمكنك الاستمرار في اتباع الجدول الزمني الأصلي المخطط له. يمكنك الحصول على جلسة إضافية كل أسبوعين إلى جانب جلستك الأصلية. إذا لم تلتزمي بجدول الجلسات، فقد لا تصلي إلى النتائج المرجوة.
	لقد أجريت جلستي بعد أكثر من 3 أيام من التاريخ المخطط له في الأصل.	
	آخر جلسة من الجلسات التي تجرى كل أسبوعين.	يمكنك اختيار إجراء جلسة إضافية كل أسبوعين إلى جانب فترة الجلسات الخاصة بك، ولكن تعتمد النتائج المطلوبة على جدول الجلسات المشار إليه، وتختلف النتائج من شخص لآخر.
	أنهيتُ مرحلة الجلسات، لكنني أشعر أنني لم أحصل على النتائج المطلوبة.	
مرحلة اللمسات الأخيرة	خارج مدة الجلسات.	يُرجى اتباع الجدول الصحيح. إذا لم تلتزمي بجدول الجلسات، فقد لا تصلي إلى النتائج المرجوة. يمكنك الاستمرار في اتباع الجدول الأصلي المخطط لك أو اختيار إعادة مرحلة الجلسات، ولكن تعتمد النتائج المطلوبة على جدول العلاج المشار إليه، وتختلف النتائج من شخص لآخر.
	لم أجر أي جلسة خلال 4 أيام من تاريخ جلستي المجدولة.	

استكشاف الأخطاء وإصلاحها

يلخص هذا الفصل المشكلات الأكثر شيوعًا التي قد تواجهينها مع جهاز Lumea الخاص بك أو إذا لم تلتزمي بجدول جلساتك. إذا لم تتمكني من حل المشكلة بالمعلومات الواردة أعلاه، ففضلي بزيارة www.philips.com/support للحصول على قائمة بالأسئلة الشائعة أو اتصلي بمركز رعاية المستهلك في دولتك كما هو موضح في نشرة الضمان الدولي المضمنة.

المشكلات الشائعة عند استخدامك لجهاز Lumea:

المشكلة

الحل

1. قد تكون بشرتك داكنة جدًا. تحقق من جدول درجات لون البشرة/لون شعر الجسم على الغلاف الداخلي لهذا الكتيب.
2. تأكدي من توصيل المرفق (بشكل صحيح) بالجزء الأمامي من جهاز Lumea.
3. تأكدي من ملائمة الجهاز لمنطقة البشرة بالكامل لمدة 5 ثوانٍ على الأقل. في بعض المناطق (المناطق الحساسة والعظمية) قد يكون من الصعب التلامس الجيد مع الجلد. يجب ألا يدخل أي ضوء بين جهاز Lumea وبشرتك. حاولي تغيير موضع جهاز Lumea.
4. تأكدي من الضغط على زر الوميض الموجود على المقبض الداخلي لاستخدام الجهاز على منطقة من جلدك.

لا يوميض جهاز Lumea / يوميض مصباح "الاستعداد للوميض" باللون البرتقالي، ثم يوميض ويفلق بشكل متكرر.

1. تأكدي من إزالة كل الشعر في المنطقة التي تريد استخدام جهاز Lumea عليها قبل استخدامه. يمكن أن يسبب حرق الشعر ظهور الرواح.
2. تأكدي من أن بشرتك نظيفة وجافة قبل بدء المعالجة.
3. تأكدي من نظافة جهاز Lumea وملحقاته وجفافه قبل بدء المعالجة. قد تتراكم بقايا الشعر، أو العرق، أو الأوساخ من عملية الوميض وقد تسبب تحسسات جلدية إذا تركت دون تنظيف.

يصدر جهاز Lumea رائحة غريبة في أثناء الوميض.

1. يعمل جهاز Lumea في وضع خطأ. اضغطي على زر التشغيل لأكثر من ثانيتين على الزر، وبعدها يعود جهاز Lumea إلى وضع الاستعداد. أو افصلي جهاز Lumea عن الكهرباء، وسيغلق.
2. انتظري 5 دقائق. سيعود جهاز Lumea تلقائيًا إلى وضع الاستعداد بعد 5 دقائق دون مساعدة.

جميع مؤشرات شدة الضوء الموجودة أعلى جهاز Lumea الخاص بي تضيء وتطفئ، كما يضيء مصباح "الاستعداد للوميض" باللون البرتقالي، ثم يضيء ويطفئ.

التخطيط لجلساتك التالية

امسحي رمز الاستجابة السريعة ضوئيًا لمشاهدة مقطع فيديو قصير.



اتبعي جدول المعالجات المحددة.

مرحلة الاختبار

قبل أول جلسة لك لكل منطقة في الجسم أو بعد تعرض بشرتك لأشعة الشمس دون وقاية، اختبري بشرتك من خلال اتباع الخطوات في قسم „اختبار بشرتك“. انتظري 24 ساعة بعد اختبار البشرة قبل البدء في إجراء جلستك.



24 ساعة

مرحلة الجلسات: 4 جلسات كل أسبوعين

خلال فترة الجلسات التي تتضمن 4 جلسات، احرصي على إجراء جلسة كل أسبوعين.

احرصي على إجراء جلسة في غضون 3 أيام قبل أو بعد تاريخ الجلسة المخطط لها خلال فترة الجلسات.



مرحلة اللمسات الأخيرة: مرة كل 4 أسابيع

بعد فترة الجلسات، استخدمِي جهاز Lumea كل 4 أسابيع للحفاظ على النتائج.


احرصي على إجراء جلسة في غضون 4 أيام قبل أو بعد تاريخ الجلسة المخطط لها خلال فترة تحسين النتائج.



ملاحظة: إننا فوّتَ جلسة، فارجعي إلى فصل „استكشاف الأخطاء وإصلاحها“ ليرشدك. إذا لم تلتزمي بجدول الجلسات، فقد لا تصلي إلى النتائج المرجوة.


5. نظفي نافذة خروج الضوء

تأكد من عدم ترك أي أوساخ، أو بقايا، أو مياه على سطح نافذة خروج الضوء.

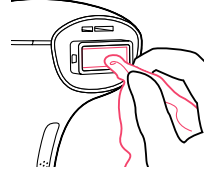


6. الضغط على الملحق

تأكد من حماية نافذة خروج الضوء لجهاز Lumea من التلف من خلال الضغط على أحد الملحقات.

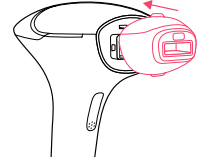

7. التخزين

خزني جميع الأجزاء في العبوة الأصلية أو في الحقيبة، أو المحفظة، أو الصندوق المرفق في مكان جاف وخالٍ من الغبار.

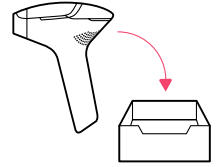

بعد العلاج

بعد الاستخدام، إذا كنت ترغبين في وضع المستحضرات أو الكريمات أو مزيلات العرق أو المرطبات أو مستحضرات التجميل على المناطق المعالجة، فتتحقي دائماً من نشرة التعليمات الخاصة بالمنتج الذي تريد تطبيقه للتأكد من أنه مناسب لحالة بشرتك في تلك اللحظة.


التنظيف والتخزين

امسحي رمز الاستجابة السريعة ضوئيًا لمشاهدة مقطع فيديو قصير.




1. إيقاف التشغيل، وفصل الكهرباء، وتبريد جهاز Lumea الخاص بك

⚠ لا تلمسي نافذة خروج الضوء خلال الاستخدام أو بعده مباشرةً، إذ يصبح هذا الجزء ساخنًا جدًا. بعد استخدام جهاز Lumea، افصليه عن الطاقة وافصلي كذلك محول الطاقة واتركيهما ليبردًا لمدة 15 دقيقة قبل تنظيفهما وتخزينهما.

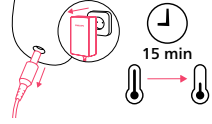

2. اغطي على الملحق


قومي بيسحيه برفق.

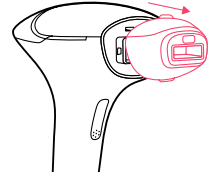

3. بللي قماشة التنظيف المرفقة مع جهاز Lumea الخاص بك قليلًا


⚠ لا تستخدمى جهاز Lumea في الأماكن الرطبة.
لا تدعي الماء يدخل إلى جهاز Lumea.

بللي قماشة التنظيف المرفقة مع جهاز Lumea الناعمة ببضع قطرات من الماء.

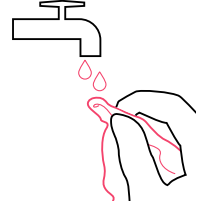

4. تنظيف الملحقات

نظفي كل جزء من المرفق الذي استخدمته، بما في ذلك العاكسات المعدنية.

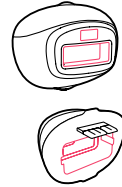

8. الوميض ثم التحريك

⚠️ **لا** تستمري في الجلسة إذا شعرتي بألم. قللي إعداد شدة الضوء أو توقفي عن استخدام جهاز Lumea.




⚠️ **لا** تعرّضي بشرتك لأشعة الشمس دون حماية عند استخدام جهاز Lumea. خلال فترة العلاج بالضوء النبضي المكثف بالكامل، وخاصة خلال الـ 48 ساعة الأولى بعد العلاج، قومي بحماية بشرتك من أشعة الشمس بالملابس أو ضعي منتجات واقية من الشمس لا تقل عن SPF50 على مناطق العلاج غير المغطاة.

اضغطي على زر الوميض. انتقلي إلى منطقة مجاورة وأومضي مرة أخرى.

9. معالجة منطقة أخرى من الجسم


⚠️ **لا** تبدئي باستخدام جهاز Lumea أو تعاودي استخدامه في غضون 4 أسابيع من تعرّض بشرتك لأشعة الشمس دون حماية (مثل الشمس الطبيعية أو الحمام الشمسي) أو في غضون أسبوعين من تعرّضك للاستمرار باستخدام وسائل صناعية (الكريمات، أو البخاخات، أو الحبوب). احرصي على إجراء اختبار البشرة قبل بدء الجلسة أو استئناؤها، اطلعي على قسم „اختبار بشرتك“.

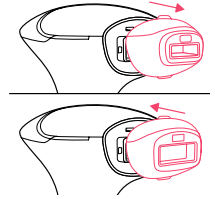

إذا كنت ترغبين في استخدامه على منطقة أخرى، فتأكدي من اختيار الملحق الصحيح واليد من جديد بالخطوة 1. وفي حال استخدامه على منطقة ما لأول مرة أو بعد التعرض لأشعة الشمس دون حماية، فتأكدي من إجراء اختبار البشرة. انتظري 24 ساعة وافحصي جلدك لمراقبة أي تحسسات.

10. الإنتهاء إيقاف التشغيل وفصل الكهرباء


⚠️ **لا** تلمسي نافذة خروج الضوء خلال الاستخدام أو بعده مباشرة، إذ يصبح هذا الجزء ساخناً جداً. بعد استخدام جهاز Lumea، افصليه عن الطاقة وافصلي كذلك محول الطاقة واتركيهما ليبردا لمدة 15 دقيقة قبل تنظيفهما وتخزينهما.

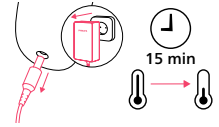

استمري في الضغط على زر التشغيل لمدة ثانيتين. افصلي جهاز Lumea ومحول الطاقة لديك من مقبس الحائط. اتركي جهاز Lumea حتى يبرد.


11. نظفي جميع أجزاء جهاز Lumea لديك بعد الاستخدام

من المهم تنظيف جميع أجزاء جهاز Lumea الخاص بك. قد تتراكم بقايا الشعر، أو العرق، أو الأوساخ من عملية الوميض وقد تسبب تحسسات جلدية إذا تراكمت دون تنظيف. اتبعي الخطوات الواردة في قسم „التنظيف والتخزين“.

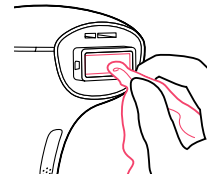

12. تحديد موعد جلستك التالية

اتبعي الخطوات الواردة في قسم „التخطيط لجلساتك التالية“.

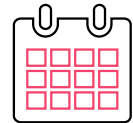

تنفيذ المعالجة


امسحي رمز الاستجابة السريعة ضوئياً لمشاهدة مقطع فيديو قصير.



اختبري تحسس بشرتك قبل استخدامه على منطقة جديدة في جسمك أو بعد تعرض بشرتك لأشعة الشمس دون حماية.


1. اتبعي الخطوات الواردة في قسم „الاستعداد للجلسة“.

2. افحصي جهاز Lumea بحثاً عن أي تلف أو أوساخ

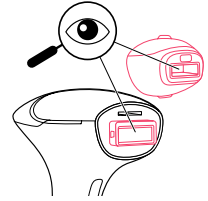

3. ركبي الملحق الصحيح

انقري فوق المرفق الخاص بالمنطقة التي تريد علاجها.

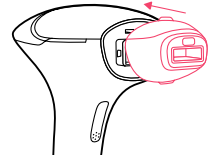

4. الشعور بالراحة

اشعري بالراحة قبل أن تبديني المعالجة.

يوصى باستخدام جهاز Lumea بعيداً عن أشعة الشمس المباشرة. يمكن أن يؤثر ضوء الشمس على مستشعر SmartSkin.

5. التوصيل بالتيار الكهربائي. اضغطي على زر التشغيل.

ستضيء مؤشرات شدة الضوء الخمسة الموجودة أعلى زر التشغيل بشكل مستمر من اليسار إلى اليمين، وستتكرر حتى تضعي جهاز Lumea على بشرتك.



6. تحديد الإعداد

اضغطي على زر التشغيل للتبديل بين الإعدادات وتحديد أكثر إعداد مريح لك لم ينتج عنه أي تحسسات جديدة.



7. بدء المعالجة

ضعي جهاز Lumea على بشرتك. سوف يضيء ضوء الاستعداد للوميض باللون الأبيض.



7. احصلي على إعداد شدة الضوء المناسب لكِ

سيتحول أحد مؤشرات شدة الضوء إلى اللون الأبيض للإشارة إلى الإعداد المقترح لكي. يتم تحديد الإعداد تلقائيًا بعد 5 ثوانٍ.



ملاحظة: إذا كان مصباح „الاستعداد للوميض“ يضيء باللون البرتقالي، فقد تكون بشرتك داكنة جدًا بحيث لا يمكن استخدام الجهاز عليها؛ فيُرجى الرجوع إلى فصل „استكشاف الأخطاء وإصلاحها“.

8. ابحتي عن الإعداد المريح لكِ

ضعي جهاز Lumea مرة أخرى على بشرتك. سوف يضيء ضوء الاستعداد للوميض باللون الأبيض.




اضغطي على زر الوميض مرة واحدة. إذا كان الوميض مريحًا، فيمكنك زيادة الإعداد. إذا شعرتي بألم، فقللي مستوى الإعداد. قومي بالتبديل بين الإعدادات عن طريق الضغط على زر الطاقة. غَيِّري الإعداد، واستخدمي جهاز Lumea على المنطقة التالية، ثم شغلي الوميض مرة أخرى. كرري ذلك حتى تصلي إلى أكثر إعداد مريح لكِ.



9. إيقاف التشغيل وفصل الكهرياء


⚠ لا تلمسي نافذة خروج الضوء خلال الاستخدام أو بعده مباشرة، إذ يصبح هذا الجزء ساخنًا جدًا. بعد استخدام جهاز Lumea، افصليه عن الطاقة وافصلي كذلك محول الطاقة واتركيهما ليبردا لمدة 15 دقيقة قبل تنظيفهما وتخزينهما.

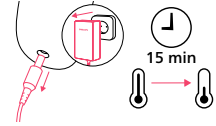

استمري في الضغط على زر التشغيل لمدة ثانيتين. افصلي جهاز Lumea ومحول الطاقة لديك من مقبس الحائط. اتركي جهاز Lumea حتى يبرد.

10. نظفي جميع أجزاء جهاز Lumea لديك بعد الاستخدام

من المهم تنظيف جميع أجزاء جهاز Lumea الخاص بك. قد تتراكم بقايا الشعر، أو العرق، أو الأوساخ من عملية الوميض وقد تسبب تحسسات جلدية إذا تركت دون تنظيف. اتبعي الخطوات الواردة في قسم „التنظيف والتخزين“.

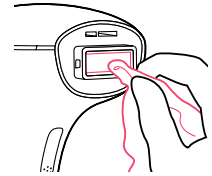

11. انتظري 24 ساعة وافحصي جلدك لمراقبة ظهور أي تحسسات

إذا لم تكن هناك تفاعلات جلدية بعد 24 ساعة، ابدأ العلاج بأعلى وضع مريح. إذا أدى الإعداد الأكثر راحة لك إلى تحسس الجلد، فاختاري إعدادًا أقل لم ينتج عنه أي تحسسات جلدية. إذا لم تختبري إعدادًا أقل، فكري اختبار الجلد.



اختبار جلدك

امسحي رمز الاستجابة السريعة ضوئياً لمشاهدة مقطع فيديو قصير.




⚠ لا تبدئي باستخدام جهاز Lumea أو تعاودي استخدامه في غضون 4 أسابيع من تعرّض بشرتك لأشعة الشمس دون حماية (مثل الشمس الطبيعية أو الحمام الشمسي) أو في غضون أسبوعين من تعرّضك للاستمرار باستخدام وسائل صناعية (الكريمات، أو البخاخات، أو الحبوب).

قبل أول جلسة لك لكل منطقة في الجسم أو بعد تعريض بشرتك لأشعة الشمس دون وقاية، اختبري تحسس جلدك، وابحثي عن الإعداد المناسب لبدء جلستك.

يوصى باستخدام جهاز Lumea بعيداً عن أشعة الشمس المباشرة. يمكن أن يؤثر ضوء الشمس على مستشعر SmartSkin.


1. اتبعي الخطوات الواردة في قسم «الاستعداد للجلسة».

2. افحصي جهاز Lumea بحثاً عن أي تلف أو أوساخ

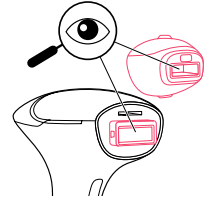


3. ركبي الملحق الصحيح

انقري فوق المرفق الخاص بالمنطقة التي تريد علاجها.

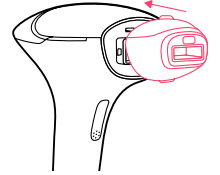

4. التوصيل بالتيار الكهربائي

ضعي مأخذ جهاز Lumea في مقبس الحائط.

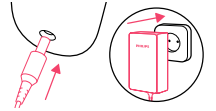


5. اضغطي على زر التشغيل.

ستضيء مؤشرات شدة الضوء الخمسة الموجودة أعلى زر التشغيل بشكل مستمر من اليسار إلى اليمين، وستتكرر حتى تضعي جهاز Lumea على بشرتك.



6. ضعِي جهاز Lumea على بشرتك.

يجب أن يتم وضع الجهاز مباشرة على بشرتك.

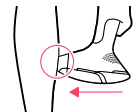


استعدي لجلستك

امسحي رمز الاستجابة السريعة ضوئيًا لمشاهدة مقطع فيديو قصير.




1. إزالة الشعر

استخدمي الشمع أو الحلاقة أو إزالة الشعر قبل استخدام جهاز Lumea. قومي بالجلسة قبل ظهور الشعر الجديد على بشرتك. إذا اخترت إزالة الشعر بالشمع، فانتظري 24 ساعة قبل استخدام Lumea.

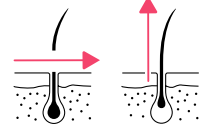

2. الجلد نظيف وجاف


يجب عدم ترك ماء أو زيت أو كريم أو حليب على جلدك.

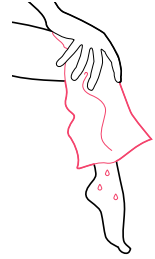

3. اختاري الملحق المناسب وافحصي جهاز Lumea بحثًا عن أي تلف أو أوساخ

⚠️ افحصي جهاز Lumea دائمًا قبل استخدامه. لا تستخدمي جهاز Lumea أو محول الطاقة في حال تلفهما. استبدلي دائمًا الجزء التالف بأحد الأنواع الأصلية.

⚠️ لا تستخدمي Lumea إذا كان هناك أي غسول أو جل على بشرتك أو على جهاز Lumea. استخدميه فقط على البشرة النظيفة والجافة. تأكدي من أن جهاز Lumea نظيف وجاف أيضًا.

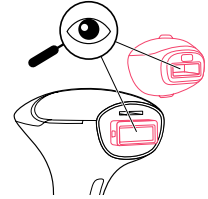

دليل جلستك

أقرني تعليمات السلامة المهمة في دليل المستخدم للتحقق مما إذا كانت جلسة جهاز Lumea مناسبة لك أو لا؛ حتى تتمكني من الحصول على الاستفادة الكاملة من جهازك.

قبل جلستك

1. استعدي لجلستك
قومي بإزالة الشعر، وتأكدي من أن بشرتك نظيفة وجافة.

2. اختار جلدك
ابحثي عن الإعداد المريح لك

3. انتظري 24 ساعة لمراقبة التحسسات الجلدية
24 ساعة

الجلسة

4. تنفيذ المعالجة

5. التنظيف والتخزين

6. التخطيط لجلساتك التالية



امسحي رمز الاستجابة السريعة ضوئيًا وشاهدي مقاطع الفيديو الإرشادية الخاصة بنا للحصول على الإرشادات خطوة بخطوة حول كيفية إجراء جلسة باستخدام جهاز Lumea، بالإضافة إلى مقاطع فيديو عامة حول جهاز Lumea.



طرق المعالجة


امسحي رمز الاستجابة السريعة ضوئيًا لمشاهدة مقطع فيديو قصير.



يحتوي جهاز Lumea على نمطين للجلسات للحصول على استخدام أكثر راحة لمختلف مناطق الجسم.


التحريك والوميض

للمناطق الأكبر (الساقان)، استمري في الضغط على زر الوميض في أثناء تحريك جهاز Lumea على بشرتك. يستمر Lumea في إطلاق الومضات حتى تتركي الزر أو في حال عدم ملائمة الجهاز للبشرة بشكل جيد.

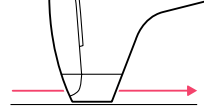

التثبيت والوميض

بالنسبة لمناطق الجسم الصغيرة (الإنبت)، ضعي جهاز Lumea على بشرتك واضغطي على زر الوميض. ارفعي جهاز Lumea عن بشرتك وضعيه بجانب المنطقة التي سبق واستخدمته عليها. اضغطي على زر الوميض مرة أخرى لتحرير ووميض آخر.

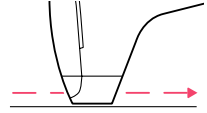

امسحي رمز الاستجابة السريعة ضوئياً لمشاهدة مقطع فيديو قصير.



يأتي جهاز Lumea مزوداً بملحقات مختلفة يمكنك استخدامها على مناطق الجسم المختلفة. حدي الملحق المخصص للمنطقة التي تريد استخدام الجهاز عليها.

نوع الملحق		BRI984	BRI983	BRI938	BRI937	BRI933	BRI932	BRI931	BRI930
ملحق الجسم - أكبر حجم للنافذة، - تصميم منحني إلى الداخل لتغطية المناطق الواقعة أسفل خط العنق ومعالجتها بشكل فعال.		✓	✓	✓	✓	✓	✓	✓	✓
ملحق الوجه - أصغر حجم نافذة لمعالجة دقيقة لأسفل عظام الوجنتين، - تصميم صغير مسطح، - فلتز مدمج إضافي لاستخدام الجهاز على الوجه أسفل عظام الوجنتين. ⚠ لا تستخدم جهاز Lumea على الوجه والرقبة. الشعر في هذه المناطق له خصائص مختلفة. نصيحة: ضعي لسانك بين شفتك العليا وأسنانك أو انفخي الخدين لتسهيل المعالجة.		✓	✓	✓	✓	✓	✓		
ملحق الدقة - حجم متوسط للنافذة؛ - تصميم منحني خارجي للاستخدام الدقيق والفعال على مناطق البكيني وتحت الإبط.				✓	✓				
ملحق الدقة لمنطقة البكيني: - حجم متوسط للنافذة؛ - تصميم منحني خارجي للاستخدام الدقيق والفعال على منطقة البكيني.		✓	✓						
ملحق الدقة للإبطين: - حجم متوسط للنافذة؛ - تصميم منحني خارجي للاستخدام الدقيق والفعال على منطقة الإبطين.		✓	✓						

مستشعر SmartSkin

يحدد مستشعر SmartSkin درجة لون بشرتك وسيساعدك على تحديد إعداد مريح لك عن طريق إضاءة أحد مؤشرات شدة الضوء أعلى جهازك.

كما أنه يمنعك تلقائياً من معالجة مناطق بشرتك الداكنة جداً للاستخدام الآمن.



ضوء الاستعداد للوميض

سيعلمك ضوء الاستعداد اللوميض الأبيض عندما يكون لديك تلامس جيد مع الجلد ويكون جهاز Lumea جاهزاً للوميض.

إنما كان مصباح «الاستعداد للوميض» برتقالياً، فإن بشرتك في هذه المنطقة من الجسم داكنة جداً بحيث لا يمكن استخدام الجهاز عليها، عالج منطقة ذات لون أفتح في جسمك.



تقنية SenseIQ 2.0

توفر تقنية SenseIQ 2.0 لك جلسات مخصصة لإزالة الشعر، إذ يجمع هذا المنتج المميز من علامتنا التجارية بين:

- أجهزة استشعار ذكية متطورة تحدد درجة لون بشرتك وتدعمك في اختيار إعداد للحصول على جلسة استخدام مريحة.
- ملحقات مميزة مصممة خصيصاً لتناسب كل منحني في جسمك.
- تطبيق Philips Lumea IPL لمساعدتك في القيام بروتينك المتعلق بالجلسات.




SenseIQ 2.0
technology



امسحي رمز الاستجابة
السريعة صوتيًا
لمشاهدة مقطع فيديو
قصير.

نظرة عامة عن جهاز Lumea

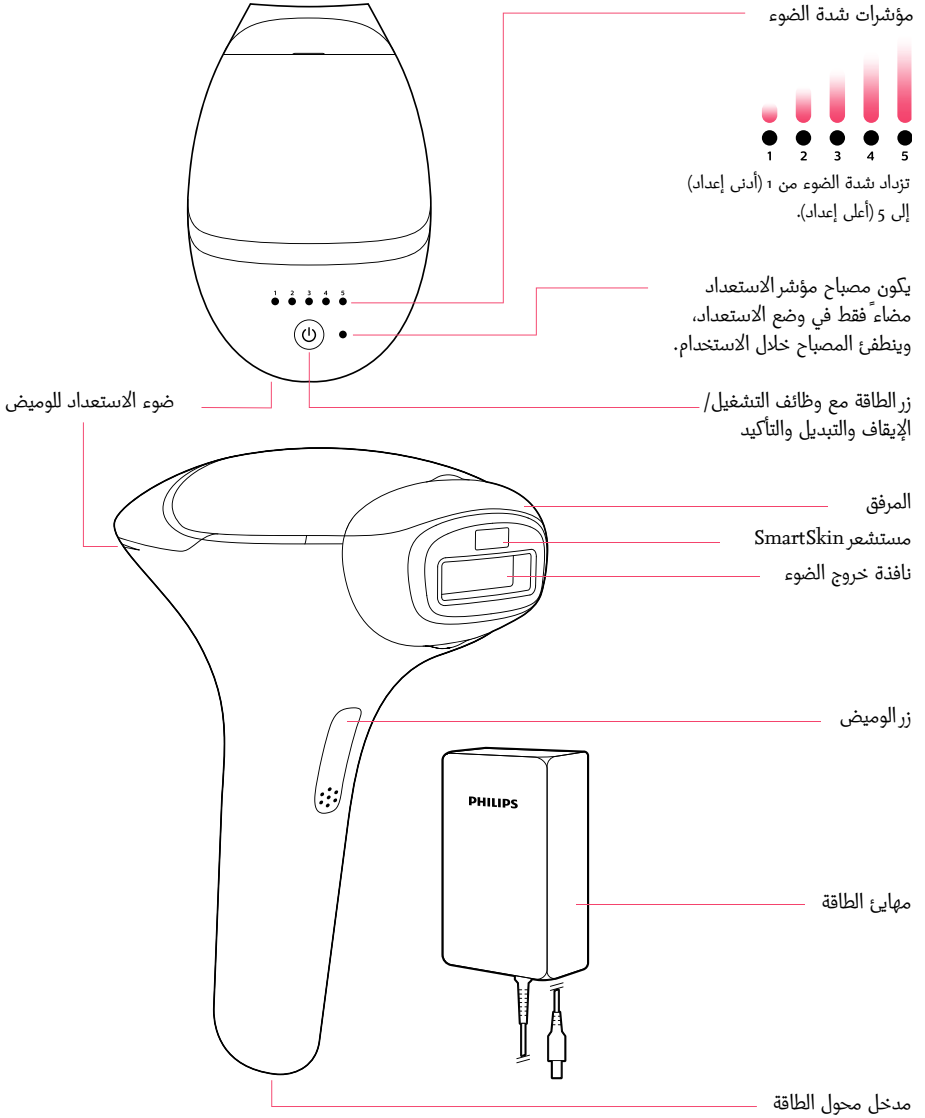


غير ظاهر:
- قماش التنظيف
- حقيبة، أو محفظة مميزة، أو صندوق فاخر لمستحضرات التجميل.

ماذا يُقصد بتقنية نبضات الضوء المكثف؟

امسحي رمز الاستجابة السريعة ضوئياً لمشاهدة مقطع فيديو قصير.



1. **IPL تعني الضوء النبضي المكثف.**
تمتص جذور شعرك نبضات خفيفة من الضوء.

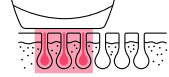

2. **تتسبب هذه النبضات من الضوء في دخول بصيالات الشعر إلى مرحلة الراحة.**
يتساقط الشعر بشكل طبيعي، ويقل نموه لبضعة أشهر.



3. **لا ينمو الشعر بشكل مستمر؛ ينمو في دورات.**
لذلك لا ينمو كل الشعر في منطقة واحدة في الوقت نفسه. نظراً لأن الضوء النبضي المكثف فعال فقط على الشعر الذي ينمو بنشاط، فمن المهم تكرار المعالجة على النحو المحدد.



أعراض جانبية غير شائعة	ما الذي تمت ملاحظته؟	ما مدة استمراره وكيف يمكن علاجه؟
التندب	قد تظهر ندبة كآثر ثانوي للحروق أو التهاب الجلد.	من المهم استشارة طبيبكِ في حال الإصابة بحرق أو التهاب الجريبات، وذلك لمنع ظهور ندبة. لا تستخدمى الجهاز على الأنسجة التي توجد بها ندبات.
الإصابة بمرضى الهريس	تعاني المستخدمات المصابات بعدوى فيروس الهريس البسيط من البثور الصغيرة أو الجروح المفتوحة، وعادةً ما تكون بالقرب من الشفاه ومنطقة الأعضاء التناسلية.	عادةً ما تستمر الإصابة لمدة تتراوح من عدة أيام إلى أسابيع ويمكن أن تقل بتناول الأدوية المضادة للفيروسات التي يصفها طبيبكِ.
الصداع والبقع الضوئية	قد يعاني الأشخاص الحساسون من الصداع أو رؤية بقع ضوئية بعد تعرض العين لوميض الضوء.	عادةً ما تستمر هذه الآثار لمدة تتراوح من بضع دقائق إلى عدة ساعات. لتجنب هذه الآثار، احرصى على استخدام جهاز Lumea في غرفة جيدة الإضاءة وعدم ترك أي مسافة بين بشرتك والجهاز قبل إطلاق وميض الضوء.

قد تؤدي إزالة الشعر بالضوء النبضي المكثف إلى زيادة نمو الشعر لدى بعض الأفراد. تتضمن مناطق الجسم المعرضة للخطر خط الفك، والرقبة، والمناطق الموجودة فوق عظام الوجنتين عند النساء، وأعلى الظهر، والكتفين عند الرجال. استنادًا إلى البيانات الحالية، فإن الفئات الأكثر خطرًا للتعرض لهذه الآثار هي الإناث من منطقة البحر الأبيض المتوسط والشرق الأوسط وجنوب آسيا.

أعراض جانبية غير شائعة	ما الذي تمت ملاحظته؟	ما مدة استمراره وكيف يمكن علاجه؟
تورم (وذمة)	يظهر التورم حول بصيالات الشعر على شكل نتوءات في المنطقة التي يُستخدم عليها الجهاز، يصاحبها احمرار أحياناً.	عادةً ما يقل هذا الأثر خلال ساعات، ولكن قد يستمر لمدة تصل إلى أسبوع واحد. يمكنك تبريد المنطقة باستخدام كمادات باردة.
جفاف الجلد	يبدو الجلد خشناً ومتقشرًا في منطقة المعالجة.	يستمر هذا الأثر عادةً ما بين يوم واحد وشهر واحد. يمكنك وضع مرطب غير معطر على المنطقة المعالجة.
الحروق والبثور	احمرار شديد بالجلد، وتورم، وربما ظهور بثور إلى جانب الشعور بالألم.	قد تستمر هذه التحسسات من يوم واحد حتى شهر واحد، وقد يتطلب الأمر الحصول على علاج طبي أو دواء. ننصحك باستشارة طبيب. انتظري حتى موعد الجلسة التالية حتى يشفى الجلد تمامًا. تأكدي من استخدام شدة ضوء أقل في الجلسات التالية.
قشرة سطحية	منطقة بنية وقشرية محددة بطريقة لافتة للنظر قد تظهر على ألوان البشرة الداكنة	
التهاب الجريبات	تورم واحمرار حول بصيالات الشعر، مصحوبًا بحكة و/أو ظهور بثور	قد يتطلب التهاب الجريبات تلقي علاج طبي أو دواء. ننصحك باستشارة طبيب. قد يستمر التهاب الجريبات لمدة تتراوح من عدة أيام وحتى عدة أسابيع. انتظري حتى موعد الجلسة التالية حتى يشفى الجلد تمامًا. تأكدي من استخدامك لشفرة حلقة نظيفة عند الحلاقة قبل الجلسة.
تغيرات صبغية/تلون البشرة	يظهر جلد المنطقة التي يُستخدم الجهاز عليها أغمق أو أفتح من الجلد المحيط بها.	عادةً ما يستمر هذا الأثر لعدة أسابيع ولمدة تصل إلى عام في بعض الحالات النادرة. لا تعالجي المناطق ذات التغيرات الصبغية/تغيرات لون البشرة حتى تستعيد بشرتك لونها الطبيعي. تأكدي من استخدام شدة ضوء أقل في الجلسات التالية.
عدوى الجلد	تظهر في شكل جلد ملتهب. يُعد هذا أثاراً جانبياً قد يحدث بعد حرق الجلد أو إصابات أخرى.	يتطلب التهاب الجلد تلقي علاج طبي أو دواء. ننصحك باستشارة طبيب، لا تستخدمى الجهاز على المناطق الملتهبة من الجلد.

الآثار الجانبية المحتملة للعلاج بتقنية الضوء النبضي المكثف (IPL)

يعرض الجدول أدناه الآثار الجانبية المعروفة التي قد تحدث بعد العلاج بأجهزة تقنية الضوء النبضي المكثف.

الآثار الجانبية الشائعة	ما الذي تمت ملاحظته؟	ما مدة استمراره وكيف يمكن علاجه؟
الإحساس بالألم/عدم الراحة	الشعور بألم موضعي أو عدم الراحة في المنطقة التي استخدم الجهاز عليها الجلسة أو بعدها	إذا شعرتي بالألم أو عدم الراحة خلال الجلسة، فاخفضي إعداد شدة الضوء أو توقف عن استخدام جهاز Lumea. قد يستمر الشعور بالألم أو عدم الراحة لمدة تصل إلى أسبوع واحد بعد الجلسة.
الإحساس بالحكة	الشعور بالحكة في المنطقة التي استخدم الجهاز عليها	عادة ما تُعالج هذه التحسسات خلال ساعات ولكن قد تستمر لمدة تصل إلى أسبوع واحد. يمكنك تبريد المنطقة باستخدام كمادات باردة.
الإحساس بالحرارة/الحرقان	الشعور بالسخونة الموضعية أو الحرقان في المنطقة التي استخدم الجهاز عليها	عادة ما يقل هذا الأثر خلال ساعات ولكن قد يستمر لمدة تصل إلى يوم واحد. يمكنك تبريد المنطقة باستخدام كمادات باردة.
احمرار الجلد (طفح وردي)	احمرار دون ظهور بقعة سوداء في المنطقة التي استخدم الجهاز عليها	ويُعد احمرار الجلد الذي يستمر أكثر من يوم واحد عرضاً جانبياً غير شائع لتقنية نبضات الضوء المكثف وعلامة على حرق الجلد. انتظري حتى موعد الجلسة التالية حتى يشفى الجلد تماماً. تأكدي من استخدام شدة ضوء أقل في الجلسات التالية.

السلامة العامة - لتجنب الإصابة

- لا تستخدم جهاز Lumea على اضطرابات الجلد الوعائية (مثل الدوالي).
- تأكد من أنه يمكن فصل محول الطاقة بسهولة عن مأخذ التيار الكهربائي في حال حدوث عطل. يُعد محول الطاقة جهاز فصل.
- لا تشاركي جهاز Lumea مع أشخاص آخرين لأسباب صحية. جهاز Lumea مخصص للاستخدام الفردي فقط.

تنويه

يجب إبلاغ Philips عن أي حادث خطير وقع فيما يتعلق بهذا الجهاز عبر www.philips.com/support وإلى السلطة المختصة وهي المركز الوطني للإبلاغ عن الأجهزة الطبية (NCMDR) على موقع الويب <https://ncmdr.sfda.gov.sa/default.aspx>

المخاطر والفوائد المتعلقة باستخدام جهاز Lumea

اختبر أمان جهاز Lumea وفعالته في دراسة سريرية شملت 94 امرأة (تتراوح أعمارهن بين 18 و64 عاماً). ولم يغادرياً من المشاركات الدراسة مبكراً بسبب مشاكل متعلقة باستخدام جهاز Lumea.

من بين 94 امرأة، عانت 31 امرأة من أثر جانبي خفيف وعابر واحد أو أكثر: ألم (26)، وحكة (8)، واحمرار (5)، والتهاب الجريبات (2)، والحرقان (1)، وتقشر سطحي، (1) وتورم (1).

النتائج بعد 4 جلسات كل أسبوعين:

- تلقت 75 امرأة الجلسات على الشفة العلوية، وقل الشعر لديهن بنسبة وصلت إلى 61.1%.
- تلقت 76 امرأة الجلسات على منطقة الإبط، وقل الشعر لديهن بنسبة وصلت إلى 60.6%.
- تلقت 77 امرأة الجلسات على منطقة البيكين، وقل الشعر لديهن بنسبة وصلت إلى 61.9%.
- تلقت 76 امرأة الجلسات على الساقين، وقل الشعر لديهن بنسبة وصلت إلى 93.1%.

استكملت 35 امرأة 8 جلسات إضافية لتحسين النتائج (مرة واحدة كل 4 أسابيع)، وقل ظهور الشعر لديهن بنسبة وصلت إلى 67.0% في المناطق التي استخدمن الجهاز عليها.

⚠ تنبيهات

الأمان الحارّي - لتجنب الألم، أو احمرار الجلد، أو الحروق

- لا تديئي أو تستأنفي استخدام Lumea في غضون 4 أسابيع بعد تعريض بشرتك لأشعة الشمس من دون حماية (أشعة الشمس الطبيعية أو جهاز التسمير) أو في غضون أسبوعين بعد التسمير الاصطناعي (باستخدام الكريمات أو البخاخات أو الحبوب). احرصي على إجراء اختبار البشرة قبل بدء الجلسة أو استئنافها، اطلعي على قسم «اختبار بشرتك».
- لا تستخدمي جهاز Lumea مع الجسم أو ملحق الدقة على الشفة العليا أو مناطق أخرى من وجهك. استخدمي ملحق الوجه فقط لأنه يحتوي على فلتر مدمج إضافي للعمل على الوجه أسفل عظام الوجنتين.
- لا تستمري في الجلسة إذا كنت تعانين من ألم أو تحسّس جلدي لأكثر من 24 ساعة. قللي إعداد شدة الضوء أو توقفي عن استخدام جهاز Lumea.
- لا تعرّضي بشرتك لأشعة الشمس دون حماية عند استخدام جهاز Lumea. خلال فترة العلاج بالضوء النبضي المكثف بالكامل، وخاصة خلال الـ 48 ساعة الأولى بعد العلاج، قومي بحماية بشرتك من أشعة الشمس بالملابس أو ضعي منتجات واقية من الشمس لا تقل عن SPF50 على مناطق العلاج غير المغطاة.
- لا تقومي بالوميض على نفس البقعة عدة مرات أثناء العلاج. ويجب أن يغطي سطح الملحق البقعة السابقة للتأكد من استخدامه على منطقة البشرة بأكملها، لكن استخدامه على المنطقة نفسها عدة مرات لن يزيد من فعاليته.
- لا تستخدمي جهاز Lumea على المنطقة الأغعمق من البشرة مقارنة بالمنطقة المحيطة بها. يمكن أن تؤدي مناطق الجلد الشديدة التصبغ هذه إلى ارتفاع درجة حرارة الجلد بشكل مفرط. قد تشمل المناطق الداكنة ما يلي:
 - البقع الداكنة مثل الشامات، أو البثور، أو النمش الكبير، أو الوحمات، أو البقع العمرية، أو الأوردة الداكنة،
 - العلامات المضافة على الجلد، مثل الوشم أو المكياج (الدائم).
- لا تستخدمي جهاز Lumea على وجه الرجل ورقبته. الشعر في هذه المناطق له خصائص مختلفة.
- لا ترسمي علامات على الجلد/المنطقة التي تستخدمين الجهاز عليها بأي نوع من أنواع الأقلام الرصاص أو الحبر.
- لا تستخدمي جهاز Lumea إذا كان هناك أي غسول، أو جل، أو شعر ظاهري على بشرتك أو على الجهاز. يُستخدم الجهاز فقط على البشرة النظيفة، والجافة، والخالية من الشعر. تأكدي من أن جهاز Lumea نظيف وجاف أيضًا.
- لا تستخدمي جهاز Lumea في حال تلف نافذة خروج الضوء أو إطارها. افحصي جهاز Lumea بحثًا عن التلف قبل استخدامه.
- لا تعرّضي جهاز Lumea لصدمات شديدة لتجنب خطر تلف المصباح.
- لا تستخدمي الجهاز على أي موضع أكثر من مرة واحدة كل أسبوعين.

سلامة العينين - لتجنب تعب عينيك

- دائماً استخدمي جهاز Lumea في غرفة مضاءة جيداً. عيناكي أكثر حساسية للضوء الشارد في الظروف المظلمة.

السلامة العامة - لتجنب الإصابة الخطيرة

- **افحصي** جهاز Lumea دائمًا قبل استخدامه. لا تستخدمِي جهاز Lumea أو محول الطاقة في حال تلفهما. استبدلي دائمًا الجزء التالف بأحد الأنواع الأصلية.
- **افصلي** دائمًا جهاز Lumea، بما في ذلك محول الطاقة، عن التيار الكهربائي بعد الاستخدام وفي حال انقطاع التيار لمنع تلفه.
- لا تستخدمِي جهاز Lumea مباشرةً فوق أي غرسات نشطة في جسمك أو بالقرب منها؛ مثل:
 - منظم ضربات القلب،
 - المحفز العصبي،
 - مضخة الأنسولين.
- تحققي من دليل الغرسات النشطة لمعرفة المسافة الآمنة بينها وبين الأجهزة الإلكترونية الأخرى.
- لا تضعي سلك محول الطاقة حول رقبتك لتجنب الاختناق.
- لا تستخدمِي جهاز Lumea إذا كنت تعانين من ضعف قدراتك الجسدية، أو الحسية، أو العقلية؛ فقد يؤدي هذا إلى الاستخدام غير الصحيح له.
- لا تستخدمِي جهاز Lumea، بما في ذلك كبل محول الطاقة، ضمن نطاق يبلغ 30 سم من الأجهزة التي تعمل على شبكة WiFi و/أو بتقنية Bluetooth (معدات اتصالات RF المحمولة) أو الهواتف الذكية لأنها قد تتسبب في تعطل جهاز Lumea.
- لا تستخدمِي جهاز Lumea بجانب جهاز آخر أو بالقرب منه. فقد يؤدي هذا إلى تشغيله بشكل غير صحيح.
- لا تتركي جهاز Lumea من دون مراقبة عند تشغيله لتجنب خطر نشوب حريق.
- **ضعي** جهاز Lumea بعيداً عن متناول أيدي الأطفال. يجب الإشراف على الأطفال للتأكد من أنهم لا يلعبون بالمعدات.
- **دائمًا** تأكدي من أن محول الطاقة متصل بمقبس حائط مثبت بشكل صحيح ويعمل لتجنب خطر نشوب حريق.

⚠ تحذيرات

الأمان الحارّي - لتجنب الألم، أو احمرار الجلد، أو الحروق

- لا تستخدم جهاز Lumea على الشفاه، أو فتحتي الأنف، أو الأنين، أو فروة الرأس، أو الحلمتين، أو الشفرين الصغيرين، أو العجان، أو حول فتحة الشرج، ويجب على الرجال عدم استخدامه حول القضيب، أو كيس الصفن. إن هذه المناطق حساسة للغاية ولها بشرة أرق وأعمق.
- لا تلمس نافذة خروج الضوء خلال الاستخدام أو بعده مباشرة، إذ يصبح هذا الجزء ساخناً جداً. بعد استخدام جهاز Lumea، افصله عن الطاقة وافصله كذلك محول الطاقة واتركيهما ليبردا لمدة 15 دقيقة قبل تنظيفهما وتخزينهما.
- لا تستخدم جهاز Lumea في حالة معاناتك أماً جلدية نشطة، فقد يؤدي ذلك إلى تفاقم الحالة أو التسبب في الشعور بألم إضافي، أو احمرار البشرة، أو الحروق. قد تشمل الأمراض الجلدية ما يلي:
 - الكدمات، أو الطفح الجلدي، أو القروح، أو الجروح المفتوحة، أو الندوب، أو حروق الشمس. انتظري حتى تلتئم الإصابة تماماً قبل استئناف العلاج.
 - تهيجات الجلد أو الالتهابات أو الإصابات النشطة. انتظري حتى يتعافى الجلد تماماً قبل استئناف العلاج.
 - تلف الجلد بعد استخدام كريمات إزالة الشعر. انتظري يومين على الأقل قبل استكمال الجلسة.
 - الجلد التالف بعد استخدام طرق تقشير الجلد أو كشط الجلد (الدقيق) أو تجديد شباب الجلد. انتظري أسبوعاً واحداً على الأقل قبل استكمال الجلسة.

سلامة العين - لتجنب الإضرار بالعين

- لا تستخدم جهاز Lumea بالقرب من أو على العينين، أو الجفون، أو الحاجبين. استخدم جهاز Lumea فقط أسفل عظام الوجنتين.

السلامة الكهربائية - لتجنب الصدمات الكهربائية

- لا تستخدم جهاز Lumea في محيط رطب (مثل: استخدامه داخل الحمام أو بالقرب من حوض الاستحمام، أو حمام السباحة).
- لا تدعي الماء يدخل إلى جهاز Lumea.
- لا تعرضي جهاز Lumea لدرجة حرارة عالية، وخاصة في بيئة ذات رطوبة عالية، من دون تركه لمدة 30 دقيقة حتى يتكيف مع درجة الحرارة المحيطة في الغرفة التي سيستخدم فيها (راجع «المواصفات الفنية»).
- لا تحاولي فتح جهاز Lumea، أو محول الطاقة، أو إصلاحه، أو تعديله، أو تعديله.
- لا تلمس أي جزء داخلي من جهاز Lumea أو محول الطاقة التالف.
- لا تستخدم أي ملحقات أو محول طاقة آخر غير المورد مع Lumea.
- دائماً تحقق مما إذا كان الجهد المشار إليه في محول الطاقة يتوافق مع جهد التيار الكهربائي المحلي أم لا قبل توصيل الجهاز.

تعليمات مهمة بشأن السلامة

⚠️ اقرئي دائماً هذه التعليمات واتبعيها بعناية بهدف الحصول على أقصى استفادة عند استخدامك لجهاز Lumea وبهدف تقليل مخاطر التعرض لإصابة، والتحسسات الجلدية، والآثار الجانبية الناتجة منه. تأكدي من حفظ هذه التعليمات للرجوع إليها في ما بعد.

⚠️ موانع الاستعمال

- لا تستخدمي جهاز Lumea إذا كنت قد تلقيت علاجاً كيميائياً أو علاجاً إشعاعياً خلال الأشهر الثلاثة الماضية.
- لا تستخدمي جهاز Lumea إذا انخفضت حساسيتك للألم بعد استخدام الأدوية (مثل مسكن الآلام أو كريم التخدير) أو بسبب حالة طبية.
- لا تستخدمي جهاز Lumea إذا كان لون بشرتك بنياً داكناً أو أسمر داكناً.
- لا تستخدمي جهاز Lumea إذا كنت حاملاً، أو مرضعة، أو ترضعين طفلك رضاعة طبيعية.
- لا تستخدمي جهاز Lumea إذا كان عمرك أقل من 18 عاماً أو أكثر من 65 عاماً.

الحساسية للضوء:

- لا تستخدمي جهاز Lumea إذا كانت بشرتك شديدة الحساسية للضوء، على سبيل المثال إذا كنت تعانين الحالات التالية:
 - حساسية من الشمس،
 - تناول أدوية أو علاجات عشبية معينة. تحقق من نشرة التعليمات الموجودة داخل عبوات أدويةك لمعرفة ما إذا كانت قد تسبب لك زيادة الحساسية لأشعة الشمس.

الحالات الطبية:

- لا تستخدمي جهاز Lumea - أو استشري طبيبك أولاً - إذا كنت تعانين الحالات التالية:
 - مرض الصرع مع الحساسية للضوء الوامض،
 - سرطان الجلد النشط أو تاريخ مرضي سابق لسرطان الجلد،
 - اضطراب تصبغ الجلد مثل البهاق،
 - أعراض سابقة من ضعف التئام الجروح،
 - حالة نقص المناعة الناجمة عن مرض مثل الإيدز أو الأدوية المثبطة للمناعة،
 - أعراض سابقة من اضطراب الكولاجين، بما في ذلك أعراض سابقة من تكوين ندبة الجدة.

مرحبًا بك في Lumea

احصلي على تقليل لطيف وسريع وفعال للشعر باستخدام جهاز Lumea. يمكنك استخدام جهاز Lumea بكل سهولة وأنت مرتاحة داخل منزلك. مستعدة للبدء؟



امسحي الرمز ضوئيًا باستخدام كاميرا الهاتف أو تطبيق QR أو تفضلي بزيارة www.philips.com/myIPL

ملاحظة: في دليل المستخدم هذا، يشير Lumea إلى أجهزة Philips Lumea IPL 7000 Series، Series 8000، و Series 9000.

الاستخدام المخصص

إن Philips Lumea هو جهاز منزلي مخصص لإزالة الشعر الزائد والاستخدام لشخص واحد. يمكن استخدام جهاز Lumea للنساء اللاتي تتراوح أعمارهن بين 18 و 65 عامًا ممن يتمتعن ببشرة بنية فاتحة إلى متوسطة، وشعر جسم طبيعي أشقر غامق، وبني، وبني غامق، وأسود. يمكن للنساء استخدام جهاز Lumea على الجسم وعلى الوجه* (أسفل عظام الوجنتين). إن Lumea عبارة عن جهاز محمول باليد مخصص للاستخدام المنزلي.

* وذلك، إذا كان جهاز Lumea لديك يتضمن ملحق الوجه (BRI932، BRI933، BRI937، BRI938، BRI983، BRI984) فقط.

هل يناسبك جهاز Lumea؟

جهاز Philips Lumea ليس مناسبًا لجميع الأشخاص. تحقق من جداول لون البشرة ولون شعر الجسم في بداية دليل المستخدم هذا واقرئي تعليمات السلامة المهمة بعناية لمعرفة ما إذا كان جهاز Lumea مناسبًا لك أم لا.

ألوان شعر الجسم غير المناسبة: لا تكون المعالجة باستخدام Lumea فعالة إذا كان لديك شعر أشقر فاتح أو رمادي أو أحمر أو أبيض لأن الشعر بهذه الألوان لا يمتص ضوءًا كافيًا.


ألوان البشرة غير المناسبة: لا يمكنك استخدام جهاز Lumea إذا كانت بشرتك ذات ألوان تتراوح من بني داكن إلى أسمر داكن، لأن ألوان البشرة هذه تمتص الضوء بسهولة أكبر؛ ما قد يؤدي إلى الشعور بالألم أو إلى احمرار الجلد أو إلى الحروق.

الفصول

69	مرحبًا بك في Lumea
69	الاستخدام المخصص
68	تعليمات مهمة بشأن السلامة
68	موانع الاستعمال
67	تحذيرات
65	تنبيهات
64	المخاطر والفوائد المتعلقة باستخدام جهاز Lumea
63	الآثار الجانبية المحتملة للعلاج بتقنية الضوء النبضي المكثف (IPL)
60	ماذا يُقصد بتقنية نبضات الضوء المكثف؟
59	نظرة عامة عن جهاز Lumea
58	تقنية SenseIQ 2.0
57	المرفقات
56	طرق المعالجة
55	دليل جلستك
54	استعدي لجلستك
53	اختبار جلدك
51	تنفيذ المعالجة
49	التنظيف والتخزين
48	بعد العلاج
47	التخطيط لجلساتك التالية
46	استكشاف الأخطاء وإصلاحها
44	المعلومات الإضافية

غير مناسبة

درجات لون البشرة - مناسبة

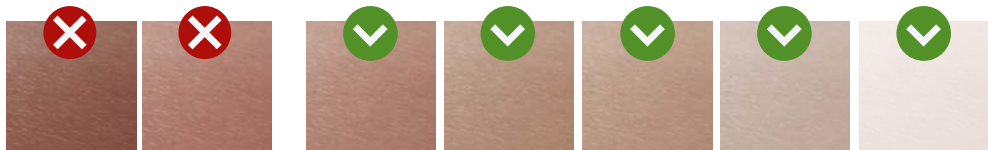

غامق

بني غامق

بني متوسط

بني فاتح


بيج

فاتح

وسط

غير مناسبة

ألوان شعر الجسم - مناسبة

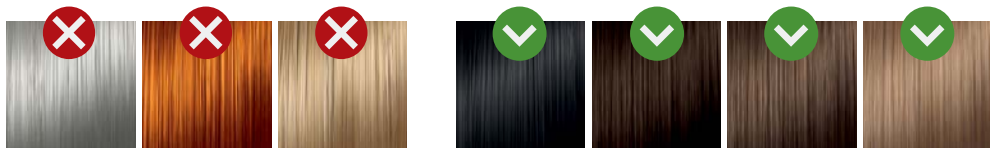

أبيض / رمادي

أحمر

أشقر فاتح

أسود

بني غامق

بني

أشقر غامق



© 2024 Koninklijke Philips N.V.
All rights reserved.
جميع الحقوق محفوظة.
3000.090.3393.1
(2024-08-08)



Philips Consumer Lifestyle B.V.
Tussendiepen 4
9206 AD Drachten
The Netherlands

فيليبس كونسيومر لايفستايل بي.في.
دراشتن، هولندا AD توسانديبين ٤، ٩٢٠٦

FSC

www.philips.com