
PHILIPS

Colgate

SonicPro for Kids

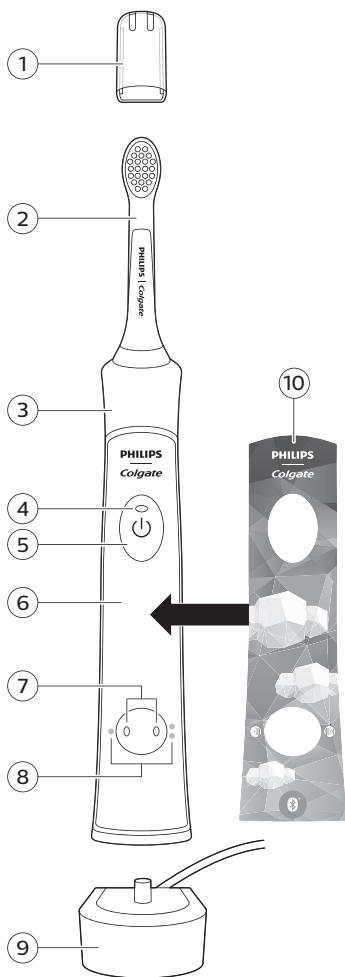English 6
Español 19
Português do Brasil 33

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **www.philips.com/welcome**.

Important safety information

Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Note: Complete firmware update, if prompted, to access the latest improvements and features.

Warnings

- Keep the chargers away from water.
- After cleaning, make sure the charger is completely dry before you connect it to the wall socket.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Only use original Philips accessories or consumables. Only use charger provided with the product.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country. If the cord is damaged, replace the charger with a new one.
- Do not use the chargers outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- This appliance has only been designed for cleaning teeth, gums and tongue.

- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

MEDICAL WARNINGS

- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables. Only use detachable supply unit HX61xx
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.

- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on recycling.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Electromagnetic fields (EMF)


This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Storage and transportation conditions

Temperature -10°C to 60°C (14 °F to 140 °F).

Explanation of symbols

The following symbols may appear on the product:
Read Operator's Manual.

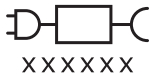Detachable Power Supply Part: The "xxxxxx" indicates the model number of the approved chargers to be used.
(HX61xx)

Introduction

Congratulations on your purchase. To fully benefit from the support that Philips Colgate SonicPro for Kids offers, to register your product at www.philips.com/welcome.

Your Philips Colgate SonicPro for Kids (Fig. 1)

- 1 Brush head cap
- 2 Brush head
- 3 Handle with soft grip
- 4 Battery charge indicator
- 5 On/off button
- 6 Panel for your favorite sticker
- 7 Brushing mode lights
- 8 Brushing modes
- 9 Charger
- 10 Sticker

Preparing for use

Applying your favorite sticker


The Philips Colgate SonicPro for Kids comes with replaceable stickers.

- 1 Apply your favorite sticker on the handle.

Note: Make sure you line up the cut-out areas on the sticker with the buttons on the front of the handle.

Attaching the brush head

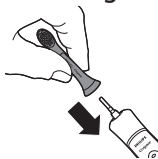- 1 Align the brush head so the bristles point in the same direction as the front of the handle.
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

Bluetooth connectivity

Before you start the connectivity procedure, make sure your mobile device runs on Bluetooth 4.0 or higher.

Note: For the latest information on device compatibility, please visit the App Store or Google Play.

Connecting your toothbrush

- 1 Download the Philips Colgate SonicPro for Kids app from the App Store or Google Play.
- 2 Confirm Bluetooth is switched on on your mobile device before you open the app. Check the settings on your mobile device to turn on Bluetooth.
- 3 Open the Philips Colgate SonicPro for Kids app and follow the app instructions displayed on your mobile device.
- 4 To connect the toothbrush to the app, turn on the toothbrush by pressing either the power on/off button or the brushing mode button.

Note: Make sure the toothbrush is not on the charger when you turn it on.

- 5 The app automatically connects with the toothbrush. If the toothbrush does not connect, follow the troubleshooting steps provided in the app.

Note: Keep your mobile device close to your toothbrush when you use your toothbrush. Make sure your mobile device is safe in a dry place.


Using your Philips Colgate SonicPro for Kids


Brushing instructions

Note: For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.

- 1 With the attached brush head, wet the bristles and apply a small amount of toothpaste.


Note: Use no more than a pea sized amount of toothpaste and spit out excess toothpaste after brushing.

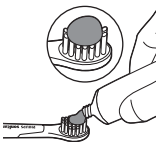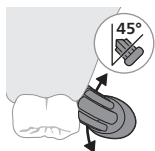
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

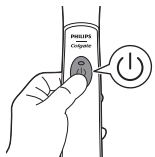- 3 Press the power on/off button to turn on the toothbrush.
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

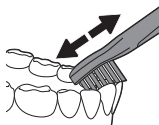
- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the KidPacer feature.

- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

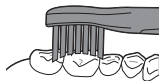
- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Note: When the Philips Colgate SonicPro for Kids is used in clinical studies, the handle should be fully changed and the EasyStart feature deactivated.

Your Philips Colgate SonicPro for Kids is safe to use on:

- Braces (brush heads wear out sooner when used on braces).
- Dental restorations (fillings, crowns, veneers).

Brushing modes

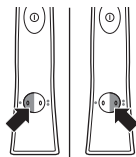The Philips Colgate SonicPro for Kids automatically starts in the default high mode.


- Press the brushing mode button to toggle between modes when the toothbrush is off.

The green mode light indicates the selected mode.

- High mode: Full-strength mode designed for older/independent brushers.
- Low mode: Reduced power mode used for training and introducing younger children to the sonic experience.

Battery status and charging

Charging your Philips Colgate SonicPro for Kids

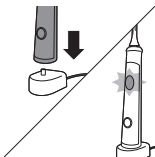- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
 - The battery charge indicator in the power on/off button flashes and the handle beeps twice to indicate that the toothbrush is charging.
 - When the toothbrush is fully charged, the battery charge indicator stops flashing and remains on.

Note: If the battery charge of the Philips Colgate SonicPro for Kids is low, you hear 3 beeps and the battery charge indicator flashes rapidly for 30 seconds after the brushing cycle.

Note: To keep the battery fully charged at all times, you may keep the Philips Colgate SonicPro for Kids on the charger when not in use. The charging indicator flashes until your toothbrush is fully charged.

Features

KidTimer

- The KidTimer plays a congratulatory sound sequence at the end of the pre-determined brushing time. Additionally, the KidTimer automatically turns off the handle when the brushing cycle is complete.
- When in high mode, the KidTimer is set at 2 minutes.
- When in low mode the KidTimer is set at 1 minute and slowly increases over time. The progression of time helps young children reach the dentist recommended 2 minutes of brushing.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

KidPacer

The Philips Colgate SonicPro for Kids uses a short series of tunes to alert the brusher to move to the next brushing quadrant of the mouth and concludes the brushing cycle with a congratulatory tune before it automatically turns off. This ensures that the child's teeth get a thorough, overall clean (see 'Brushing instructions').

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

EasyStart

The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing power of the Philips Colgate SonicPro for Kids. The toothbrush comes with the EasyStart feature activated.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the EasyStart ramp-up cycle.

Activating or deactivating the EasyStart feature

- Place the handle on the charger.
- To activate EasyStart: Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

- To deactivate EasyStart: Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the EasyStart feature beyond the initial ramp-up period is not recommended and reduces Philips Colgate SonicPro for Kids effectiveness in removing plaque.

Cleaning

You should clean your Philips Colgate SonicPro for Kids regularly to remove toothpaste and other residue. Failure to clean your toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean brush heads, the handle or the charger in the dishwasher.

Warning: Do not clean with boiling hot water.

Toothbrush handle


- 1 Remove the brush head from the handle and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal at the top of the handle with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.


Note: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean handle as this may cause discoloration.

Brush head

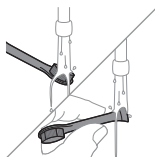- 1 Remove the brush head from the handle and rinse it thoroughly.
- 2 Rinse the brush head connection with warm water at least once a week.

Charger

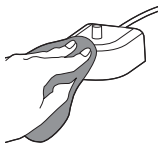- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the Philips Colgate SonicPro for Kids for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

- Replace Philips Colgate SonicPro for Kids brush heads every 3 months to achieve optimal results.
- Use only Philips Colgate SonicPro for Kids replacement brush heads.
- Brush head can be thrown away with the normal household waste.

Warranty and support

If you need information or support, please visit **www.philips.com/support** or read the international warranty leaflet.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver.

Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Colgate SonicPro for Kids and let it run until it stops. Repeat this step until you can no longer turn on the Philips Colgate SonicPro for Kids.
- 2 Remove and discard the brush head.


- 3 Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle.


- 4 Insert the screwdriver into the gap and pry the bottom cap out of the handle.

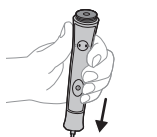- 5 Hold the handle upside down and push down on the shaft to release the internal components of the handle.

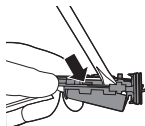- 6 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier.

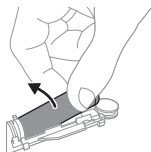- 7 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an

official collection point for recycling. By doing this, you help to preserve the environment.

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©2019 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV.
Colgate and the Colgate Logo are the registered trademarks of Colgate-Palmolive Company.