PHILIPS

5100 Series/系列

PHILIPS
sonicare

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300014409641 (2024-12-19)

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
© 2024 Koninklijke Philips N.V. (KPNV). All
rights reserved. Philips and the Philips shield are
trademarks of KPNV. Sonicare and the Sonicare
together are trademarks of Philips Oral Healthcare,
LLC and/or KPNV.


www.philips.com/Sonicare


1

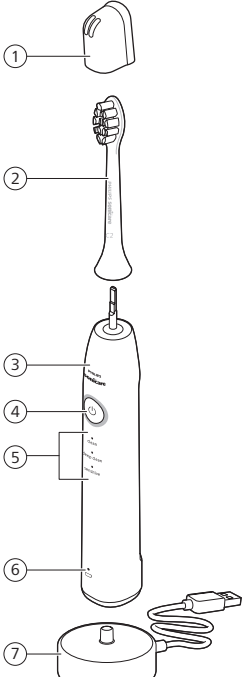2

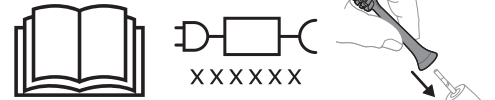
5

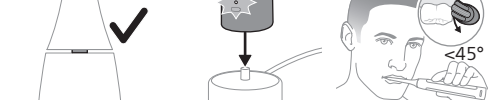8

11

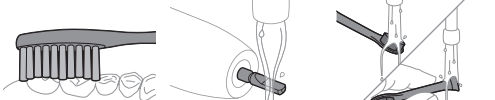14

English

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support').
- If the USB cable or USB power supply unit is damaged, stop using it.
- Keep the charger away from water.
- After cleaning, make sure the charger is completely dry before you connect it.
- Do not clean any part of the product in the dishwasher.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Make sure the handle is completely dry before you connect the USB charger.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not use adapters and chargers outdoors or near heated surfaces.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.

Medical warnings

- This appliance is a personal care device and is not intended for use on multiple persons.
- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5 V, $\geq 2V$. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country (see the international warranty leaflet for contact details). Using a non-certified supply unit may cause hazards or serious injuries.
- Charge, use and store the product at a temperature between 0 °C and 40 °C.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.

Storage and transportation conditions

Temperature -10°C to +60°C (14 °F to 140 °F).

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Your Philips Sonicare (Fig. 1)

- Brush head cap
- Brush head
- Handle
- Power on/off and mode button
- Brushing modes
- Battery indicator
- Charger

Note: Appearance may slightly differ, depending on the model.

Note: Wall adapter not included. You can use any compatible IPX4 splash-proof adapter (identifiable by the "IPX4" text or symbols on the adapter).

The following symbols may appear on the product. This symbol means read the instructions for use, before using the appliance (Fig. 2). This symbol means to use the indicated power supply (Fig. 3).

Getting started

Attaching the brush head

- Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 4).
- Press the brush head down onto the metal shaft until there is a small 1-2 mm gap (Fig. 5).

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

Charging your Philips Sonicare

- Place the power toothbrush handle on the charger (Fig. 6).
- Connect the charger to a 5V USB wall adapter.
- Your toothbrush continues to charge until it is fully charged.
 - When the toothbrush is charging, the battery indicator pulsates green.
 - When fully charged, the battery indicator will light up solid green for 30 seconds before turning off.

Note: It can take up to 10 hours to fully charge the battery, but you can use the toothbrush before it is fully charged.

Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.

Brushing instructions

- Wet the bristles and apply a small amount of toothpaste.
 - Place the toothbrush bristles against the teeth at a slight angle (45 degrees) (Fig. 7). The handle has two 45° degree indication bars to help guide your Bass brushing technique. Press the brush head gently against the gum to make the bristles reach the gum line or slightly beneath the gum line (Fig. 8).
 - Press the power on/off button to turn on the Philips Sonicare.
- Note: Keep the center of the brush in contact with the teeth at all times.
- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. Do not scrub. It is not recommended to scrub as you would with a manual toothbrush.

- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 9).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the QuadPacer feature (see chapter 'Features').

- Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth) (Fig. 10).
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer (Fig. 11).

Your Philips Sonicare toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Battery status

This Philips Sonicare toothbrush is designed to provide at least 80 brushing sessions after fully charged for 10 hours, each session lasting 2 minutes. When the battery is low, the battery indicator starts to flash amber.

Note: When the battery is completely empty, the toothbrush switches off. Place the handle on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your Philips Sonicare toothbrush on the charger when not in use.

Note: Battery usage time is expected to reduce over product life time.

Features

Brushing modes

Your power toothbrush is equipped with three different modes:

- Clean**
For exceptional everyday cleaning (2 minutes).
- Deep Clean**
For invigorating and thorough cleaning (3 minutes).
- Sensitive**
For gentle cleaning of gums (2 minutes).

Note: When you use the toothbrush for the first time, the default mode is 'Clean'.

To toggle between the brushing modes:

- Press the power on/off button once to turn on the toothbrush.
- Soft Start, 3 seconds.
- Press a second time **within 2 seconds** to change the brushing mode.

After the toothbrush has been on for more than 2 seconds, if you want to change modes, you'll need to pause the toothbrush first (click the mode button), then press the mode button **within 2 seconds**.

EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare. The EasyStart feature works with Clean mode only. To deactivate EasyStart see 'Activating or deactivating features'.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

QuadPacer

The Quadpacer is an interval timer that beeps every 30 seconds to remind you to brush all teeth evenly by dividing your mouth into 4 sections (Fig. 10).

Pressure sensor

Your Sonicare toothbrush measures the pressure you apply while brushing to protect your gums and teeth from damage. If you apply excess pressure, the handle changes its vibration until you reduce the pressure.

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
 - Pressure Sensor Feedback
- To activate or deactivate these features, follow the instructions below:

EasyStart

- Put the handle on the plugged-in charger.
- Keep the power on/off button pressed until you hear a single short beep.
- Release the power on/off button.
 - Triple tone of low-medium-high means the EasyStart feature has been activated. The battery indicator blinking green twice shows the EasyStart is activated.
 - Triple tone of high-medium-low means the EasyStart feature has been deactivated. The battery indicator blinking amber once shows the EasyStart is deactivated.

Note: To achieve optimal efficacy, EasyStart needs to be deactivated.

Pressure Sensor Feedback

- Put the handle on the plugged-in charger.
- Keep the power on/off button pressed until you hear three short beeps.
- Release the power on/off button.
 - Triple tone of low-medium-high means the Pressure Sensor Feedback feature has been activated. The battery indicator blinking green twice shows the Pressure Sensor Feedback is activated.
 - Triple tone of high-medium-low means the Pressure Sensor Feedback feature has been deactivated. The battery indicator blinking amber once shows the Pressure Sensor Feedback is deactivated.

Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean brush heads, the handle or the charger in the dishwasher.

Warning: Do not clean with boiling hot water.

Toothbrush handle

- Remove the brush head and rinse the metal shaft area with warm water (Fig. 12). Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

Wipe the entire surface of the handle with a damp cloth.

Caution: Do not tap the handle on the sink to remove excess water.

Brush head

- Rinse the brush head and bristles after each use (Fig. 13).
- Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Charger

- Unplug the charger before you clean it.
- Wipe the surface of the charger with a damp cloth.

Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter 'Cleaning') and store it in a cool and dry place away from direct sunlight.

Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number HX507x.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.	This is normal. A small gap is necessary to allow movement in order to deliver the vibrations needed to clean your teeth.	
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the handle.	Remove the brush head from the handle and re-set it on the metal shaft leaving small gap. (1-2mm)
I have charged my toothbrush, but it only works for a short time.	The toothbrush handle is not placed vertically on the charger, resulting in no charge.	Place handle vertically on charger base and make sure the charger is plugged in.

Problem	Possible cause	Solution
The toothbrush stops working.	The battery is empty.	Charge the toothbrush handle on the charger.

Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 14).
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

CONTAINS LITHIUM ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the appliance is disconnected from the wall socket and that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

Only remove the rechargeable battery when you discard the product. Before you remove the battery, make sure that the product is disconnected from the wall socket and that the battery is completely empty.

- Remove the inner body from the handle.
- Remove the rechargeable battery from the inner body.

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