

PHILIPS

sonicare

*DiamondClean
Essential*
7000 Series/系列

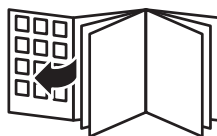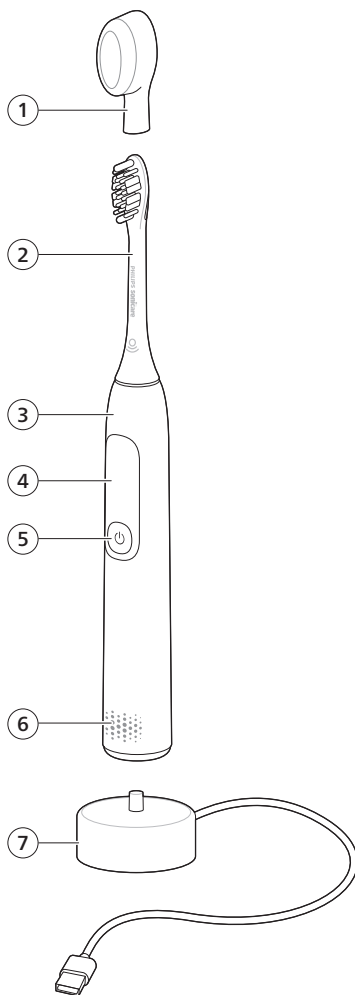English 6

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Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **www.philips.com/welcome**.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

Warnings

- Only use original Philips accessories or consumables. Only use the USB charger and USB cable provided with the product.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support'). If the cord is damaged, replace the charger.
- Keep the charger away from water.
- Make sure the handle is completely dry before you connect the USB charger.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Do not clean any part of the product in the dishwasher.
- Do not use adapters and chargers outdoors or near heated surfaces.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.

- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.

Medical warnings

- Stop using the appliance and consult your dentist/doctor if excessive bleeding occurs after use, if bleeding continues to occur after 1 week of use or if you experience discomfort or pain.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use this appliance.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use this appliance.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Battery safety instructions

- This appliance contains batteries that are only to be removed by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.

- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- Charge only with Philips WAA1001, WAA2001, UL or equivalent safety approved IPX4 adapters with input ratings (100-240 V~; 50/60 Hz; 3.5 W) and output ratings (5 V dc; 2.5 W). For the list of adapters, please see **www.philips.com/support**.
- Charge, use and store the product at a temperature between 0 °C and 40 °C.
- This appliance contains a Lithium-ion rechargeable cell. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

Storage and transportation conditions


Temperature -20°C to +60°C (-4 °F to 140 °F).

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

The following symbols may appear on the product:
Read Operator's Manual.

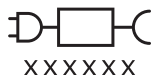Detachable Power Supply Part: The "xxxxxx" indicates the model number of the approved chargers to be used.

Intended use

Sonicare power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare power toothbrushes are intended for consumer home use. Use by children should be with adult supervision.

Your Philips Sonicare toothbrush

- 1 Brush head cap(s)
- 2 Smart brush head(s)
- 3 Handle
- 4 Smart display
- 5 Power on/off and mode button
- 6 Pressure sensor indicator
- 7 USB charger

Note: Wall adapter not included. You can use any compatible IPX4 splash-proof adapter (identifiable by the "IPX4" text or symbols on the adapter).

Getting started

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it stops.


Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.


To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gradually increases the power over the first 14 brushing sessions to help you adapt to brushing with your Sonicare toothbrush (see 'EasyStart').

Brushing instructions

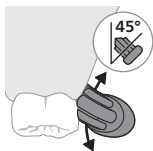
- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Press the power on/off button to turn on the Philips Sonicare. Then select the brushing mode.
- 3 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

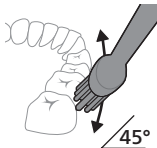- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

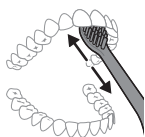- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the QuadPacer feature (see chapter 'Features').

- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.


Your Philips Sonicare toothbrush should be safe to use on braced (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Note: When the Philips Sonicare toothbrush is used in clinical studies, it must be used in Clean mode with the handle fully charged and the EasyStart feature deactivated.

Charging your Philips Sonicare

- 1 Attach the USB charger to a wall adapter and plug the wall adapter into an electrical outlet.
- 2 Place the handle on the charger.
 - The battery icon with lightning symbol is shown on the display and you hear two beeps. The light ring lights up white for 5 seconds.

Note: It can take up to 16 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.

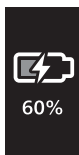Note: This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day).

Features

Brushing modes

Based on the smart brush head you use; the optimal mode will be selected automatically by the BrushSync technology. To change modes, press the mode button while the handle is turned on. Your power toothbrush is equipped with five different modes:

Clean

For exceptional everyday cleaning (2 minutes).

Deep Clean

For invigorating and thorough cleaning (3 minutes).

Gum Pro

For exceptional cleaning and gentle gum massaging (3 minutes).

Sensitive

For gentle cleaning of gums (2 minutes).

White

For removing stains from teeth (2 minutes 30 seconds). After brushing your whole mouth (2 minutes), spend the remaining time brushing the front teeth.

Note: When you use the toothbrush with the smart brush head for the first time, the default mode is 'Gum Pro'.

To toggle between the brushing modes:

- Press the power on/off button once to turn on the toothbrush.
- Press a second time **within 5 seconds** to change the brushing mode.
- After the toothbrush has been on 5 seconds, if you want to change modes, you'll need to turn the toothbrush off, then on and press the power on/off button **within 5 seconds**.

Post Brushing Feedback

Your power toothbrush comes with four different types of feedback:

Duration Feedback

It records if the last brushing cycle reached 2 minutes.

Interval Feedback

It indicates if the time interval between the last two brushing cycles is controlled in the proper duration to keep plaque growth low.

Pressure Feedback

It tells you if the pressure was applied properly during the last brushing cycle.

Cleanliness Feedback

It assesses the overall cleaning effectiveness based on consistent adherence to the criteria of duration, interval, and pressure during brushing.

EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

To deactivate EasyStart see 'Activating or deactivating features'.

Brush head replacement reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When a new smart brush head is attached for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and the display shows the brush head symbol.

- 2 Over time, based on the pressure you apply and how long the brush head has been used, the handle will track brush head wear in order to determine the optimal time to change your brush head. When the life of the brush head ends, the display shows the brush head replacement reminder symbol and indicates that you should replace the brush head as soon as possible.

Pressure Sensor

Your Philips Sonicare is equipped with an advanced sensor that measures the pressure you apply while brushing. The light ring at the bottom of the handle provides instant feedback. Adjust your pressure for a gentle and effective clean:

- Blinking purple with vibration change means you're pressing too hard.
- Steady green means your pressure is optimal.
- Breathing purple means you're pressing too lightly.

QuadPacer

The QuadPacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. This Philips Sonicare toothbrush comes with the QuadPacer feature activated.

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush will pause. After a pause of 15 seconds, the SmarTimer resets.

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder

To activate or deactivate these features, follow the instructions below:

EasyStart

- 1 Put the handle on the plugged-in charger.

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- 2 Keep the power on/off button pressed until you hear a single short beep (after 3 seconds).
 - 3 Release the power on/off button.
- Triple tone of low-medium-high means the EasyStart feature has been activated. The icon on the display shows that the EasyStart has been activated.

- Triple tone of high-medium-low means the EasyStart feature has been deactivated. The icon on the display shows that the EasyStart has been deactivated.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

Brush head replacement reminder

Note: You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

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- 1 Put the handle on the plugged-in charger.
 - 2 Press and hold the power on/off button while the handle remains on the charger.
 - 3 Keep the power on/off button pressed until you hear two short beeps.
 - 4 Release the power on/off button.
- Triple tone of low-medium-high means the brush head replacement reminder feature has been activated. The icon on the display shows that the brush head replacement reminder has been activated.
 - Triple tone of high-medium-low means the brush head replacement reminder feature has been deactivated. The icon on the display shows that the brush head replacement reminder has been deactivated.

Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean the product or accessories with sharp objects or in the dishwasher, microwave, with chemicals or in boiling hot water

Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

Brush head

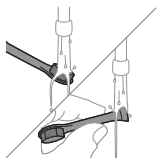- 1 Rinse the brush head and bristles after each use.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Storage

If you are not going to use the product for a long time, unplug it, and clean it (see chapter 'Cleaning'). Store it in a cool and dry place away from direct sunlight.

Charge the toothbrush regularly to prevent the battery from discharging completely. A completely discharged battery may fail to charge again.

Brush head replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX389x).

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This is normal. A small gap is necessary to allow movement in order to deliver the vibrations needed to clean your teeth.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the handle.	Remove the brush head from the handle and re-set it on the metal shaft leaving small gap. (1-2 mm)
	The battery is low.	Charge the toothbrush handle on the charger.
I have charged my toothbrush, but it only works for a short time.	The toothbrush handle is not placed vertically on the charger, resulting in no charge.	Place handle vertically on charger base and make sure the charger is plugged in.

Problem	Possible cause	Solution
The toothbrush stops working.	The battery is empty.	Charge the toothbrush handle on the charger.

Warranty and support

If you need information or support, please visit **www.philips.com/support** or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

Li-ion rechargeable battery

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the appliance is disconnected from the wall socket and that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

When you handle batteries, make sure that your hands, the product and the batteries are dry.


If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.


Before removing the battery, make sure that the battery is completely empty.


To remove the rechargeable battery, you need a flat-head (standard) screwdriver.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Remove the film at the bottom of handle, then release the screw.

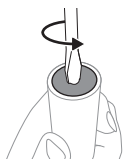
- 3 Remove the end cap from the toothbrush handle.

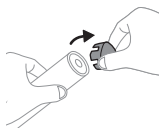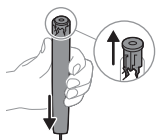
- 4 Holding the handle upside down, press the shaft down on a hard surface.

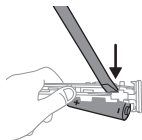- 5 Wedge the screwdriver between the battery and the white frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the white frame.

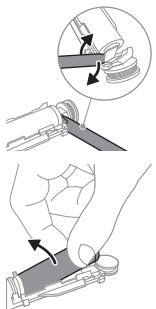- 6 Insert the screwdriver between the bottom of the battery and the white frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

- 7 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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©2025 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.
Lifestyle B.V.