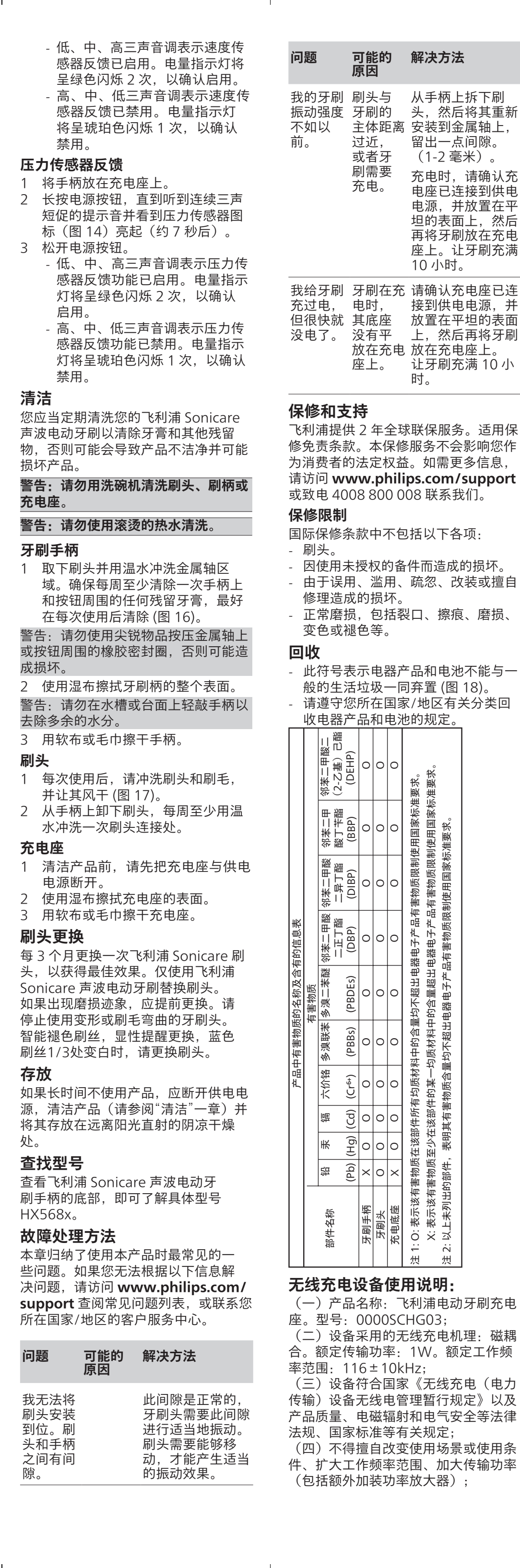
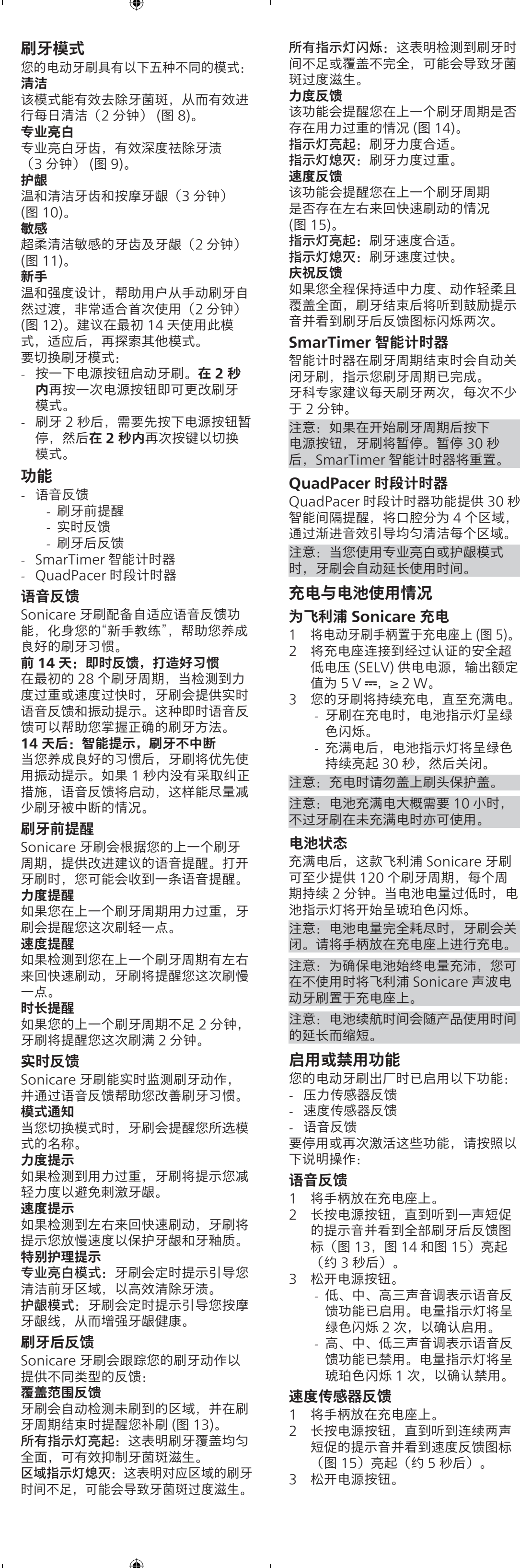
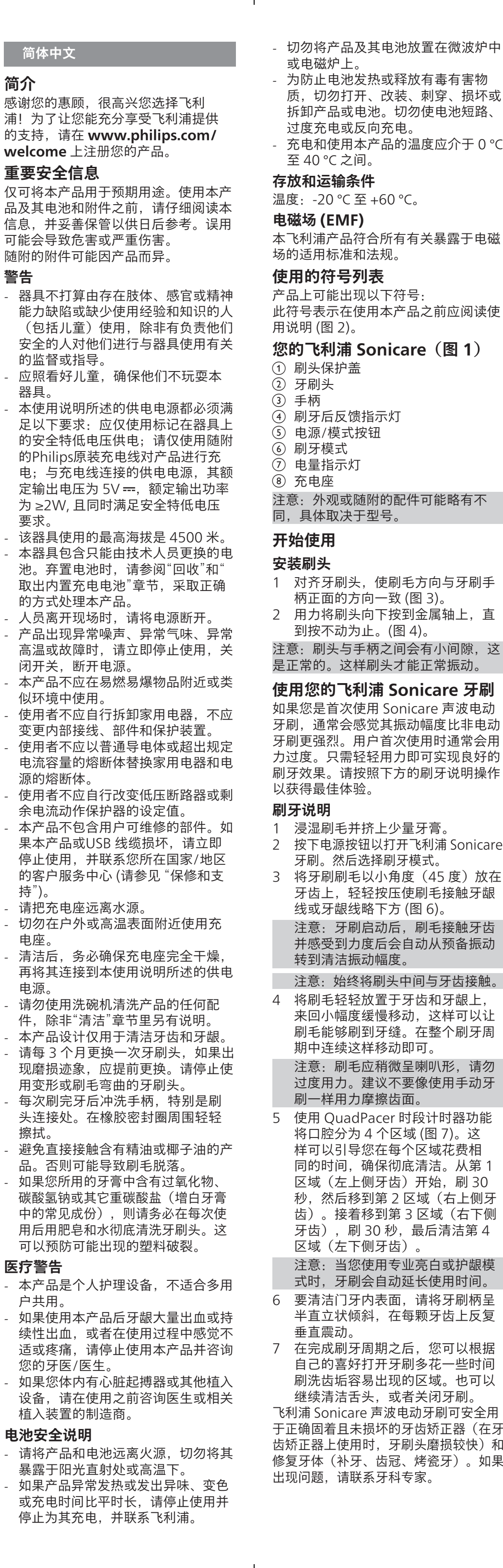
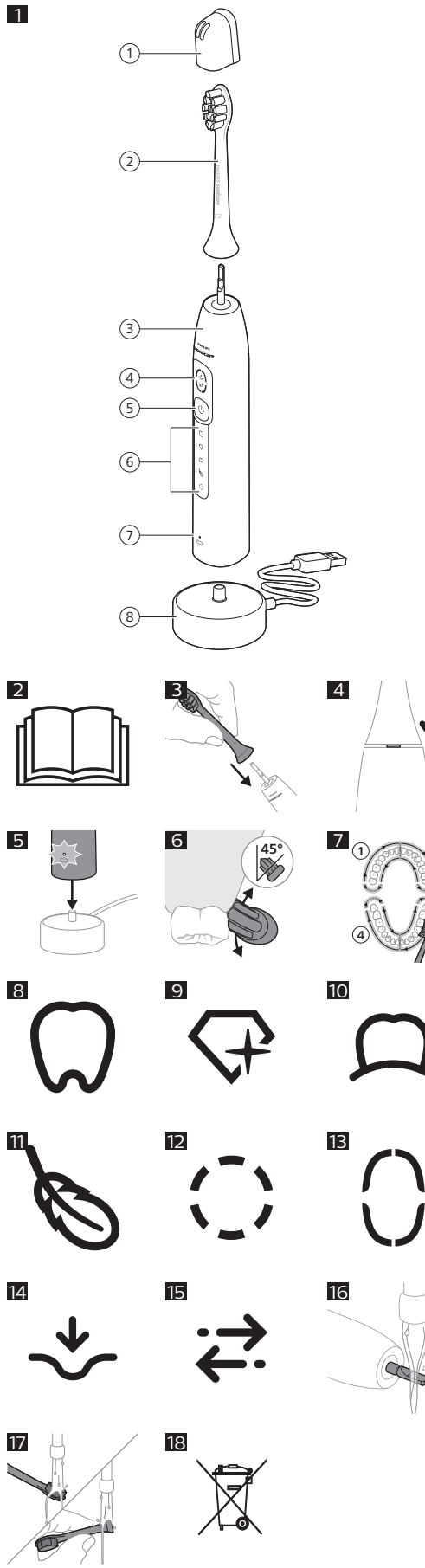
English

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.

Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see "Warranty and support").
- If the USB cable or USB power supply unit is damaged, stop using it.
- Keep the charger away from water.
- After cleaning, make sure the charger is completely dry before you connect it.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher.
- Do not use the charger outdoors or near heated surfaces.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
 - The maximum altitude of use is 4500 meters.
- Medical warnings**
- This appliance is a personal care device and is not intended to be shared by multiple users.
 - Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
 - If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

- Battery safety instructions**
- This appliance contains batteries that are only removable by skilled persons.
 - Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
 - If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
 - Do not place products and their batteries in microwave ovens or on induction cookers.
 - Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
 - This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
 - To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5 V $\overline{\text{rms}}$ $\pm$ 2 V. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country (see the international warranty leaflet for contact details). Using a non-certified supply unit may cause hazards or serious injuries.
 - Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).
- Storage and transportation conditions**
- Temperature: -20°C to +60°C (14 °F to 140 °F).
- Electromagnetic fields (EMF)**
- This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

- LIST OF SYMBOLS USED**
- The following symbols may appear on the product. This symbol means read the instructions for use, before using the appliance (Fig. 2).
- Your Philips Sonicare (Fig. 1)**
- Brush head cap
 - Brush head
 - Handle
 - Post brushing feedback indicator
 - Power/mode button
 - Brushing modes
 - Battery indicator
 - Charger
- Note: Wall adapter not included. Appearance may slightly differ, depending on the model.

- Getting started**
- Attaching the brush head**
- Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 3).
 - Firmly press the brush head down onto the metal shaft until it will not go any further (Fig. 4).
- Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.
- Using your Philips Sonicare toothbrush**
- If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.

- Brushing instructions**
- Wet the bristles and apply a small amount of toothpaste.
 - Press the power button to turn on the Philips Sonicare. Then select the brushing mode.
 - Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline (Fig. 6).
- Note: The vibration activates automatically once the toothbrush is turned on and the bristles touch your teeth.
- Note: Keep the center of the brush head in contact with the teeth at all times.
- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.
- Note: The bristles should slightly flare. Do not scrub. It is not recommended to scrub as you would with a manual toothbrush.
- Divide your mouth into 4 sections with the QuadPacer feature (Fig. 7). This ensures thorough cleaning by guiding you to spend equal time in each section. Begin brushing in section 1 (left upper teeth) and brush for 30 seconds before moving to section 2 (right upper teeth). Continue brushing in section 3 (right lower teeth) and brush for 30 seconds before moving to section 4 (left lower teeth).
- Note: Additional time is added when you use the Expert White or Gum Care mode..

- To clean the inside surfaces of the front teeth, tilt the brush handle up and make several vertical overlapping brushing strokes on each tooth.
 - After you have completed the brushing cycle, you can spend additional time brushing the areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.
- Your Philips Sonicare toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered to and not compromised. If a problem occurs, please follow up with a dental professional.

- Brushing modes**
- Your power toothbrush is equipped with five different modes:
- Clean**
- It removes plaque effectively for daily cleaning (2 minutes) (Fig. 8).
- Expert White**
- Professionally whitens and brightens teeth, helps to deeply remove stains on teeth (3 minutes) (Fig. 9).
- Gum Care**
- Gentle cleaning for teeth and gum massage (3 minutes) (Fig. 10).
- Sensitive**
- Extra gentle cleaning for sensitive gums and teeth (2 minutes) (Fig. 11).
- Easy**
- Designed with gentle intensity to ease the transition from manual brushing, ideal for first-time user (2 minutes) (Fig. 12). It is recommended to use this mode for the first 14 days to build comfort before exploring other modes.

- To toggle between the brushing modes:
- Press the power button once to turn on the toothbrush. Press a second time **within 2 seconds** to change the brushing mode.
 - After 2 seconds of brushing, you will need to press the power button to pause first, and then press it again **within 2 seconds** to switch modes.

- Features**
- Voice Feedback
 - Pre-brushing reminder
 - Instant feedback
 - Post-brushing feedback
 - SmartTimer
 - QuadPacer

- Voice Feedback**
- Your Sonicare toothbrush comes with adaptive voice feedback designed to be your "Newbie Coach" for better brushing habits.
- First 14 days: Instant habit-building support**
- During the initial 28 brushing sessions, it provides real-time voice feedback and vibration alert when excess pressure or scrubbing is detected. This instant voice feedback helps you master the right technique.
- After 14 days: Uninterrupted brushing with smart alerts**
- After you have built good habits, it prioritizes vibration alerts. If no corrective action is made within 1 second, the voice feedback activates to minimize interruptions.

- Pre-brushing reminder**
- Your Sonicare toothbrush provides voice reminders that suggest improvements based on your last brushing session. You might hear one when you turn on the toothbrush.
- Pressure reminder**
- It reminds you to apply moderate pressure if you brushed too hard during your last brushing session.
- Scrubbing reminder**
- It reminds you to slow scrubbing speed if aggressive side-to-side motions were detected in your last brushing session.
- Duration reminder**
- It reminds you to brush for the full 2 minutes if your last brushing session ended early.

- Instant feedback**
- Your Sonicare toothbrush measures your brushing motions and provides voice feedback in real-time to help you improve your brushing habits.
- Mode notification**
- It notifies you the name of the selected mode when you switch modes.
- Pressure alert**
- It alerts you with vibration intensity reduction when excessive pressure is detected to prevent gum irritation.
- Scrubbing alert**
- It alerts you when aggressive side-to-side brushing motions are detected to protect gums and enamel.
- Special Care reminder**
- Expert white mode:** It navigates your brushing with time cues to polish front teeth for effective stain removal.
- Gum care mode:** It navigates your brushing with gum massage intervals to strengthen gum health.

- Post-brushing feedback**
- Your Sonicare toothbrush tracks your brushing motions to provide different types of feedback:
- Coverage Feedback**
- It detects missed areas and reminds you to brush them again at the end of the brushing session. (Fig. 13)
- All lights on:** You brushed evenly and effectively to prevent plaque growth.
- Area light off:** This area received insufficient brushing time, which might lead to excessive plaque growth.
- All lights blink:** Insufficient brushing time or incomplete coverage detected, this may lead to excessive plaque growth.
- Pressure Feedback**
- It indicates if too much pressure was detected during the last brushing session. (Fig. 14)
- Light on:** Pressure was moderate.
- Light off:** Pressure was too heavy.
- Scrubbing Feedback**
- It indicates if aggressive side-to-side motions were detected in the last brushing session. (Fig. 15)
- Light on:** Gentle, controlled motions maintained.
- Light off:** Aggressive scrubbing detected.

- Celebration Feedback**
- If you brush consistently with moderate pressure, gentle motions, and complete coverage, you will hear a celebration comment and see the post-brushing feedback icons blink twice.
- SmartTimer**
- The SmartTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.
- Note: If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 30 seconds, the SmartTimer resets.

- QuadPacer**
- The QuadPacer is a 30-second interval timer that divides your mouth into 4 sections and uses progressive tone to guide you brushing equally in each section.
- Note: Additional time is added when you use the Expert White or Gum Care mode.
- Charging and battery status**
- Charging your Philips Sonicare**
- Place the power toothbrush handle on the charger (Fig. 5).
 - Connect the charger to a certified safety extra low voltage (SELV) supply unit with output rating 5V $\overline{\text{rms}}$ $\pm$ 2V.
 - Your toothbrush continues to charge until it is fully charged.
 - When the toothbrush is charging, the battery indicator pulsates green.
 - When fully charged, the battery indicator will light up solid green for 30 seconds before turning off.
- Note: It takes up to 10 hours to fully charge the battery, but you can use the toothbrush before it is fully charged.

- Battery status**
- This Philips Sonicare toothbrush is designed to provide at least 120 brushing sessions after fully charged, each session lasting 2 minutes. When the battery is low, the battery indicator starts to flash amber.
- Note: When the battery is completely empty, the toothbrush switches off. Place the handle on the charger to charge it.
- Note: To keep the battery fully charged at all times, you may keep your Philips Sonicare toothbrush on the charger when not in use.
- Note: Battery usage time is expected to reduce over product life time.

- Activating or deactivating features**
- Your toothbrush comes with the following features activated:
- Pressure Sensor Feedback
 - Scrubbing Sensor Feedback
 - Voice Feedback
- To deactivate or reactivate these features, follow the instructions below:
- Voice Feedback**
- Put the handle on the plugged-in charger.
 - Keep the power button pressed until you hear a single short beep and see all post-brushing feedback icons (Fig. 13, Fig. 14, and Fig. 15) light up (after 3 seconds).
 - Release the power button.
 - Triple tone of low-med-high means the Voice Feedback feature has been activated.
 - Triple tone of high-med-low means the Voice Feedback feature has been deactivated.

- Scrubbing Sensor Feedback**
- Put the handle on the plugged-in charger.
 - Keep the power button pressed until you hear a series of two short beeps and see the Scrubbing Feedback icon (Fig. 15) light up (after 5 seconds).
 - Release the power button.
 - Triple tone of low-med-high means the Voice Feedback feature has been activated.
 - Triple tone of high-med-low means the Voice Feedback feature has been deactivated.
- Pressure Sensor Feedback**
- Put the handle on the plugged-in charger.
 - Keep the power button pressed until you hear a series of three short beeps and see the Pressure Sensor icon (Fig. 14) light up (after 7 seconds).
 - Release the power button.
 - Triple tone of low-med-high means the Pressure Sensor Feedback feature has been activated.
 - Triple tone of high-med-low means the Pressure Sensor Feedback feature has been deactivated.

- Warranty and support**
- Philips offers a 2-year warranty. Warranty exceptions apply. This warranty does not affect your rights as a consumer under law. For more information, visit www.philips.com/support or call us.
- Warranty restrictions**
- The terms of the international warranty do not cover the following:
- Brush heads.
 - Damage caused by use of unauthorized replacement parts.
 - Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
 - Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- Recycling**
- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 18).
- Follow your country's rules for the separate collection of electrical products and batteries.

- Remove of built-in rechargeable battery**
- CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.**
- Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.**
- When you handle batteries, make sure that your hands, the product and the batteries are dry.**

- Wipe the entire surface of the handle with a damp cloth.
- Caution: Do not tap the handle on the sink or counter to remove excess water.**
- Dry the handle with a soft cloth or towel.
- Brush head**
- Rinse the brush head and bristles after each use and allow to air dry (Fig. 17).
 - Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

- Charger**
- Unplug the charger before you clean it.
 - Wipe the surface of the charger with a damp cloth.
 - Dry the charger with a soft cloth or towel.

- Brush head replacement**
- Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.
- Replace the brush head if signs of wear appear. Stop using a brush head with crushed or bent bristles. Smart fading bristles provide a visible replacement reminder: when the blue bristles fade to white at one-third of their length, replace the brush head.

- Storage**
- If you are not going to use the product for a long time, unplug it, clean it (see chapter "Cleaning") and store it in a cool and dry place away from direct sunlight.
- Locating the model number**
- Look on the bottom of the Philips Sonicare toothbrush handle for the model number HX568x.
- Troubleshooting**
- This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.	This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.	
My toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm).
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

- Warranty and support**
- Philips offers a 2-year warranty. Warranty exceptions apply. This warranty does not affect your rights as a consumer under law. For more information, visit www.philips.com/support or call us.
- Warranty restrictions**
- The terms of the international warranty do not cover the following:
- Brush heads.
 - Damage caused by use of unauthorized replacement parts.
 - Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
 - Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.
- Recycling**
- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 18).
- Follow your country's rules for the separate collection of electrical products and batteries.

- Remove of built-in rechargeable battery**
- CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.**
- Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.**
- When you handle batteries, make sure that your hands, the product and the batteries are dry.**

- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come in contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

- The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the appliance is disconnected from the wall socket and that the battery is completely empty.
- Remove the inner body from the handle.
 - Remove the rechargeable battery from the inner body.

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