

PHILIPS












sonicare

Compact Flosser
1000

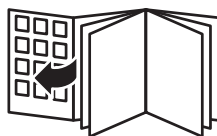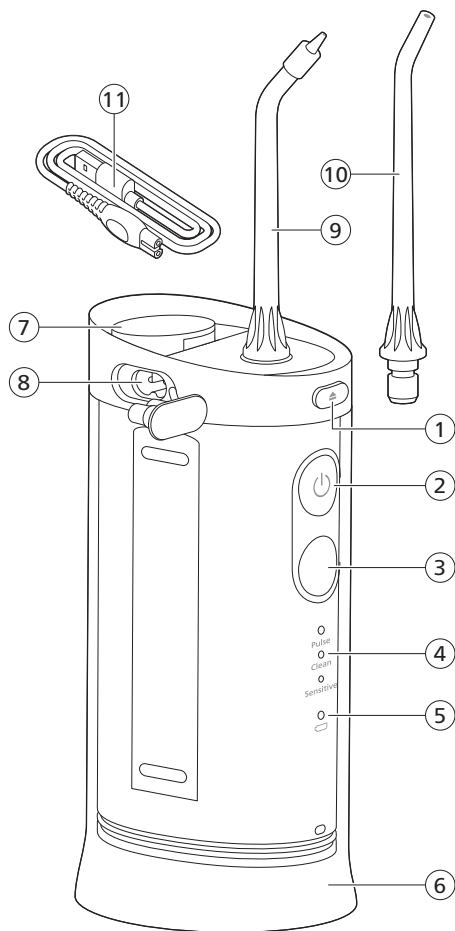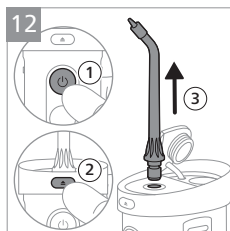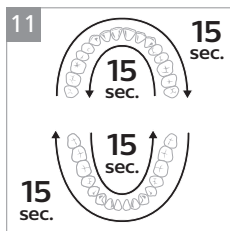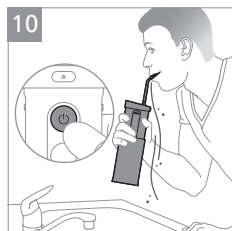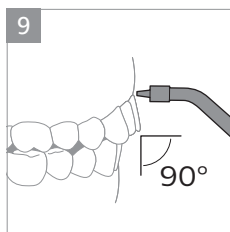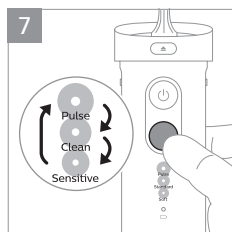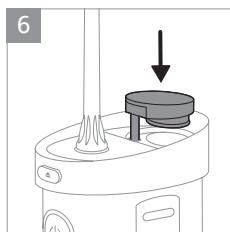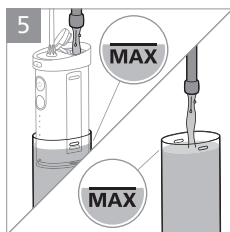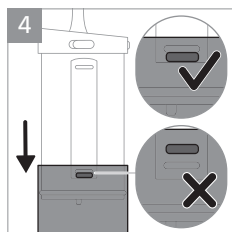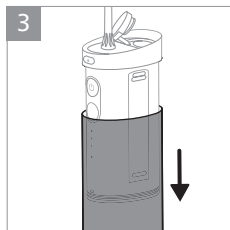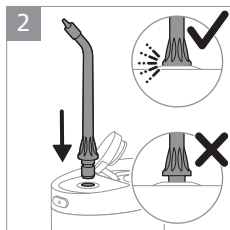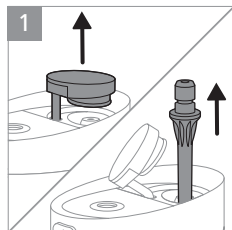English 6
Azərbaycanca 23
Čeština 41
Dansk 58
Deutsch 76
Eesti 97
Español 114
Français 132
Hrvatski 150
Italiano 168
Latviešu 186
Lietuviškai 203
Magyar 220
Nederlands 238
Norsk 256
O'zbek 273
Polski 291
Português 311
Română 329
Shqip 347
Slovenščina 365
Slovensky 383
Srpski 401
Suomi 419
Svenska 437
Türkçe 454
Türkmençe 471
Ελληνικά 489
Български 509
Кыргызча 529
Македонски 547
Русский 565
Тоҷикӣ 585

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Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

WARNINGS

- Keep the charging cable away from water.
- After cleaning, make sure the charging cable is completely dry before you connect it to the wall socket.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not charge the product outdoors or near heated surfaces.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country. If the USB cable or USB power supply unit is damaged stop using it.
- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- To avoid physical injury, do not fill the reservoir with hot water. Do not put foreign objects in the reservoir.
- This appliance has only been designed for cleaning teeth, gums.
- Do not use to clean eyes, ears or nose.
- The maximum altitude of use is 4500 meters.
- Water may accumulate on the counter or floor and present a slipping hazard if the device is not used over the sink.
- Clean the product (especially inside the reservoir) on regular basis.
- Charge only with Philips WAA1001, WAA2001, UL or equivalent safety approved adapters with input ratings (100 - 240V~; 50/60Hz; 3.5W) and output ratings (5Vdc; 2.5W). For the list of adapters, please see www.philips.com/support.

MEDICAL WARNINGS

- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- This appliance is a personal care device and is not intended to be shared by multiple users.
- Consult your dental professional before you use this product if you have had oral or gum surgery in the previous 2 months.

BATTERY SAFETY INSTRUCTIONS

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- Charge and use the product at a temperature between 3°C and 40°C (37 °F and 104 °F).
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- Power consumption in:
 - Off mode (no attachments) < 0.1 W

8 English

- Low energy mode < 0.3 W
- Charging: USB Charging (5VDC, 1A)

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Storage and transportation conditions

Temperature (Transport and Storage) between -20°C to 60°C (-4°F and 140 °F)

SAVE THESE INSTRUCTIONS

The following symbols may appear on the product:

This symbol means read the instructions for use, before using the appliance.

XXXXXX

This symbol means to use the indicated power supply. Detachable Power Supply Part: The "xxxxxx" indicates the model number of the approved chargers to be used.

Introduction

Welcome to Philips! Congratulations on your new Sonicare Compact Flosser. To fully benefit from the support that Philips offers, register your product at philips.com/welcome.

Intended use

The Philips Sonicare Compact Flosser is an oral irrigator intended to remove plaque and food debris between teeth and on the tooth surface in order to reduce tooth decay and improve and maintain gum health. This

product is intended for consumer home use, normally used at least once a day.

Your Philips Sonicare Compact Flosser

- 1 Nozzle release button
- 2 Power on/off button
- 3 Mode button
- 4 Mode indicator light
- 5 Battery indicator light
- 6 Reservoir
- 7 Filling port/Nozzle storage
- 8 Charging port


Accessories

- 9 Comfort nozzle
- 10 Standard nozzle
- 11 USB-A to small plug (wall adapter not included)


Note: Accessories included may vary based on the model purchased.

Getting started


Set up your Compact Flosser

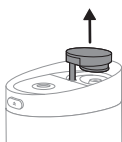- 1 Gently pull on the cover of the filling port.

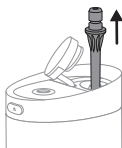
- 2 Remove the nozzle from the storage compartment.

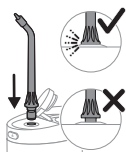
- 3 Insert your nozzle into the device.

Note: Ensure the nozzle is fully inserted and flush with the top of the device.

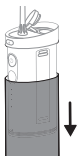- 4 Ensure the cover of the filling port is open and pull down the reservoir.


Note: When the filling port cover is open, it is easier to pull down the reservoir.

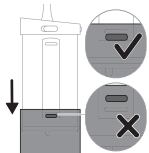
- 5 Pull down the reservoir until the device's stopper is above the stopper inside the reservoir.

Filling the reservoir


Warning: Remove the charger from the Compact Flosser prior to filling the reservoir and using the device.

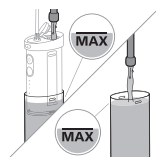- 1 Open the filling port and fill the reservoir with lukewarm water.
Or

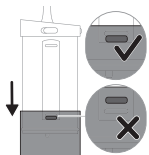
Water can also be filled directly into the reservoir. Gently pull the reservoir off the Compact Flosser and fill with water to the max fill line. Reattach the reservoir and make sure the stopper of the device and reservoir are correctly in place.

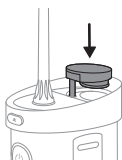2 Gently close the filling port.

Note: When filling reservoir, do not exceed the max fill line (200ml).

Modes

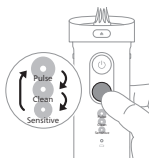Your Compact Flosser is equipped with three different modes. To toggle between modes, press the mode button.

The modes have different intensities. Listed from low to high:

Sensitive

This mode, with pacer, provides gentle water pressure. The pacer indicates when you are supposed to move to the next section of your mouth via a short pause in the water flow every 15 seconds.

Pulse

This mode, alternates strong and gentle pulses, and is designed to massage and stimulate the gums. The pacer indicates when you are supposed to move to the next section of your mouth via a short pause in the water flow every 15 seconds.

Clean

This mode, with pacer, provides constant water pressure for a great clean. The pacer indicates when you are supposed to move to the next section of your mouth via a short pause in the water flow every 15 seconds.

Note: For first time use, start the device in Sensitive mode for the lowest intensity. It is recommended to increase the intensity you use over the course of a few weeks. This can give better cleaning benefits as you are able to tolerate higher intensities.

Note: This Compact Flosser comes with the pacer activated. To activate or deactivate the pacer with the Compact Flosser off, press and hold the mode selection button for 3 seconds. When deactivated, the battery indicator light will flash amber once. When activated, the battery indicator light will flash green twice.


Using your Compact Flosser

- 1 Before turning on the device, point the nozzle tip at the sink. Allow water to spray into the sink. Turn off the device once water dispenses from the device.


- 2 Place the nozzle tip in your mouth. Lean over the sink and partially close your lips over the nozzle tip to avoid splashing. Leave enough space to allow water to flow from your mouth into the sink.

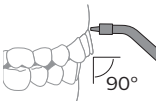- 3 Position the nozzle tip just above the gum line at a 90° angle.


- 4 Press the power on/off button on the handle to turn on the Compact Flosser and to start flossing.

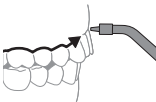
- a Starting from your back teeth in your upper or lower jaw, slide the nozzle tip along the gum line and pause briefly between each tooth.

Note: Rotate the nozzle as desired for easier positioning along the gumline.

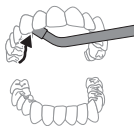- b Follow the pacer to clean both the outside (15 seconds) and the inside of your teeth (15 seconds) for a total of 30 seconds in one arch of your mouth.

Note: Remember to clean the insides of your teeth.

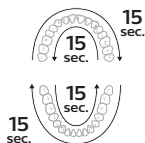
- c Follow the pacer and repeat in the other arch of your mouth for another 15 seconds on the inside and 15 seconds on the outside until all spaces around and between teeth have been cleaned for a total of at least one minute.

Note: You can stop the water flow during use by pressing the on/off button on the Compact Flosser.


Note: If you choose to refill the reservoir during use, turn the device off while refilling the reservoir and then pick up where you left off until all spaces around and between your teeth have been cleaned for a total of at least one minute.

After each use

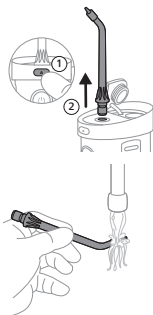- 1 Remove the nozzle by pulling the nozzle and holding the nozzle release button at the same time.
- 2 Rinse the tip under warm water.

14 English

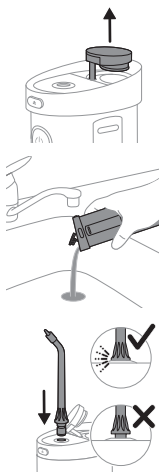3 Gently open the filling port on top of the Compact Flosser.

4 Pour out any excess water. Dry the nozzle, the outside of the handle and the reservoir.

5 Insert the nozzle back into the handle.

Note: To facilitate drying, keep the filling port open or leave the Compact Flosser disassembled between uses (the reservoir can be pulled off completely).

Note: Do not leave the Compact Flosser unassembled longer than necessary.

Charging and battery status


Your Sonicare Compact Flosser is designed to provide at least 24 flossing sessions, each session lasting 1 minute on a fully charged battery. The device will shut off after 2 minutes of use.

Note: It can take up to 4 hours to fully charge the battery before first use, but you can use the Compact Flosser before it is fully charged.

Charging


1 Gently pull the charging port cover on the side of the Compact Flosser to expose the charging port.

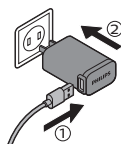- 2** Connect the USB-A to small plug charging cable into a 5V USB wall adapter.

Note: Wall adapter is not included in all models. Only use a 5V adapter (Type A).

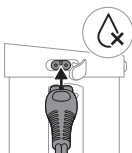- 3** Plug the USB-A to small plug charging cable into the charging port.

Note: Make sure the handle is completely dry before you connect the USB-A to small plug charging cable.

- 4** The battery indicator flashes while charging.
- 5** Close the charging port cover on the Compact Flosser when not actively charging.

Battery status

Battery Status

The battery indicator shows the battery status when you complete the flossing session, or when you pause the Compact Flosser.

When the Compact Flosser is charging

Battery status	Battery indicator
Full	Solid green
Charging	Flashing green

Note: The battery indicator will flash until fully charged. When it becomes fully charged while plugged in, the light will turn solid green for 30 seconds and then it will turn off.

When the Compact Flosser is not charging

Battery status	Battery indicator
Full	Solid green
Medium	Flashing green
Low	Flashing amber
Empty	Flashing amber for 5 seconds and the Compact Flosser turns off automatically


Note: When the low battery indicator appears, the Compact Flosser is designed to provide at least 3 flossing sessions, each session lasting 1 minute.

Cleaning


The Compact Flosser needs to be cleaned once a week. Clean the device and the accessories as follows. Accessories may vary by model.

Note: Do not run or clean device with high alcohol content cleaning solutions.


Warning: Do not clean the Compact Flosser while the USB-A to small plug charging cable is attached. Unplug the device and close the charging port cover before cleaning.

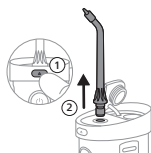
- 1 Remove the nozzle by pulling the nozzle and holding the nozzle release button at the same time.

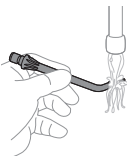
- 2 Rinse the tip under warm water.

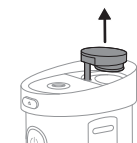- 3 Gently open the filling port on top of the Compact Flosser.

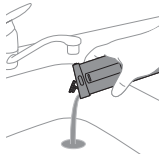- 4 Pour leftover water in the sink.

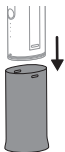
- 5 Gently remove the reservoir from the Compact Flosser by pulling it off. The reservoir can be hand washed with warm water and mild soap. When necessary, you can clean the handle by wiping with a soft cloth and mild soap.


Storage

Your Compact Flosser is designed such that you can store one nozzle in the handle. To use this functionality:


- 1 Remove the nozzle by pulling the nozzle and holding the nozzle release button at the same time.

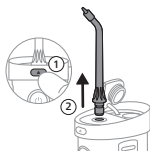
- 2 Gently open the filling port on top of the Compact Flosser.

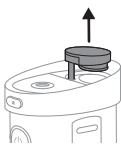- 3 Insert the nozzle into the filling port.

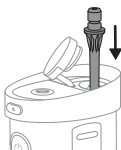
- 4 Collapse the Compact Flosser by pushing the reservoir up.

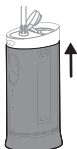Note: When the filling port cover is open, it is easier to collapse the reservoir.

- 5 Gently close the filling port.

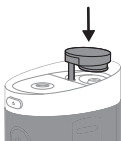Replacement

Replace nozzles every 6 months or immediately when damaged. To buy accessories or spare parts, visit **philips.com** or go to your local Philips Sonicare dealer.

Note: The Compact Flosser is only compatible with Philips Sonicare Compact Flosser Nozzles. It is not compatible with other Philips Sonicare nozzles such as those for the Power Flosser.

Troubleshooting

Question	Answer
Can I use mouthwash or essential oils in my reservoir?	It is not recommended to use essential oils in the reservoir. Mouthwash can be mixed with water and used in the reservoir. After using mouthwash, rinse the device to prevent clogging by partially filling reservoir with warm water and running product with tip pointed into sink until reservoir is empty.
Will other nozzles fit on my Philips Sonicare Compact Flosser?	Only use Philips N1 and N2 nozzles for this device. Other Philips or non-Philips nozzles are not compatible with this product.
How often do I need to replace my nozzle?	We advise you to replace your nozzle after 6 months of use. You can find more information about buying replacement parts at philips.com
I'm having difficulty extending the reservoir	Open the filling port and remove the nozzle before attempting to extend the reservoir. This should allow the reservoir to easily extend.

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in

the qualifying period is provided. Philips recommends use of genuine Sonicare nozzles for optimal performance. Use of other nozzles may damage your Compact Flosser and limit your warranty protection. Contact our Consumer Care Center at 1-800-682-7664 or **www.philips.com/sonicare-warranty** (North America), outside North America contact your local Philips Consumer Care Center. For product information and online user instructions, visit: **<https://www.usa.philips.com/c-w/support-home/sonicare>** (North America) or **www.philips.com/support** (outside North America).

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Nozzles
- Damage caused by use of unauthorized replacement parts or unauthorized nozzles
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE

EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.
- Information on REACH substances, contained in Philips products, can be found on www.philips.com/REACH


Removing the rechargeable battery

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Instructions for removal of built-in rechargeable batteries can be found on **www.philips.com/support**.


- This appliance contains Lithium-ion rechargeable cell. **MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

To remove the rechargeable battery, you need a screwdriver and a wire cutter. Observe basic safety precautions when you follow the procedure outlined below.


- 1 To deplete the rechargeable battery of any charge, let the Compact Flosser run until it stops.

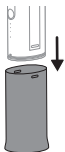
- 2 Remove the reservoir by pulling it off.

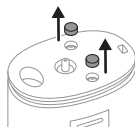
- 3 Remove the two screw caps.

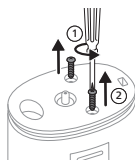
- 4 Remove the screws.

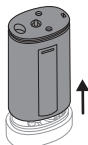- 5 Pull the internal components from the body.

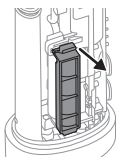- 6 Remove the battery cover.

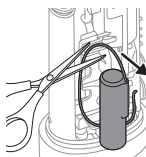- 7 Pull the battery from the cradle. Cut wires to remove the battery.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

- 9 The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

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©2025 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.
Philips Oral Healthcare, LLC
Bothell, WA 98021 U.S.A
Philips Electronics UK Limited