PHILIPS **sonicare**

4100 - 5100

— |

| —

— |

| —

English 6
Azərbaycanca 9
Čeština 12
Dansk 15
Deutsch 18
Eesti 21
Español 24
Español 27
Français 30
Français (Canada) 33
Hrvatski 36
Italiano 39
Latviešu 42
Lietuviškai 45
Magyar 48
Nederlands 51
Norsk 54
O'zbek 57
Polski 60
Português 63
Română 66
Shqip 69
Slovenščina 72
Slovensky 75
Srpski 78
Suomi 81
Svenska 84
Tiếng Việt 87
Türkçe 90
Türkmençe 93
Ελληνικά 96
Български 99
Кыргызча 102

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Removing the rechargeable battery

The built-in rechargeable battery can be removed for disposal. If you are not comfortable disassembling electronics or if otherwise instructed in the product manual, please seek professional assistance at your local electronics recycling or repair shop.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

Be sure to protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the toothbrush, and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

Do not reinstall battery once it has been removed for recycling. Consumer installation of battery may increase risk of fire or other potential safety hazards.


If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Removing the rechargeable battery


Before removing the battery, make sure that the battery is completely empty.

To remove the rechargeable battery, you will need these tools: A T8 or T10 torx screwdriver, tweezers, a plastic pry tool (or guitar pick), and a hard surface such as a table.


- 1 Deplete the rechargeable battery of any charge. To do this, remove the toothbrush from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare toothbrush.
- 2 Remove the brush head from the toothbrush.

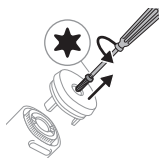
- 3 Using a torx screwdriver, unscrew and remove (counterclockwise) the screw from the base of the toothbrush handle.

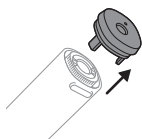- 4 Remove the endcap from the base of the toothbrush. If the endcap is tight and does not release easily from the housing tube, use a plastic pry tool to help separate the endcap from the housing tube.

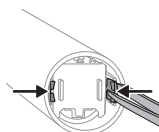
- 5 Use tweezers or a plastic pry tool to pry the two retention tabs away from the inside walls of the handle.

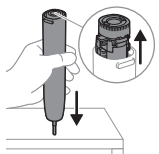- 6 Once the tabs are away from the wall, release the internal assembly from the housing tube by inverting the handle and firmly pressing the metal shaft onto a solid surface (do not slam the handle down with force). If the internal components do not easily release from the housing tube, double check that you have completed step 5 above.
- 7 Gently pull the internal assembly to remove it from the base of the housing tube.

Caution: To avoid injury or damage to the product, be aware of the sharp edges and exposed electronic components.

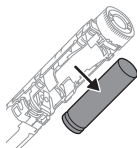- 8 Use a plastic pry tool to push or pry the battery out of the battery tray.

Caution: Do not use screwdriver or metal object to pry the battery out of the battery tray.

8 English

- 9 For battery disposal, cover the battery contacts of the used battery with tape to prevent any electrical short from residual battery charge. The used rechargeable battery can now be recycled through your community battery or electronics disposal process.

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©2025 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.