

PHILIPS

sonicare

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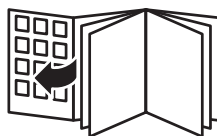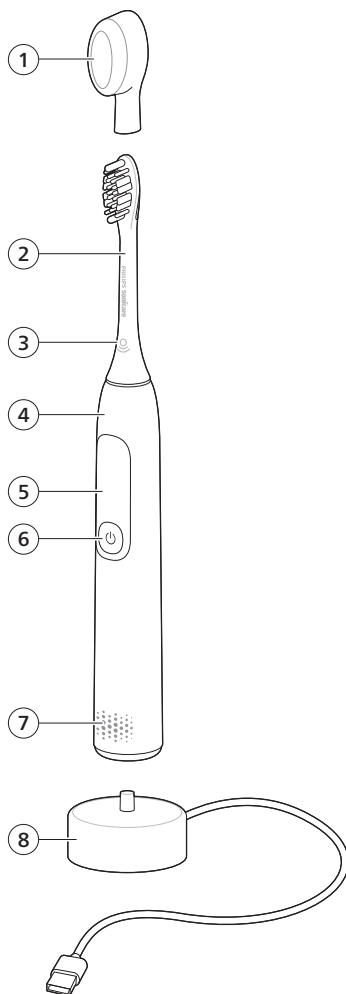English 6
Español 25
Français (Canada) 46

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **www.philips.com/welcome**.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.

Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country.
- If the USB cable or USB power supply unit is damaged, stop using it.
- Keep the charger away from water.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not use the charger outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- After cleaning, make sure the charger is completely dry before connecting to power.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.

- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.

Medical warnings

- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended to be shared by multiple users.

Battery safety instructions

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Disposal.

- To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5V, $\geq 1.5\text{W}$. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country (see the international warranty leaflet for contact details). Using a non-certified supply unit may cause hazards or serious injuries.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).

Storage and transportation conditions

Temperature -20°C to +60°C (-4 °F to 140 °F).

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC Compliance

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).
- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1 This device may not cause interference.
- 2 This device must accept any interference, including interference that may cause undesired operation of the device.


SAVE THESE INSTRUCTIONS

LIST OF SYMBOLS USED

The following symbols may appear on the product:
This symbol means read the instructions for use, before using the appliance.

This symbol means to use the indicated power supply.

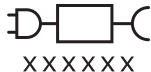Intended use

This Sonicare power toothbrush is intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. This Sonicare power toothbrush is intended for consumer home use. Use by children should be with adult supervision.


Your Philips Sonicare (Fig 1.)

- 1 Brush head cap
- 2 Smart brush head
- 3 BrushSync symbol
- 4 Handle
- 5 Smart display
- 6 Power/mode button
- 7 Light ring
- 8 USB charger (wall adapter not included)

Note: Accessories included may vary based on the model purchased.

Getting started

Attaching the brush head

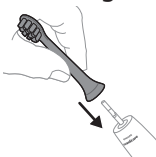- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using

a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gradually increases the power over the first 14 brushing sessions to help you adapt to brushing with your Sonicare toothbrush.

Brushing instructions


- 1 Wet the bristles and apply a small amount of toothpaste.


- 2 Press the power button to turn on your Philips Sonicare. Then select a brushing mode.

Note: The first 3 seconds of all modes start with a gentle vibration. This prevents toothpaste from splattering and lets you cycle through the modes to select your preferred brushing mode.

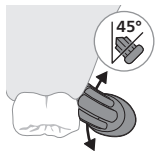
- 3 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

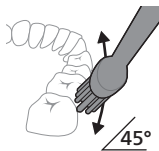- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

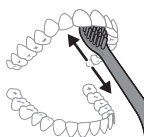
- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

Note: Additional time is added when you use the White mode to focus on the front section of your teeth.

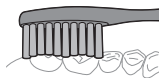- 7 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

Note: When using the Philips Sonicare toothbrush in clinical studies, ensure the handle is fully charged and the EasyStart feature is deactivated.


Modes

Your power toothbrush is equipped with two different modes: Clean and White.

Based on the smart brush head you use; the optimal mode will be selected automatically by the BrushSync technology.

To manually change the mode, press the power button when the toothbrush is on. All Philips Sonicare brush heads work with each mode.

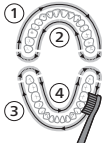Mode	Clean	White
Mode image		
Benefit	Plaque removal	Plaque and surface stains removal
Total brushing time	2 minutes	2 minutes and 30 seconds
How to brush	Brush each segment for 30 seconds 	Brush each segment for 30 seconds,  then brush the upper and lower front teeth for 15 seconds each.

BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The  symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- Brush head mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to www.philips.com/toothbrush-heads for more information.

Features

- SmarTimer
- Brushing reports
- QuadPacer
- EasyStart
- Brush head replacement reminder
- Brush head mode pairing

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 15 seconds, the SmarTimer resets.

Brushing reports

Your Philips Sonicare regularly provides a report on your brushing behavior.

After brushing the display shows a report with the brushing time.


- Brushing time of 2 minutes or more

- Brushing time of less than 2 minutes
- During brushing the display shows the current brushing mode and brushing time

QuadPacer

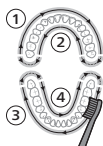The QuadPacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. This Philips Sonicare toothbrush comes with the QuadPacer feature activated.

EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

To deactivate EasyStart (see 'Activating or deactivating features').

Brush head replacement reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When a new smart brush head is attached for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and the display shows the brush head symbol.
- 2 Over time, based on the pressure you apply and how long the brush head has been used, the handle will track brush head wear in order to determine the optimal time to change your brush head. When the life of the brush head ends, the display shows the brush head replacement reminder symbol and indicates that you should replace the brush head as soon as possible.

Note: Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Brush head mode pairing


Brush heads with BrushSync technology are equipped with a microchip which communicates with the handle and automatically pairs it with the recommended mode and intensity. Most brush heads with BrushSync technology will auto pair to Clean Mode. Brush heads with specific benefits will pair with the corresponding mode (when available), examples include:
W Optimal White: White Mode

Note: If you choose to change mode from the recommended setting, the brush head will remember the chosen setting for future brushing sessions.

Charging and battery status

Charging your Philips Sonicare

- 1 Attach the USB charger to a wall adapter and plug the wall adapter into an electrical outlet.
- 2 Place the handle on the plugged-in charger.

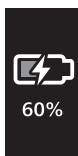- The battery icon with lightning symbol is shown on the display and you hear two beeps. The light ring at the bottom of the handle will also light up white for 5 seconds.

Note: It can take up to 16 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.

Note: This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day).

Battery status

The LCD screen will show the battery level after each brushing session, when placing the handle on the charger, and when removing it from the charger. While charging, press the power button to check the battery level.

The battery icon changes color with different battery levels: it turns yellow if the battery level is below 40%; while low battery is indicated by a red icon flashing with two short beeps.

Battery status	Battery icon color	Battery percentage
Low	Red	0% - 20%
Medium	Yellow	21% - 40%
Full	Blue	41% - 100%

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder

To activate or deactivate these features, follow the instructions below:

EasyStart

- 1 Place the handle on the plugged-in charger.
- 2 Keep the power button pressed until you hear a single short beep (after 3 seconds).

3 Release the power button.

- Triple tone of low-medium-high means the EasyStart feature has been activated. The icon on the display shows that the EasyStart has been activated.

- Triple tone of high-medium-low means the EasyStart feature has been deactivated. The icon on the display shows that the EasyStart has been deactivated.

Note: To achieve clinical efficacy, EasyStart needs to be deactivated.

Brush head replacement reminder

Note: You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

- 1 Place the handle on the plugged-in charger.
 - 2 Press and hold the power button while the handle remains on the charger.
 - 3 Keep the power button pressed until you hear two short beeps (after 6 seconds).
 - 4 Release the power button.
- Triple tone of low-medium-high means the brush head replacement reminder feature has been activated. The icon on the display shows that the brush head replacement reminder has been activated.
 - Triple tone of high-medium-low means the brush head replacement reminder feature has been deactivated. The icon on the display shows that the brush head replacement reminder has been deactivated.


Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure

to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

Warning: Do not clean the product or accessories with sharp objects or in the dishwasher, microwave, with chemicals or in boiling hot water

Toothbrush handle

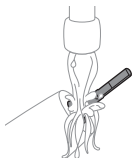- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste on the handle and around the buttons at least once a week, but preferably after each use.


Caution: Do not push on the rubber seal on the metal shaft or around the buttons with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Caution: Do not tap the handle on the sink or counter to remove excess water.

- 3 Dry the handle with a soft cloth or towel.

Brush head

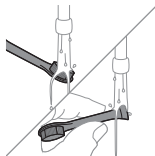- 1 Rinse the brush head and bristles after each use and allow to air dry.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Storage

If you are not going to use the product for a long time, unplug it, and clean it (see chapter 'Cleaning'). Store it in a cool and dry place away from direct sunlight. Charge the toothbrush regularly to prevent the battery from discharging completely. A completely discharged battery may fail to charge again.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (e.g. HX123A).

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm). When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

Problem	Possible Cause	Solution
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided.

Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection.

Contact our Consumer Care Center at 1-888-744-5477 (North America), outside North America contact your local Philips Consumer Care Center. Internet information: **www.philips.com/support**.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH

ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Disposal

CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

- Contact your local town or city officials for electronics recycling and battery disposal information. You can also call 1-800-822-8837 or visit **www.call2recycle.org** for drop-off locations.
- For assistance, visit our website **www.philips.com/support** or call 1-888-744-5477 toll free.
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit: **www.recycle.philips.com**.
- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

**CONTAINS LITHIUM-ION RECHARGEABLE CELL.
MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

When you handle batteries, make sure that your hands, the product and the batteries are dry.


If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.


Before removing the battery, make sure that the battery is completely empty.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Remove the film at the bottom of handle, then release the screw.

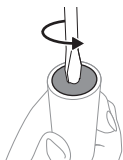
- 3 Remove the end cap from the toothbrush handle.

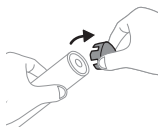
- 4 Holding the handle upside down, press the shaft down on a hard surface.

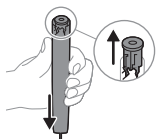
- 5 Wedge the screwdriver between the battery and the white frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the white frame.

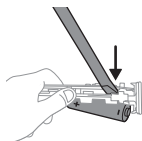- 6 Insert the screwdriver between the bottom of the battery and the white frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

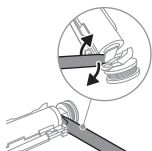- 7 Grab the battery and pull it away from the internal components to break the second metal battery tab.

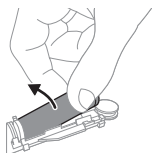Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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©2025 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.
Philips Oral Healthcare, LLC.