
# PHILIPS

## LUMEA

IPL HAIR REMOVAL 7000, 8000, 9000

BRI931, BRI932, BRI937, BRI984

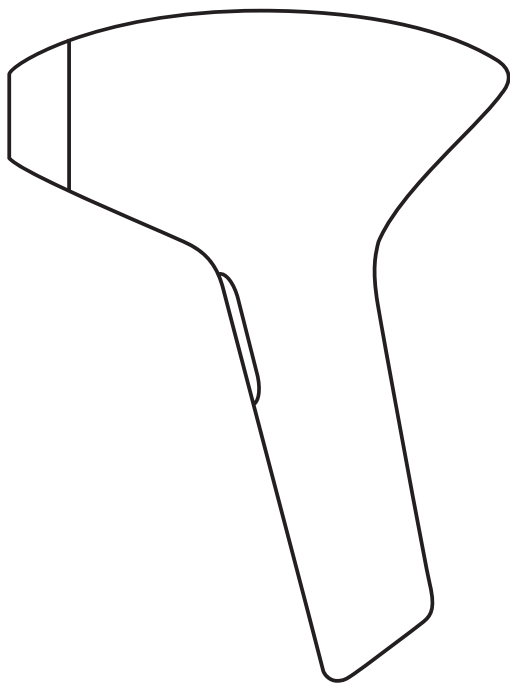**Skin tones – Suitable**

**Teints de peau – Compatibles**

Fair  
Pâle

Light  
Clair

Beige

Light  
brown  
Brun  
clair


Medium  
brown  
Brun  
moyen

Dark  
brown  
Brun  
foncé

Dark  
Foncé

**Body hair colors – Suitable**

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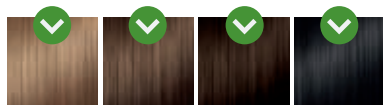Dark  
blond  
Blond  
foncé

Brown  
Brun

Dark  
brown  
Brun  
foncé

Black  
Noir

**Unsuitable**

**Non compatible**

Light  
blond  
Blond  
clair

Red  
Roux

White/  
Grey  
Blanc/  
Gris

English	5
Français (Canada)	39

# Chapters

Welcome to Lumea	6
Indications for use	6
Important safety instructions	7
Contraindications	7
Warnings	8
Cautions	10
Risks and benefits of using Lumea	11
Potential side effects of IPL treatment	12
What is IPL?	15
Lumea overview	16
SenseIQ technology	17
Attachments	18
Treatment modes	19
Your treatment guide	20
Prepare for treatment	21
Test your skin	22
Perform a treatment	24
Clean and store	26
After treatment	27
Plan your next treatments	28
Troubleshooting	29
Additional information	31

# Welcome to Lumea

Achieve gentle, fast and effective hair reduction with Lumea. You can conveniently use Lumea in the comfort of your home. Are you ready to get started?

Scan with your phone camera, QR app, or visit [www.philips.com/myIPL](http://www.philips.com/myIPL)

Note: In this user manual, Lumea refers to the Philips Lumea IPL 7000 Series, 8000 Series and 9000 Series.

## Indications for use

Philips Lumea IPL is indicated for the removal of unwanted hair for a single user. It is also indicated for the permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

### Is Lumea suitable for you?

Philips Lumea is not suitable for everyone. Check the skin tone and hair color tables at the start of this user manual and read the important safety instructions carefully to see if Lumea is suitable for you.

Philips Lumea can be used by persons in the age range of 18-65 years with fair up to medium brown skin tones, and naturally dark blond, brown, dark brown and black hair. Unsuitable body hair colors: Treatment with Lumea is not effective if you have light blond, grey, red or white hairs as hairs with these colors do not absorb enough light.

Unsuitable skin tones: You cannot use Lumea if you have dark brown to dark skin tones as these skin tones absorb light more easily which could result in pain, skin redness or burns.

Philips Lumea, attached with a face attachment, can be used on the female face (below the cheekbone). Philips Lumea, attached with a body or precision attachment, can be used on the body (below the neckline). Philips Lumea is a handheld device intended for home use.

# Important safety instructions

**⚠ ALWAYS** read and follow the instructions carefully to fully benefit from the use of your Lumea and to minimize the risk of injury, skin reactions and side effects. Make sure to save these instructions for future reference.

## **⚠ Contraindications**

- **DO NOT** use Lumea if you have had chemotherapy or radiation therapy in the last 3 months.
- **DO NOT** use Lumea if your sensitivity to pain is reduced after using medication (e.g. painkiller, numbing cream) or due to a medical condition.
- **DO NOT** use Lumea if you have a dark brown or darker skin tone.
- **DO NOT** use Lumea if you are pregnant, lactating or breastfeeding.
- **DO NOT** use Lumea if you are under the age of 18 years or above 65 years.

### **Light sensitivity:**

- **DO NOT** use Lumea if your skin is extra sensitive to light, for example if you:
  - have a sun allergy;
  - are taking certain medication or herbal remedies. Check the instruction leaflet inside your medication package to see if what you are taking can cause increased sensitivity to sunlight.

### **Medical conditions:**

- **DO NOT** use Lumea – or consult your doctor first – if you have:
  - epilepsy with flashlight sensitivity;
  - an active skin cancer or a history of skin cancer;
  - a skin pigmentation disorder such as vitiligo;
  - a history of poor wound healing;
  - an immunocompromised state caused by a disease such as AIDS or by immunosuppressive drugs;
  - a history of collagen disorders, including a history of keloid scar formation.

## Warnings

### Thermal safety - To avoid pain, skin redness or burns

- **DO NOT** use Lumea on lips, nostrils, ears, scalp, nipples, labia minora, perineum, around the anus, penis, or scrotum. These areas are highly sensitive and have thinner and darker skin.
- **DO NOT** touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.
- **DO NOT** use Lumea on an active skin condition, as this may worsen the condition or cause additional pain, skin redness or burns. Skin conditions may include:
  - Active bruises, rashes, sores, open wounds, scars or sunburns. Wait until the injury is completely healed before resuming treatment.
  - Active skin irritations, inflammations or infections. Wait until the skin is completely recovered before resuming treatment.
  - Damaged skin after using depilatory (hair removal) creams. Wait at least 2 days before resuming treatment.
  - Damaged skin after using skin peeling, (micro)dermabrasion or skin rejuvenation methods. Wait at least 1 week before resuming treatment.

### Eye safety - To avoid eye damage

- **DO NOT** use Lumea close to or on the eye, eyelid or eyebrow. Use Lumea only below the cheekbones.

### Electrical safety - To avoid electrical shock

- **DO NOT** use Lumea in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).
- **DO NOT** let water run into Lumea.
- **DO NOT** expose Lumea to large temperature changes, especially in a high humidity environment, without giving Lumea time to acclimatize for 30 minutes to the ambient temperature of the room in which it will be used (see 'Technical specifications').
- **DO NOT** attempt to repair or modify Lumea. Any repairs or alterations made by the user may void the warranty and compliance of Lumea. Changes or modifications made to Lumea not expressly approved by Philips may void the user's authority to operate this equipment.
- **DO NOT** touch any inner part of a damaged Lumea or the power adapter.
- **DO NOT** use any other attachments or power adapter than supplied with Lumea.
- **ALWAYS** check if the voltage indicated on the power adapter corresponds to the local mains voltage before you connect the equipment.

### General safety - To avoid serious injury

- **ALWAYS** check Lumea before you use it. Do not use Lumea or the power adapter if it is damaged. Always replace a damaged part with one of the original type.
- **ALWAYS** unplug Lumea, including the power adapter, after use and in case of a power failure to prevent damage to the equipment.
- **DO NOT** use Lumea directly over or close to any active implants such as a:
  - pacemaker;
  - neurostimulator;
  - insulin pump.Check the active implants' manual for the safe distance to other electronic devices.
- **DO NOT** wear the power adapter cord around your neck to avoid strangulation.
- **DO NOT** use Lumea if you have reduced physical, sensory or mental capabilities as this could result in incorrect usage.
- **DO NOT** use Lumea, including the power adapter cable, within 12 inches (30 cm) of Wi-Fi and/or Bluetooth devices (portable RF communication equipment) or smartphones as it may cause Lumea to malfunction.
- **DO NOT** use Lumea next to or stacked with other equipment. This could result in improper operation.
- **DO NOT** leave Lumea unattended when it is turned on to avoid the risk of fire.
- **KEEP** Lumea out of the reach of children. Children should be supervised to ensure that they do not play with the equipment.
- **ALWAYS** check that the power adapter is connected to a properly installed and working power outlet to avoid the risk of fire.

## Cautions

### Thermal safety - To avoid pain, skin redness or burns

- **DO NOT** start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills). Do a skin test before starting/resuming your treatment, see section 'Test your skin'.
- **DO NOT** use Lumea with the body or precision attachment on the upper lip or other areas of the face. Only use the face attachment as it has an extra integrated filter for treatment on face below the cheekbones.
- **DO NOT** continue the treatment if you experience pain or a skin reaction that lasts longer than 24 hours. Reduce the light intensity setting or stop using Lumea.
- **DO NOT** expose your skin to sunlight unprotected when using Lumea. During the entire IPL treatment period, and especially in the first 48 hours after treatment, protect your skin from sunlight with clothing or apply sun-blocking products of not less than SPF50 on the non-covered treatment areas.
- **DO NOT** flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase performance.
- **DO NOT** use Lumea on skin that is darker than the surrounding treatment area. These highly pigmented skin areas can lead to excessive heating of the skin. Darker areas may include:
  - dark spots such as moles, warts, large freckles, birthmarks, age spots or dark veins;
  - marks added to the skin, such as tattoos or (permanent) makeup.
- **DO NOT** use Lumea on male neck and face. The hairs in these areas have different characteristics.
- **DO NOT** mark your skin/treatment area with any type of pencil or pen.
- **DO NOT** use Lumea if there is any lotion, gel or visible hair on your skin or on the device. Only treat skin that is clean, dry and hair-free. Make sure Lumea is clean and dry as well.
- **DO NOT** use Lumea if the light exit window or window frame is damaged. Check Lumea for damage before use.
- **DO NOT** subject Lumea to heavy shocks to avoid risk of damage to the lamp.
- **DO NOT** treat more often than once every 2 weeks.

### Eye safety - To avoid eye discomfort

- **ALWAYS** use Lumea in a well-lit room. Your eyes are more sensitive to stray light in darker conditions.

### General safety - To avoid injury

- **DO NOT** use Lumea over vascular skin disorders (e.g. varicose veins).
- **ENSURE** that the power adapter can be easily disconnected from the electrical outlet in case of a malfunction. The power adapter is considered a disconnection device.
- **DO NOT** share Lumea with other persons for hygienic reasons. Lumea is intended for single person use only.

## Risks and benefits of using Lumea

The safety and effectiveness of Lumea was tested in a clinical study involving 94 women (aged 18-64 years). None of the subjects left the study early due to problems associated with the use of Lumea.

Out of the 94 subjects, 31 experienced one or more mild and transient side effects: pain (26), itchiness (8), redness (5), folliculitis (2), burning sensation (1), superficial crusting (1) and swelling (1).

Results after 4 bi-weekly treatments:

- 75 women received treatments on their upper lip and achieved up to 61.1 % hair reduction.
- 76 women received treatments on their armpit area and achieved up to 60.6 % hair reduction.
- 77 women received treatments on their bikini area and achieved up to 61.9 % hair reduction.
- 76 women received treatments on their legs and achieved up to 93.1 % hair reduction.

35 women completed 8 additional touch-up treatments (once every 4 weeks) and achieved up to 67.0 % hair reduction over all treated areas.

# Potential side effects of IPL treatment

The table below lists known side effects that may occur following treatment with IPL devices.

<b>Common side effects</b>	<b>What is observed?</b>	<b>How long does it last and how can it be helped?</b>
Pain/ uncomfortable sensation	Local pain or discomfort on the treated area during or after treatment	If your treatment feels painful/uncomfortable, reduce the light intensity setting or stop using Lumea. Pain or discomfort after treatment can last up to 1 week.
Itching sensation	Perceived itching in the treated area	These reactions typically resolve within hours but can last up to 1 week. You can cool the area with a cold pack.
Warm/burning sensation	Local warm or burning feeling on the treated area	
Skin redness (erythema)	Redness without a clear edge in the treated area	This reaction typically resolves within hours but can last up to 1 day. You can cool the area with a cold pack. Skin redness lasting longer than 1 day is an uncommon side effect of IPL and a sign of a skin burn. Wait with the next treatment until the skin has healed completely. Make sure you use a lower light intensity for next treatments.

<b>Uncommon side effects</b>	<b>What is observed?</b>	<b>How long does it last and how can it be helped?</b>
Swelling (edema)	Swelling around hair follicles appears as bumps in the treated area, sometimes accompanied by redness.	This reaction typically resolves within hours but can last up to 1 week. You can cool the area with a cold pack.
Skin dryness	Skin appears rough and flaky in the treatment area.	This reaction typically lasts between 1 day and 1 month. You can apply a non-scented moisturizer to the treated area.
Burns and blisters	Extensive skin redness, swelling and possibly blisters, accompanied by pain	These reactions can last between 1 day and 1 month and may require medical treatment or medication. We advise you to consult a doctor.
Superficial crust	A sharply defined brownish and crusty area that can occur on darker skin tones	Wait with the next treatment until the skin has healed completely. Make sure you use a lower light intensity for next treatments.
Folliculitis	Swelling and redness around hair follicles, accompanied by itching and/or pustule formation	Folliculitis may require medical treatment or medication. We advise you to consult a doctor. Folliculitis may last several days up to several weeks. Wait with the next treatment until the skin has healed completely. Make sure you use a clean razorblade when shaving before treatment.
Pigmentary changes/skin discoloration	Skin of the treated area appears darker or lighter than the surrounding skin.	This reaction typically lasts several weeks and, in rare cases, up to a year. Do not treat areas with pigmentary changes/skin discolorations until your skin has regained its normal skin tone. Make sure you use a lower light intensity for next treatments.


<b>Uncommon side effects</b>	<b>What is observed?</b>	<b>How long does it last and how can it be helped?</b>
Skin infection	Appears as inflamed skin. This is a secondary effect which may occur after a skin burn or other injuries.	Skin infection requires medical treatment or medication. We advise you to consult your doctor. Do not treat over skin infections.
Scarring	A scar may occur as a secondary effect after a burn or skin infection.	It is important to consult your doctor if you have a burn or folliculitis to prevent formation of a scar. Do not treat over scar tissue.
Herpes outbreak	Small blisters or open sores, typically near the lip or genitals, experienced by users infected with Herpes Simplex virus	An outbreak typically lasts for several days to weeks and can be shortened by taking anti-viral medications prescribed by your doctor.
Headache and afterimages	Sensitive persons may experience headache or afterimages when the light flash reaches the eye.	These reactions typically last between minutes to several hours. To prevent them, make sure you use Lumea in a well-lit room and press the device firmly onto your skin before flashing.

Hair removal by Intense Pulsed Light can cause increased hair growth in some individuals. Body areas at risk include women's jawline, neck and areas above the cheekbone and men's upper back and shoulders. Based upon current data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage.

# What is IPL?

Scan the QR code to watch a short video.
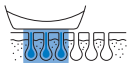
## 1. IPL stands for Intense Pulsed Light.

Gentle pulses of light are absorbed by the roots of your hair.

## 2. These pulses of light cause the hair follicles to go into a resting phase.

The hair sheds naturally, and growth is reduced for a certain number of months.


## 3. Hair does not grow continuously: it grows in cycles.

Not all the hairs in one area are growing at the same time. Since IPL is only effective on hair that is actively growing, it is important to repeat treatments as indicated.

# Lumea overview


Scan the QR code to watch a short video.
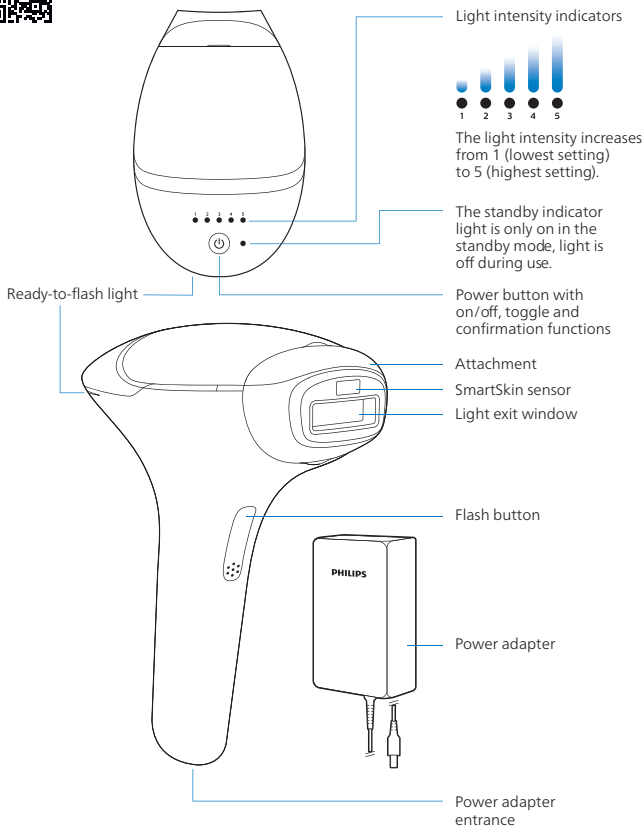
## SmartSkin sensor

SmartSkin sensor determines your skin tone and will help you select a comfortable setting by lighting up one of the light intensity indicators on top of your device.

It also automatically helps you prevent treating areas of your skin that are too dark for safe use.
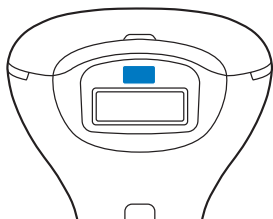
## Ready-to-flash light

The white ready-to-flash light will let you know when you have good skin contact and your Lumea is ready to flash.

If the ready-to-flash light blinks orange, your skin in this body area is too dark to treat. Treat a lighter colored body area.

## SenseIQ technology

**SenseIQ**  
Lumea app  
for optimal results

**SenseIQ** technology for a personalized hair removal treatment. Our signature combination of:

- SmartSkin sensors that measure your skin tone and support you in selecting a setting for a comfortable treatment.
- Unique attachments designed for each curve of your body.
- Philips Lumea IPL app to support your treatment routine.

# Attachments

Scan the QR code to watch a short video.

Your Lumea comes with different attachments you can use on different body areas. Select the attachment that is intended for the area you want to treat.


	Attachment type	BRI931	BRI932	BRI937	BRI984
	<b>Body attachment</b> <ul style="list-style-type: none"> <li>- largest window size;</li> <li>- a curved-in design to effectively cover and treat areas below the neckline.</li> </ul>	✓	✓	✓	✓
	<b>Face attachment</b> <ul style="list-style-type: none"> <li>- smallest window size for precise treatment below the cheekbones;</li> <li>- small, flat design;</li> <li>- extra integrated filter for treatment on the face below the cheekbones.</li> </ul> <p><b>⚠ DO NOT</b> use Lumea on male neck and face. The hairs in these areas have different characteristics.</p> <p><b>Tip:</b> Place your tongue between your upper lip and your teeth or puff out your cheeks to make the treatment easier.</p>		✓	✓	✓
	<b>Precision attachment</b> <ul style="list-style-type: none"> <li>- medium window size;</li> <li>- curved-out design for precise and effective use on the bikini and armpit areas.</li> </ul>			✓	
	<b>Precision attachment for bikini</b> <ul style="list-style-type: none"> <li>- medium window size;</li> <li>- curved-out design for precise and effective use on the bikini area.</li> </ul>				✓
	<b>Precision attachment for armpits</b> <ul style="list-style-type: none"> <li>- medium window size;</li> <li>- curved-out design for precise and effective use on the armpit area.</li> </ul>				✓

# Treatment modes

Scan the QR code to watch a short video.

Your Lumea has two treatment modes for more convenient use on different body areas.
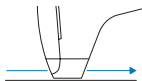
## Slide & flash

For larger areas (legs), keep the flash button continuously pressed as you slide Lumea over your skin. Lumea keeps releasing flashes until you release the button or if you lose good skin contact.

## Stamp & flash

For small body areas (bikini, armpit), place Lumea on your skin and press the flash button. Lift Lumea from your skin and place it beside the spot previously treated. Press the flash button again to release another flash.

# Your treatment guide

Please read the important safety instructions in the user manual to check if a Lumea treatment is suitable for you, so you can fully benefit from your device.

## Before your treatment

### 1. Prepare for treatment

Remove hair, ensure that your skin is clean and dry

### 2. Test your skin

Find your comfortable setting

### 3. Wait 24h for skin reactions

## Treatment

### 4. Perform a treatment

### 5. Clean and store

### 6. Plan your next treatments


Scan the QR code and watch our how-to videos for step-by-step guidance on performing a Lumea treatment, along with general videos about your Lumea device.

## Prepare for treatment

Scan the QR code to watch a short video.
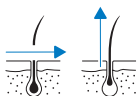
### 1. Remove hair

Wax, shave or epilate before using Lumea.

Do a treatment before new hair is visible on your skin.

If you choose to wax, wait 24 hours before you use Lumea.

### 2. Clean and dry skin

No water, oil, cream or body milk should be left on your skin.

### 3. Choose the right attachment and check your Lumea for damage and dirt

**⚠ ALWAYS** check Lumea before you use it. Do not use Lumea or the power adapter if it is damaged. Always replace a damaged part with one of the original type.

**⚠ DO NOT** use Lumea if there is any lotion or gel on your skin or on Lumea. Only treat skin that is clean and dry. Make sure Lumea is clean and dry as well.

## Test your skin


Scan the QR code to watch a short video.


**!** **DO NOT** start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills).

Before your first treatment on each body area or after exposing your skin to sunlight unprotected, test your skin's reaction and find your suitable setting for treatment. It is recommended to use Lumea out of direct sunlight. The sunlight can influence the SmartSkin sensor.

**1. Follow the steps in section 'Prepare for treatment'**

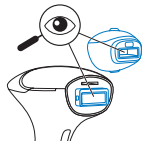
**2. Check your Lumea for damage and dirt**

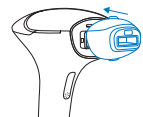**3. Attach the right attachment**

Click on the attachment for the area you want to treat.

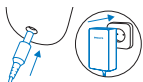**4. Plug in**

Plug in Lumea to the power outlet.

**5. Press the power button**

The 5 light intensity indicators above the power button will light up continuously from left to right and back, until you place Lumea on your skin.

**6. Place Lumea on your skin**

It must sit flush against your skin.

### 7. Get your light intensity setting

One of light intensity indicators will turn white to indicate the proposed setting for you. The setting is selected automatically after 5 seconds.

Note: If the ready-to-flash light blinks orange your skin may be too dark to treat, please consult chapter 'Troubleshooting'.

### 8. Find your comfortable setting

Place Lumea back on your skin. The ready-to-flash light will light up white.

Press the flash button 1 time. If the flash is comfortable, you can increase the setting. If you experience pain, lower the setting. Toggle through the settings by pressing the power button.

Change the setting, move Lumea to the next spot and flash again. Repeat until you find the highest comfortable setting.


### 9. Turn off, unplug

**DO NOT** touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.

Hold the power button for 2 seconds. Unplug your Lumea and the power adapter from the power outlet. Let your Lumea cool down.
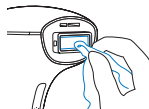
### 10. Clean all parts of your Lumea after use

It is important to clean all parts of your Lumea. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned. Follow the steps in section 'Clean and store'.

### 11. Wait 24 hours and check your skin for any reactions

If there are no skin reactions after 24 hours, start the treatment with the highest comfortable setting. If the highest comfortable setting resulted in a skin reaction, choose a lower setting that did not result in any skin reactions. If you did not test a lower setting, repeat the skin test.

## Perform a treatment

Scan the QR code to watch a short video.

Test your skin's reaction before treating a new body area or after exposing your skin to sunlight unprotected.


### 1. Follow the steps in section 'Prepare for treatment'

### 2. Check your Lumea for damage and dirt

### 3. Attach the right attachment

Click on the attachment for the area you want to treat.
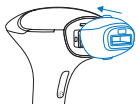
### 4. Be comfortable

Make yourself comfortable before you start the treatment.

It is recommended to use Lumea out of direct sunlight. The sunlight can influence the SmartSkin sensor.

### 5. Plug in. Press the power button

The 5 light intensity indicators above the power button will light up continuously from left to right and back, until you place Lumea on your skin.

### 6. Select setting

Press the power button to toggle through the settings and select your highest comfortable setting that did not result in any skin reactions.


### 7. Begin treatment

Place Lumea on your skin.

The ready-to-flash light will light up white.

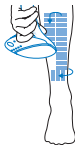
## 8. Flash. Then move

**⚠ DO NOT** continue the treatment if you experience pain. Reduce the light intensity setting or stop using Lumea.

**DO NOT** flash the same spot multiple times during a treatment. The surface of the attachment should overlap the previous spot to make sure you flash the complete skin area, but treating the same area multiple times does not increase performance.

Press the flash button. Move to an adjacent spot and flash again.
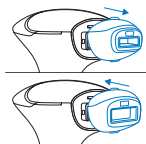
## 9. Treat another body area

**⚠ DO NOT** start or resume using Lumea within 4 weeks after exposing your skin to sunlight unprotected (natural sun or tanning bed) or within 2 weeks after artificial tanning (creams, sprays or pills). Do a skin test before starting/resuming your treatment, see section 'Test your skin'.


If you want to treat another body area, make sure you select the right attachment and start again with step 1. If you are treating a body area for the first time or after exposure to sunlight unprotected, make sure you perform a skin test. Wait 24 hours and check your skin for any reactions.

## 10. Finished. Turn off and unplug

**⚠ DO NOT** touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.

Hold the power button for 2 seconds. Unplug your Lumea and the power adapter from the power outlet. Let your Lumea cool down.
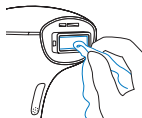
## 11. Clean all parts of your Lumea after use

It is important to clean all parts of your Lumea. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned. Follow the steps in section 'Clean and store'.

## 12. Schedule your next treatment

Follow the steps in section 'Plan your next treatments'.

## Clean and store


Scan the QR code to watch a short video.

### 1. Turn off, unplug and cool your Lumea

**⚠ DO NOT** touch the light exit window during or directly after use, this part becomes very hot. After use, unplug Lumea and the power adapter and let them cool down for 15 minutes before cleaning and storing.
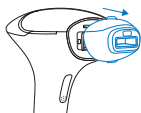
### 2. Click off the attachment

Gently pull it off.


### 3. Lightly moisten your Lumea cleaning cloth

**⚠ DO NOT** use Lumea in wet surroundings.


**DO NOT** let water run into Lumea.

Dampen the soft Lumea cleaning cloth with a few drops of water.

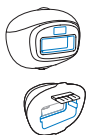
### 4. Clean the attachments

Clean each part of the attachment you used, including the metal reflectors.

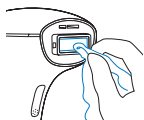
### 5. Clean the light exit window

Make sure there is no dirt, residue or moisture left on the surface of the light exit window.
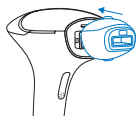
### 6. Click on an attachment

Make sure the light exit window of your Lumea is protected from damage by clicking on one of the attachments.
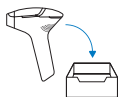
### 7. Store

Store all parts in the original packaging or in the provided bag, pouch or box in a dust-free and dry place.

## After treatment

After use, if you want to apply lotions, creams, deodorants, moisturizers or cosmetics to the treated areas, always check the instruction leaflet of the product you want to apply to make sure it is suitable for the condition of your skin at that moment.

# Plan your next treatments

Scan the QR code to watch a short video.

Follow the indicated treatment schedule.

24 hours

## Test phase

Before your first treatment on each body area or after exposing your skin to sunlight unprotected, do a skin test by following the steps in section 'Test your skin'. Wait 24 hours after the skin test before you start your treatment phase.

## Treatment phase: 4 bi-weekly treatments

During the treatment phase which lasts 4 treatments, perform a treatment every 2 weeks.

Perform a treatment within 3 days before or after the planned treatment date during the treatment phase.

## Touch-up phase: once every 4 weeks

After the treatment phase, use your Lumea once every 4 weeks to maintain results.

Perform a treatment within 4 days before or after the planned treatment date during the touch-up phase.

Note: If you miss a treatment, consult chapter 'Troubleshooting' for guidance. If you deviate from the treatment schedule, you may not achieve the intended results.

# Troubleshooting

This chapter summarizes the most common problems you could encounter with your Lumea or if you deviate from your treatment schedule. If you are unable to solve the problem with the information below, visit [www.philips.com/support](http://www.philips.com/support), or call 1-888-744-5477 for assistance.

Common problems while using your Lumea:

## Problem

---

My Lumea is not flashing / the ready-to-flash light turns orange and starts to blink.

---

My Lumea produces a strange smell while flashing.

---

The light intensity indicators on top of my Lumea are blinking and the ready-to-flash light is orange and blinking.

## Solution

---

1. Your skin may be too dark. Check the skin tone/hair color tables inside the cover of this booklet.
2. Make sure your attachment is (properly) connected to the front of your Lumea.
3. Make sure you have full skin contact for at least 5 seconds. In some areas (sensitive and bony area) it might be difficult to make good skin contact. No light should be able to enter between Lumea and your skin. Try repositioning your Lumea.
4. Make sure you press the flash button on the inside handle to treat an area of your skin.

- 
1. Make sure you have removed all hairs in the area you want to treat, before using your Lumea. Burned hairs can cause smells.
  2. Make sure your skin is clean and dry before starting a treatment.
  3. Make sure your Lumea and its attachments are clean and dry before starting a treatment. Residue like hair, sweat, or dirt can build up from the flashing process and may cause skin reactions if left uncleaned.

- 
1. Lumea is in error mode. Press the power button for more than 2 seconds, Lumea should go back to standby. Or unplug Lumea from the mains, it should turn off.
  2. Wait 5 minutes. Lumea will automatically go back to standby after 5 minutes without help.

What to do if you deviate from your treatment schedule:

## Phase

---

Treatment phase

## Problem

---

Missed treatment.

*I did my treatment more than 3 days after the originally planned date.*

## Solution

---

Try to stay closer to your planned treatment date. You can continue with the original planned schedule. You may add one additional bi-weekly treatment to your treatment phase. If you deviate from the treatment schedule, you may not achieve the intended results.

---

Last bi-weekly treatment.

*I finished my treatment phase, but I feel I did not reach the claimed results.*

---

You may choose to add an additional bi-weekly treatment to your treatment phase, but claimed results are based on the indicated treatment schedule and individual results vary.

---

Touch-up phase

---

Out of treatment window.

*I did not treat within 4 days of my scheduled treatment date.*

---

Please follow the correct schedule. If you deviate from the treatment schedule, you may not achieve the intended results. You can continue with the original planned schedule or you may choose to restart the treatment phase, but claimed results are based on the indicated treatment schedule and individual results vary.

# Additional information

## Full Two-Year Warranty

Philips warrants its products for a specific period of time after the first date of purchase. If the product is purchased in the USA, Philips North America LLC (USA) or, if the product is purchased in Canada, Philips Electronics Ltd., warrants this new product against defects in materials or workmanship for a period of two years from the date of purchase (as evidenced by a valid proof of purchase), and agrees to replace any defective product without charge. Philips will refund the purchase price if it is unable to provide replacement.

Parts that are subject to normal wear and tear, accessories, replacement batteries and consumables such as toothbrush heads, are not covered by the terms of this warranty.

**IMPORTANT:** This warranty does not cover damage resulting from accident, misuse or abuse, lack of reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified voltage.\*

**NO RESPONSIBILITY IS ASSUMED FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES.** For more information and in order to obtain warranty service, simply go to [www.philips.com/support](http://www.philips.com/support).

NOTE: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state, province to province or country to country. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above exclusion and limitations may not apply to you. This warranty is transferable.

\*Read enclosed instructions carefully.

### THE WARRANTY DOES NOT COVER THE FOLLOWING:

- Indirect or consequential damages, such as loss of time, loss of data, or loss of income. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.
- Compensation for activities done by yourself, such as regular maintenance, installing firmware updates, or saving data.

### THE WARRANTY DOES NOT APPLY:

- If the proof of purchase, model number, serial number, or production date code on the product has been made illegible.
- If the product was purchased as faulty stock.
- If the defect is the result of abuse or excessive use outside the intended purpose.
- If the defect is caused by environmental conditions, connecting peripherals, additional equipment, or accessories (including software) that are not in accordance with the user manual.
- If the unit has been damaged, including but not limited to damage by animals, lightning, abnormal voltage, fire, natural disaster, transportation, or water (unless the user manual expressly states that the product may be rinsed).
- For non-electrical parts, normal wear and tear, accessories, replacement batteries or replacement of consumable parts such as toothbrush heads.
- If the product does not function properly because it was not originally designed for use in the country where the product is used, which may occur if you import the product.
- If the product does not function properly due to problems with connection with service providers, such as for the internet.

## 90-Day Money-Back-Guarantee

If you are not fully satisfied with your product, send the product back and Philips will refund you the full purchase price. The product must be shipped prepaid by insured mail, insurance prepaid, have the original sales receipt, indicating purchase price and date of purchase. Philips cannot be responsible for lost mail. The product must be postmarked no later than 90 days after the date of purchase. Philips reserves the right to verify the purchase price of the product and limit refunds not to exceed suggested retail price.

To obtain the money-back guarantee service, please go to [www.philipsmoneyback.com](http://www.philipsmoneyback.com), and follow the simple 3-step process. Delivery of refund check will occur 6-8 weeks after receipt of the returned product.

## Disposal

Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit [www.recycle.philips.com](http://www.recycle.philips.com).

## Electromagnetic compatibility

This product has been tested and found to comply with IEC 60601-1-2 for electromagnetic compatibility (EMC) as Class B according to CISPR 11. Tests levels are listed in the accompanying tables. Only the power adapter that is supplied with Lumea should be used with your device.

### Guidance and manufacturer's declaration on electromagnetic emissions and immunity

Lumea is intended for use in the electromagnetic environment listed here. The user of the device should assure that it is used in such an environment to prevent improper operation.

Emissions	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	Lumea does not use RF energy for its intended function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	Lumea is suitable for use in locations in residential environments and in establishments directly connected to a low voltage power supply network which supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

Immunity test	Immunity test level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	Floors can be a source of static electricity. When the air is very dry (<20 % RH) the user is advised to use Lumea in a room with concrete or wooden floor (instead of using it in a room with a synthetic floor).
Radiated RF EM fields IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	Lumea is immune against any RF wireless communication device such as smart phone, Wi-Fi or Bluetooth device at a distance of 30 cm or more. Avoid a closer distance as this may result in improper operation.
Proximity fields from RF wireless Communications Equipment IEC 60601-1-2	385 (27), 450 (28), 710 (9), 745 (9), 780 (9), 810 (28), 870 (28), 930 (28), 1720 (28), 1845 (28), 1970 (28), 2450 (28), 5240 (9), 5500 (9), 5785 (9) MHz (V/m)	
RATED power frequency magnetic fields IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	
Proximity magnetic fields IEC 61000-4-39	30 kHz, CW, 8 A/m 134.2 kHz, Pulse modulation 2.1 kHz, 65 A/m 13.56 MHz, Pulse modulation 50 kHz, 7.5 A/m	
Electrical fast transient/burst IEC 61000-4-4	±2 kV 100 kHz repetition frequency	Mains power quality expected of a typical commercial or residential environment.
Surges IEC 61000-4-5	±0.5 kV, ±1 kV Line-to-line	
Conducted disturbances induced by RF fields IEC 61000-4-6	3 V r.m.s., before modulation is applied 0.15 MHz – 80 MHz 6 V r.m.s., before modulation is applied in ISM and amateur radio bands between 0.15 MHz and 80 MHz 80 % AM at 1 kHz	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1cycle and 70 % UT; 25/30 cycles Single phase at 0° 0 % UT; 250/300 cycle	Mains power quality expected of a typical commercial or residential environment. At power interruptions Lumea will switch off.

**Note:** UT is the AC mains voltage prior to application of the test level.

If in the unlikely event Lumea is locked inoperative in a function by disturbances from other devices beyond the levels stated in the table above, Lumea will shut-off completely and must be restarted by the power button. In that case make sure the distance to the other device causing this disturbance is enlarged.

## Electromagnetic fields (EMF)

This device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Supplier's Declaration of Conformity 47 CFR § 2.1077 Compliance Information

Unique Identifier: BRC332, BRC344

### Responsible Party – U.S. Contact Information

Philips Personal Health  
a division of Philips North America LLC  
P.O. Box 10313  
Stamford, CT 06904  
United States of America

### FCC Compliance Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

No changes shall be made to the equipment without the manufacturer's permission as this may void the user's authority to operate the equipment.

## Technical specifications

Electrical specification	The combination of the power adapter and Lumea is considered: Class II ME EQUIPMENT
Attachment	Type BF applied part The applied part (attachment) can become hot during use <140 °F (<60 °C). If this occurs, let the device cool down and limit the contact duration to less than 1 minute.
Mode of operation	Continuous operation
Intended operator	The user is the intended operator
IP rating power adapter	IP22
Light output	
Risk group classification	Exempt
Wavelength range of emitted optical radiation	Body and precision attachments: 520 nm – 1200 nm Face attachment: 590 nm – 1200 nm
Spot size	Body attachment: 4.1 cm <sup>2</sup> (32.7 mm x 12.6 mm) Precision attachment: 3 cm <sup>2</sup> (30.0 mm x 10.0 mm) Face attachment: 2 cm <sup>2</sup> (25.0 mm x 8.0 mm)
Optical exposure	Setting 1: 2.5 ± 0.80 J/cm <sup>2</sup> Setting 2: 3.1 ± 0.85 J/cm <sup>2</sup> Setting 3: 3.8 ± 0.91 J/cm <sup>2</sup> Setting 4: 4.4 ± 0.96 J/cm <sup>2</sup> Setting 5: 5.0 ± 1.00 J/cm <sup>2</sup>
Maximum optical energy	Body attachment: 23.0 J Precision attachment: 16.8 J Face attachment: 11.2 J
Optical homogeneity	Max. +/-30 % deviation from average optical exposure in treatment area
Pulse sequence	Single pulse, for all settings
Pulse interval	1.0 – 2.4 s, for all settings
Pulse width	1.25 +/-0.4 ms (FWHM)
Transport and storage conditions	Lumea remains operational in NORMAL USE within its specification after transport or storage in the following environmental range: - a temperature range of -13 °F to 158 °F (-25 °C to 70 °C); and - a relative humidity up to 95 % RH, noncondensing; and - an atmospheric pressure range of 700 hPa to 1060 hPa. If the device has been stored in a hot or cold environment, place it in an environment with a temperature of 68 °F (20 °C) for 30 minutes to let it reach a temperature within the usage conditions (41 °F to 104 °F/5 °C to 40 °C) before you use it.
Operating conditions	Lumea complies with its specifications when operated in NORMAL USE under the following environmental operating conditions: - a temperature range of 41 °F to 104 °F (5 °C to 40 °C); - a relative humidity range of 15 % to 90 % RH, noncondensing; and - an atmospheric pressure range of 700 hPa to 1060 hPa.

<b>Power adapter</b>	
Power adapter input	100-240 V ~ 50/60 Hz 1.5 A
Power adapter output	19.5 V $\equiv$ 3.33 A 65.0 W
Average active efficiency	$\geq$ 88.0 %
Efficiency at low load (10 %)	$\geq$ 79.0 %
No-load power consumption	<0.21 W
<b>Device weight</b>	Power adapter ~ 0.7 lb (~ 310 g) Handpiece ~ 0.9 lb (~ 430 g)
<b>Device size</b>	Power adapter ~ 4.3 x 2.6 x 2.7 inch (~ 10.8 x 6.5 x 6.8 cm) Handpiece ~ 78.7 x 5.7 x 3.7 inch (~ 20.0 x 14.5 x 9.5 cm)
<b>Service life</b>	5 years, based on Lumea treatment schedule with 14 full body treatments per year.
<b>Materials</b>	
Housing	PC (Polycarbonate), ABS (Acrylonitrile-butadiene-styrene), PA (Polyamide)
Cable	TPE (Thermoplastic elastomer)

## Warranty and support

Philips offers a 2-year warranty. Warranty exceptions apply. This warranty does not affect your rights as a consumer under law. For more information, visit [www.philips.com/support](http://www.philips.com/support) or call us.

## Explanation of symbols

The safety signs and symbols are essential to ensure that you use this product safely and correctly and to protect you and others from injury. Below you find the meaning of the safety signs and symbols on the label and in the user manual.

**IEC 60417-5333 Type BF applied part**  
This symbol identifies the type BF applied part complying with IEC 60601-1.

**ISO 7000-2620 Humidity limitation**  
This symbol indicates the acceptable upper and lower limits of relative humidity for transport and storage.

**ISO 7000-2621 Atmospheric pressure limitation**  
This symbol indicates the acceptable upper and lower limits of atmospheric pressure for transport and storage.

**ISO 7000-0632 Temperature limit**  
This symbol indicates the maximum and minimum temperature limits at which the item shall be stored, transported or used.

**ISO 7000-0626 Keep dry**  
This symbol indicates that the device shall be kept dry.

**IP22**

This symbol indicates the protection against ingress of solid foreign objects  $\geq 12.5$  mm and against harmful effects due to the ingress of dripping water when tilted at 15 degrees.

**IEC 60417-5172 Class II equipment**  
This symbol indicates 'Class II Equipment'. The power adapter is double-insulated.

**ISO 7000-3082 Manufacturer**  
This symbol indicates the manufacturer of the product.

**IEC 60417-6049 Country of manufacture**  
This symbol indicates the country of manufacture and the manufacturing date. 'CC' is the two-letter country code.

**LOT**

**ISO 7000-2492 Batch number**  
This symbol indicates the manufacturer's batch number.

**REF**

**ISO 7000-2493 Catalogue number**  
This symbol indicates the manufacturer's catalogue number.

**IEC 60417-5032 Alternating current**  
This symbol indicates alternating current.

**IEC 60417-5031 Direct current**  
This symbol indicates direct current.

**ISO 7000-1641 Consult instructions for use**  
This symbol indicates that operating instructions should be considered when operating the device.

**ISO 7010-W017 Warning; Hot surface**  
This symbol warns about a hot surface. The background color is yellow, the color of the symbol is black.

***t<sub>a</sub>***

**IEC 61558-1 Rated maximum ambient temperature**  
This symbol indicates the maximum ambient temperature at which the transformer may be operated continuously under normal conditions of use.

**IEC 60417-5957 For indoor use only**  
This symbol means that the device is designed for indoor use.

**IEC 60417-5947 Safety isolating transformer, short-circuit proof**  
This symbol means that the transformer is designed to be short-circuit proof.

**IEC 60417-6190 Power supply unit, switch mode**  
This symbol means that the transformer is considered a Switch Mode Power Supply (SMPS) unit.

#### IEC 60417-6343 Maximum altitude

This symbol indicates that the appliance is intended to be usable up to the maximum altitude 3000 m (9842 ft.).

#### ISO 7000-3706 Single patient multiple use

This symbol indicates that the device may be used multiple times on a single user.

#### SenseIQ technology

This symbol indicates that the device is equipped with SenseIQ technology.

**ISO/IEC 16022** Information technology – Automatic identification and data capture technologies – Data Matrix barcode symbology specification. The data matrix (QR-code) is used to capture information in a machine-readable form.

#### ISO 15223-1 Unique Device Identifier

Indicates the Unique Device Identifier that is used to identify a specific product.

#### ISO 7010-W001 General warning sign

This symbol means that the user should consult the instructions for use for important safety information such as contraindications and warnings. On the device, the background color is yellow, the color of the symbol is black.

#### ISO 7000-0434A Caution

This symbol means that the user should consult the instructions for use for important safety information on cautions.

A recycling symbol to identify the material of which an item is made, to facilitate recycling or other reprocessing. The symbol includes a number to identify the material code.

#### IEC 60417-5009 Stand-by

This symbol identifies the switch or switch position by means of which part of the equipment is switched on in order to bring it into the stand-by condition. The symbol also identifies the control to shift to or to indicate the state of low power consumption.

This symbol indicates Forest Stewardship Council®. The FSC trademarks enable consumers to choose products that support forest conservation, offer social benefits, and enable the market to provide an incentive for better forest management.

#### ISO 15223-1 Medical device

This symbol indicates that this is a medical device.

The FCC mark is voluntary mark, which indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission.

MEDICAL — GENERAL MEDICAL EQUIPMENT PHILIPS IPL HAIR REMOVAL AS TO ELECTRICAL SHOCK, FIRE AND MECHANICAL HAZARDS ONLY IN ACCORDANCE WITH AAMI ES60601-1 (2005) + AMD 1 (2012) and the applicable Amendments and/or Particular and/or related Collateral Standards for which the product has been found to comply by UL.

MEDICAL — GENERAL MEDICAL EQUIPMENT PHILIPS IPL HAIR REMOVAL AS TO ELECTRICAL SHOCK, FIRE AND MECHANICAL HAZARDS ONLY IN ACCORDANCE WITH CSA-C22.2 No. 60601-1:2014 and the applicable Amendments and/or Particular and/or related Collateral Standards for which the product has been found to comply by UL.

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All rights reserved.
Philips Consumer Lifestyle B.V.