
**PHILIPS**


**sonicare**

For Kids

Design a Pet Edition

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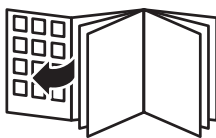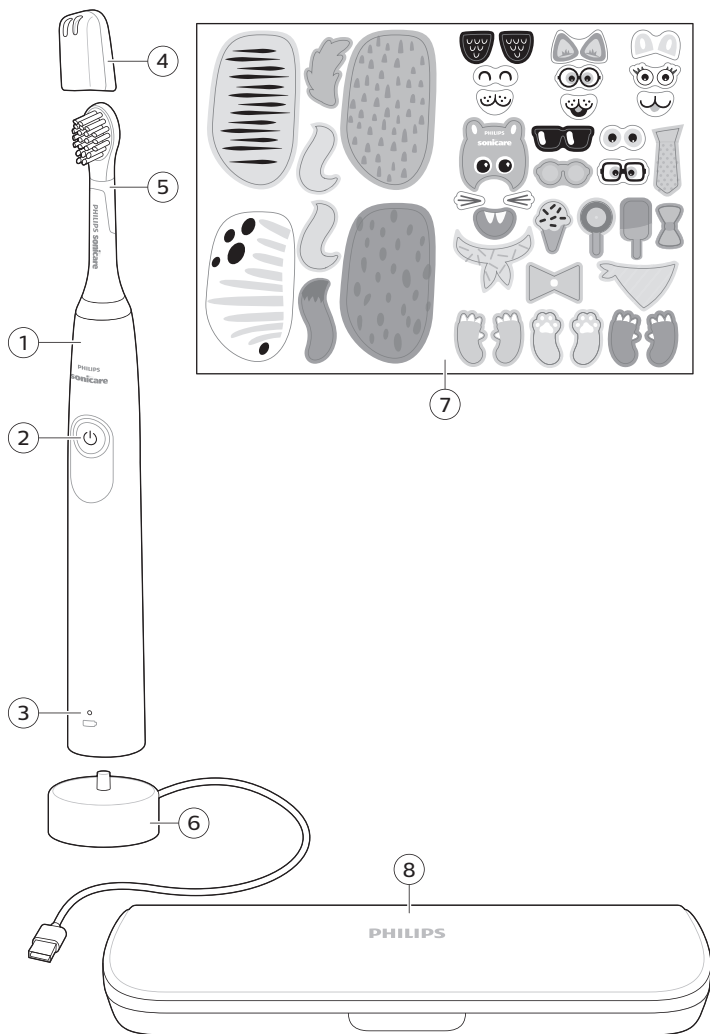## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

## Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products. NOT A TOY. Check brush heads regularly, cracks from dropping can lead to a choking hazard.

## Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country.
- Keep charging cables and wall adapters away from water.
- After cleaning, make sure the wall adapters and charging cables are completely dry before connecting to power.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not use charging cables and wall adapters outdoors or near heated surfaces.

- If the charging cable or wall adapter is damaged, stop using it.
- This appliance has only been designed for cleaning teeth and gums.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Make sure brush head and toothbrush are dry before storing in travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.
- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

### **Battery safety instructions**

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.

- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Disposal.
- Charge only with Philips WAA1001, WAA2001 safety approved wall adapters with input ratings (100 -240V~; 50/60Hz; 3.5W) and output ratings (5V ; 2.5W). For the list of approved adapters, please see [www.philips.com/support](http://www.philips.com/support).
- Charge and use the product at a temperature between 0 °C and 40 °C (32 °F and 104 °F).

### Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

### Storage and transportation conditions


Temperature -20 °C to +60 °C (-4 °F to 140 °F).

### LIST OF SYMBOLS USED

The following symbols may appear on the product:  
This symbol means read the instructions for use, before using the appliance.

This symbol means to use the indicated power supply.
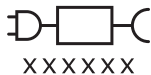
### Intended use

Sonicare for Kids - Design a Pet Edition power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare for Kids -

Design a Pet Edition power toothbrushes are intended for personal use. Use by children should be with adult supervision.

## Your Philips Sonicare (Fig 1.)

- 1 Handle
- 2 Power on/off button
- 3 Battery indicator

### Accessories

- 4 Brush head cap
- 5 Brush head
- 6 USB charger (wall adapter not included)
- 7 Stickers
- 8 Travel case

Note: Accessories included may vary based on the model purchased.

## Getting started

### Sonicare for Kids App

Download the Sonicare for Kids App to support healthy habit building. The App is compatible with wide range of iOS and Android (™) smart devices.

Note: If your toothbrush is not Bluetooth capable you can still use the Sonicare for Kids app on iOS and Android (™) smart devices to enhance your brushing experience.

### Applying your favorite sticker


The Sonicare for Kids - Design a Pet Edition comes with several stickers to choose from. The stickers can be moved around to create the design you choose.

- 1 Apply your favorite sticker(s) to the handle


Note: Make sure the handle is clean and dry, and that stickers are fully pressed onto the handle with no air gaps to ensure full adhesion.

Note: Kids should not be left unattended when applying the stickers.

## Attaching the brush head

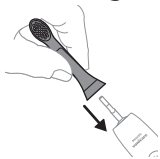- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

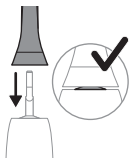- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

Note: It is normal to see a slight gap between the brush head and the handle.

## Using your Philips Sonicare toothbrush


If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience. To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gently increases the power over the first 14 brushings to help you adapt to brushing with a Philips Sonicare toothbrush.

## Brushing instructions

Note: For children under age 8, it is recommended that a parent or adult caregiver perform toothbrushing.


- 1 Wet the bristles and apply a small amount of toothpaste.

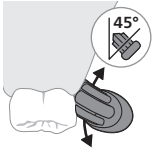
- Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

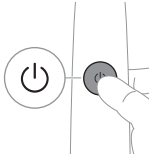
- Press the power on/off button to switch on the Philips Sonicare.

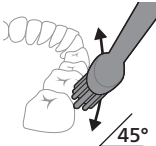- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.


Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

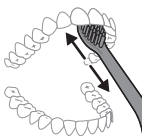- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

## KidPacer

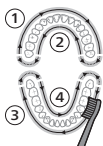To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the KidPacer feature. The KidPacer divides the brushing time into 4 equal segments and indicates when you should move to the next area. The segments are indicated with a short tune. The toothbrush automatically stops at the end of a brushing session with a congratulatory tune.

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush switches off.

## Intensities

Your power toothbrush is equipped with 2 different intensities: gentle and extra gentle.

Note: When you use the toothbrush for the first time, the default setting is extra gentle.

To change intensities:

- Press the power on/off button once to turn on the toothbrush
- Press a second time within 2 seconds to change the intensity

Press a third time within 2 seconds to pause. After 2 seconds of brushing, pressing the power on/off button again will pause the toothbrush

	Extra gentle	Gentle
<b>Benefit</b>	Reduced power intensity used for training and introducing younger children to the sonic experience.	Full-strength intensity designed for older/independent brushers.
<b>Recommended brush heads</b>	Sonicare for Kids Compact or Standard brush head	Sonicare for Kids Compact or Standard brush head
<b>Total brushing time</b>	2 minutes	2 minutes

## How to brush

## Features

- KidPacer
- SmarTimer
- EasyStart

### KidPacer

The KidPacer is an interval timer that helps your child to brush evenly through the mouth by dividing your mouth into 4 sections.

### SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. After a pause of 30 seconds, the SmarTimer resets.

### EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare.


To deactivate EasyStart (see 'Activating or deactivating features').

## Activating or deactivating features


You can activate or deactivate the following features of your toothbrush:

- EasyStart

Step 1: Place the handle on the charger.

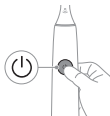
Step 2: Press and hold power on/off button for up to 3 seconds.

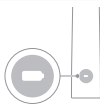Step 3: Let go of the power on/off button when you hear 1 beep.

If you see the battery indicator flash green 2 times and hear 3 tones from low to high, then the feature has been activated

OR

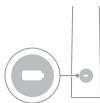If you see the battery indicator flash amber 1 time and hear 3 tones high to low, then the feature has been deactivated.

Note: Using the EasyStart feature beyond the initial ramp-up period is not recommended and reduces Sonicare's effectiveness in removing plaque.

## Charging and battery status


This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

## Charging

### USB charging

Note: Accessories included may vary based on the model purchased.

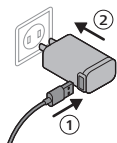- 1 Plug the USB cord of the charging stand into a wall adapter that is connected to an electrical outlet.
- 2 Place the handle on the charger.
  - The battery indicator lights up and the handle beeps twice. This indicates that the toothbrush is charging.

### Battery status

When the handle is placed on the charger, the battery indicator shows the battery level.

Battery status	Battery indicator
Charging	Flashes green
Full	Solid green for 30 seconds then switches off

### Battery status (when handle is not placed on charger)

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	Battery indicator	Sound
Full	Solid green	-
Partially full	Flashes green	-
Low	Flashes amber	The handle emits 3 beeps after brushing
Empty	Flashes amber	The handle emits 2 sets of 5 beeps after brushing


Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush on the charger when not in use.

## Cleaning

### Toothbrush handle

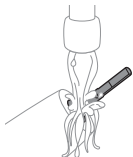
- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

**Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.**

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

### Brush head

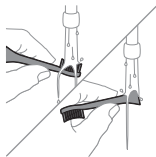- 1 Rinse the brush head and bristles after each use and allow to air dry.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

### Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.
- 3 Dry the charger with a soft cloth or towel.

### Travel case

Rinse with warm water and use a damp cloth to clean the travel case.

## Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

## Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Only use Sonicare for Kids Compact or Standard replacement brush heads.

## Locating the model number

The model number begins with "HX" (e.g. HX123A) and can be found on the bottom of your toothbrush.

## Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit [www.philips.com/support](http://www.philips.com/support) and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.

Problem	Possible Cause	Solution
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2 mm). When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

## Warranty and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the international warranty leaflet.

## Warranty exclusions

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

## Removing the rechargeable battery

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Before removing the battery, make sure that the battery is completely empty.

**CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.**

**If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.**


**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.


- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

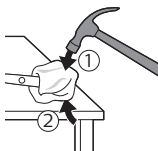
Note: You may have to hit on the end several times to break the internal snap connections.

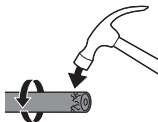
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

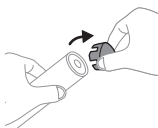
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

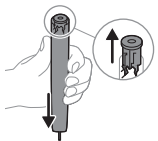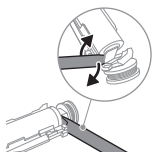- 6 Wedge the screwdriver between the battery and the frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the frame.
- 7 Insert the screwdriver between the bottom of the battery and the frame to break the metal tab connecting the battery to the printed circuit board. This will release the bottom end of the battery from the frame.

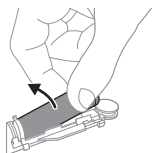- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

**Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.**

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

**Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.**

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©2026 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.