

PHILIPS sonicare

For Kids 3+ Smart

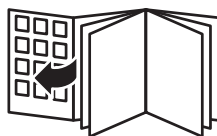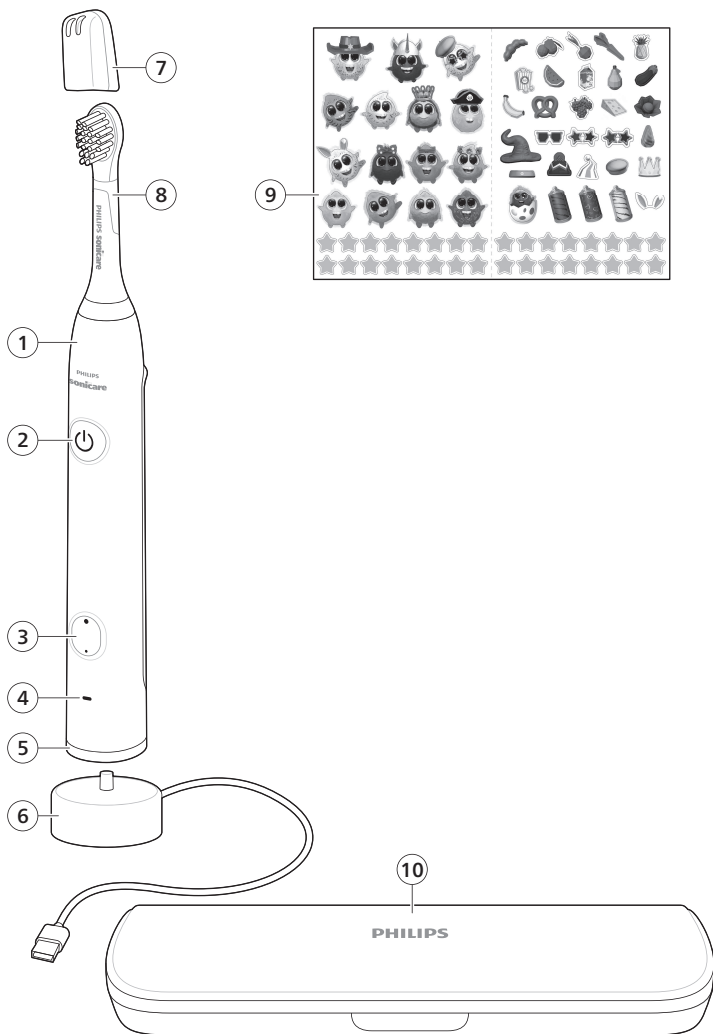English 6
Español 25
Français (Canada) 45

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **www.philips.com/welcome**.

Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products. NOT A TOY. Check brush heads regularly, cracks from dropping can lead to a choking hazard.

WARNINGS

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country.
- Keep charging cables and wall adapters away from water.
- After cleaning, make sure the wall adapters and charging cables are completely dry before connecting to power.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not use charging cables and wall adapters outdoors or near heated surfaces.

- If the charging cable or wall adapter is damaged, stop using it.
- This appliance has only been designed for cleaning teeth and gums.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Make sure brush head and toothbrush are dry before storing in travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.

MEDICAL WARNINGS

- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

Battery safety instructions

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.

- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Disposal.
- Charge only with Philips WAA1001, WAA2001 safety approved wall adapters with input ratings (100 -240V~; 50/60Hz; 3.5W) and output ratings (5V ; 2.5W). For the list of approved adapters, please see **www.philips.com/support**.
- Charge and use the product at a temperature between 0 °C and 40 °C (32 °F and 104 °F).

Storage and transportation conditions

Temperature -20 °C to 60 °C (-4 °F to 140 °F).

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC Compliance

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that

interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).
- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:


- 1 This device may not cause interference.
- 2 This device must accept any interference, including interference that may cause undesired operation of the device.

SAVE THESE INSTRUCTIONS

LIST OF SYMBOLS USED

The following symbols may appear on the product:
This symbol means read the instructions for use, before using the appliance.

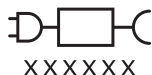This symbol means to use the indicated power supply.

Intended Use

Sonicare for Kids power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare for Kids power toothbrushes are intended for personal use. Use by children should be with adult supervision.

Your Philips Sonicare toothbrush

- 1 Handle
- 2 Power on/off button
- 3 Mode button
- 4 Battery indicator
- 5 Light ring
- 6 USB charger (wall adapter not included)

Accessories:

- 7 Brush head cap
- 8 Brush head
- 9 Stickers
- 10 Travel case

Note: Accessories included may vary based on the model purchased.

Getting started

Sonicare for Kids App

The Sonicare for Kids app connects with your toothbrush to deliver an interactive brushing experience. The app collects information about your kids' brushing habits and provides tips to enhance their oral healthcare routine. The included Parent Dashboard feature helps you understand your child's brushing progress.

The app is compatible with a wide range of iOS and Android™ mobile devices. To start using the app:

- 1 Download the Sonicare For Kids app to your mobile device.
- 2 Enable Bluetooth®/MD on your mobile device.
- 3 Open the app and follow the on-screen instructions.
- 4 To connect the toothbrush with the app, press the power button or mode button to turn it on.

Note: Do not place the toothbrush on the charger when turning it on.

- 5 Register a Philips account in the app to access the Parent Dashboard and to enable back up and restoration of your child's data.
- 6 Complete any updates, if prompted, to access the latest improvements and features.
- 7 The app will automatically connect to the toothbrush. If it does not connect, follow the troubleshooting steps in the app.

Tips: Always place your mobile device in a safe, dry location. Brush regularly with the toothbrush connected to the app for the best interactive experience.


Tips: Keep your mobile device close to the toothbrush during use for optimal data transfer. For the latest improvements and features, keep updating to the latest version of the app.

Note: You can remove the toothbrush by managing Bluetooth connections in your mobile device settings.

Note: To keep your brushing data available in the app, make sure the toothbrush is charged regularly. If the toothbrush battery is fully depleted for an extended period of time, or if the battery is replaced, previously stored brushing data may no longer be available.

Applying your favorite sticker

This toothbrush comes with several stickers to choose from.

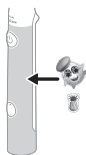- 1 Apply sticker(s) on the handle to create your personalized toothbrush.

Note: Apply stickers only on the handle and travel case. Do not apply stickers on the brush head.


Note: Make sure the handle or travel case is clean and dry, and that stickers are fully pressed onto the handle with no air gaps to ensure full adhesion.

Note: Kids should not be left unattended when applying the stickers.

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

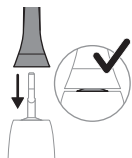- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

Note: It is normal to see a slight gap between the brush head and the handle.


Using your Sonicare for Kids toothbrush

If you are using your Sonicare for Kids toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience.


Your Sonicare for Kids toothbrush comes with the EasyStart feature activated to help you transition to using a power toothbrush. This feature gently increases the power over the first 14 brushing sessions to help you adapt to brushing with a Philips Sonicare toothbrush.

Note: For children under age 8, it is recommended that a parent or adult caregiver perform or supervise toothbrushing.

Brushing instructions

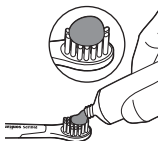
- 1 Wet the bristles and apply a small amount of toothpaste.

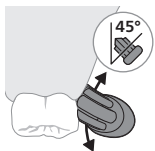
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times to ensure the whole tooth surface is brushed.

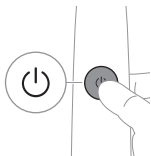- 3 Press the power on/off button to switch on the toothbrush.

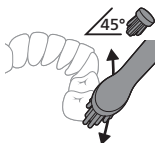
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the Sound Pacer feature (see chapter 'Features').

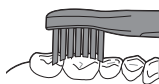
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If you experience any issues, consult your dental professional.

***Disclaimer**

If the Philips Sonicare for Kids toothbrush is used in clinical studies, it must be used in Gentle mode with the handle fully charged and the EasyStart feature deactivated.

Modes

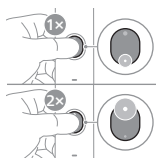
Your power toothbrush is equipped with two different modes: Extra Gentle, and Gentle.

You can press the mode button to change the mode.


Mode	Extra Gentle	Gentle
Benefit	Reduced power mode used for training and introducing younger children to the sonic experience.	Full-strength mode designed for older/independent brushers.
Total brushing time	2 minutes	2 minutes

How to brush

Brush each segment for 30 seconds

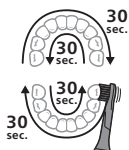Brush each segment for 30 seconds

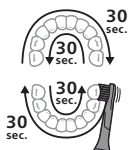Features

- 2-minute Light and Sound Timer
- Light and Sound Pacer
- EasyStart

2-minute Light and Sound Timer

The 2-minute Light and Sound Timer plays a celebratory tune and displays a light ring effect at the end of the predetermined brushing time. Additionally, the timer automatically turns off the handle when the brushing cycle is complete.


- In Extra-Gentle Mode, the timer is set at 1 minute and slowly increases over time. This gradual increase helps young children build up to the dentist-recommended 2 minutes of brushing.

- In Gentle Mode, the timer is set at 2 minutes.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: Instruct your child to brush until they hear the celebratory tune and the toothbrush turns off.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. If you press the on/off button within 30 seconds, the timer resumes to allow you to finish your brushing cycle. After a pause of 30 seconds, the timer resets.

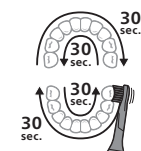Light and Sound Pacer

To make sure your child brushes evenly throughout the mouth, divide the mouth into 4 sections using the Light and Sound Pacer feature.

The pacer splits the brushing time into 4 equal segments and signals with a short tune with light effect when it's time to move to the next area.


EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushing sessions to help you get used to brushing with the Philips Sonicare.


Note: Each of the first 14 brushing sessions must be at least 1 minute in length to properly advance through the EasyStart ramp-up cycle.

Activating or deactivating the EasyStart feature

Step 1: Place the handle on a plugged in charger.

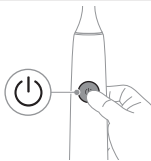
Step 2: Press and hold the power on/off button until you hear a single short beep (for 3 seconds).

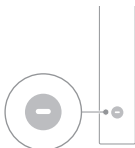If you see the battery indicator flash green 2 times and hear 3 tones low to high, then the feature has been activated.

OR

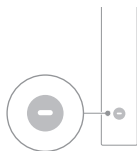If you see the battery indicator flash amber 1 time and hear 3 tones high to low, then the feature has been deactivated.


Step 3: Release the power on/off button.

Charging and battery status

This Philips Sonicare toothbrush is designed to provide on average 42 brushing sessions, each session lasting 2 minutes (21 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

Charging

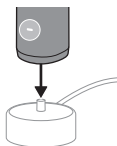- 1 Attach the USB charging cord to a wall adapter, plug the wall adapter into an electrical wall outlet.
- 2 Place the handle on the charger. The handle performs a bottom-to-top swoosh animation, sequentially fades the light ring, battery, and mode indicator lights, and emits two short beeps to confirm proper placement.
- 3 The flashing green light of the battery indicator shows that the toothbrush is charging.
- 4 Once the battery is fully charged, the green light remains solid for 30 seconds and then turns off.

Battery status (when handle is on plugged-in charger)

When the handle is plugged into the charger, the battery indicator at the bottom of the handle shows the battery level.

Battery status	LED color
Charging	Flashing green
Full	Solid green for 30 seconds and then switches off

Battery status (when handle is not on plugged-in charger)

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	LED color
Full	Solid green
Medium low	Flashing green
Low	Flashing amber and 3 beeps after brushing
Empty	Flashing amber and 2 sets of 5 beeps


Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush plugged into the charger when not in use.

Cleaning

Toothbrush handle

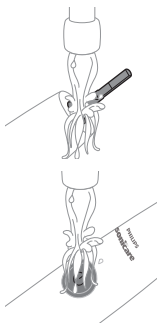- 1 Remove the brush head and rinse the metal shaft area with warm water.


Caution: Do not push on the rubber seal on the metal shaft or around the buttons with sharp objects, as this may cause damage.

- 2 Make sure you remove any residual toothpaste on the handle and around the buttons at least once a week, but preferably after each use.
- 3 Wipe the entire surface of the handle with a damp cloth.

Caution: Do not tap the handle on the sink or counter to remove excess water.

- 4 Dry the handle with a soft cloth or towel.

Brush head

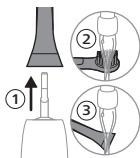- 1 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.
- 2 Rinse the brush head and bristles after each use and allow to air dry.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.
- 3 Dry the charger with a soft cloth or towel.

Travel case

Rinse with warm water and use a damp cloth to clean the travel case and allow to air dry.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean

it and store it in a cool and dry place away from direct sunlight.

Brush head replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Locating the model number

The model number begins with "HX" (e.g. HX123A) and can be found on the bottom of your toothbrush.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.

Problem	Possible Cause	Solution
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2 mm). When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

Full Two-Year Warranty

Philips warrants its products for a specific period of time after the first date of purchase. If the product is purchased in the USA, Philips North America LLC (USA) or, if the product is purchased in Canada, Philips Electronics Ltd., warrants this new product against defects in materials or workmanship for a period of two years from the date of purchase (as evidenced by a valid proof of purchase), and agrees to replace any defective product without charge. Philips will refund the purchase price if it is unable to provide replacement.

Parts that are subject to normal wear and tear, accessories, replacement batteries and consumables such as toothbrush heads, are not covered by the terms of this warranty.

IMPORTANT: This warranty does not cover damage resulting from accident, misuse or abuse, lack of

reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified voltage.*

NO RESPONSIBILITY IS ASSUMED FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES.

For more information and in order to obtain warranty service, simply go to **www.philips.com/support**.

NOTE: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state, province to province or country to country. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above exclusion and limitations may not apply to you. This warranty is transferable.

*Read enclosed instructions carefully.

THE WARRANTY DOES NOT COVER THE FOLLOWING:

- Indirect or consequential damages, such as loss of time, loss of data, or loss of income. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.
- Compensation for activities done by yourself, such as regular maintenance, installing firmware updates, or saving data.

THE WARRANTY DOES NOT APPLY:

- If the proof of purchase, model number, serial number, or production date code on the product has been made illegible.
- If the product was purchased as faulty stock.
- If the defect is the result of abuse or excessive use outside the intended purpose.
- If the defect is caused by environmental conditions, connecting peripherals, additional equipment, or accessories (including software) that are not in accordance with the user manual.
- If the unit has been damaged, including but not limited to damage by animals, lightning, abnormal voltage, fire, natural disaster, transportation, or water (unless the user manual expressly states that the product may be rinsed).

- For non-electrical parts, normal wear and tear, accessories, replacement batteries or replacement of consumable parts such as toothbrush heads.
- If the product does not function properly because it was not originally designed for use in the country where the product is used, which may occur if you import the product.
- If the product does not function properly due to problems with connection with service providers, such as for the internet.

Disposal

CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

- Contact your local town or city officials for electronics recycling and battery disposal information. You can also call 1-800-822-8837 or visit **www.call2recycle.org** for drop-off locations.
- For assistance, visit our website **www.philips.com/support** or call 1-888-744-5477 toll free.
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit: **www.recycle.philips.com**.
- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removal of built-in rechargeable battery

- The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Before removing the battery, make sure that the battery is completely empty.

- Instructions for removal of built-in rechargeable batteries can be found on **www.philips.com/support**. Use the Philips support home search bar to enter the toothbrush model number and find your product. The model number begins with "HX" (e.g. HX123A) and can be found on the bottom of your toothbrush. The battery removal instructions are in the Troubleshooting and Repair section of the product's support page.
- When handling batteries, make sure that hands, the product and the batteries are dry.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

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— |

| —

Philips Oral Healthcare, LLC.
22100 Bothell Everett Highway
Bothell, WA 98021
U.S.A.
1-888-744-5477
www.philips.com/Sonicare
www.philips.ca/sonicare

© 2026 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC
and/or KPNV.

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The Bluetooth®/MD word mark and logos are registered trademarks owned by
Bluetooth SIG, Inc. and any use of such marks by Koninklijke Philips N.V. is under
license.

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Healthcare, LLC et (ou) KPNV.
Bluetooth SIG, Inc., et toute utilisation de ces marques par Koninklijke Philips