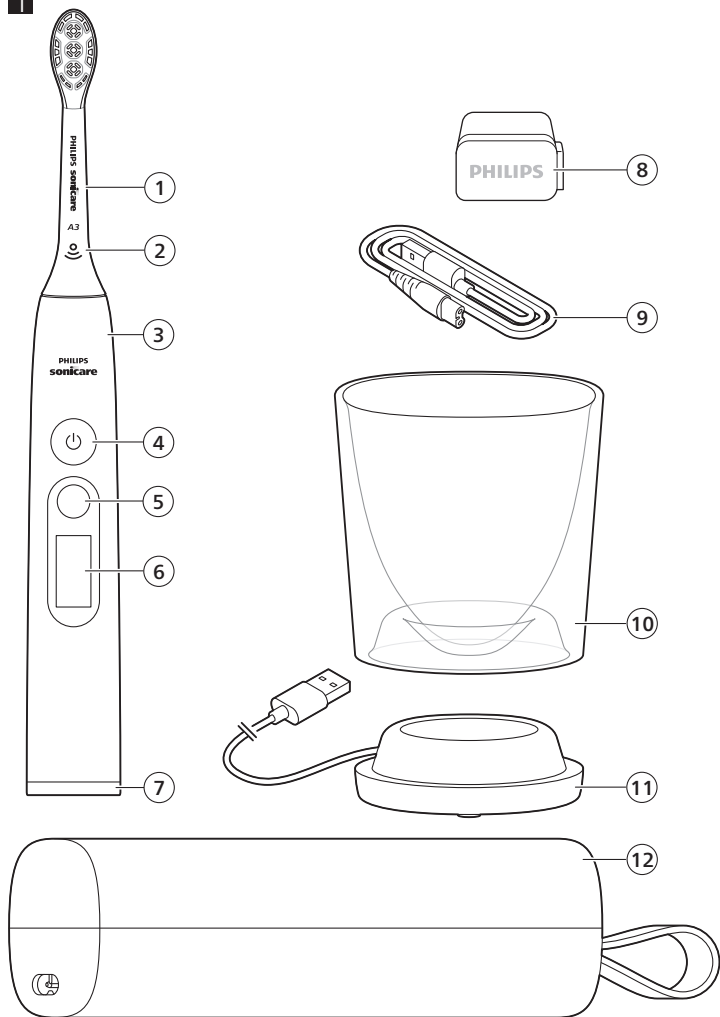
PHILIPS sonicare

9900

1

— |

| —

— |

| —

English 5
Español 43
Français (Canada) 85

Congratulations on your new Next-Generation Philips Sonicare power toothbrush! Superior plaque removal, whiter teeth, and healthier gums are at your fingertips. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time.

Welcome to the Philips Sonicare family!

IMPORTANT SAFETY INSTRUCTIONS

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

WARNINGS

- Keep the charging cables, charging travel cases, and wall adapters away from water.
- After cleaning, make sure charging cables, charging travel cases, and wall adapters are completely dry before connecting to power.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not use charging cables, charging travel cases, and wall adapters outdoors or near heated surfaces.
- Only use original Philips accessories or consumables, if supplied, only use the original USB cable, USB power supply unit, or non-USB power supply unit.
- If the charging cable, charging travel cases or wall adapters is damaged, stop using it.
- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- The maximum altitude of use is 4500 meters.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Make sure brush head and toothbrush are dry before storing in travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.

6 English

- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country.
- The Real Time Coverage feature is intended to support your full-mouth brushing. Users are advised to refrain from brushing beyond the touch-up session to avoid overbrushing. Always follow the advice of a Dental Professional.
- This appliance has only been designed for cleaning teeth and gums.

MEDICAL WARNINGS

- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

BATTERY SAFETY INSTRUCTIONS

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- Charge only with Philips WAA1001 or WAA2001, and HYTCxx or CBA62xx. For the list of approved adapters, please see philips.com/support.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- Charge and use the product at a temperature between 0 °C and 40 °C (32 °F and 104 °F).

- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- This appliance contains a Lithium-ion rechargeable cell. **MUST BE RECYCLED OR DISPOSED OF PROPERLY.**
- The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Instructions for removal of built-in rechargeable batteries can be found on www.philips.com/support.

Storage and transportation conditions

Temperature (Transport and Storage) between -20°C to 60°C (-4°F and 140 °F)

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC Compliance

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).

8 English

- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1 This device may not cause interference.
- 2 This device must accept any interference, including interference that may cause undesired operation of the device.

LIST OF SYMBOLS USED

The following symbols may appear on the product:

This symbol means read the instructions for use, before using the appliance.

XXXXXX

Detachable Power Supply Part, if supplied. The "xxxxxx" indicates the model number of the approved power supply unit to be used.

Intended use

This Sonicare power toothbrush is intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. This Sonicare power toothbrush is intended for consumer home use. Use by children should be with adult supervision.

Your Philips Sonicare (Fig. 1)

- 1 Premium All-in-One brush head (A3)
- 2 BrushSync symbol
- 3 Handle
- 4 Power on/off button
- 5 Confirm Button
- 6 Display
- 7 Brushing Feedback Light
- 8 Wall adapter
- 9 USB-A charging cable
- 10 Charging glass
- 11 Charging base
- 12 Charging travel case

Note: The content of the box may vary based on the model purchased.

App Overview

App Features

The Sonicare app collects information about your brushing habits from smart sensors on your toothbrush, then provides feedback and tips to improve your oral care routine, including:

- Customize your Sonicare toothbrush settings
 - Choose your preferred instant feedback
 - Adjust your brushing intensity
 - Change your brushing mode
- Receive live feedback on your brushing technique during your brushing session
- Track all your oral healthcare activities, including brushing, flossing, mouth rinse, and tongue cleaning
- Receive personalized tips and recommendations to improve your oral health
- View your brushing and oral healthcare habits over time
- Receive brush head replacement reminders
- Stay up to date on the latest features and upgrades

To enhance your oral care routine, visit the Sonicare app for a detailed review and your complete brushing history.

The Sonicare app is compatible with a wide range of mobile phones. To start using the Sonicare app:

Sonicare app - Getting started

- 1 Download the Sonicare app to your phone.
- 2 Ensure that your phone's Bluetooth is turned on.
- 3 Open the Sonicare app and follow the instructions on your phone screen to pair your toothbrush with the Sonicare app.
- 4 Log in to your account via the Sonicare app.

Note: When you first download the Sonicare app from the app store, you will have the latest version, so there's no need for any updates at that time. However, if any updates are available for the toothbrush or the Sonicare app, you will receive a prompt to complete those updates.

For an optimal app experience:

- Sync regularly. When you regularly sync the toothbrush with the app, you will receive real-time tips to help you improve your oral care routine.
- Keep your phone close to your toothbrush for optimal data transfer.
- Make sure your phone's Bluetooth is turned on when using the Sonicare app, allowing your toothbrush to transfer and update your brushing data in the app.

Note: The Sonicare app does not need to be open to collect your brushing data. Your toothbrush stores data and syncs to the app.

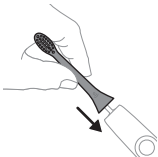
To sync automatically: Allow location permission while setting up the app. By allowing location permission, your phone knows when it is in the connection range of your toothbrush and can refresh your brushing data in the app to provide the latest insights and recommendations.

To sync manually: Connect your toothbrush with the Sonicare app at least every two weeks to keep your brushing habits and personalized insights up to date.

If you have questions about why your brushing data is collected, be sure to review the Philips Dental Privacy Notice, available throughout the app setup process.

Getting Started

Attaching the Brush Head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.
- 3 When a Sonicare Smart brush head is attached to your toothbrush, the toothbrush displays the remaining brush-head sessions and automatically applies the assigned mode and intensity.

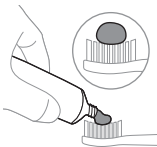
Note: It is normal to see a slight gap between the brush head and the handle.

Using Your Philips Sonicare Toothbrush


If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, use a lower intensity level for the first 1-2 weeks then transition to a higher intensity level.

Brushing Instructions

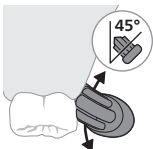
- 1 Wet the bristles and apply a small amount of toothpaste.

12 English

- Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

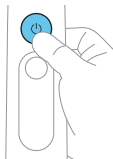Note: Keep the center of the brush in contact with the teeth at all times.


- Press the power on/off button to switch on the Philips Sonicare toothbrush.

- Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout your brushing cycle. Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

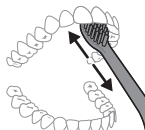Note: When Instant Feedback Pressure is selected, if you apply too much pressure, you will be alerted by a change in vibration of the handle and a flashing brushing feedback light.

- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

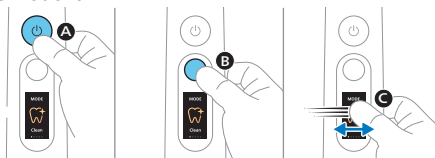
Toothbrush Interface

Your Sonicare toothbrush has two buttons and a display with swipe detection to navigate menu options.

A- Power Button: Turns on and pauses brushing

B- Confirm Button: Confirms your selection on the display with the Confirm button

C- Swipe Display: The display on your toothbrush can react to left and right swipe motions

When you use your Sonicare toothbrush for the first time after purchase or after Factory Reset, you will see the following options on the screen.

First Time Setup

The display on your Philips Sonicare toothbrush wakes up (with enough battery) when you first pick it up. The first screen you will see displays the language options, with English selected by default. You can swipe left or right to browse through other language choices.

Note: You may need to press the On/Pause or Select Button to wake up the toothbrush.

Note: We recommend to fully charge your device before starting to use it. If your toothbrush's battery is completely empty (at 0%), it will show a "battery empty" screen and turn off (refer to Battery Status). Pressing the On/Pause or Select button to wake up the toothbrush will produce two sets of five beeps. You won't be able to proceed with the first-time setup; you will only see the empty battery screen before the device powers off again.

Select Preferred Language

Swipe left through the Language options and press the Confirm Button when your preferred language is displayed. Swipe right to go back through the previous options.

Navigating Your Display

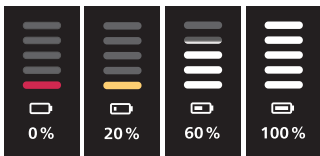
Your Sonicare toothbrush features a swipe display with a menu system that helps you keep track of all its available features and benefits.

To navigate the display, follow these steps:


- 1 Select language.

Note: Additional languages may be available for selection.

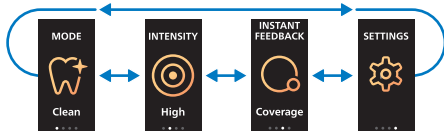- 2 A welcome animation appears on the screen, followed by the Battery Status screen.

- 3 The display will then show the Main Menu selection.

Note: If a smart brush head is attached to the handle during the First-Time Setup, the device will first display the Battery Status screen. After that, it will show the remaining brush head sessions screen before proceeding to the Main Menu. Depending on the status of the brush head, the screen may also display a "Brand New Brush Head" celebration animation.

Main Menu

When the Main Menu is referenced, this refers to the top level of menu options available whenever your toothbrush is first awakened. The Main Menu will contain options such as Mode, Intensity, Instant Feedback, and Settings.

Selection Menu

The Selection Menu will appear after you confirm your initial choice on the Main Menu. The choices presented are specific to the category you have chosen on the Main Menu.

Note: The Selection Menu has options such as different modes and intensities, as well as Instant Feedback options like Pressure + SenseIQ, Coverage, and Scrubbing Feedback.

An example of navigating the Selection Menu is shown in the following steps:

16 English

- 1 From the Main Menu, swipe left or right to go to the Mode screen.

- 2 On the Mode screen, press the Confirm Button to access the Selection Menu of modes.

- 3 A check mark in the top right corner of the screen indicates the currently selected mode.

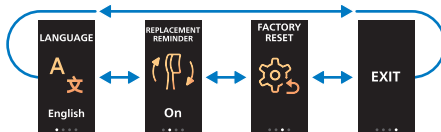
- 4 Swipe to browse available modes. All Sonicare brush heads work with each mode.
- 5 Press the Confirm Button to select a mode. A blue check mark confirms your choice and returns you to the Main Menu.
- 6 To exit without changes, select Exit in the Selection Menu.
- 7 The selected mode appears on the Mode screen.

Note: Sub-menu selections are available in Settings. To access those options, follow the steps below:

- 1 On the Main menu screen, select Settings.

- 2 On the Settings Selection Menu swipe to choose options such as Language, Replacement Reminder, or Factory Reset.

Modes

Your Sonicare toothbrush features up to eight different modes. After you wake it up, your Sonicare toothbrush will default to the Modes screen. To manually select your desired mode, follow the steps below:

- 1 From the Main Menu, swipe left or right to go to the Mode screen.

18 English

- 2 On the Mode screen, press the Confirm Button to access the Selection Menu of modes.
- 3 A check mark in the top right corner of the screen indicates the currently selected mode.

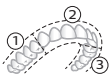
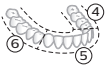
- 4 Swipe to browse available modes. All Sonicare brush heads work with each mode.

- 5 Press the Confirm Button to select a mode. A blue check mark confirms your choice and returns you to the Main Menu.
- 6 The selected mode appears on the Mode screen.
- 7 To exit without changing modes, select Exit in the Selection Menu.

Note: Always use original Sonicare brush heads with BrushSync technology as they are designed to work with every mode in a Sonicare toothbrush.

Modes	Clean	Deep Clean	Gum Health	Sensitive
Mode Image				
Benefits	Plaque removal	Superior plaque removal	Plaque removal and gum massage	Extra gentle mode for sensitive teeth and gums
Recommended intensity level	High	High	High	Low
Total brushing time	2 minutes	3 minutes	3 minutes and 20 seconds	2 minutes
How to brush	Brush each segment for 20 seconds.	Brush each segment for 30 seconds.	Brush each segment for 20 seconds, then brush segments 1,3,4, and 6 for 20 seconds each.	Brush each segment for 20 seconds.
			 	

Modes	White	Comfort	Max Care	Polish
Mode image				
Benefits	Plaque and surface stains removal	Gentle gum stimulation	For a deep clean	For extra polish
Recommended intensity level	High	High	High	High
Total brushing time	2 minutes and 40 seconds	2 minutes	3 minutes	2 minutes
How to brush	Brush each segment for 20 seconds, then brush segments 2 and 5 for 20 seconds each.  	Brush each segment for 20 seconds.	Brush one segment for 20 seconds, then brush the same segment for 10 seconds, repeat for all segments.	Brush each segment for 20 seconds.

Intensities

Our Sonicare toothbrush has three intensity settings: High, Medium, and Low. To change the intensity manually:

- 1 From the Main Menu, swipe left or right to go to the Intensity screen.
- 2 On the Intensity screen, press the Confirm Button to access the Selection Menu of Intensity.
- 3 A check mark in the top right corner of the screen indicates the currently selected intensity.
- 4 Swipe to browse available Intensities.
- 5 Press the Confirm Button to select an intensity. A blue check mark confirms your choice and returns you to the Main Menu.
- 6 The selected intensity appears on the Intensity screen.
- 7 To exit without changing Intensities, select Exit in the Selection Menu.

Note: The intensity can be adjusted either prior to brushing or during the brushing by pressing the Confirm Button.

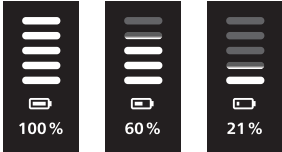
Battery Status

This Philips Sonicare toothbrush is designed to provide up to 14 days of regular brushing on a full charge. The battery level screen will briefly appear on the display:

- When the toothbrush awakens.
- Before it goes to sleep, after 30 seconds of inactivity.
- When the toothbrush is placed on the charger and begins charging.


Note: The battery can take up to 24 hours to fully charge before first use, but you can use the Sonicare toothbrush before it is fully charged.

Battery Status


Battery Status Display	Battery Level	Battery Status' Color
	Battery Empty (Two sets of five beeps will indicate an empty battery)	Red (flashing)
	Low	Yellow
		
	Medium to Full	White
		
		

Charging

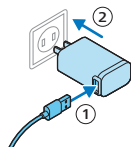
Charging on the charger

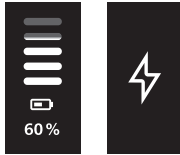- 1 Plug the USB cord of the charging base into a power supply unit.

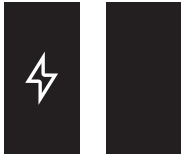
- 2 Place the toothbrush in the plugged-in charging base with the glass on top. The battery status screen appears briefly, then switches to a charging icon with a blinking symbol. The toothbrush emits two beeps to indicate that charging has started.


- 3 Leave the toothbrush on the charger until it is fully charged. The display will turn off when the toothbrush has finished charging.

Note: It may take up to 24 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.

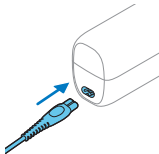
Charging with the travel case

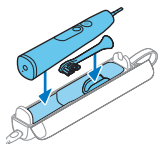
- 1 Plug the charging cord into the travel case and into a power supply unit.

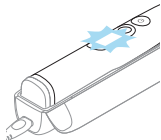- 2 Place toothbrush in the travel case.

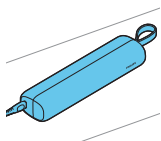- 3 The battery status will appear on the display with a charging icon while the toothbrush emits two beeps to indicate that charging has started.

- 4 Leave the toothbrush in the travel case until fully charged. The display will turn off when the toothbrush has finished charging.

Note: Make sure the travel case is placed on its side for better stability.

Note: Not all models include a charging travel case.

Select Your Instant Feedback

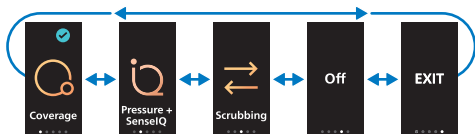
Your Sonicare toothbrush provides instant feedback during brushing that helps you enhance your brushing technique, including motions and habits. Only one type of feedback can be activated at a time on the handle, and the light ring only indicates the selected feedback.

Note: Within the Sonicare app, you can view all your feedback results, including those not selected on the handle.

To activate these features on your toothbrush:

- 1 From the Main Menu, swipe left or right to the Instant Feedback screen and press the Confirm Button.

- 2 On the Selection Menu, swipe left or right between the different types of feedback and press the Confirm Button for the instant feedback you want.
- 3 A check mark appears on the screen to confirm that your selected feedback is active.
- 4 To exit without changes to instant feedback, select Exit in the Selection Menu.

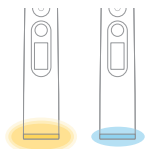
Note: If you select an alternative instant feedback, that feedback will become your new default.

Instant Feedback: Coverage (only 9900/9000)


Your Sonicare toothbrush provides guidance on how long to brush each segment of your mouth, supporting you in achieving complete and even coverage. The light ring gives you instant feedback throughout your brushing session. While you brush, the toothbrush detects the location in your mouth, and the light ring glows yellow to indicate any segments that haven't been brushed long enough.

Once you've spent sufficient time cleaning a particular segment, the toothbrush beeps and the light ring glows solid blue to signal that segment has been thoroughly brushed, alerting you to move on to the next segment. As you move the toothbrush around your mouth, the feedback updates in real time, reflecting whether you have spent long enough time on that segment.


26 English

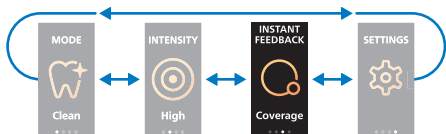To access or change this Instant Feedback selection:

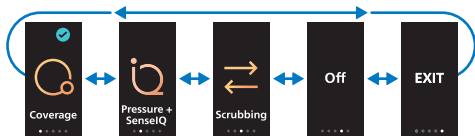- 1 From the Main Menu, swipe left or right to the Instant Feedback screen and press the Confirm Button. Upon first-time use, the default instant feedback for the 9000/9900 Series is Coverage.


- 2 If you previously selected a different instant feedback, that selection will be shown on the display. To change it, swipe to your desired instant feedback and press the Confirm Button.

Best Practices to Increase the Accuracy of Instant Feedback Coverage

- **Follow these "Brushing Instructions"**
 - Keep the brush head at the correct angle to the gum line as stated in the Brushing Instructions.
 - Finish brushing the current segment before you move on to the next. Your toothbrush will alert you once you have finished brushing the current segment.
- **Minimize Body Movement**

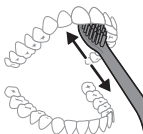- Maintain a steady posture while brushing. Avoid excessive head or body movement, including bending or rotating, which can interfere with sensor accuracy.
- **Minimize Interruptions During Brushing**
 - Frequently removing and reinserting the toothbrush during a session can disrupt the toothbrush from tracking its location.


Note: To manually improve the accuracy of the Instant Feedback Coverage, brush the inside of the upper or lower front teeth for two seconds to reset the position.

Instant Feedback: Pressure + SenseIQ (only 9900)


Your Sonicare toothbrush is designed with the Pressure + SenseIQ feature to protect your teeth and gums. If you apply excess pressure for an extended period of time, your toothbrush will automatically lower the intensity by 1 level. You will experience a brief pause in the pressure feedback and feel the intensity adjust to the next lower level. At the end of your session, the handle's display shows a mouth map that highlights only the segments with excessive pressure for an extended time.

Note: Every new brushing cycle the intensity will go back to your pre-selected setting.

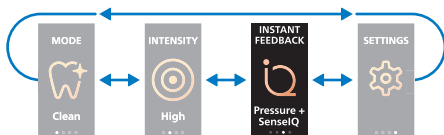
By default, Pressure Sensor + SenseIQ is turned off. To turn this feature on:

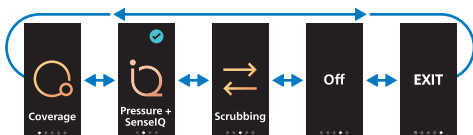- 1 From the Main Menu, swipe to the Instant Feedback screen and press the Confirm Button.

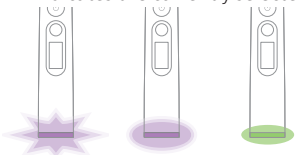- 2 Swipe to the Pressure + SenseIQ and press the Confirm Button to select this instant feedback. A check mark in the top right corner of the screen indicates the currently selected Instant Feedback.

During brushing, the light ring gives you instant feedback on how much pressure you are using. It glows green when using the optimal amount of pressure, pulses purple with a gentle vibration if pressing too hard, and slowly fades to purple in an up-and-down motion when you're not applying enough pressure.

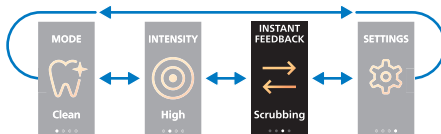
Instant Feedback: Scrubbing

Your Sonicare toothbrush measures the motions while brushing to detect scrubbing behavior (see brushing instructions for optimal technique).


When Scrubbing Feedback is activated, the light ring at the bottom of the handle glows pink, and the toothbrush changes its vibration as a reminder to stop scrubbing. The light ring turns off when you stop scrubbing.

To turn this feature on:

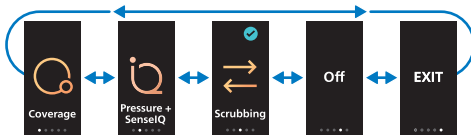
- 1 From the Main Menu, swipe to the Instant Feedback screen and press the Confirm Button.

- 2 Swipe to the Scrubbing and press the Confirm Button to select this instant feedback. A check mark in the top right corner of the screen indicates the currently selected Instant Feedback.

Note: When Instant Scrubbing Feedback is selected, the displayed Mouth Map is the Coverage Mouth Map.

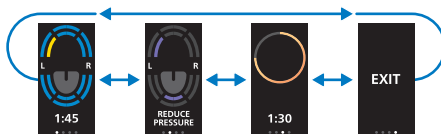Post Brushing Feedback


When you finish or pause your brushing session, you can swipe the screen to see additional brushing information. Depending on the model of your

toothbrush, this includes separate screens for the Coverage Mouth Map, Pressure Mouth Map, Brushing Duration, and the option to exit.

Note: You must brush your teeth for at least 30 seconds before the Mouth Map appears.

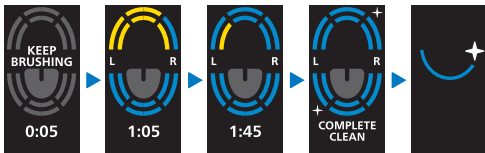Understanding the Mouth Map (only for 9900/9000)

The Coverage Mouth Map is the first thing you will see in post-brushing feedback. After each brushing session, the display will show a mouth map consisting of 12 segments. The outer segments represent the outside of your teeth, and the inside segments represent the inside of your teeth. The letters L and R on the side help show your mouth's left and right sides. The tongue icon helps you distinguish between your top and bottom teeth.

Note: You can also view a more detailed mouth map within the Sonicare app.

Post-Brushing Feedback: Coverage (only 9900/9000)

To help you achieve complete and even coverage, when you finish or pause your brushing, the display will show a mouth map consisting of 12 segments of your teeth. The segments that have been thoroughly brushed will appear blue, and segments that need more brushing time appear yellow. You need to brush for at least 30 seconds before this coverage map is available.

Once all areas have been thoroughly brushed, the toothbrush will celebrate with an all-clean sound, and all areas will show solid blue. You can turn the toothbrush back on to touch up any segments you have not completely brushed. The Mouth Map will be updated to reflect the retouched areas based on your brushing.

Note: Make sure to touch up any missed segments before your toothbrush goes to sleep; otherwise, the toothbrush will treat the brushing as a brand new session. We recommend performing only one touch-up session to avoid overbrushing.

Post-Brushing Feedback: Pressure

If you apply too much pressure in any area of your mouth for an extended period (5 consecutive seconds), that segment on the mouth map will be shown in purple. This is a reminder to be gentler when pressing the brush head against your teeth and gums. Pressing too hard reduces efficacy and wears out your brush head faster. The gray area on the Mouth Map represents either optimal or under pressure.

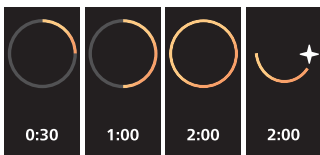
Note: The Post-Brushing Pressure Mouth Map will only appear if you apply too much pressure for an extended time on at least one segment. If no overpressure (for an extended period) occurs, no Post-Brushing Pressure Feedback will be displayed.

Note: You can also view your results in the Sonicare app.

Post-Brushing Feedback: Duration

The display shows how long you have been brushing, and the circle shows the progress of your brushing session. The display will show a full circle when the brushing session is completed.

Note: Brushing duration may vary based on the mode you selected.

BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The  symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- BrushSync mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to

www.philips.com/toothbrush-heads for more information.

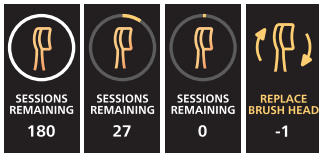
BrushSync Mode Pairing

Sonicare brush heads with BrushSync technology are equipped with a microchip that communicates with the handle and automatically pairs it with the recommended mode and intensity. Most brush heads with BrushSync technology will automatically pair to Clean Mode. Brush heads with specific benefits will pair with the corresponding mode (when available), examples include:

- G3 Premium Gum Care: Gum Health Mode
- W3 Premium White or W DiamondClean (NAM) / W OptimalClean (ROW): White Mode
- S2 Sensitive: Sensitive Mode

Note: If you choose to change the mode/intensity from the recommended setting, the brush head will remember your chosen setting for future brushing sessions.

Brush head life


When you attach a genuine Sonicare smart brush head to your Sonicare toothbrush for the first time, the display will indicate the number of remaining brushing sessions, which is 180, until a brush head replacement is recommended. The “Replace Brush Head -1” message appears when the number of brush head sessions remaining drops below 0. This message reminds the user to change the brush head.

Note: Replace Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Brush Head Replacement Reminder

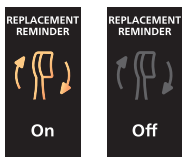Your Sonicare toothbrush is equipped with BrushSync technology, which tracks the wear of your brush head. The Brush Head Replacement Reminder will notify you when the number of brush head sessions remaining drops below 0. This feature guarantees the best cleaning and care of your teeth.


By default, the Brush Head Replacement Reminder feature is enabled. To manually turn this feature on or off on your Sonicare toothbrush, follow the steps below:


- 1 From the Main Menu, swipe left or right to go to the Settings screen.
- 2 On the Settings screen, press the Confirm Button to access the Selection Menu.

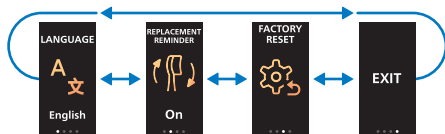34 English

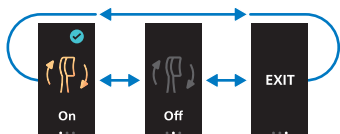
- 3 Swipe to the Replacement Reminder screen and press the Confirm Button. A check mark in the top right corner of the screen indicates the current selection.

- 4 Swipe to browse On/Off or Exit.
- 5 Press the Confirm Button to make your selection. A check mark confirms your choice and returns you to the Settings sub menu. To return to the main Menu, select Exit.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

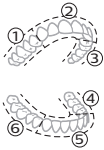BrushPacer

The BrushPacer divides the brushing time into 6 mouth segments and indicates when to move to the next area via a brief pause in vibration. The toothbrush automatically stops at the end of the brushing session.

When Instant Feedback Coverage is activated, the BrushPacer continues to operate but adapts to the Coverage guidance to guide you between segments. Instant Feedback Coverage divides the mouth into 12 segments and provides guidance based on the time spent in each segment.

Note: If any other type of Instant Feedback is selected, the BrushPacer remains active.

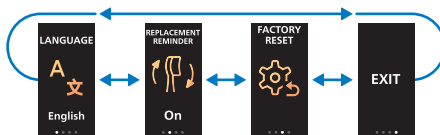
Settings

Factory Reset


To reset all settings and restore your Sonicare toothbrush's original factory settings, follow the steps below:

- 1 From the Main Menu, swipe to go to the settings screen.
- 2 On the Settings screen, press the Confirm Button.

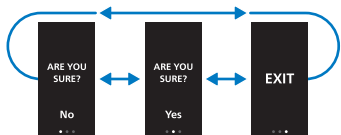- 3 Swipe left or right until the display is on the Factory Reset screen. Press the Confirm Button to confirm your selection.


- 4 A screen appears with the message "Are you sure?"
- 5 Swipe left or right to select "Yes" or "No" (the Exit screen is also available to cancel the Factory Reset). If "Yes" is selected, a Resetting status screen will appear.
- 6 You will know that your Factory Reset was successful if the display prompts you to select your desired language.

Language

The default language of your toothbrush is English. To change the language of your display:

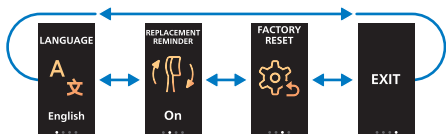- 1 From the Main Menu, swipe to the Settings screen.

- 2 On the Settings screen, press the Confirm Button to access the Selection Menu of languages.

- 3 A check mark in the top right corner of the screen indicates the currently selected language.

- 4 Swipe to browse available languages.
- 5 Press the Confirm Button to select the language of your choice. A check mark confirms your choice and returns you to the Main Menu. To exit without changes, select Exit in the Selection Menu.
- 6 The selected language appears on the Language screen.

Additional Features

SmarTimer

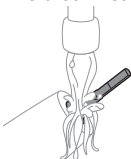
The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing session. Dental professionals recommend brushing for at least 2 minutes twice a day.

Note: If you press the power button after you start the brushing session, the toothbrush will pause. After a 30-second pause, the SmarTimer resets.

Cleaning

Brush Head and Toothbrush

The brush heads and toothbrush can be cleaned by rinsing it with warm water.

- 1 Remove the brush head from the toothbrush and rinse it thoroughly.
- 2 Rinse the entire toothbrush, especially the brush head connection. Gently clean around the rubber seal. Repeat at least once a week.

Note: Do not push on the rubber seal at the top of the toothbrush.

Travel Case and Charger

- 1 Unplug the charger and travel case before you clean them.
- 2 Use a damp cloth to wipe the surface of the charger and travel case.

Danger: Unplug the USB wall adapter and chargers before cleaning them.

Cautions

- Do not clean the product or accessories in the dishwasher.
- Do not use isopropyl alcohol, vinegar, bleach, or any other household cleaning products to clean the product or accessories as this may cause discoloration.
- Make sure brush head and toothbrush are dry before storing in travel case.
- Do not use essential oils to clean the brush head, the product or accessories as this can cause damage.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Locating the Model Number


Look on the bottom of the Philips Sonicare toothbrush handle for the model number.

Frequently asked questions

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		A small gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the handle.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm).
The toothbrush stops working.	The battery is empty.	Charge the toothbrush handle on the charger.
I charged my toothbrush, but it only works for a short time.	The toothbrush handle is not placed vertically on the charger (with plastic cover or glass on top), resulting in no charge.	Place the handle vertically on the charger (with plastic cover or glass on top) and make sure the charger is plugged in.

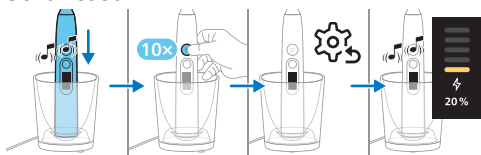Soft Reset

Follow the below steps to fix any screen issues or lights that are not lining up:

- 1 Place the handle on a valid and plugged-in charger base to start charging. You can confirm this by the charging UI and the sound of two beeps.
- 2 Press the on/pause button 10 times consecutively.

- 3 The handle will then perform the reset process.
- 4 Once the reset is complete, the handle will wake up and return to the charging state, this will be confirmed again by the charging UI and two beeps.

Warranty and support

Philips warrants its products for two years after the date of purchase. Defects due to faulty materials and workmanship will be repaired or replaced at Philips' expense provided that convincing proof of purchase in the qualifying period is provided.

Philips recommends use of genuine Sonicare brush heads for optimal performance. Use of other brush heads may damage your Sonicare toothbrush and limit your warranty protection.

Contact our Consumer Care Center at 1-888-744-5477 (North America), outside North America contact your local Philips Consumer Care Center.

Internet information: www.philips.com/support.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY

SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Software Licenses and Disclaimers

Apache 2.0 open-source software license disclaimer

Copyright [yyyy] [name of copyright owner]

Licensed under the Apache License, Version 2.0 (the "License"); you may not use this file except in compliance with the License. You may obtain a copy of the License at

<http://www.apache.org/licenses/LICENSE-2.0> Unless required by applicable law or agreed to in writing, software distributed under the License is distributed on an "AS IS" BASIS, WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied. See the License for the specific language governing permissions and limitations under the License.

BSD 2-Clause open-source software license disclaimer

Copyright (c) 2013, Mike Ryan,

Based on micro-ecc, Copyright (c) 2013, Kenneth MacKay

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

* Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.

* Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE

DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT HOLDER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

MIT open-source software license disclaimer

Copyright (C) 2020 Amazon.com, Inc. or its affiliates. All Rights Reserved. Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions:

The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT, TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

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©2026 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.
de Philips Oral Healthcare, LLC o KPNV.