
**PHILIPS**  
**sonicare**

**8000 - 9900**

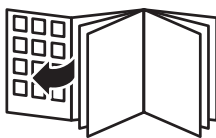— |

| —

— |

| —

English 6  
Azərbaycanca 9  
Bahasa Melayu 12  
Čeština 15  
Dansk 18  
Deutsch 21  
Eesti 24  
Español 27  
Español 30  
Français 33  
Français (Canada) 36  
Hrvatski 39  
Italiano 42  
Latviešu 45  
Lietuviškai 48  
Magyar 51  
Nederlands 54  
Norsk 57  
O'zbek 60  
Polski 63  
Português 66  
Português do Brasil 69  
Română 72  
Shqip 75  
Slovenščina 78  
Slovensky 81  
Srpski 84  
Suomi 87  
Svenska 90  
Tiếng Việt 93  
Türkçe 96  
Türkmençe 99  
Ελληνικά 102

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## Removing the rechargeable battery

The built-in rechargeable battery can be removed for disposal. If you are not comfortable disassembling electronics, please seek professional assistance at your local electronics recycling or repair shop.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

**Be sure to protect your eyes, hands, fingers, and the surface on which you work.**

When you handle batteries, make sure that your hands, the toothbrush, and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

Do not reinstall battery once it has been removed for recycling. Consumer installation of battery may increase risk of fire or other potential safety hazards.


**If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.**

## Removing the rechargeable battery

- Tools: T8 TORX screwdriver, #1 Phillips-head screwdriver with a shaft diameter of max. 3.5 mm and plastic pry tool (or guitar pick).
- 1 When possible, deplete the rechargeable battery of any charge. To do this, remove the toothbrush from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare toothbrush.
- 2 Remove the brush head from the toothbrush.


- 3 Using a T8 TORX screwdriver, unscrew and remove (counterclockwise) the screw from the base of the toothbrush handle.

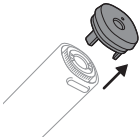- 4 Remove the endcap from the base of the toothbrush. If the endcap is tight and does not release easily from the housing tube, use a plastic pry tool to help separate the endcap from the housing tube.


- 5 Use a Phillips screwdriver to rotate the metal locking tab inside the bottom center of the toothbrush. Turn counterclockwise until rotation stops (less than a 1/4 of a turn).

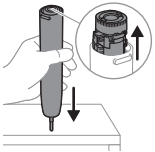- 6 To release the internal components from the housing tube, firmly press the metal shaft onto a solid surface (do not slam the handle down with force). If the internal components do not easily release from the housing tube, double check that you have completed step 5 above.
- 7 Gently pull the internal assembly to remove it from the base of the housing tube.

**Caution:** To avoid injury or damage to the product, be aware of the sharp edges and exposed electronic components.

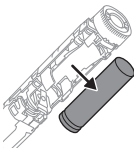- 8 Pinch the battery and pull it away from the battery tray.

**Caution:** Do not use screwdriver or metal object to pry the battery out of the battery tray.

- 9 For battery disposal, cover the battery contacts of the used battery with tape to prevent any electrical short from residual battery charge. The used rechargeable battery can now be recycled through your community battery or electronics disposal process.

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©2026 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.