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












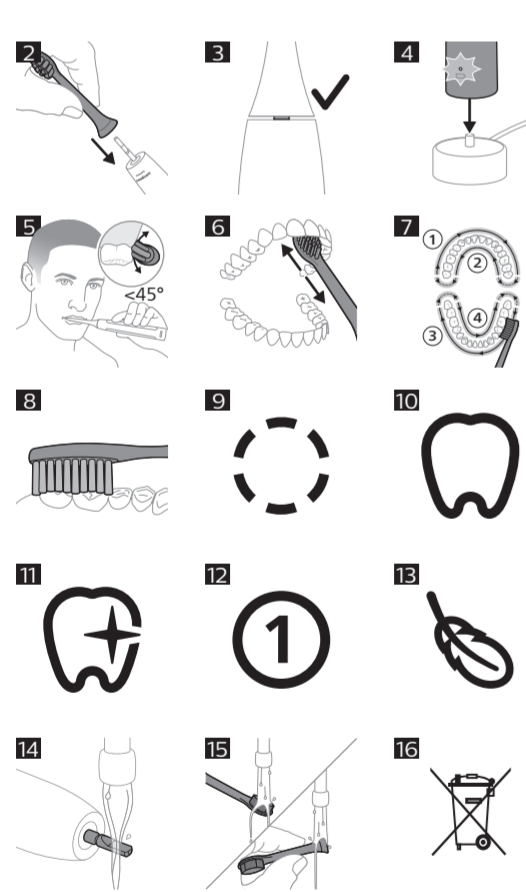
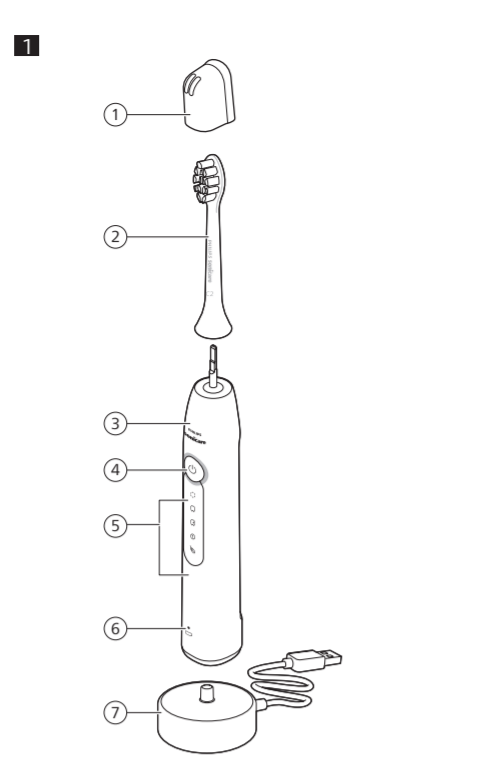
## 1
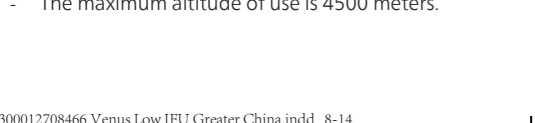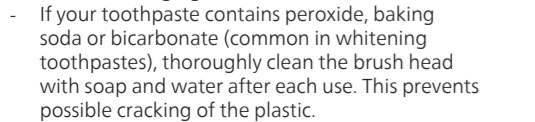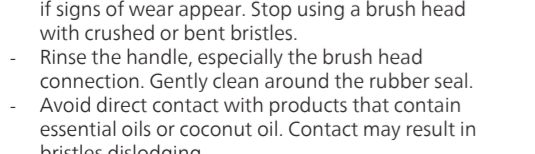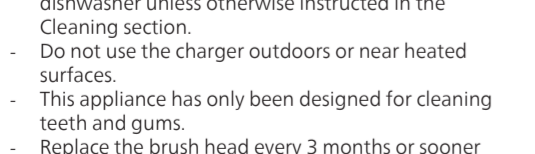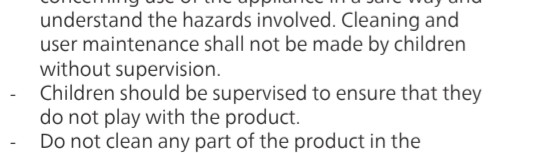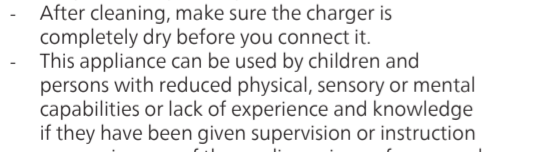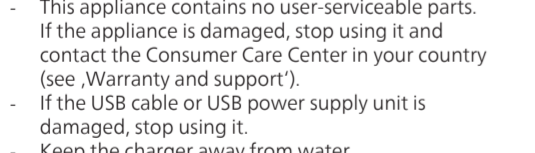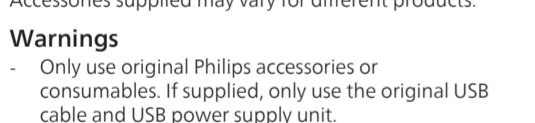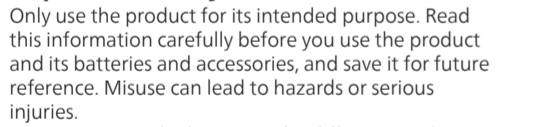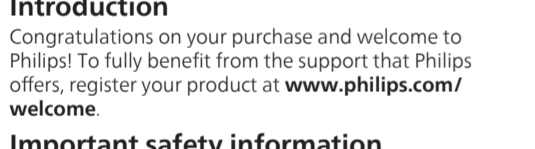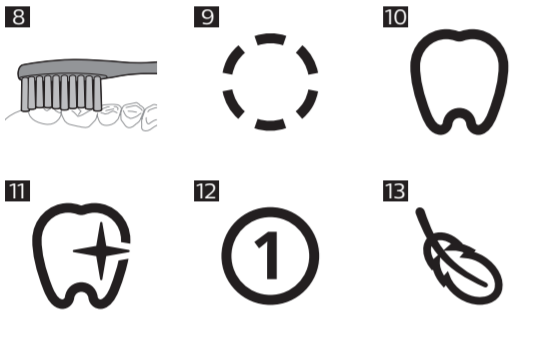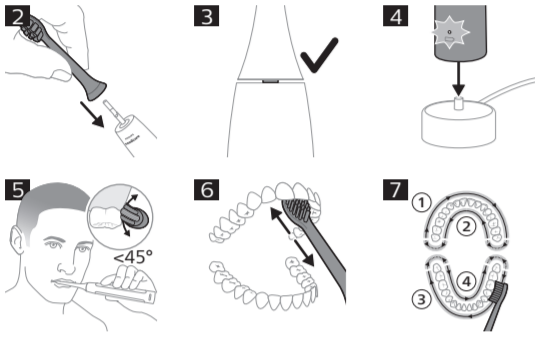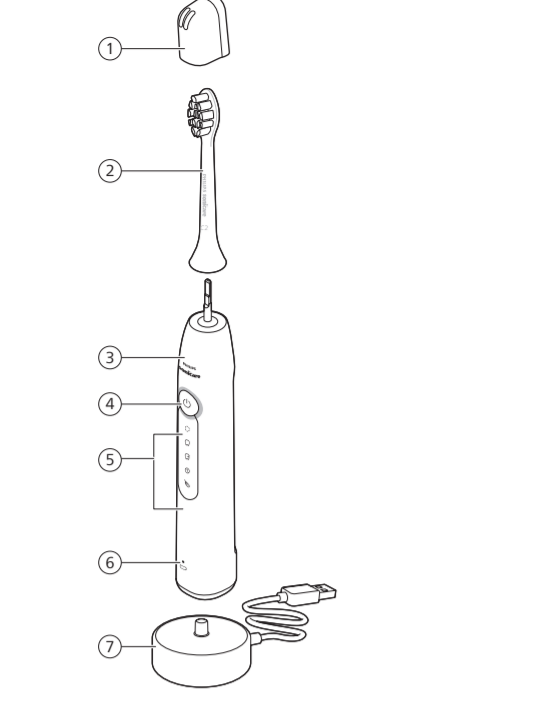
## Medical warnings

- Stop using the appliance and consult your dentist/ doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

## Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5 V ± 2V. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country. Using a non-certified supply unit may cause hazards or serious injuries.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).

## Storage and transportation conditions

Temperature -20°C to +60°C (-4 °F to 140 °F).

## Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Your Philips Sonicare (Fig. 1)

① Brush head cap

② Brush head

③ Handle

④ Power and mode button

⑤ Brushing modes

⑥ Battery indicator

⑦ Charger

Note: Wall adapter not included. Appearance may slightly differ, depending on the model.

## Getting started

## Attaching the brush head

- Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 2).

- Firmly press the brush head down onto the metal shaft until it will not go any further (Fig. 3).

## Charging your Philips Sonicare

- Place the power toothbrush handle on the charger (Fig. 4).
- Connect the charger to a certified safety extra low voltage (SELV) supply unit with output rating 5V ± 2V.
- Your toothbrush continues to charge until it is fully charged.
  - When the toothbrush is charging, the battery indicator pulsates green.
  - When fully charged, the battery indicator will light up solid green for 30 seconds before turning off.

Note: It can take up to 10 hours to fully charge the battery, but you can use the toothbrush before it is fully charged.

## Battery status

This Philips Sonicare toothbrush is designed to provide at least 10 hours brushing sessions after fully charged for 10 hours, each session lasting 2 minutes. When the battery is low, the battery indicator starts to flash amber.

Note: When the battery is completely empty, the toothbrush switches off. Place the handle on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your Philips Sonicare toothbrush on the charger when not in use.

Note: Battery usage time is expected to reduce over product life time.

## Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow instructions below for the best experience.

## Operation instructions

**Power on/Pause**  
Press the power button briefly to turn the toothbrush on or pause operation.

**Switch brushing modes**  
Press and hold the power button to switch between brushing modes. Release the button to enter the selected mode.

**Soft start**  
All models start with a gentle vibration for the first 3 seconds to prevent toothpaste splatter and to help you gradually get used to the vibration intensity.

## Brushing modes

Your power toothbrush is equipped with five different modes:

**Easy**  
Designed with gentle intensity to ease the transition from manual brushing, ideal for first-time user (2 minutes) (Fig. 9). It is recommended to use this mode for the first 14 days to build comfort before exploring other modes.

**Clean**  
It removes plaque effectively for daily cleaning (2 minutes) (Fig. 10).

**Deep Clean**  
It helps to remove plaque (3 minutes) (Fig. 11).

**Refresh**  
For removing stains from teeth (1 minute) (Fig. 12).

**Sensitive**  
Gentle cleaning for sensitive gums and teeth (2 minutes) (Fig. 13).

Note: When you use the toothbrush for the first time, the default mode is ‚Easy‘.

## Brushing instructions

- Wet the bristles and apply a small amount of toothpaste.
- Place the toothbrush bristles against the teeth at a slight angle (45 degrees) (Fig. 5). The handle has two 45 - degree indication bars to help guide your Bass brushing technique. Press the brush head gently against the gum to make the bristles reach the gum line or slightly beneath the gum line.
- Press the power button to turn on the Philips Sonicare.
  - Note: Keep the center of the brush in contact with the teeth at all times.
- Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.
  - Note: The bristles should slightly flare. Do not scrub. It is not recommended to scrub as you would with a manual toothbrush.
- To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 6).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections using the QuadPacer feature (see chapter ‘Features‘).

- Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth) (Fig. 7).
- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs (Fig. 8).

Your Philips Sonicare toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

## Features

### Smarterim

The Smarterim indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle. Dental professionals recommend brushing no less than 2 minutes twice a day.

### QuadPacer

The Quadpacer is an interval timer that beeps every 30 seconds to remind you to brush all teeth evenly by dividing your mouth into 4 sections (Fig. 7). Refresh mode (1 min) beeps every 15 seconds.

### Pressure sensor

Your Sonicare toothbrush measures the pressure you apply while brushing to protect your gums and teeth from damage. If you apply excess pressure, the handle changes its vibration until you reduce the pressure.

### Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

**Warning: Do not clean brush heads, the handle or the charger in the dishwasher.**

**Warning: Do not clean with boiling hot water.**

### Toothbrush handle

- Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste on the handle and around the buttons at least once a week, but preferably after each use (Fig. 14).

**Caution: Do not push on the rubber seal on the metal shaft or around the buttons with sharp objects, as this may cause damage.**

- Wipe the entire surface of the handle with a damp cloth.

**Caution: Do not tap the handle on the sink or counter to remove excess water.**

- Dry the handle with a soft cloth or towel.

### Brush head

- Rinse the brush head and bristles after each use and allow to air dry (Fig. 15).
- Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

### Charger

- Unplug the charger before you clean it.
- Wipe the surface of the charger with a damp cloth.
- Dry the charger with a soft cloth or towel.

## Brush head replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

- Replace the brush head if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- Smart fading bristles provide a visible replacement reminder: when the blue bristles fade to white at one-third of their length, replace the brush head.

## Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter ‘Cleaning’) and store it in a cool and dry place away from direct sunlight.

Charge the toothbrush regularly to prevent the battery from discharging completely. A completely discharged battery may fail to charge again.

## Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (e.g HX123A).

## Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance: If you are unable to solve the problem with the information below, visit **www.philips.com/support** for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm). <p>When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 10 hours.</p>
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 10 hours.

I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 10 hours.
--	---	--

## Warranty and support

Philips offers a 2-year warranty. Warranty exceptions apply. This warranty does not affect your rights as a consumer under law. For more information, visit **www.philips.com/support** or call us.

### Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 16).
- Follow your country’s rules for the separate collection of electrical products and batteries.

## Removing the rechargeable battery

**CONTAINS LITHIUM-ION RECHARGEABLE CELL. MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.**

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the appliance is disconnected from the wall socket and that the battery is completely empty.

- Remove the inner body from the handle.
- Remove the rechargeable battery from the inner body.

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©2026 Koninklijke Philips N.V. (KPNV). 保留所有权利。 Philips Sonicare 是 Philips Oral Healthcare, LLC 的注册商标。 Philips Sonicare 5300 Series 是 Philips Oral Healthcare, LLC 的注册商标。