
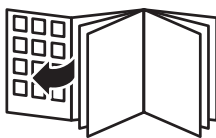
PHILIPS

sonicare

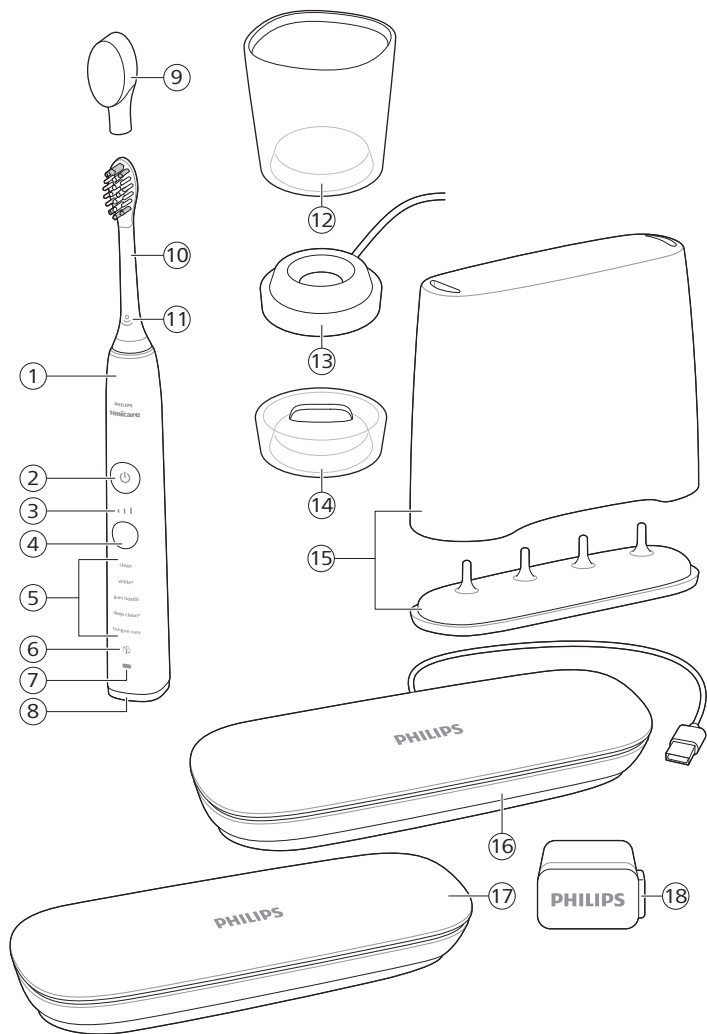*DiamondClean
Smart*

English 6
Español 30
Français (Canada) 57

Introduction

Welcome to the Philips Sonicare family!
This toothbrush enables you to obtain superior plaque removal, whiter teeth and healthier gums. Using Sonicare's combination of gentle sonic technology and clinically developed and proven features, you can be confident that you're getting the very best clean, every time.

Further support and product registration are available to you at: **www.philips.com/support**

IMPORTANT SAFETY INSTRUCTIONS

- Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.

WARNINGS

- Keep the charging cables and charging travel cases away from water.
- After cleaning, make sure the charging cables and charging travel cases are completely dry before connecting to power.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not clean any part of the product in the dishwasher unless otherwise instructed in the Cleaning section.
- Do not use charging cables or charging travel cases outdoors or near heated surfaces.
- Only use original Philips accessories or consumables. Only use the power supply unit provided with the product.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support').
- If charging cables or charging travel cases are damaged, stop using it.
- This appliance has only been designed for cleaning teeth and gums.

- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- The maximum altitude of use is 4500 meters.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.
- Make sure brush head and toothbrush are dry before storing in travel case.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.

MEDICAL WARNINGS

- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.
- This appliance is a personal care device and is not intended for use on multiple patients in a dental practice or institution.

BATTERY SAFETY INSTRUCTIONS

- This appliance contains non-replaceable batteries that are only to be removed by skilled persons.
- Accessories supplied may vary for different products.
- Only use original Philips accessories and consumables.
- Charge only with Philips CBAxxxx (input ratings: 100 - 240V~; 50/60Hz; 1.3W), TCAxxxx (input ratings: 100 - 240V~; 50/60Hz; 1.5W), WAAxxxx (input ratings: 100 - 240V~; 50/60Hz; 3.5W). For the list of approved adapters, please see philips.com/support.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.

- Do not modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- Charge and use the product at a temperature between 0 °C and 40 °C (32 °F and 104 °F).
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- When you handle batteries, make sure that your hands, the product and the batteries are dry.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal object (e.g. coins, hairpins, rings). Do not wrap batteries in aluminium foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- This appliance contains a Lithium-ion rechargeable cell. **MUST BE RECYCLED OR DISPOSED OF PROPERLY.**
- The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Instructions for removal of built-in rechargeable batteries can be found on **www.philips.com/support**.

Storage and transportation conditions

Temperature (Transport and Storage) between -20°C to 60°C (-4°F and 140 °F)

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

FCC Compliance

Changes or modifications not expressly approved by Philips could void the user's FCC authorization to operate the toothbrush.

This device complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment (handle) off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna (of the radio or television).
- Increase the separation between the equipment (handle) and receiver.
- Connect the equipment (charger) into an outlet on a circuit different from that to which the receiver is connected.

ISED Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1 This device may not cause interference.

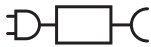
- This device must accept any interference, including interference that may cause undesired operation of the device.

SAVE THESE INSTRUCTIONS

The following symbols may appear on the product:

This symbol means read the instructions for use, before using the appliance.

XXXXXX

Detachable Power Supply Part, if supplied. The "xxxxxx" indicates the model number of the approved power supply unit to be used.

Intended use

Sonicare DiamondClean Smart power toothbrushes are intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. Sonicare DiamondClean Smart power toothbrushes are intended for personal use. Use by children under the age of 8 should be with adult supervision.

Your Philips Sonicare (Fig 1.)

- Handle
- Power on/off button
- Intensity indicator
- Mode/intensity button
- Mode indicator
- Brush head replacement reminder
- Battery indicator
- Brushing feedback light

Accessories:

- Brush head cap

- 10 Smart brush head(s)
- 11 BrushSync symbol
- 12 Charging glass
- 13 Charging base
- 14 Charging stand
- 15 Brush head holder
- 16 Charging travel case with integrated USB cord
- 17 Travel case
- 18 USB-A wall adapter

Note: Accessories, features, and benefits included may vary based on the model purchased.

Sonicare app - Getting started

The Sonicare app pairs with your toothbrush to provide you a connected experience. By connecting your toothbrush to your Sonicare app account, you will be able to:

- Track your brushing progress.
- Receive personalized tips and actionable recommendations to improve your oral health.
- Access the full range of benefits and receive ongoing upgrades.

The app is compatible with a wide range of iPhone and Android™ smartphones.

To start using the Sonicare app:

- 1 Download the Sonicare app to your phone.
- 2 Ensure your phone's Bluetooth is turned on.
- 3 Pick up your toothbrush to ensure it is active (lights on).
- 4 Open the app and follow the steps.
- 5 Pair your toothbrush with the app.
- 6 Create your account via the app. Complete firmware update, if prompted, to access the latest improvements and features.
- 7 Brush regularly with your toothbrush connected to the app. When you regularly sync the toothbrush with the app, you will receive updates to help improve your oral healthcare.
- 8 Sync regularly.

9 Keep your phone close to your toothbrush for optimal data transfer.

To sync manually: Pair/connect your toothbrush with the app every two weeks to benefit from the app features/capabilities.

To sync automatically: Allow location permission while setting up the app. By allowing location permission, your phone knows when it is in the connection range of your toothbrush and can refresh your brushing data in the app to provide the latest insights and recommendations.

Note: Make sure your phone's Bluetooth is switched on when using the app so that your toothbrush can transfer and update your brushing data in the app. If you have questions about why your brushing data is collected, be sure to review the Sonicare Privacy Statement, available throughout the app setup process.

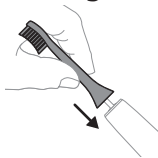
App features

Your Sonicare power toothbrush is equipped with smart sensors to provide you with feedback about your brushing habits, including:

- Daily brushing sessions
- Time spent during each brushing session
- Applied brushing pressure
- Personalized tips and actionable recommendations
- Brush head replacement reminder (based on actual use)

The Sonicare app collects the data provided by the smart sensors and tracks your brushing habits over time.

Getting started

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it will not go any further.

Note: It is normal to see a slight gap between the brush head and the handle.

Using your Philips Sonicare toothbrush

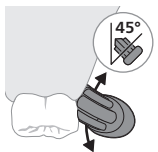
If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first-time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow the brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, use a lower intensity setting for the first 1-2 weeks then transition to a higher setting.

Brushing instructions


- 1 Wet the bristles and apply a small amount of toothpaste.


- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees). Pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

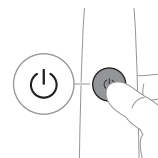Note: Keep the center of the brush in contact with the teeth at all times.

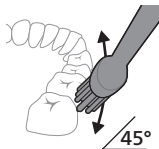
- 3 Press the power on/off button to switch on the Philips Sonicare.

- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

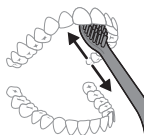
Note: If you apply too much pressure, you will be alerted by a change in vibration of the handle and a flashing purple brushing feedback light.

- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.
- 6 After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush switched on or off, as you prefer.

Note: Your Philips Sonicare toothbrush is safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised.

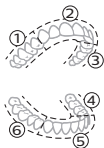BrushPacer

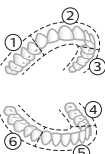
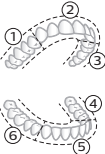
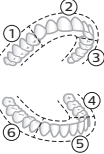
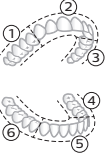
To make sure you brush evenly throughout the mouth, divide your mouth into 6 sections using the BrushPacer feature. The BrushPacer divides the brushing time into 6 equal segments and indicates when you should move to the next area. The segments are indicated with a brief pause in vibration. The toothbrush automatically stops at the end of the brushing session.

Note: Additional segments may be added when you use the White+ or Gum Health Modes.

Modes


Your power toothbrush is equipped with up to 5 different modes: Clean, White+, Gum Health, Deep Clean+ and Tongue Care (specific models only). Based on the smart brush head you use, the optimal mode and intensity will be selected automatically by the BrushSync

technology. To change modes, press the mode/intensity button while the handle is turned off.

	Clean and Deep Clean+	White+	Gum health	Tongue Care (specific models only)
Benefit	Plaque removal	Plaque and surface stains removal	Plaque removal and gum massage	Tongue cleaning
Recommended intensity level	High	High	High	Low
Total brushing time	Clean: 2 minutes Deep Clean+: 3 minutes*	2 minutes and 40 seconds	3 minutes and 20 seconds	20 seconds
How to brush	<p>Clean: Brush each segment for 20 seconds</p>  <p>Deep Clean+: Brush each segment for 30 seconds*</p> 	<p>Brush each segment for 20 seconds,</p>  <p>then brush segments 7 and 8 for 20 seconds each.</p> 	<p>Brush each segment for 20 seconds,</p>  <p>then brush segments 7, 8, 9 and 10 for 20 seconds each.</p> 	<p>1. Wet the micro bristles and brush your tongue with a back and forth motion. Clean as far back as comfortable. 2. Rinse your mouth and repeat the procedure two more times.</p>

*Note: When brushing with the Philips Sonicare app in 'Deep Clean+' mode, the brushing time will be 2 minutes when connected to the app and 3 minutes when not connected to the app.

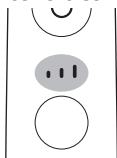Intensities

Your power toothbrush comes with 3 different intensity settings:

- low intensity (one bar)
- medium intensity (two bars)
- high intensity (three bars)

When attached to the handle, the smart brush head will automatically select the recommended intensity. If you wish to change to a different intensity, you can do so by pressing the mode/intensity button while brushing. The intensity cannot be changed as long as the handle is powered off or paused.

BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The  symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- BrushSync mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to

www.philips.com/toothbrush-heads for more information.

Features

- Pressure Sensor Feedback
- BrushPacer
- SmartTimer
- Brush Head Mode Pairing
- Easy Start
- Brush Head Replacement Reminder

Pressure Sensor Feedback

Your Sonicare toothbrush measures the pressure you apply while brushing, to protect your gums and teeth from damage.

If you apply excess pressure, the handle changes its vibration and the brushing feedback light flashes purple at the bottom of the handle until you reduce the pressure.

Note: Your product comes with the Pressure Sensor activated. To deactivate this feature (see 'Activating or deactivating features').

BrushPacer

The BrushPacer is an interval timer that helps you to brush evenly through the mouth by dividing your mouth into 6 sections. (see 'Brushing instructions')

SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power on/off button after you started the brushing cycle, the toothbrush pauses. After a pause of 30 seconds, the SmarTimer resets.

Brush Head Mode Pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates with the handle and automatically pairs it with the recommended mode and intensity.

Most brush heads with BrushSync technology will auto pair to Clean Mode. Brush heads with specific benefits will pair with the corresponding mode (when available), examples include:

- (G3) 2-in-1: Plaque Removal + Gum Care: Gum Health mode
- TongueCare+: Tongue Care mode

Note: If you choose to change mode/intensity from the recommended setting, the brush head will remember the your chosen setting for future brushing sessions.

EasyStart

This Philips Sonicare model comes with the EasyStart feature deactivated. To activate EasyStart (see 'Activating or deactivating features'). The EasyStart feature gently increases the power over the first 14 brushings to help you get used to brushing with the Philips Sonicare. The EasyStart feature works with all modes.

Disclaimer

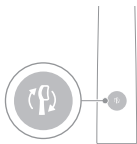
When the Philips Sonicare toothbrush is used in clinical studies, it must be used in clean mode at high intensity with the handle fully charged and the EasyStart feature deactivated.

Brush Head Replacement Reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head.

Note: This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When you attach a new smart brush head for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and begins tracking brush head wear.
- 2 Over time, based on the pressure you apply and the total time the brush head has been used, the handle tracks brush head wear in order to determine the optimal time to replace your brush head. This feature guarantees the best cleaning and care of your teeth.
- 3 When the brush head replacement reminder indicator lights up amber, you should replace your brush head.

All smart brush heads come with the brush head replacement reminder feature activated. To deactivate the brush head replacement reminder see 'Activating or deactivating features'.

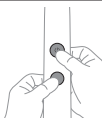
Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush Head Replacement Reminder
- Pressure Sensor Feedback

Note: You can only activate or deactivate the Brush Head Replacement Reminder when a smart brush head is on the handle.


Step 1: Place the handle on the charging stand or in the charging glass.

Step 2: Press and hold the mode/intensity button.


EasyStart	Pressure Sensor Feedback	Brush Head Replacement Reminder
Press the power on/off button 1 time.	Press the power on/off button 2 times.	Press the power on/off button 3 times.
↓	↓	↓

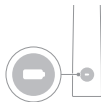Step 3: Let go of the mode/intensity button.

If you see the battery indicator flash white 2 times and hear 2 beeps then the feature has been activated.

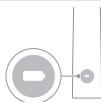OR

If you see the battery indicator flash amber 1 time and hear 1 beep, then the feature has been deactivated.

Charging and battery status

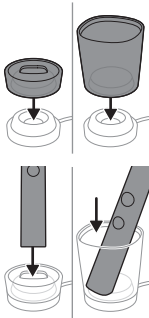
This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day) on a fully charged battery. The battery indicator shows the battery status when you complete the 2-minute brushing session, when you pause the toothbrush, or while it is charging.

Note: It can take up to 24 hours to fully charge the battery before first use, but you can use the Philips Sonicare toothbrush before it is fully charged.

Charging

Charging on charging base

- 1 Plug the cord of the charging base into an electrical outlet.
- 2 Place the charging stand or the charging glass on the charging base.
- 3 Place the handle on the charging stand or in the charging glass.
 - The battery indicator and light ring illuminate white and the handle beeps twice. This indicates the toothbrush has started charging.


- 4 The charger will charge your toothbrush until it is fully charged.
 - The handle is charging when the battery indicator is flashing white. When the handle is fully charged, the battery indicator lights up white for 30 seconds and then switches off.


Charging with the travel case

Note: Accessories, features, and benefits included may vary based on the model purchased.

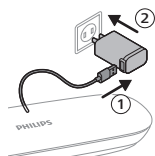- 1 Take the USB cable from the bottom of the charging travel case.
- 2 Connect the charging travel case via the USB cable to 5V DC input power either from:

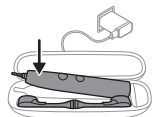- The USB adapter provided
- a UL-listed or certified to ANSI/UL 60950-1 personal computer's standard USB connector or,
- a UL-listed or certified to ANSI/UL 60950-1 LPS power supply with USB output port or,
- a UL-listed or certified 1310 Class 2 power supply.

- 3 Place the handle in the travel case. The travel case charges your toothbrush until your toothbrush is fully charged.

Battery status

Battery status (handle on plugged in charger/charging travel case)

Handle on plugged in charger/charging travel case:

When the handle is placed on the charger or in the charging travel case, the battery indicator shows the battery level.

Battery status	Battery indicator
Charging	Flashes white
Full	Solid white for 30 seconds and then switches off

Battery status (when handle is not placed on charger/charging travel case)

Handle is not placed on charger/charging travel case:

When the toothbrush is awake, the battery indicator at the bottom of the handle shows the battery level.

Battery status	Battery indicator	Sound
Full	Solid white	-
Partially full	Flashes white	-
Low	Flashes amber	The handle emits 3 beeps after brushing

Empty

Flashes amber

The handle emits 2 sets of 5 beeps after brushing

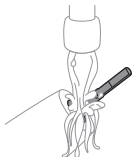
Note: To save energy, the battery indicator will switch off when not in use.

Note: When the battery is completely empty, the toothbrush switches off. Place the toothbrush on the charger to charge it.

Note: To keep the battery fully charged at all times, you may keep your toothbrush on the charger when not in use.

Cleaning

Toothbrush handle


- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste at least once a week, but preferably after each use.

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Note: Do not tap the handle on the sink to remove excess water.

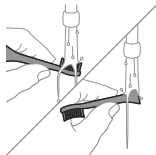Brush head

- 1 Rinse the brush head and bristles after each use and allow to air dry.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week and allow to air dry. Rinse the travel cap as often as needed and air dry. Use only on a dry brush head.

Brush head holder

Rinse with warm water and use a damp cloth to clean the cover and the base and allow to air dry.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

Travel case

Rinse with warm water and use a damp cloth to clean inside of travel case and allow to air dry.

Charging travel case

- 1 Unplug the charging travel case before cleaning.
- 2 Use damp cloth to wipe the surface of the travel case and allow to air dry.

Storage

If you are not going to use the product for an extended period of time, unplug it from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Replacement

Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX99XX).

Frequently asked questions

Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit www.philips.com/support and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Possible Cause	Solution
I cannot fully attach the brush head. There is a gap between the brush head and handle.		This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations.
My Philips Sonicare toothbrush vibration is less powerful than before.	The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged.	Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2 mm). When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.
I charged my toothbrush, but it only works for a short time.	The base of the toothbrush is not placed flat on the charger when attempting to charge.	Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.

Warranty and support

Philips warrants its products for a specific period of time after the first date of purchase. If the product is purchased in the USA, Philips North America LLC (USA) or, if the product is purchased in Canada, Philips Electronics Ltd., warrants this new product against defects in materials or workmanship for a period of two years from the date of purchase (as evidenced by a valid proof of purchase), and agrees to replace any defective product

without charge. Philips will refund the purchase price if it is unable to provide replacement.

Parts that are subject to normal wear and tear, accessories, replacement batteries and consumables such as toothbrush heads, are not covered by the terms of this warranty.

IMPORTANT: This warranty does not cover damage resulting from accident, misuse or abuse, lack of reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified voltage.*

NO RESPONSIBILITY IS ASSUMED FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES.

For more information and in order to obtain warranty service, simply go to **www.philips.com/support**.

NOTE: This warranty gives you specific legal rights, and you may also have other rights which vary from state to state, province to province or country to country. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above exclusion and limitations may not apply to you. This warranty is transferable.

*Read enclosed instructions carefully.

THE WARRANTY DOES NOT COVER THE FOLLOWING:

- Indirect or consequential damages, such as loss of time, loss of data, or loss of income. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.
- Compensation for activities done by yourself, such as regular maintenance, installing firmware updates, or saving data.

THE WARRANTY DOES NOT APPLY:

- If the proof of purchase, model number, serial number, or production date code on the product has been made illegible.
- If the product was purchased as faulty stock.
- If the defect is the result of abuse or excessive use outside the intended purpose.

- If the defect is caused by environmental conditions, connecting peripherals, additional equipment, or accessories (including software) that are not in accordance with the user manual.
- If the unit has been damaged, including but not limited to damage by animals, lightning, abnormal voltage, fire, natural disaster, transportation, or water (unless the user manual expressly states that the product may be rinsed).
- For non-electrical parts, normal wear and tear, accessories, replacement batteries or replacement of consumable parts such as toothbrush heads.
- If the product does not function properly because it was not originally designed for use in the country where the product is used, which may occur if you import the product.
- If the product does not function properly due to problems with connection with service providers, such as for the internet.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- In mouth attachments including brush heads and nozzles.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

IMPLIED WARRANTIES

ANY IMPLIED WARRANTIES, INCLUDING IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THE EXPRESS WARRANTIES SET FORTH ABOVE. IN SOME STATES LIMITATIONS ON DURATION OF IMPLIED WARRANTIES DO NOT APPLY.

LIMITATION OF REMEDIES

IN NO EVENT SHALL PHILIPS OR ANY OF ITS AFFILIATED OR SUBSIDIARY COMPANIES BE LIABLE FOR ANY

SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, TORT, OR ANY OTHER LEGAL THEORY. SUCH DAMAGES INCLUDE, WITHOUT LIMITATION, LOSS OF SAVINGS OR REVENUE; LOSS OF PROFIT; LOSS OF USE; THE CLAIMS OF THIRD PARTIES INCLUDING, WITHOUT LIMITATION, DENTISTS AND DENTAL HYGIENISTS; AND COST OF ANY SUBSTITUTE EQUIPMENT OR SERVICES. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES.

Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

Removing the rechargeable battery

The built-in rechargeable battery **must only be removed by a qualified professional** when the appliance is discarded. Before removing the battery, make sure that the battery is completely empty.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the product and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in

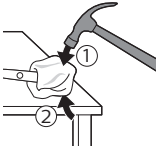
aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.


- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

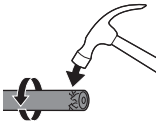
Note: You may have to hit on the end several times to break the internal snap connections.

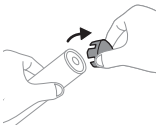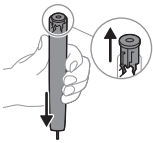
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.


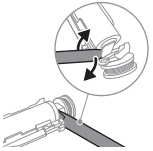
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

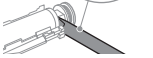- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

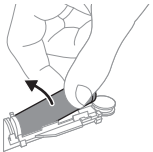- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

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©2026 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.
de Philips Oral Healthcare, LLC et (ou) KPNV.
de Philips Oral Healthcare, LLC o KPNV.