PHILIPS

ReAura

LED MASK LTM7100

|          |    |            |     |
|----------|----|------------|-----|
| English  | 3  | Italiano   | 59  |
| Deutsch  | 16 | Nederlands | 72  |
| Español  | 31 | Polski     | 86  |
| Français | 45 | Svenska    | 100 |

## WELCOME TO THE WORLD OF LIGHT AND BEAUTY!

Light has always been a symbol of renewal — and now, it's literally illuminating the path to better skin. Once confined to high-tech dermatology clinics, LED (Light Emitting Diode) light sessions have quietly become one of the most effective, non-invasive tools in modern beauty routines. You've likely seen the futuristic masks all over social media — but behind those glowing selfies lies real science, decades in the making.

## DISCOVERING THE PHILIPS REAURA LED MASK

We are proud to call Philips a born innovator. For over 130 years, we have been improving people's lives with a steady flow of ground-breaking innovations.

Now we are proud to bring this spirit of innovation to the 3-in-1 ReAura face mask backed by decades of Philips LED innovation. Just use it for 10 minutes every other day and enjoy healthy-looking glow.

## GENERAL DESCRIPTION (Fig. 1)

- 1 Top part of LED mask for eye area
- 2 Main mask strap with hook and loop fastening
- 3 Bottom part of LED mask for chin or neck
- 4 Neck straps with hook and loop fastening
- 5 USB-C cable

Note: Make sure you use a suitable USB-C supply unit to charge the controller (5 V, max. 2 A).

- 6 Controller
- 7 On/off and mode selection button
- 8 Battery status light

- 9 Session progress lights
- 10 Blemish care mode light
- 11 Age care mode light
- 12 Daily glow mode light
- 13 Pouch

For all figures referred to in this user manual, see the foldout page and the inside of the cover at the front of this booklet.


## HOW IT WORKS

Think of a LED light session as a gentle conversation with your skin. Instead of using heat or chemicals, the light delivers energy in the form of different colors — or wavelengths — that penetrate the skin at varying depths. Each color has a unique role:

- Red + infrared light: to visibly reduce the appearance of fine lines and wrinkles.
- Blue + red light: to visibly reduce the appearance of blemishes and redness.
- Blue + red + infrared light: a combination of different light colors to maintain your healthy-looking glow.

Usage of the LED mask is entirely non-invasive, painless, and suitable for most skin types. It doesn't burn, peel, or require downtime — which explains why it has become a go-to for those seeking results without disruption.

It has been designed to minimize the exposure of your eyes to direct LED light, allowing you to comfortably multi-task during your skincare routine without straining your eyes.

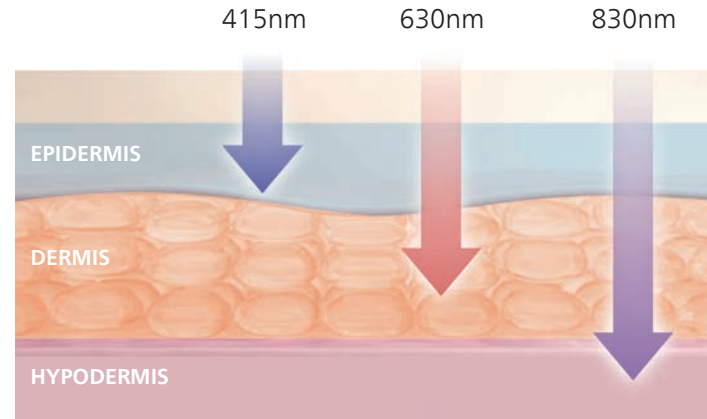SKIN LAYERS

## WHAT TO EXPECT

A session with the LED mask typically lasts 10 minutes and feels like a spa ritual — no discomfort, just quiet meditation. With consistent use, many people experience healthy-looking glow. This clinically proven technology enables visible reduction in appearance of fine lines, wrinkles, redness and blemishes. The key word here is consistency: much like exercise, results build gradually with regular sessions.

## PRECAUTIONS, CONTRAINDICATIONS AND POTENTIAL SIDE EFFECTS

**This device is intended for use by people in good health on intact skin. Before you use the Philips LED mask for the first time, carefully read all safety information (contraindications, precautions, side effects and warnings) in this user manual and in the separate safety information leaflet. Do not look at the LED light directly. Before you use the device for the first time, you must test your skin for light sensitivity.**

These recommendations are derived from published research regarding precautions and contraindications as of the printing date.

For up-to-date information, please visit us online at [www.philips.com](http://www.philips.com)

**STOP using the device immediately at the first sign of discomfort.**

## PRECAUTIONS

Due care is required in these circumstances and device use may need to be modified. Consult a medical professional if you currently have or suspect you may have any of the following conditions or if you have any questions.

- Recent injury, surgery, or facial treatment including neurotoxin, dermal filler, microneedling, laser and chemical peel until the skin has fully healed.
- Current herpes simplex virus breakout
- Broken skin.
- Retinol application before use of red LED light.

## CONTRAINDICATIONS

Do not use if any of the following contraindications applies to you:

- You have vitiligo, melasma or hyperpigmentation (especially if exacerbated by mild warmth).

- You have skin rash (e.g. psoriasis or eczema), open wounds, blisters, local tissue inflammation, infection, bruises, or tumors on the face, hairline, ears or neck.
- You are pregnant, nursing, maybe pregnant, or if you are under the age of 18.
- You have any abnormal sensations (e.g. numbness).
- You have cancer, tumors or suspicious lesions.
- You have epilepsy.
- You have cardiopathy (heart disease), a pacemaker or other implanted electronic device.
- You have implants in the face.
- You have light sensitive skin disorder (e.g. lupus, porphyria), photosensitive eczema or albinism.
- You have a photosensitive disorder of the eyes (e.g. glaucoma or genetic/congenital retinal conditions), have had recent laser eye surgery, or are extremely sensitive to light.

- You have undergone cataract surgery, refractive lens exchange, or any other surgical procedure involving removal of the natural lens of the eye (with or without an intraocular lens implant).
- You are taking:
  - medications that cause light sensitivity of the eyes and/or skin
  - medications for severe inflammatory acne
  - retinol, other sun-sensitive medications or products, or benzoyl peroxide
- You have light-induced headaches.
- You have tattoos or permanent make-up in the area that you wish to use the device on.
- You have an allergy to silicone or Spandex (Elastane).
- You are using perfume, anti-bacterial soaps, deodorants or petroleum-based products in the area you wish to use the device on. These products can cause chemical burns under LED light.
- You have skin that is burned or has been recently burned from sun or tanning equipment.

- You have varicose veins in the area that you wish to use the device on.
- You have dark brown or black spots such as birthmarks, moles or warts in the area you wish to use the device on.

POTENTIAL SIDE EFFECTS

Some users may experience the following reactions, which are common when using LED masks:

- Skin sensations: temporary redness, swelling, local warmth, itchiness, dryness, flakiness or minor inflammation.
- Vision and headaches: temporary light-induced headaches, eye dryness, redness, blurred vision or general irritation.
- Prolonged discoloration: changes in skin tone or pigment.

Most of these effects are usually short-lived and resolve within hours (generally less than 24 hours). If symptoms persist, discontinue use and consult a healthcare professional.

OPERATING MODES

Only use the mask for one light session every other day.

The LED mask is equipped with 3 modes:

- Age care mode (Fig. 2): The age care mode is a combination of red and near-infrared light (630 nm and 830 nm) clinically proven to enable visible reduction in appearance of fine lines and wrinkles. The mode duration is 10 minutes. The recommended usage is every other day, for a minimum of 8 weeks.
- Blemish care mode (Fig. 3): The blemish care mode is a combination of blue and red light (415 nm and 630 nm) clinically proven to enable visible reduction in appearance of blemishes and redness. The mode duration is 10 minutes. The recommended usage is every other day, for a minimum of 4 weeks.

- Daily glow mode (Fig. 4): The daily glow mode is a combination of blue, red, and near-infrared light (415 nm, 630 nm and 830 nm) in shorter sessions to maintain the glow. Mode duration is 5 minutes. You can use the daily glow mode every other day for as long as you like.

| Mode                                                                                                                         | Session    | Frequency       | Period              |
|------------------------------------------------------------------------------------------------------------------------------|------------|-----------------|---------------------|
| Age care mode (Fig. 2)                                                                                                       | 10 minutes | Every other day | At least 8 weeks    |
| Blemish care mode (Fig. 3)                                                                                                   | 10 minutes | Every other day | At least 4 weeks    |
| Once you have achieved the desired effects with the age care mode and/or blemish care mode, you can use the daily glow mode. |            |                 |                     |
| Daily glow mode (Fig. 4)                                                                                                     | 5 minutes  | Every other day | As long as you like |

## USAGE AREAS

This two-piece mask can be used in 3 different ways to give you either a full-face session, a targeted session for the eye area or a targeted session for the neck.

- Full face (Fig. 5): Use both parts of the mask attached to each other. Align the bottom part of the mask with the top part. The bottom part will attach itself to the top part with magnets. Use the mask in combination with the main mask strap.
- Eye area (Fig. 6): Only use the top part of the mask. Use this part in combination with the main mask strap.
- Neck (Fig. 7): Only use the bottom part of the mask. Use this part in combination with the neck straps.

## GETTING STARTED

**Before you use the Philips LED mask for the first time, carefully read the separate leaflet with important safety information as well as chapter 'Precautions and contraindications' in this booklet.**

### SKIN TEST

Before you use the device on your face or neck for the first time, you must test your skin for light sensitivity in the following way:

- 1 Make sure that your skin is clean and dry.
- 2 Place the mask in contact with the skin, ideally on your arm, the back of your hand or inner thigh, with the LED light shining towards your skin (Fig. 8).
- 3 Switch on the device by pressing the on/off button for at least 2 seconds (Fig. 9).

- 4 Let the device complete a full session at any mode. The LED mask switches off automatically after the session.

**It is normal that the mask feels warm on your skin. If your skin feels too hot or uncomfortable at any time during the session, stop using the device immediately.**

- 5 After this full usage session, wait at least 6 hours. If you experience any side effects (such as redness or itching), then your skin is not suitable for sessions with the LED mask. In this case, please contact the Philips Consumer Care Center in your country ([www.philips.com/support](http://www.philips.com/support)).

### CHARGING

Charge the controller before first use and when the battery is empty.

- 1 Insert one end of the USB-C cable into the controller (Fig. 10).
- 2 Insert the other end of the USB-C cable in a universal USB-C supply unit (not provided) (Fig. 11).

- 3 Put the USB-C supply unit in a compatible wall socket (100-240V) (Fig. 12).
- 4 During charging, the battery status light flashes white (Fig. 13) to indicate that the controller is charging.
- 5 It takes up to 6 hours to fully charge the controller. The battery is fully charged when the battery status light lights up green continuously (Fig. 14).
- 6 After charging, disconnect the controller from the USB-C cable and remove the supply unit from the wall socket (Fig. 15).

### BATTERY STATUS INDICATIONS

The controller shows the battery status during mode selection and after the device has switched off. During the session, the battery status will only be shown if the battery is running low.

- Solid orange when you switch on the device (Fig. 16): the battery is almost empty and the device will switch off within 30 seconds.

- Flashing orange during use (Fig. 17): the battery is running low. The mask may switch off before completion of the session.
- Solid green (Fig. 14): the battery contains sufficient charge to operate the device.

## HOW TO USE THE PHILIPS REAURA LED MASK

**Warning: Keep the LED mask, the controller, the USB cable and the USB supply unit (not provided) dry (Fig. 18).**

**Only use the mask for one light session every other day.**

**If you notice any prolonged or severe side effects, stop using the mask.**

- 1 Unpack your device. The device may come partly charged, but we recommend charging it fully before first use.

- 2 Perform the skin test before use the LED mask on your face or neck for the first time.
- 3 Begin with a clean dry face.

**Warning: Make sure that the skin you are using the LED mask on is free from perfume, anti-bacterial soaps, deodorants, petroleum-based products and products containing retinol.**

- 4 Ensure the most comfortable fit of the mask by taking the following things into account:
  - The main strap can be threaded through 2 different positions on the mask to cater for different face shapes and sizes. The high slot on the top part of the mask is recommended for wider faces (Fig. 19). The lower slot on the top part of the mask is recommended for smaller faces (Fig. 20).
  - Move the strap lower or higher at the back of your head (Fig. 21) to find the angle that makes the mask sit most comfortably on your face.

- Adjust the strap on both sides of the mask (Fig. 22), not just on one side.
- If the mask feels uncomfortable, remove the strap from the mask completely and attach it to the alternative slots.
- The mask should feel snug but not too tight on your face.
- Wearing your hair in a low pony tail makes it easier to position and adjust the strap properly.

Note: If you only want to use the bottom part of the mask, use the two small straps (Fig. 23).

- 5 Connect the charged controller to the mask with the USB-C cable provided.
  - When you use both parts of the mask or only the bottom part, insert the USB-C plug into the socket on the bottom part of the mask (Fig. 24).
  - When you only use the top part of the mask, insert the USB-C plug into the socket on top part of the mask (Fig. 25).

- 6 Place the mask on your face.

Note: This device has built-in eyeshields. Normal use of the eyeshields may result in redness around your eye area. Some redness is normal and should subside within 5-10 minutes after using the device.

- 7 Switch on the device by pressing the on/off button for at least 2 seconds (Fig. 9). The battery status light goes on to show the battery status.

**Warning: Only switch on the mask when it is placed properly and securely on your face. Do not stare directly into the LEDs when the mask is switched on. Only use the mask if the LEDs are not damaged. Re-adjust the position of the mask if the LED light enters your eyes during use.**

**8** Press the on/off button briefly one or more times to select the required mode (Fig. 26):

- Age care mode (Fig. 2)
- Blemish care mode (Fig. 3)
- Daily glow mode (Fig. 4)

**9** If you do not press the on/off button again, the last mode you selected will start automatically after 2 seconds.

**10** The light of the selected mode flashes 5 times (Fig. 27).

Note: As long as the light flashes, you can still change the selected mode by pressing the on/off button.

**11** When the light of the selected mode stops flashing, the LEDs in the mask go on and the session starts.

Note: If the battery is running low during the session, the battery light starts flashing orange. The mask may switch off before completion of the session.

**12** The lights on the controller show the progress of the session (Fig. 28). When the session starts, all 4 progress lights are on and the top light flashes slowly. After ¼ of the session, the top light goes out. The light below it starts flashing and goes out after ¼ of the session. And so on until the bottom light goes out to mark the end of the session.

Note: If you want to pause the session, press the on/off button once. The mode lights start pulsating in the color of the selected mode (Fig. 29). To continue the session, press the on/off button once more.

**13** The LED mask switches off automatically after the session. The battery status light goes on to show the battery status.

**14** When the session has stopped, take off the mask and continue with your daily skincare routine.

**15** Clean the mask after each use.

**16** Store the device in the pouch supplied in a cool, dry and safe place. For long-term storage, store the device with a fully charged battery.

## CLEANING THE PHILIPS REAURA LED MASK

**Warning: Keep the LED mask, the controller, the USB cable and the USB supply unit (not provided) dry (Fig. 18). Do not use chemical cleaners or alcohol-based cleaners on any parts of the mask and controller. Never use compressed air, scouring pads, abrasive cleaning agents or aggressive liquids to clean the device.**

- 1** Make sure that the mask is switched off and unplugged from any power source.
- 2** After every use, clean the inside and outside of the device with a damp or dry microfibre cloth or with a non-alcohol based anti-bacterial wipe:
  - Wipe the inside of the mask with a damp or dry cloth (Fig. 30).
  - Wipe the outside of the mask with a damp or dry cloth (Fig. 31).
- 3** Let the device air-dry completely before you use or charge it again.

## SUPPORT

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support).

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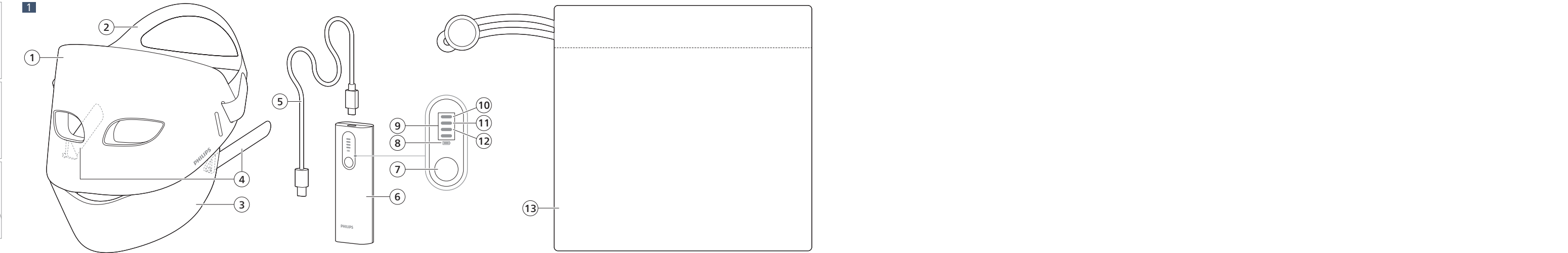
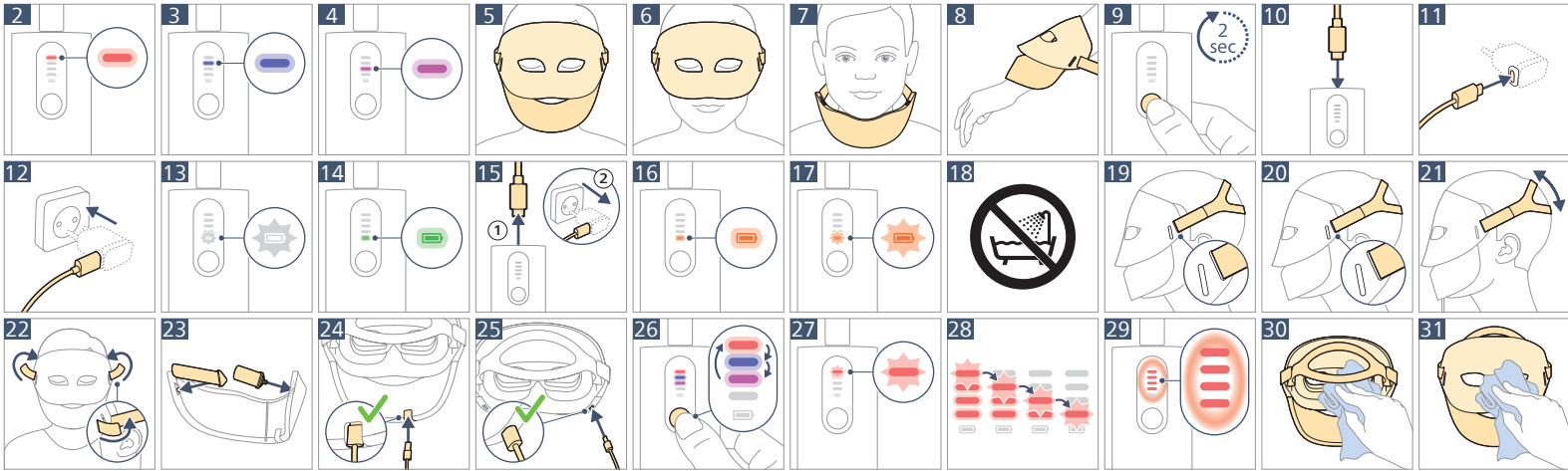
**Important safety information**  
Only use the product for its intended household purpose. Read this important information carefully (Fig. 1) before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries.

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- You have vitiligo, melasma or hyperpigmentation (especially if exacerbated by mild warmth).
- You have skin rash (e.g. psoriasis or eczema), open wounds, blisters, local tissue inflammation, infection, bruises, or tumors on the face, hairline, ears or neck.
- You are pregnant, nursing, maybe pregnant, or if you are under the age of 18.
- You have any abnormal sensations (e.g. numbness).
- You have cancer, tumors or suspicious lesions.
- You have epilepsy.
- You have cardiopathy (heart disease), a pacemaker or other implanted electronic device.
- You have implants in the face.
- You have light sensitive skin disorder (e.g. lupus, porphyria), photosensitive eczema or albinism.
- You have a photosensitive disorder of the eyes (e.g. glaucoma or genetic/ congenital retinal conditions), have had recent laser eye surgery, or are extremely sensitive to light.
- You have undergone cataract surgery, refractive lens exchange, or any other surgical procedure involving removal of the natural lens of the eye (with or without an intraocular lens implant).
- You are taking:
  - medications that cause light sensitivity of the eyes and/or skin
  - medications for severe inflammatory acne
  - retinol, other sun-sensitive medications or products, or benzoyl peroxide
- You have light-induced headaches.
- You have tattoos or permanent make-up in the area that you wish to use the device on.
- You have an allergy to silicone or Spandex (Elastane).
- You are using perfume, anti-bacterial soaps, deodorants or petroleum-based products in the area you wish to use the device on. These products can cause chemical burns under LED light.
- You have skin that is burned or has been recently burned from sun or tanning equipment.
- You have varicose veins in the area that you wish to use the device on.
- You have dark brown or black spots such as birthmarks, moles or warts in the area you wish to use the device on.

#### WARNINGS

- Before using the device, you must carry out a skin test. Please follow the instructions in the the user manual.
- To charge the product, only use a certified safety extra low voltage (SELV) supply unit with output rating 5 V, max. 2 A. If you need support finding the correct USB supply unit, contact the Consumer Care Center in your country. Using a non-certified supply unit may cause hazards or serious injuries. Philips is not responsible for damages that may occur due to the use of third-party chargers.
- Keep the LED mask, the controller, the USB cable and the USB supply unit (not provided) dry (Fig. 2).
- The device is not waterproof. Only charge and use the device indoors in a well-ventilated, dry location. Do not charge or use the device outdoors, in a bathroom, or within 3 meters of a bathtub, shower or pool or in another moist environment. Do not reach for an appliance that has fallen into or come into contact with water. Unplug immediately. Clean the device according to the instructions in the user manual.
- This device is intended for use by people in good health on intact skin.
- This device is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- This device and its pouch are not intended for use by children under 18 years old. Children should be supervised to ensure they do not play with the device and its pouch.
- Because of hygiene, only one person should use the appliance.
- Never use or charge the device if the USB supply unit (not provided), USB cable, mask or controller is damaged in any way, if the device is not working properly, if it has been dropped or dropped into water, if the housing appears to be degraded or if the transparent silicone layer of the mask is damaged. Inspect all parts for damage before use. Philips is not responsible for damages that may occur due to the use of third-party chargers. Do not attempt to repair or disassemble the device, as this may result in an electric shock or fire.
- Do not use the supply unit in or near wall sockets that contain an electric air freshener to prevent irreparable damage to the supply unit.
- Do not use the device for more than one light session per day.
- Only switch on the mask when it is placed properly and securely on your face. Do not stare directly into the LEDs when the mask is switched on. Only use the mask if the LEDs are not damaged. Re-adjust the position of the mask if the LED light enters your eyes during use.
- Stop using the device immediately if you experience pain or discomfort beyond what is expected from an ice pack or a hot shower. Discontinue use of the device if it overheats or becomes uncomfortably hot.
- If you notice any prolonged or severe side effects, stop using the mask.
- Do not modify the supply unit.
- Do not use the device or the USB supply unit in the presence of explosive atmospheres (gaseous fumes, dust or flammable materials). Sparks may be generated, possibly causing a fire.
- The device should be kept away from radio frequency (RF) equipment and other devices that emit electromagnetic waves to avoid interference.
- Unplug the USB cable and the USB supply unit when you are not using or charging the device. Store away when not in use. Do not stretch the USB cable or place the cable under strain. Do not handle the cable with wet hands.
- Store the device in a cool, dry and safe place. For long-term storage, store the device with a fully charged battery.
- Charge and use the product at a temperature between 0 °C and 30 °C.
- Never use compressed air, scouring pads, abrasive cleaning agents or aggressive liquids to clean the appliance.
- Do not use chemical cleaners or alcohol-based cleaners on any parts of the mask and controller.
- Only use original Philips accessories or consumables.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.

- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries. Do not place a battery near anything that may cause a short circuit, such as coins, keys, or screws in your pocket. A short circuited battery may cause fire and personal injury.
- If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

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#### Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

#### Warranty and support

For all product support, power consumption details and the EU Declaration of Conformity, please visit [www.philips.com/support](http://www.philips.com/support). Philips offers a 2-year warranty. Warranty exceptions apply. This warranty does not affect your rights as a consumer under law. For more information, visit [www.philips.com/support](http://www.philips.com/support) or call us.

#### Replacement of battery unit

The battery unit is replaceable. For information about how to replace it, visit [www.philips.com/support](http://www.philips.com/support).

#### Recycling

This symbol means that electrical products and batteries shall not be disposed of with normal household waste (Fig. 3). Follow your country's rules for the separate collection of electrical products and batteries.

#### Removal of built-in rechargeable battery

The built-in rechargeable battery must only be removed by a qualified professional when the appliance is discarded. Before removing the battery, make sure that the appliance is disconnected from the wall socket and that the battery is completely empty.

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**

**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**

- 1 Remove the screw in the back of the controller with a screwdriver.
- 2 Slide the battery compartment cover off the controller housing.
- 3 Remove the battery unit from the controller.

#### Explanation of symbols

|  |                                                                                                                                                                                                     |
|--|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | This symbol identifies the manufacturer of the product: Philips Consumer Lifestyle B.V., Tussendiepen 4, 9206 AD Drachten, The Netherlands.                                                         |
|  | This symbol means that you have to read the user manual carefully before use and to save the user manual for future reference (Fig. 1).                                                             |
|  | This symbol means that the manufacturer has taken all necessary measures to ensure that the product complies with the applicable safety legislation for distribution within the European Community. |
|  | This symbol means that the user should read this leaflet before use because it contains important safety information such as contraindications, warnings and precautions.                           |
|  | This symbol means that the device is in accordance with UK directives on health, safety and environmental protection standards.                                                                     |
|  | This symbol indicates the unique device identifier.                                                                                                                                                 |
|  | This symbol indicates the serial number of the device.                                                                                                                                              |
|  | Indicates manufacturer's catalog number of the product.                                                                                                                                             |
|  | This symbol means that you have to keep the LED mask, the controller, the USB cable and the USB supply unit (not provided) dry.                                                                     |
|  | This symbol means that the device has been manufactured in China.                                                                                                                                   |
|  | This symbol means that the device is protected against solid objects over 12.5 mm in diameter and vertical water droplets when the device is tilted 15°.                                            |
|  | This symbol indicates the manufacturing date of the device.                                                                                                                                         |
|  | This symbol means the separate collection of waste electrical and electronic equipment.                                                                                                             |

#### Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit [www.philips.com/support](http://www.philips.com/support) for a list of frequently asked questions or contact the Consumer Care Center in your country.

| Problem                                              | Possible cause                      | Solution                                             |
|------------------------------------------------------|-------------------------------------|------------------------------------------------------|
| The charging light on the controller does not go on. | The battery level is low.           | Charge the controller with the USB-C cable provided. |
|                                                      | The device is faulty.               | Contact the Consumer Care Center in your country.    |
| The battery does not charge.                         | The USB-C supply unit is defective. | Replace the USB-C supply unit.                       |
|                                                      | The USB-C cable is defective.       | Replace the USB-C cable.                             |
|                                                      | The battery is defective.           | Contact the Consumer Care Center in your country.    |

| Problem                                                                   | Possible cause                                               | Solution                                                                 |
|---------------------------------------------------------------------------|--------------------------------------------------------------|--------------------------------------------------------------------------|
| The LED lights of the mask do not switch off after the session has ended. | The device is faulty.                                        | Contact the Consumer Care Center in your country.                        |
| The mask, controller, packaging or USB-C cable is damaged in any way.     |                                                              | Do not use the device. Contact the Consumer Care Center in your country. |
| Only part of the mask lights up.                                          | The connection between the components of the mask is faulty. | Contact the Consumer Care Center in your country.                        |

#### Technical specifications

|                                   |                                                                                          |
|-----------------------------------|------------------------------------------------------------------------------------------|
| Model number                      | LTM7100                                                                                  |
| Operating temperature             | 0 °C - 30 °C                                                                             |
| Operating humidity                | 0 % - 85 % relative humidity                                                             |
| Storage and transport temperature | -10 °C - 50 °C                                                                           |
| Maximum storage humidity          | 0 % - 85 % relative humidity                                                             |
| Electrical output                 | 3.6 V                                                                                    |
| LEDs                              | Blue: 415 nm ± 10 nm<br>Red: 630 nm ± 10 nm<br>Near-infrared: 830 nm ± 10 nm             |
| Mode duration                     | Age care mode: 10 minutes<br>Blemish care mode: 10 minutes<br>Daily glow mode: 5 minutes |

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