
PHILIPS
sonicare

For Kids 3+

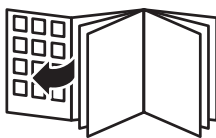— |

| —

— |

| —

English	6	Русский	105
Azərbaycanca	9	Тоҷикӣ	108
Čeština	12	Українська	111
Dansk	15	Қазақша	114
Deutsch	18	Հայերեն	117
Eesti	21	ქართული ენა	120
Español	24	עברית	123
Español	27		
Français	30		
Français (Canada)	33		
Hrvatski	36		
Italiano	39		
Latviešu	42		
Lietuviškai	45		
Magyar	48		
Nederlands	51		
Norsk	54		
O'zbek	57		
Polski	60		
Português	63		
Română	66		
Shqip	69		
Slovenščina	72		
Slovensky	75		
Srpski	78		
Suomi	81		
Svenska	84		
Türkçe	87		
Türkmençe	90		
Ελληνικά	93		
Български	96		
Кыргызча	99		
Македонски	102		

Removing the rechargeable battery

The built-in rechargeable battery can be removed for disposal. If you are not comfortable disassembling electronics or if otherwise instructed in the product manual, please seek professional assistance at your local electronics recycling or repair shop.

Battery removal should only be performed by adults. Keep batteries and tools out of reach of children to prevent hazards.

Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.

Be sure to protect your eyes, hands, fingers, and the surface on which you work.

When you handle batteries, make sure that your hands, the toothbrush, and the batteries are dry.

To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.

Do not reinstall battery once it has been removed for recycling. Consumer installation of battery may increase risk of fire or other potential safety hazards.


If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.

Removing the rechargeable battery


Before removing the battery, make sure that the battery is completely empty.

To remove the rechargeable battery, you will need these tools: A T8 or T10 torx screwdriver, a small cloth, a 5.5mm Allen wrench, tweezers, a plastic pry tool (or guitar pick), and a hard surface such as a table.


- 1 Deplete the rechargeable battery of any charge. To do this, remove the toothbrush from the charger, turn on the toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare toothbrush.
- 2 Remove the brush head from the toothbrush.
- 3 Using a torx screwdriver, unscrew and remove (counterclockwise) the screw from the base of the toothbrush handle.

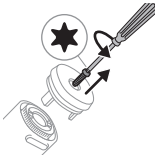
- 4 Wrap a small, thin cloth around the end of your Allen wrench, insert it into the end cap. Then, carefully pry the end cap out of the housing tube.

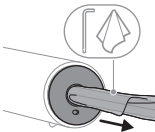- 5 Use tweezers or a plastic pry tool to pry the two retention tabs away from the inside walls of the handle.

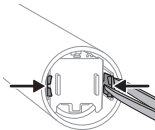
- 6 Once the tabs are away from the wall, release the internal assembly from the housing tube by inverting the handle and firmly pressing the metal shaft onto a solid surface (do not slam the handle down with force). If the internal components do not easily release from the housing tube, double check that you have completed step 5 above.
- 7 Gently pull the internal assembly to remove it from the base of the housing tube.

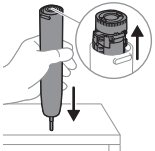Caution: To avoid injury or damage to the product, be aware of the sharp edges and exposed electronic components.

8 English

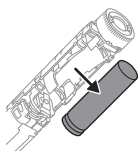- 8 Use a plastic pry tool to push or pry the battery out of the battery tray.

Caution: Do not use screwdriver or metal object to pry the battery out of the battery tray.

- 9 For battery disposal, cover the battery contacts of the used battery with tape to prevent any electrical short from residual battery charge. The used rechargeable battery can now be recycled through your community battery or electronics disposal process.

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©2026 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.