




PHILIPS

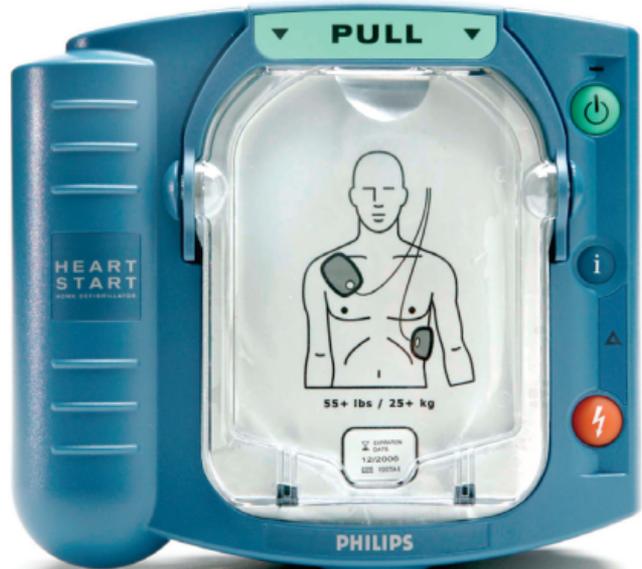

HeartStart Home Owner's Manual

Your guide to set up, operation,
maintenance, and accessories

IN AN EMERGENCY, CALL EMS FIRST!

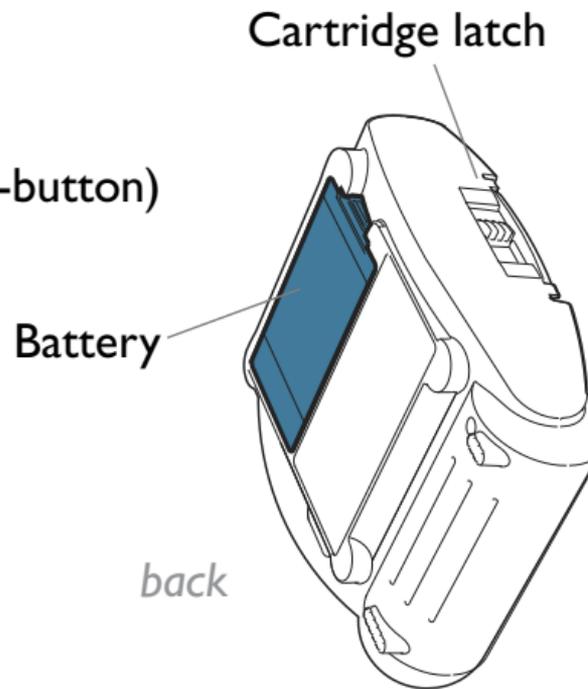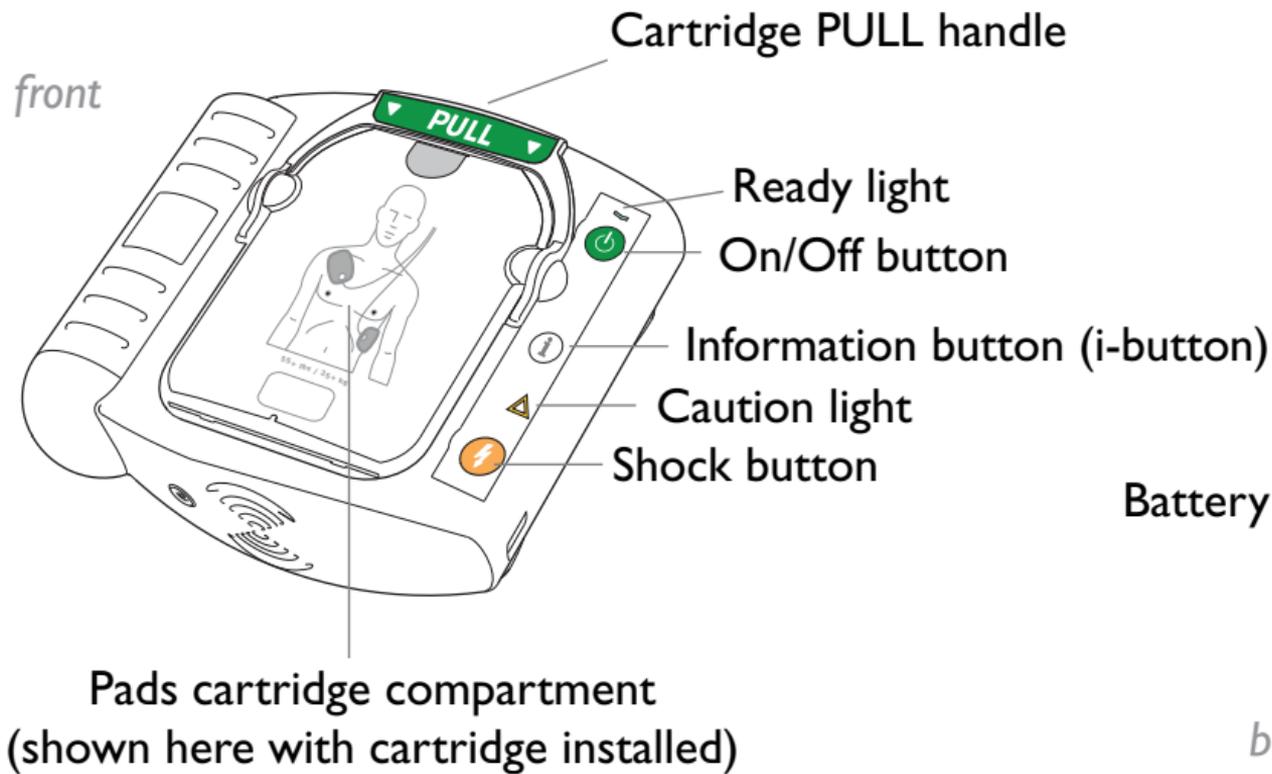

IMPORTANT NOTE

Time is critical in treating sudden cardiac arrest! Survival rates are directly related to how soon victims are defibrillated. For every minute of delay, the chance of survival drops by about 10%. However, it is important to understand that defibrillation cannot always assure survival. In some cases, the cause of the cardiac arrest simply cannot be survived no matter how quickly the victim is treated.

HeartStart Home Owner's Manual

Edition 10

PHILIPS

ABOUT THIS EDITION

The information in this guide applies to the model M5068A Philips HeartStart Home Defibrillator. This information is subject to change. Please contact us at 1.866.333.4246 or www.philips.com/AEDSupport for information on revisions.

EDITION HISTORY

Edition 10

Publication date: March 2015

Publication number: 453564544051

LANGUAGES

Please contact Philips at 1.866.333.4246 with questions about the HeartStart Defibrillator in other languages.

NOTICES

© 2015 Koninklijke Philips N.V.

All rights reserved.

Specifications are subject to change without notice.

Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

CONTENTS

Introduction	1	When to replace the battery	21
The HeartStart Defibrillator	1	When to replace the pads	22
Sudden cardiac arrest	1	How to clean it	23
Training and practice	3	How to have it repaired	24
State and local requirements	4	Common Questions	25
Setting up the HeartStart	5	Is it safe for the user and bystanders?	25
How to set it up the first time	5	Can I use it on an infant or child?	26
Recommended accessories	9	Can I use it on wet or metal surfaces?	27
Using the HeartStart	10	Can I travel with it?	27
When to use it	10	What if I donate the HeartStart Home to someone else?	28
How to use it	11	Ordering Accessories	29
How to know if CPR is needed	16	For More Information	30
What to do next	16	Warnings and Precautions	31
When EMS arrives	17	Glossary	35
Contacting Philips after a use	18	Warranty	37
Maintaining the HeartStart	19	Index	39
How to check that it's ready for use	19		
What to do when it chirps	19		

NOTES

Use this space to record information about the HeartStart Home Defibrillator.

HeartStart Product Line: HSI

One of the Philips HSI family of defibrillators.

Model Number: M5068A

Serial Number: _____

(on back of device)

Date of Purchase: _____

Place of Purchase: _____

Date of Registration: _____

Customer Support: 1.866.333.4246

SMART Pads Cartridge

installed on: _____

expiration date: _____

installed on: _____

expiration date: _____

HeartStart Battery

installed on: _____

expiration date: _____

installed on: _____

expiration date: _____

INTRODUCTION

THE HEARTSTART DEFIBRILLATOR

The Philips HeartStart Home Defibrillator is designed to be easy to use. The HeartStart talks the user through each step of treating someone who is in sudden cardiac arrest.

If you have concerns about your health or an existing medical condition, talk to your doctor. A defibrillator is not a replacement for seeking medical care.

SUDDEN CARDIAC ARREST

Sudden cardiac arrest (SCA) occurs when the heart unexpectedly stops pumping blood through the body. This is usually caused by a disorganized heartbeat called VF, for “ventricular fibrillation.”

Many victims of SCA do not have warning signs or symptoms. Some people may have a higher risk for SCA than others. Causes vary and may be different for infants and

children than for adults. Talk to your doctor about your risk factors.

The Philips HeartStart Home Defibrillator is used to treat VF by sending a shock across the heart, so it can start beating regularly again. Unless this is successful within the first few minutes after the heart stops beating, the victim is not likely to survive.

SCA is not the same as a heart attack, even though sometimes a heart attack can lead to SCA. If you feel the signs of a heart

attack — such as pain, pressure, or a squeezing feeling in your chest or anywhere in your upper body, shortness of breath, and stomach upset — you should call Emergency Medical Services (EMS) right away and also call your doctor.

For people with certain medical conditions who are likely to suffer SCA, an implantable cardioverter-defibrillator (ICD) may be advised. The HeartStart is not intended as a substitute for an ICD. Your doctor can tell you if an ICD might be right for you.

TRAINING AND PRACTICE

The HeartStart Home Defibrillator comes with training materials to help you become familiar with the signs of sudden cardiac arrest and how to use the HeartStart to help a victim of SCA.

Anyone who might use the HeartStart should watch the training video that comes with it and carefully read both this Owner's Manual and the Guide. Remember, [the HeartStart cannot be used by the victim of](#)

[SCA](#) — another person must use it to treat the victim.

It is important to review these materials periodically. They are intended to help make sure that, in an emergency, you will be able to use the HeartStart quickly and confidently.

In addition, you can buy reusable supplies from Philips for practicing how to use your HeartStart. See “Ordering Accessories” for ordering information.

Anyone who might use the HeartStart should also have training in CPR (cardio-pulmonary resuscitation). Many training organizations offer a combined course in CPR and defibrillator use. For a list of organizations offering this training in your area, contact us by phone at 1.866.DEFIBHOME (1.866.333.4246) or online at www.HeartStartHome.com.

STATE AND LOCAL REQUIREMENTS

Check with your state health department to see if there are any local or state requirements about owning and using a defibrillator. You can find contact information for learning about your state's requirements online, at www.philips.com/AEDSupport.



SETTING UP THE HEARTSTART


HOW TO SET IT UP THE FIRST TIME

It is easy to set up your HeartStart. Just follow these steps.

- 1 **Check the contents** of the HeartStart Home Defibrillator box.
- 2 **Take the HeartStart out** of its red case.

NOTE: Leave the pads cartridge installed. To avoid drying out the adhesive pads, do not open the cartridge until you need to use the pads in an emergency.

- 3 **Pull out the red tab** to start the automatic battery self-test.

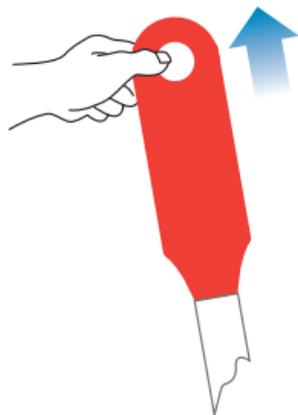

4 Press the flashing orange Shock button when the voice tells you to. (No shock will be given.) Wait about one minute. When the HeartStart passes the self-test, it turns off and the green Ready light starts blinking.

Your HeartStart came wrapped with a set of cardboard setup instructions. Tear the cardboard along the line at the bottom and remove the wrapping. Then reinstall the defibrillator in its carry case.

If there was no voice during the self-test, if the HeartStart tells you the Shock button was not verified when you pressed it, or if the green Ready light is not blinking at the end of the test, call 1.866.DEFIBHOME (1.866.333.4246) for customer service.

NOTE: Leave the battery installed. The HeartStart runs a self-test every day to help be sure that it is ready for use.

5 Watch the short training video that comes with the HeartStart. Plan to review it at least once a year.

Review the Quick Reference Guide, found inside the carry case. It provides step-by-step illustrated emergency directions for users who are hearing impaired or who are using the HeartStart in a noisy area. Then replace the Quick Reference Guide, face up, in the clear plastic window inside the carry case. (This Owner's Manual is also designed to be stored in the case, in the

pocket behind the Quick Reference Guide window.)

Write your address on the Call 911 card on the *back* of the carry case. (If 911 is not your local EMS number, use the back of the card and write in your EMS number.) A smaller Call 911 card is on the *front* of the case. (Use the back of the card if 911 is not your EMS number.)



6 Decide where to keep the HeartStart.

Some typical storage places include the family room, kitchen, or living room.

- It should be in a convenient, central area, close to a telephone, so that in an emergency the user can call Emergency Medical Services (usually 911) and get the HeartStart without wasting time.
- It should be a place where you can see and hear the HeartStart, so you can easily check that its green Ready light is blinking and hear it chirping if a self-test reveals a problem.

- It should be somewhere the temperature stays between 50° F and 109° F. (It will work at temperatures between 32° to 122° F.)

Be sure everyone who may use it knows where the HeartStart is kept.

Help us stay in touch with you. Fill out the HeartStart Support Program enrollment card in the box, or call us and we will enroll you. This allows Philips or the FDA to send you important notifications about the product.

RECOMMENDED ACCESSORIES

It is always a good idea to have a spare battery and a spare pads cartridge.

There are also some other things that are useful to keep with the HeartStart:

- a pair of scissors — to cut off clothing (comes with your HeartStart)
- disposable gloves — to protect the user
- a disposable razor — to shave the chest if hair prevents good pads contact
- a pocket mask or face shield — to protect the user

- a towel or absorbent wipes — to dry the victim's skin for good pads contact

Philips has a Fast Response Kit with all these items.

You can buy special Infant/Child SMART Pads,* for use on children less than 8 years or less than 55 pounds. Practice Kits with reusable adult or training pads and pad placement guides are also available. See “Ordering Accessories” for ordering information.

* Available by prescription only. Talk to your doctor.

USING THE HEARTSTART

WHEN TO USE IT


The HeartStart Home Defibrillator should *only* be used to treat someone you think may be a victim of SCA.

How can you tell if a person is in SCA? A person in SCA does not respond when shaken, and is not breathing normally.* Try to wake the person. Gently shake or tap

* Infants and children more often become unresponsive due to breathing or other issues than to heart problems. See page 26 for directions on using the HeartStart on an infant or child.

the person's shoulder and shout, "Are you OK?" Someone who is just sleeping will respond. Someone who is not responsive but not in sudden cardiac arrest will be breathing normally.

To check for breathing, put your ear near the person's mouth and look at the chest to see if it rises and falls.

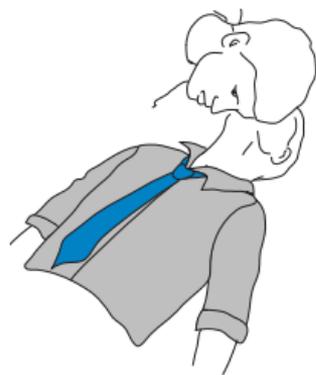

The HeartStart should not be used on a person who responds when shaken, or is breathing normally.

If you are not sure the person is in SCA, apply the HeartStart and follow its instructions. It will not let you shock someone unless a shock is needed.

HOW TO USE IT


If you find someone who may be in SCA, act calmly and quickly. Time is critical!

▶ **Call Emergency Medical Services** (usually 911) immediately or, if possible, have someone else call while you use the HeartStart. Clearly tell the emergency operator:


“This is (give your name and location).
Someone here has collapsed.
I have a HeartStart Home Defibrillator.
Please send help immediately.”

▶ **Quickly get the defibrillator.** Put it on the same side of the person that you are on.


Unzip the case and open it flat. To turn on the HeartStart, **lift the cartridge PULL handle**. Set the cover aside and follow the HeartStart's voice instructions.

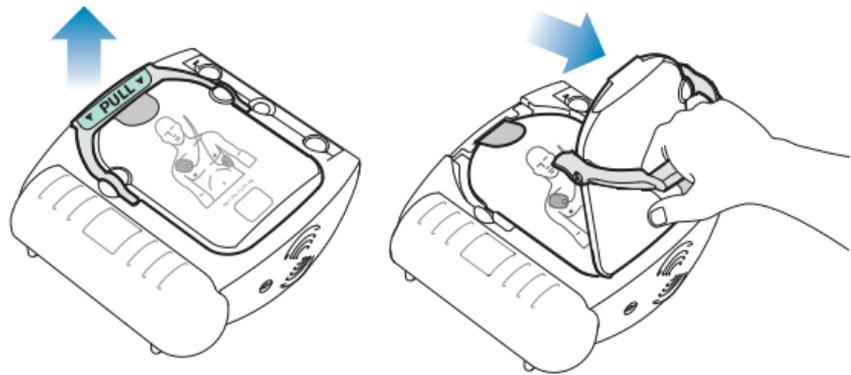

▶ Quickly **remove clothing from the person's chest**. If necessary, cut off the clothing with the scissors. (If the victim is an infant or young child, remove all clothing from the upper body, to bare both the chest and the back.)

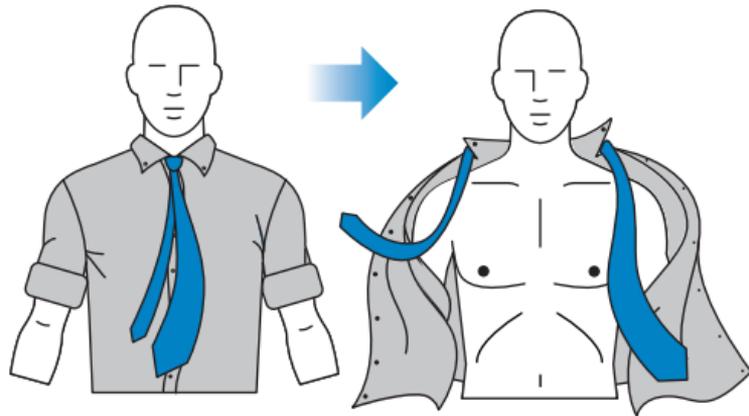

▶ Pull the tab at the top of the pads cartridge to **peel off the film seal.**

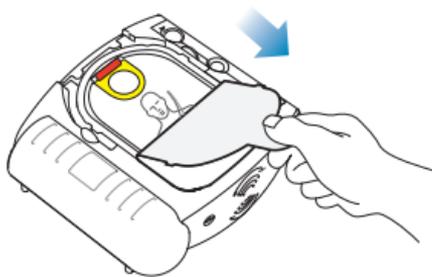

Inside are two pads on a yellow plastic liner. Lift the pads out of the cartridge. You will notice a picture on each pad showing where to place it on the patient.

NOTES: Before applying the pads to the patient, remove any medicine patches from the person's chest, any headphones, and other electrical devices with a wire or a cable.

It is okay to use the HeartStart on someone who has an implanted device such as a pacemaker or ICD, but do not place the pads directly over the implant. You can tell where an implant is located by a lump in the skin and a scar.

It is safe to use the HeartStart on someone wearing an oxygen mask.



▶ **Peel one pad off the liner.** Look at the picture on the pad. Place the pad with the sticky side down, on the person's bare skin,




as shown in its picture.

▶ Press the pad against the skin firmly. Then **peel and place the other pad** as shown in its picture.

Where to place pads on adults and children 8 years old or 55 pounds or more.

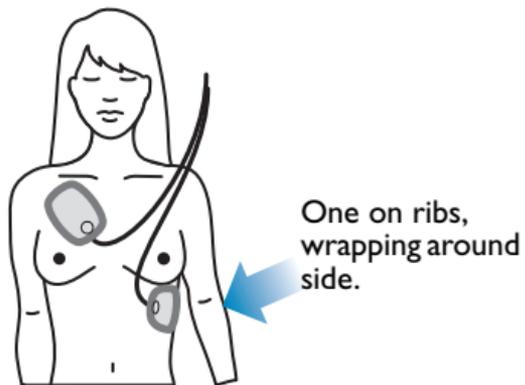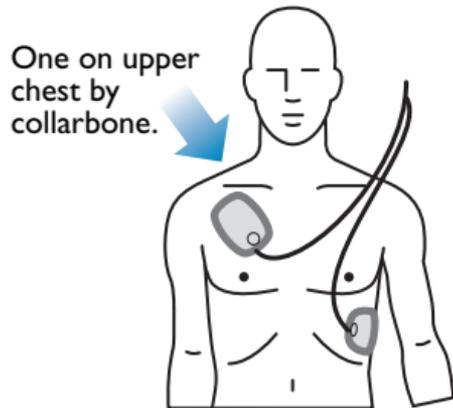

Where to place pads on children less than 8 years old or 55 pounds.

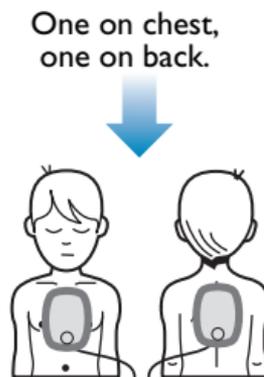

NOTE: Be sure to take each pad off the liner before applying it.

The pads have to stick well to the skin in order to work. If the person's skin is wet, dry it before applying the pads. If hair on the chest prevents the pads from sticking to the skin, clip or shave off the hair from the area where the pad will go.

As soon as both pads are in place, the HeartStart begins to analyze the patient's heart rhythm. It tells you that no one should be touching the patient and the Caution  light starts flashing.

▶ **If a shock is needed**, the HeartStart tells you to **press the flashing orange button** .

When you press the Shock button, the HeartStart gives a shock. It tells you the shock has been delivered, then pauses so you can perform CPR if needed. The HeartStart then analyzes the heart rhythm again to see if another shock is needed.



NOTE: No shock is given unless you press the button when it is flashing.

▶ **If a shock is not needed**, the HeartStart tells you it is safe to touch the patient (the Caution light  turns off). It reminds you to check the patient and **perform CPR if needed**. When a shock is not needed, it may be because the person is not in cardiac arrest and may need other medical attention. Make sure EMS has been called.

HOW TO KNOW IF CPR IS NEEDED

If the person is breathing, you do not need to do CPR. If the person is not breathing,

you should perform CPR. The HeartStart tells you to press the i-button  for coaching on how to perform CPR.

NOTE: For CPR coaching, you must press the i-button while it is flashing.

WHAT TO DO NEXT

Leave the HeartStart turned on. Leave the pads on. Follow the HeartStart's instructions and tend to the patient until EMS arrives. If another shock is needed, the HeartStart will give you step-by-step instructions.

Even if you did not use the HeartStart, be sure the patient gets immediate medical care. EMS should already have been called. If possible, contact the patient's personal doctor. *Shocking the heart can save a life, but it will not cure the underlying cause of SCA.*

WHEN EMS ARRIVES

When the EMS team arrives, they will take over caring for the patient. Explain as clearly as you can what has happened. Then tell them that if they hold down the i-button  until it beeps, the HeartStart

Home Defibrillator will report how many shocks it gave and how long ago it was turned on. (More detailed information can be transferred from the HeartStart's memory to a computer running special software.)

After care is transferred to EMS, *turn off the HeartStart*: press and hold down the green On/Off button  until it turns off (about one second).

NOTE: The HeartStart pads should be taken off before another defibrillator is used.

CONTACTING PHILIPS AFTER A USE

If you have used or seen use of the HeartStart Home Defibrillator, we urge you to contact us by phone at 1.866.DEFIBHOME (1.866.333.4246) or email at www.HeartStartHome.com. You will be helping us learn more about sudden cardiac arrest.

We are here to answer any questions you may have about the experience.

Also, please let us know if you ever have any problems using or maintaining your HeartStart Home Defibrillator.

MAINTAINING THE HEARTSTART

HOW TO CHECK THAT IT'S READY FOR USE

Check the HeartStart regularly to be sure the green Ready light is blinking. You can see the Ready light through the window in the red carry case. With the battery installed, the HeartStart tests itself every day. As long as the green Ready light is blinking, you can be confident that the defibrillator has passed its latest self-test and is ready for use.

WHAT TO DO WHEN IT CHIRPS

Your HeartStart tests itself daily, weekly, and monthly so you can be sure it is ready to use in an emergency. Any time the HeartStart chirps, it is telling you something is wrong. The green Ready light will be off, which means you need to take action so that the HeartStart will be ready to save a life.

- If you hear a pattern of single chirps (♪ ... ♪ ... ♪...) coming from your HeartStart, the i-button  will be flashing. Press it, and the HeartStart will tell you what to do next.

NOTE: If the HeartStart says to take out the battery then put it back in, wait for at least five seconds before putting it back in.

- If you hear a pattern of triple chirps (♪♪♪... ♪♪♪... ♪♪♪...) coming from your HeartStart, it could mean something serious is wrong with your AED. It may not be ready to use in an emergency, so

it's important that you call technical support.

- As long as you're not in an emergency, call Philips for technical support right away at 1.866.333.4246 if you hear triple chirps.
- If you're in an emergency and you hear triple chirps, try to use the device. You can take out the battery and put it back in, which will sometimes allow the HeartStart to work. As soon as the emergency is over, call Philips even if the device isn't chirping anymore.

WARNING: If you take out the battery and put it back in more than once when you hear a pattern of triple chirps, your HeartStart might seem ready to use, even though there's a chance it won't be ready to use in an emergency. If you hear a pattern of triple chirps from your HeartStart, *please call Philips immediately*, unless you're in an emergency and need to try to use your AED.

WHEN TO REPLACE THE BATTERY

The lithium battery used in the HeartStart normally lasts for about four years in standby mode. (“Standby mode” means the HeartStart is ready for use, with the battery installed and the green Ready light

blinking.) When the HeartStart is turned on for use, a new battery normally lasts for 4 hours, or can give up to 200 shocks.

After you put the battery in the HeartStart, the self-test checks it every day. When battery power gets low, the HeartStart will start chirping to alert you.

The battery is not rechargeable. Do not try to open, crush, or burn it. After it is used up, throw it away in household waste or according to local regulations.

NOTE: Practicing or showing how to use the HeartStart with training pads uses battery power. For example, one hour of training will use about 10% of the battery power.

See “Setting Up the HeartStart” for how to put in the battery. See “Ordering Accessories” for ordering information.


WHEN TO REPLACE THE PADS

The SMART Pads are designed to last for two years in an unopened, undamaged cartridge. Each cartridge has an expiration date printed on the film seal and the package. The HeartStart checks the pads

every day as part of its self-test. If for any reason the pads are no longer good, the HeartStart will start chirping and the green Ready light will turn off.

The pads are designed for only one use. Replace the pads when they expire or are used.

To replace used pads or swap adult pads and infant/child pads cartridges, slide the latch at the top of the HeartStart over to the side. Then lift out the cartridge, and put in the new one.

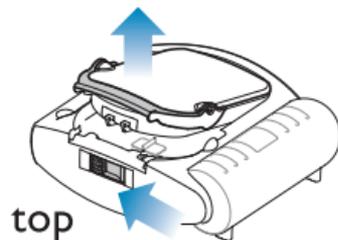

If used, throw away the pads and cartridge in household waste or according to local regulations.

See “Setting Up the HeartStart” for directions on how to put in a new cartridge. See “Ordering Accessories” for ordering information.

HOW TO CLEAN IT

You can clean the outside of the HeartStart and its red carry case with a soft cloth dampened with:

- soapy water
- chlorine bleach
(2 Tbsp/quart or liter of water)
- ammonia-based cleaners

Do not put any part of the HeartStart into water or other fluids. Do not use alcohol, rough materials, cleaners, or strong solvents.

HOW TO HAVE IT REPAIRED

The HeartStart is not designed for repair by the user. Do not open the HeartStart, remove its covers, or try to repair it. You could be injured by dangerous high voltages and currents.

If the HeartStart is damaged, please call 1.866.DEFIBHOME (1.866.333.4246) for service.

COMMON QUESTIONS

Some of the most common questions about using the HeartStart are answered here.

See “Warnings and Precautions” for more information.

Is it safe for the user and bystanders?

There is no risk of harm to the user or a bystander as long as:

- the HeartStart is used according to directions, and

- no one is touching the patient when the Shock button is pressed.

The HeartStart cannot deliver a shock unless it is used on someone whose heart needs a shock.

Even if you press the orange Shock button, the HeartStart will not deliver a shock unless the defibrillator tells you to push it.

Can I use it on an infant or a child?

Yes. For children younger than 8 years or those who weigh less than 55 pounds, special HeartStart Infant/Child SMART Pads are recommended. When used with these special pads, the HeartStart delivers less energy, and the optional CPR instructions are tailored for use on a child.

Infants and children more often become unresponsive due to breathing or other issues than to heart problems. Perform CPR first before calling for 911/EMS or using the HeartStart.

For children 8 years or older, or who weigh 55 pounds or more, use the adult pads that come with the HeartStart.

WARNING: *If you are not sure about the child's exact age or weight, or if infant/child pads are not available, DO NOT DELAY TREATMENT.*

Use the adult pads, but place them on the child's chest and back, as shown on page 14. (With the adult pads on a child, the optional CPR coaching will be for an adult.) **To treat an infant or child:**

- Perform infant/child CPR while a bystander calls 911/EMS and brings the HeartStart.
- If no bystander is available, perform 1-2 minutes of CPR before calling 911/EMS, then apply the HeartStart.
- If you witnessed the child's collapse, call EMS *immediately* and *then* get the HeartStart.

Can I use it on wet or metal surfaces?

Yes. It is safe to use the HeartStart on someone lying on a wet surface. It is also safe to use it on someone lying on a metal surface. As always, no one should touch the person when the Shock button is pressed. Naturally, the user should get the victim out of standing water, such as a pool or bathtub, before applying the HeartStart. It is important to dry the chest completely, so the pads will stick well to the skin.

Can I travel with it?

Yes. The HeartStart is designed to be easily

portable. Just remember that if you take it along when you travel, you should keep it where it can stay within its temperature range (see page 8). When you get back, put the HeartStart in its normal place, so everyone will know where to find it in an emergency.

There may be rules against taking a defibrillator into some countries. Before you leave, find out if there are restrictions.

Also, you may want to check with the airlines before boarding the aircraft with scissors in the case.

What if I donate the HeartStart Home to someone else?

The training materials included with this HeartStart Home defibrillator are designed to help users to use the device in a home setting, or on family outings. **If this device will be given or donated to another individual for use in the home**, be sure to provide the training materials included with the HeartStart Home and visit www.HeartStartHome.com to find additional training resources.

HeartStart defibrillators that are used by an organization, or in a facility that is not a home, require different considerations. These may include additional training, a review of the organization's medical emergency response policies, a site review to ensure proper placement, and medical oversight, among other things. **If this device will be given or donated to an organization where it will be used in a facility**, be sure that the recipient organization seeks medical oversight.

ORDERING ACCESSORIES

Accessories for the HeartStart Home Defibrillator include:

- M5070A battery pack
- M5071A Adult SMART Pads cartridge
- M5072A Infant/Child SMART Pads cartridge*
- M5073A Adult Practice Kit with training pads cartridge and pads placement guide
- M5074A Infant/Child Practice Kit with training pads cartridge and pads placement guide

* Available by prescription only. Talk to your doctor.

- M5090A adult pads placement guide
- M5075A standard carry case (for extra storage room)
- M5076A slim carry case (comes with the HeartStart Home)
- 68-PCHAT Fast Response Kit (pouch with a pair of scissors, 2 pairs of disposable gloves, a disposable razor, a pocket mask, and an absorbent wipe)

Call 1.866.DEFIBHOME (1.866.333.4246)
or visit www.philips.com/AEDSupport to order.

FOR MORE INFORMATION

If you have any questions about the HeartStart Home Defibrillator after you watch the video and read this *Owner's Manual*, please contact us at 1.866.DEFIBHOME (1.866.333.4246).

We will be happy to answer your questions and to provide you with more information.

You can find detailed *Instructions for Use* for the HeartStart Home Defibrillator online at www.philips.com/AEDSupport. Clinical summaries of research studies about the technology used in the HeartStart are also available online.

WARNINGS AND PRECAUTIONS

It is important to understand how to use your HeartStart Home Defibrillator safely. Please read these warnings and precautions carefully.

A warning describes something that could cause serious personal injury or death. A precaution describes something that could cause minor personal injury, damage to the HeartStart, loss of data stored in the

HeartStart, or less chance of successful defibrillation.

NOTE: The HeartStart Home Defibrillator is designed to be used only with Philips-approved accessories. The HeartStart may perform improperly if non-approved accessories are used.

WARNINGS

- **Flammable gases.** If the HeartStart is used to give a shock in the presence of flammable gases such as in an oxygen

tent, there is a risk of explosion.
(However, it is safe to use the HeartStart on someone wearing an oxygen mask.)

- **Battery.** The HeartStart M5070A battery is not rechargeable. Do not try to recharge, open, crush, or burn the battery, or it may explode or catch fire.
- **Fluids.** Do not let fluids to get into the HeartStart. Avoid spilling any fluids on the HeartStart or its accessories. Spilling fluids into the HeartStart may damage it

or cause a fire or shock hazard. Do not sterilize the HeartStart or its accessories.

- **Accessories.** Using damaged or expired equipment or accessories may cause the HeartStart Home Defibrillator to perform improperly, and/or injure the victim or the user.
- **Patient Handling.** Performing CPR or otherwise handling or moving the patient while the HeartStart is analyzing heart rhythm can cause an incorrect or delayed analysis. If the HeartStart tells you a

shock is advised while you are handling or moving the patient, stop the vehicle or CPR and keep the patient as still as possible for at least 15 seconds. This will give the HeartStart time to reconfirm the analysis before telling you to press the Shock button.

- **Children.** Keep the HeartStart out of reach of children to avoid the potential risk of inhalation or swallowing of small parts or strangulation by the pads cables.

- **Cell Phones.** The HeartStart can work correctly when it is fairly close to equipment like emergency two-way radios and cell phones. Normally, using a cell phone near the patient should not cause a problem for the HeartStart. However, it is best to keep such equipment only as close as necessary to the patient and the HeartStart.

PRECAUTIONS

- **Device Handling.** The HeartStart was designed to be sturdy and reliable for many different use conditions. However, handling the HeartStart too roughly can damage it or its accessories and will invalidate the warranty. Check the HeartStart and accessories regularly for damage, according to directions.
- **Maintenance.** Improper maintenance may damage the HeartStart or cause it to function improperly. Maintain the HeartStart according to directions.
- **Skin Burns.** Do not let the pads touch each other or other electrodes, lead wires, dressings, medicine patches, etc. Such contact can cause electrical arcing and skin burns during a shock and may also divert the electrical current away from the patient's heart. During a shock, air pockets between the skin and pads can cause skin burns. To help prevent air pockets, make sure pads stick well to the skin. Do not use dried out pads because they will not provide good contact with the skin.

GLOSSARY

AED. Automated external defibrillator. The HeartStart Home Defibrillator is an AED.

CPR. Cardiopulmonary resuscitation, a way to provide artificial breathing and heartbeats to keep a victim of SCA alive until a shock can be delivered.

i-button. An information button on the front of the HeartStart. Press the i-button while it flashes during a pause to get CPR coaching. Press the i-button if it flashes at other times to get troubleshooting information. At any time, you can press and hold the i-button (for at least

two seconds) to get summary information about defibrillator's last use and condition.

On/Off button. A green button on the front of the HeartStart Home Defibrillator. Press the On/Off button to turn the defibrillator on or off, disarm the HeartStart, or stop the self-test that runs when a battery is put in.

Ready light. A small, bright green light that blinks to show that the HeartStart is ready for use. When the HeartStart is being used, the Ready light stays on without blinking.

Shock button. An orange button with a lightning bolt symbol, on the front of the HeartStart. The Shock button flashes when a shock is needed. You must press the button for the shock to be given.

SMART Pads. The pads, packed in a cartridge, used with the HeartStart. Lifting the cartridge PULL handle turns on the HeartStart and opens the cartridge. The pads stick to the patient's bare skin. The HeartStart uses the pads to sense the patient's heart rhythm and to transfer the shock, if needed.

Standby mode. The status of the HeartStart when the battery is in place, and the HeartStart

is ready for use when needed. A blinking green Ready light means the HeartStart is in standby mode.

SCA. Sudden cardiac arrest. The sudden stopping of the heart's ability to pump blood. A person in SCA does not respond when shaken and is not breathing normally.

Ventricular fibrillation. A disturbance of the normal heart rhythm that causes a disorganized heartbeat that cannot pump blood through the body and brain. It is usually the cause of sudden cardiac arrest.

WARRANTY

HeartStart Home Defibrillator

Philips will replace or repair the HeartStart Home Defibrillator, model M5068A, within eight years of the original retail purchase, if it is found to be defective in materials or workmanship. This warranty, applicable only to residential installations, is valid for the original retail purchaser from the date of initial retail purchase. Save your receipt as proof of date of purchase. This warranty gives you specific legal rights and you may also have other rights which vary from state to state. Replacement or repair will be at Philips' sole discretion.

HeartStart HSI Defibrillator Family Accessories

Philips will replace or repair accessories for the HeartStart HSI family of defibrillators, within the specified

warranty period if they are found to be defective in materials or workmanship or if they do not meet the specifications published by Philips. The HSI battery is warranted for four years from the date of installation. The SMART Pads cartridge is warranted until the expiration date listed on the cartridge. All other accessories are warranted for one year from the date of purchase. This warranty is valid for the original purchaser from the date of initial retail purchase. Save your receipt as proof of date of purchase. This warranty gives you specific legal rights and you may also have other rights which vary from state to state.

For All HeartStart Products

What is not covered under warranty: normal wear and tear; products that are modified or subjected to unusual

physical or electrical stress; misuse; neglect; alteration; improper installation, repair or testing. All HeartStart products must be used in accordance with instructions and properly maintained as required in the Owner's Manual in order to be eligible for warranty protection.

How to Obtain Warranty Service

For prompt warranty service, call 1.866.DEFIBHOME (1.866.333.4246). A customer service representative will give you detailed instructions on how to fulfill your HeartStart Home Defibrillator warranty.

Implied Warranties

Any implied warranties, including implied warranties of merchantability and fitness for a particular purpose, shall be limited to the duration of the express warranties set forth in this limited warranty. Some states do not allow

limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

Limitation of Remedies

In no event shall Philips be liable for any special, incidental or consequential damages based upon breach of warranty, breach of contract, negligence, tort, or any other legal theory. Such damages include, without limitation, loss of savings or revenue, loss of profit, loss of use, the claims of third parties, and cost of any substitute equipment or services. Some states do not allow the exclusion or limitations of incidental or consequential damages so the above limitation or exclusion may not apply to you.

Replacement

HeartStart Home Defibrillator replacement accessories can be purchased from Philips at 1.866.DEFIBHOME (1.866.333.4246) or www.philips.com/AEDSupport.

INDEX

A

accessories 9, 29

B

battery replacement 21

C

checking the HeartStart 19

chirping 19

cleaning the HeartStart 23

CPR coaching 16

I

infants and children 14, 26

P

pads position 14

pads replacement..... 22

pressing the Shock button 15

R

repairing the HeartStart 24

S

safety 25, 31

set up 5

standby mode 21, 36

storing the HeartStart 8

sudden cardiac arrest (SCA) 1, 36

symptoms of SCA 10

T

training 3

traveling with the HeartStart 27

turning off the HeartStart 17

turning on the HeartStart 12

NOTES



HAVE YOU USED YOUR HEARTSTART?

If you have been in a situation that called for using the HeartStart Home Defibrillator, we urge you to contact us by phone at 1.866.DEFIBHOME (1.866.333.4246) or email at www.HeartStartHome.com to discuss the experience.

PHILIPS

Philips Medical Systems
22100 Bothell Everett Highway
Bothell, WA 98021- 8431, USA



Authorized EU Representative

Philips Medizin Systeme Boeblingen GmbH
Hewlett-Packard Strasse 2
71034 Boeblingen, Germany



453564544051