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English

1 Important

Read this user manual carefully before you use the appliance, and save the user manual for future reference. You can also find this user manual at Philips website: www.philips.com/support

- Before each use, check the manual breast pump. If any damage or crack is found, stop using it immediately.
- Use the supplied components only.
- Use the manual breast pump as instructed in this user manual. Improper usage could harm the breast.
- Never use feeding teats as dummies, pacifiers, or soothers.
- Keep all the components away from your baby.
- This manual breast pump is intended for personal use only for hygiene reasons.
- Food coloring can discolor components, but it will not affect your usage.
- Make sure that the teething baby does not bite the teat too long. The continuous and prolonged bite will damage the teat.
- For hygiene reasons, replace teats after three months of use. Use only Philips AVENT teats and bottles together.


2 Your manual breast pump

Introduction

The silicone diaphragm of Philips AVENT manual breast pump provides 100% reliable vacuum, whilst the Let-down Massage Cushion naturally imitates your baby's suckling action to provide fast milk flow – silently, comfortably, and gently. Doctors recommend that breast milk is the best nutrition for babies during the first year, combined with solid food after the first six months. Your milk is specially adapted to your baby's needs and contains antibodies which help protect your baby against infection and allergies.

Your manual breast pump will help you breast-feed longer. You can express and store your milk so that – whether you are out for a few hours, taking a well deserved rest, or back at work – your baby can still enjoy the benefits of your breast milk, even if you can not be there to give it yourself. As the pump is easy to assemble, silent and discreet you can take it with you anywhere, allowing you to express milk at your own convenience and maintain your milk supply.

Overview

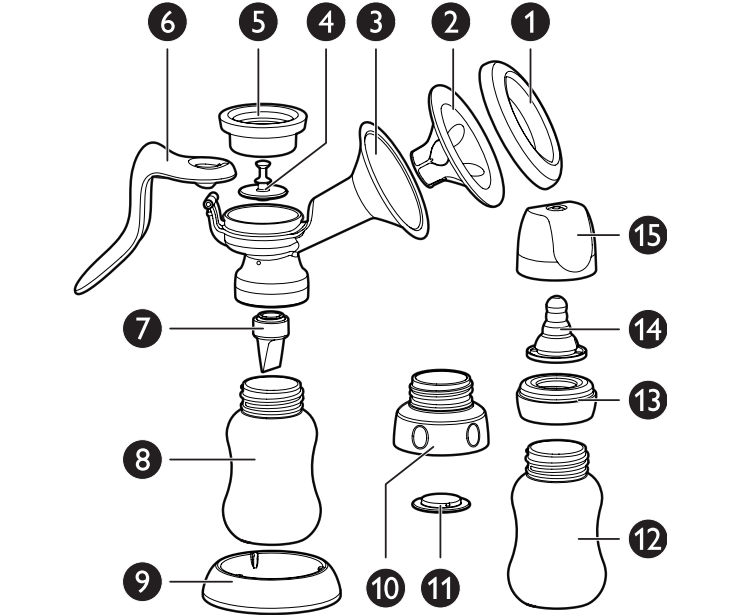Manual breast pump	Bottle (Standard)
1 Pump cover	10 Adaptor ring for wide neck bottles
2 Let-down Massage Cushion	11 Sealing disc
3 Pump body	12 Bottle (120 ml) – only available for SCF900/00
4 Stem	13 Screw ring
5 Silicone diaphragm	14 Teat (0m+)
6 Handle	15 Dormal cap
7 Valve	
8 Bottle (120 ml)	
9 Stand	

- ✱ **Tip**
 - Extra Philips AVENT travel pouches/ bags are sold separately.

3 Use your manual breast pump

Clean and sterilize

Before each use, clean and sterilize the pump.


- **To clean:**
 - Disassemble all the components, wash them in warm soapy water, and then rinse them thoroughly.
 - Never use abrasive, anti-bacterial cleaning agents, chemical solvents, or excessive concentration of detergents.
 - If necessary, clean the components with a Philips AVENT cleaning brush sold separately.
 - Take care when cleaning the valve. Rub it gently between your fingers in warm soapy water. Do not insert objects as this may cause damage. If it gets damaged, your breast pump does not function properly.
- **To sterilize:**
 - Sterilize all the components in a Philips AVENT electric or microwave steam sterilizer. The pump remains sterile in the sterilizer for a minimum of six hours, provided the lid is not removed. For details, see the sterilizer user manual.
 - Sterilize all the components in boiling water for five minutes. Make sure that the components do not contact the container too long, to avoid heat deformation.
 - If you use a sterilizing solution, do not leave the components in the solution for longer than recommended.
 - Do not store the components in drying/sterilizing cabinets with heating elements. (The maximum temperature of durance is 110°C.)

After each use, disassemble all the components, clean, and sterilize them as instructed in this user manual, to ensure hygiene.

Assemble


- ✳ **Note**
 - Make sure that you clean and sterilize the pump as instructed in this user manual.
 - Make sure that you assemble the pump correctly as instructed in this user manual.
 - To avoid leakage: Remove the residue around the rim before assembling. Do not overtighten.
 - Wash your hands thoroughly before you take sterilized components.

- 1 Insert the valve into the pump body from underneath.

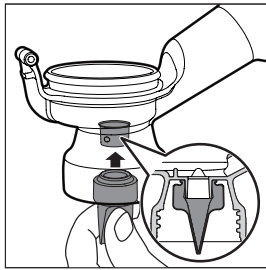
- 2 Put the bottle onto the pump body and turn it anticlockwise to fix it.


- ✳ **Note**
 - Do not overtighten when assembling the bottle to the pump body.

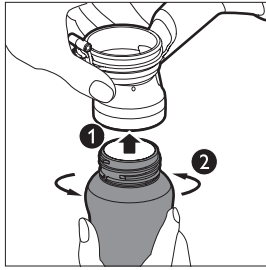- If you use a Philips AVENT wide neck bottle, put on the adaptor ring.

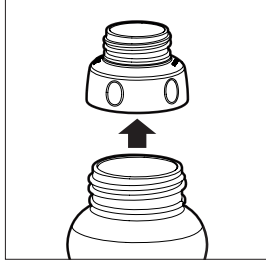
- 3 Put the stem into the silicone diaphragm.

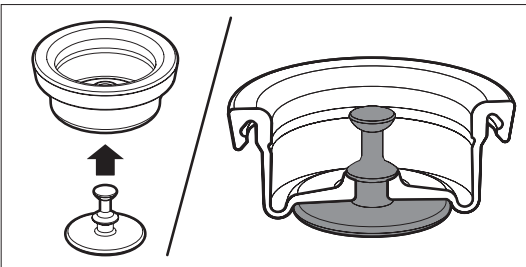
- 4 Put the silicone diaphragm onto the pump body and press the rim to ensure a perfect seal.


- ✱ **Tip**
 - It is easier if assembled wet.

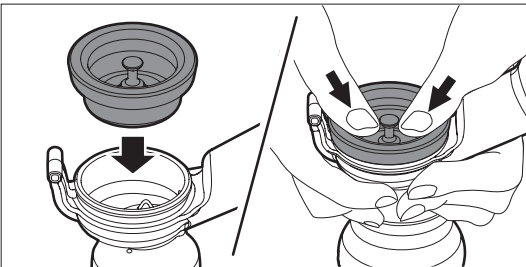- 5 Fasten the circle end of the handle on the stem and press the handle gently until it clicks into place.

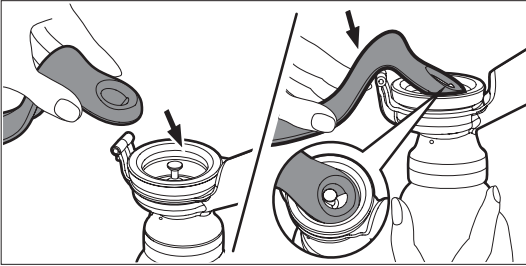- 6 Insert the Let-down Massage Cushion into the pump funnel and press the rim to ensure a perfect seal. (This is easier if assembled wet.)

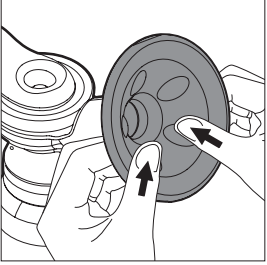- 7 To make the bottle stands stable, insert it into the stand.

- ✳ **Note**
 - To prevent toppling, insert the bottle into the stand after each use.

When to express milk

If possible, wait until your milk supply and breast-feeding schedule are well established, normally at least 2-4 weeks after birth, unless told otherwise by your doctor or breast-feeding advisor. The exceptions are:

- If you are expressing milk for your baby to be given in hospital.
- If your breasts are engorged (painful or swollen). You can express a small amount of milk before or between feeds to ease the pain and to help your baby latch on more easily.
- If you have sore or cracked nipples and you may wish to express your milk until they are healed.
- If you are separated from your baby and wish to continue to breast-feed when you are reunited. You should express your milk regularly to stimulate your milk supply.

Right times of the day to express your milk:

- Early in the morning when your breasts are full, just before or after your baby's first feed.
- After a feed, if your baby has not emptied both breasts.
- Between feeds or in a break if you have returned to work.

Hints to help you succeed

The pump is easy to use. Take some practice and you can get used to expressing your milk.

- Choose a time when you are not rushed and will not be interrupted.
- Having your baby or a photograph of your baby nearby can help encourage 'let-down'.
- Placing a warm cloth on the breasts for a few minutes before you start pumping can encourage milk flow and soothe painful breasts.
- Warmth and relaxation can encourage milk flow. Try expressing after a bath or shower.
- Try expressing from one breast while your baby is feeding from the other, or continue expressing just after a feed.
- Repositioning the pump on the breast from time to time can help stimulate your milk ducts.

- ✱ **Tip**
 - Practice with your pump to find the best technique for you. However, if the process becomes painful, stop and consult your breast-feeding advisor

Express breast milk

- ✱ **Tip**
 - Make sure that the Let-down Massage Cushion is clean.
 - If you find that you regularly express more than 120 ml per session, you can use a 240 ml Philips AVENT bottle, sold separately. Do not overflow the Philips AVENT bottle but keep upright, to avoid leakage.
 - If the process becomes painful or uncomfortable, stop using the pump and consult your health professional or breast-feeding advisor.

- 1 Wash your hands thoroughly and make sure that your breasts are clean. Gently squeeze a little milk from each nipple to make sure that milk ducts are not blocked.
- 2 Make sure that you sterilize and correctly assemble your pump as instructed in the "Clean and sterilize" and "Assemble" sections.
- 3 Express milk at the right time. For details, see the "Hints to help you succeed" section.
- 4 Relax in a comfortable chair, leaning slightly forward (using cushions to support your back). Make sure that you have a glass of water nearby.

Correction position	Maximum angle	Incorrect position

- 5 Press the funnel and the let-down massage cushion of the pump firmly against your breast. Make sure that no air escapes – otherwise there is no suction.
- 6 As you gently press the handle, you will feel the suction on your breast. You do not need to depress the handle fully to create a vacuum – only as much as it is comfortable. Your milk will soon start flowing, even though you do not use all the suction that the pump can generate.
- 7 Begin by pumping 5-6 times rapidly to initiate let-down. Then, hold the handle down for 2-3 seconds, and allow it to return to its resting place. These 2-3 second cycles imitate your baby's natural suckling pattern and allow the milk to flow between strokes.
- 8 Milk begins to flow after the first few strokes. Do not worry if your milk does not flow immediately. Relax and continue pumping. Do not continue pumping for more than 5 minutes at a time if no result is achieved. Try expressing at another time during the day.
- 9 In average, you need to pump for 10 minutes to express 60-120 ml. This depends from person to person. If you notice that you can rhythmically express more than 120ml each time, use the respective size of Philips AVENT milk bottle for storage. Never tilt the breast pump when the milk has reached the maximum line on the bottle. Replace with a new bottle before re-pumping to avoid spillage.

Store your breast milk

- ✱ **Caution**
 - Never refreeze breast milk.
 - Do not add fresh breast milk to already frozen milk.
 - Discard any milk that is left over at the end of a feed.


- Store expressed milk in Philips AVENT bottles fitted with sealing discs and refrigerated in the fridge or freezer immediately.
- You can store expressed breast milk in the refrigerator for up to 48 hours (not in the door) or in the freezer for up to 3 months. Label the bottles with the date and time of expression and use up older milk first.
- Store only milk collected with a sterile pump to feed your baby.

Feed your baby your expressed milk

- If you store breast milk in a Philips AVENT bottle, replace the sealing disc with a sterilized Philips AVENT teat before feeding.
- Defrost frozen breast milk overnight in the refrigerator. However, in an emergency, you can defrost it in a jug of hot water. Once it is thawed, use within 24 hours.
- To heat breast milk:
 - Stand the Philips AVENT bottle in a bowl of hot water.
 - Use the Philips AVENT bottle and baby food warmer. Always check milk or food temperature before feeding.
 - Never heat breast milk in a microwave as this may destroy valuable nutrients and antibodies. In addition, the uneven heating may cause hot spots in the liquid which could cause burns to the throat.
 - Never immerse a frozen bottle in boiling water to avoid crack or split.

Choose a right teat for your baby

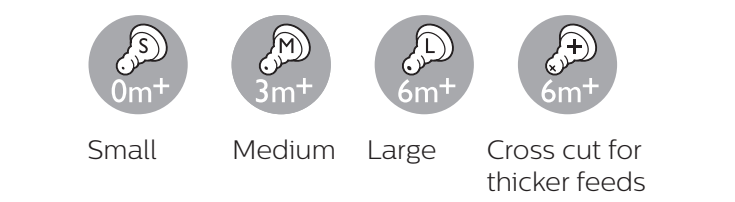Philips AVENT teats in four flow rates are designed for babies at different ages. The replacement teats are sold separately.

The four flow rates are indicated on the teats as S (small), M (medium), L (large), and + (cross cut for thicker feeds).

Philips AVENT standard teats are designed for the Philips AVENT standard bottles.

Choose a teat proper for your baby. The age proposed is for your reference only.

Store the manual breast pump

- Keep the small components away from children.
- Keep all the components away from the source of heat or direct sunlight.
- For hygiene maintenance, disassemble the manual breast pump and store all the components in a dry and covered container.

4 Troubleshooting

Problem	Solution
Lack of suction	• Make sure that all the components are assembled correctly and that the valve is the right way up. • Make sure that the diaphragm and Let-down Massage Cushion are fitted securely to the pump body and a perfect seal is created.
No milk being expressed	• Make sure that suction is created and the pump is correctly assembled. Relax and try again, expression may improve after practice. Please refer to the "Hints to help you succeed" section.
Pain in the breast area when expressing	• You may pump too hard. You do not need to use all the suction that the pump can generate. Try pressing the handle halfway and remember the 2-3 seconds pumping cycles. • Consult your health professional/ breast-feeding advisor (see section "Hints to help you succeed").
Milk accumulates in the Let-down Massage Cushion	• Remove and refit the Let-down Massage Cushion. Make sure that it is installed correctly and firmly. • Try leaning lightly forward.
Cracking/ Discoloring	• Avoid contact antibacterial or abrasive detergents/cleaners. • Combinations of detergents, cleaning products, sterilization solution, softened water, and temperature fluctuations may cause the plastic to crack. If any crack is found, stop using the pump immediately. Contact Philips AVENT consumer care center. • Food coloring can discolor components, but it will not affect your usage. • Keep all the components away from the source of heat or direct sunlight as it can discolor components.
Components broken or lost	• Contact Philips AVENT consumer care center.

If you are unable to solve the problem with the provided information, contact the Philips AVENT Consumer Care staffs. You can also consult your health professional or breast-feeding advisor.

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All rights reserved.