PHILIPS

AVENT

Single/Double
Electric Breast
Pump

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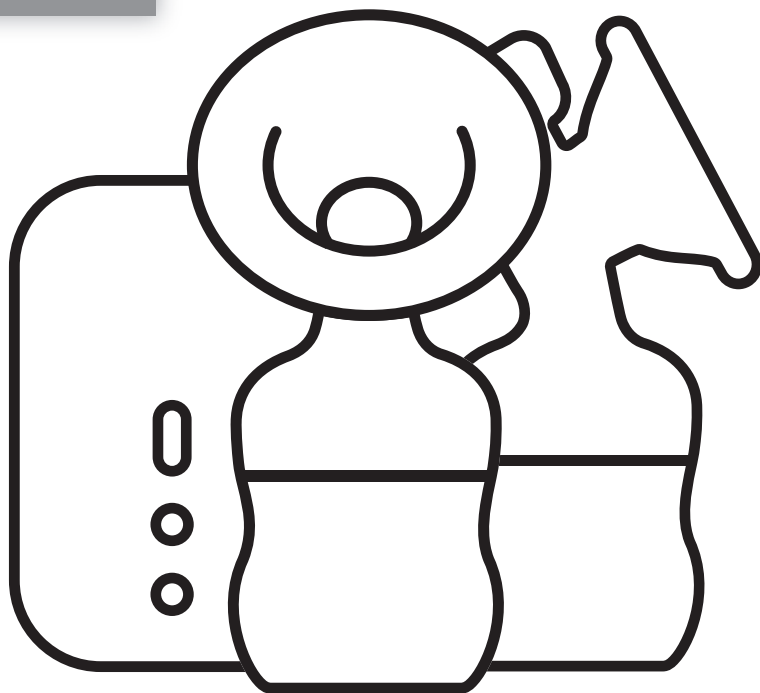
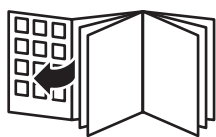
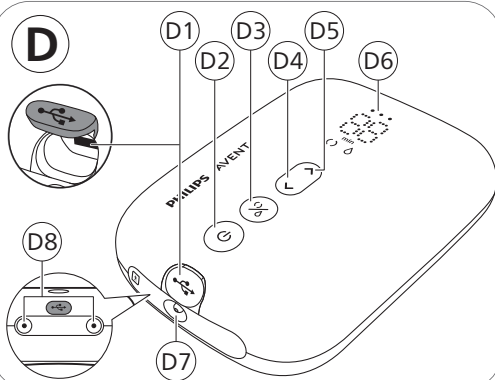
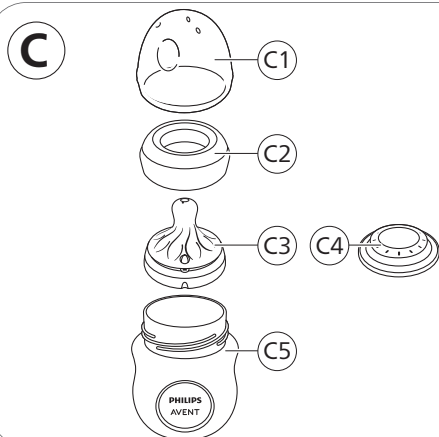
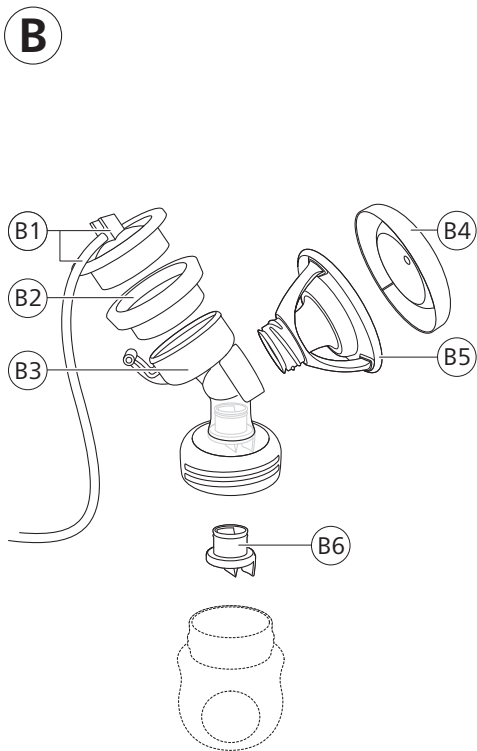
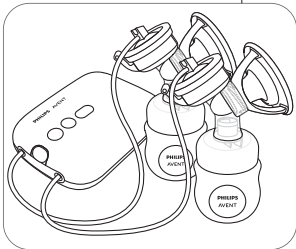
English 6
Español 38
Français (Canada) 70

Contents

Introduction	6
General description	7
Indications for use	8
IMPORTANT SAFETY INSTRUCTIONS	8
Before first use	10
Cleaning and disinfecting	11
Using the breast pump	15
Using the breast pump when you have symptoms of engorgement	22
After use	22
Feeding breast milk with the bottle	24
Compatibility	27
Ordering accessories and spare parts	27
Disposal	27
Full Two-Year Warranty	28
Support	28
Troubleshooting	28
Undesirable side effects	30
Storage conditions	30
Electromagnetic fields (EMF)	31
Technical information	31
Explanation of symbols	34

Introduction

Congratulations on your purchase and welcome to Philips Avent! To fully benefit from the support that Philips Avent offers, register your breast pump at www.philips.com/mybreastpump.

The Philips Avent single/double electric breast pump Advanced is inspired by babies and their natural suction movement. The soft silicone flange stimulates the nipple with gentle massaging movements, helping you to express milk comfortably and effectively. The combination of suction and nipple stimulation triggers fast milk flow. The no-lean-forward design enables you to sit in a relaxed, upright position while expressing.

To breastfeed longer, you can express and store your milk so that your baby can still enjoy the benefits of it, even if you cannot be there to provide it yourself. As the pump is compact, quiet and discrete to use, you can take it with you anywhere, allowing you to express milk at your own convenience and to maintain your milk supply.

If breastfeeding goes well, it is advisable (unless recommended otherwise by your healthcare professional) to wait until your milk supply and breastfeeding schedule are established (normally at least 2 to 4 weeks after giving birth) before you start expressing.

General description

Note: The numbers mentioned below refer to the numbers in the figure on the front foldout page of this user manual.

Note: The single breast pump (SCF392/62) includes one set and the double breast pump (SCF394/62) includes two sets of the breast pump kit and Philips Avent Natural bottle shown in Figures B and C.

Product description (Fig. A)

- | | |
|--------------------|--|
| A1 Motor unit | A5 Double breast pump |
| A2 Breast pump kit | * The adapter supplied may vary for different countries. |
| A3 Bottle | |
| A4 Adapter* | |

Breast pump kit (Fig. B)

- | | |
|--------------------------|----------------|
| B1 Silicone tube and cap | B4 Cover |
| B2 Silicone diaphragm | B5 Flange |
| B3 Pump body | B6 White valve |

Philips Avent Natural bottle (Fig. C)

- | | |
|-------------------|-----------------|
| C1 Bottle cap | C4 Sealing disc |
| C2 Screw ring | C5 Bottle |
| C3 Feeding nipple | |

Motor unit part description (Fig. D)

- | | |
|---|--|
| D1 USB power inlet | D5 'Level up' button |
| D2 On/off button with pause/play function | D6 LED display |
| D3 Mode selection button | D7 Tubing port for single breast pump |
| D4 'Level down' button | D8 Tubing ports for double breast pump |

Accessories

There are multiple versions of this breast pump package, which all come with a different combination of the following Philips Avent accessories:

- | | |
|--------------------|------------------------------|
| - Breast pads | - Storage cups with adapter* |
| - Travel bag | - Breast milk storage bags* |
| - Insulation pouch | - Natural bottles |
| - Cleaning brush | - Breast pump belt |

* For these accessories, a separate user manual is provided.

Indications for use

The Philips Avent single/double electric breast pump Advanced is intended to express and collect milk from the breast of a lactating woman, and to alleviate engorgement of the breast. The device is intended for a single user.

IMPORTANT SAFETY INSTRUCTIONS

Read this user manual carefully before you use the breast pump and save it for future reference.

This user manual can also be found online at: www.philips.com/support.

This product is operated by the lactating woman herself.

You can use this breast pump in an everyday environment.

CONTRAINDICATION

- **Never use the breast pump while you are pregnant, as pumping can induce labor.**

WARNINGS

Warnings to avoid choking, strangulation and injury:

- **Do not allow children or pets to play with the motor unit, the adapter, packaging materials or accessories.**

- **Always unplug the breast pump after use and when fully charged. Store the breast pump when not in use.**

- This breast pump is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge. Such persons can only use this breast pump if they are supervised by or have received proper instruction concerning use of the appliance by a person responsible for their safety.

Warning to prevent transmitting infectious diseases:

- If you are a mother who is infected with Hepatitis B, Hepatitis C or Human Immunodeficiency Virus (HIV), expressing breast milk will not reduce or remove the risk of transmitting the virus to your baby through your breast milk.

Warnings to avoid electric shock:

- Inspect the breast pump, including the adapter, for signs of damage before each use. Do not use the breast pump if the adapter or plug is damaged, if it does not work properly or if it has been dropped or submerged in water.
- If the power adapter outer casings become loose or are separating, contact the device manufacturer immediately. Separation of the casing can result in electric shock.
- Only use the adapter (S009AHz050yyyy) provided with this breast pump.

- Always ensure that it is easy to disconnect the device from the mains in case of malfunction. The adapter is considered the disconnection device.

Warnings to avoid contamination and to ensure hygiene:

- For hygienic reasons, the breast pump is only intended for repeated use by a single user.
- Clean and disinfect all parts that come into contact with the breast and breast milk before first use and after each use.
- Do not use antibacterial or abrasive cleaning agents to clean the breast pump parts as this may cause damage.
- Wash your hands thoroughly with soap and water before you touch breast pump parts and breasts to prevent contamination. Avoid touching the inside of containers or lids.

Warnings to avoid breast and nipple problems including pain:

- Do not use the breast pump if the silicone diaphragm appears to be damaged or broken. See 'Ordering accessories and spare parts' for information on how to obtain replacement parts.
- Only use accessories and parts recommended by Philips Avent.
- No modification of the breast pump is allowed. If you do so, your warranty becomes invalid.
- There are no user serviceable parts inside the motor unit of the breast pump.
- Never use the breast pump while you are sleepy or drowsy to avoid lack of attention during use.
- Always switch off the breast pump before you remove the pump body from your breast, in order to release the vacuum.
- If the vacuum is uncomfortable or causes pain, switch off and remove the breast pump from your breast.
- Do not continue pumping for more than five minutes at a time if you do not succeed in expressing any milk. Try to express at another time during the day.
- If the process becomes very uncomfortable or painful, stop using the breast pump and consult your healthcare professional.

Warnings to prevent damage and malfunction of the breast pump:

- Portable radio frequency (RF) communications equipment (including mobile phones and peripherals such as antenna cables and external antennas) should not be used closer than 30 cm (12 in) to any part of the breast pump, including adapter. This could have a negative effect on the performance of the breast pump.
- Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation.

Warnings to avoid increased electromagnetic emissions or decreased electromagnetic immunity of this device, resulting in improper operation:

- Do not use any other adapter or charging cable than supplied with the device.
- Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Battery safety instructions

- Only use the product for its intended purpose. Read this important information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.
- To charge the battery, only use the detachable adapter (S009AHZ050yyyy) provided with the product. Do not use an extension cord.

- Charge, use and store the product at the temperature ranges indicated in the technical specifications table in this user manual.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.

CAUTIONS

Cautions to prevent damage to and malfunctioning of the breast pump:

- Prevent the adapter and motor unit from coming into contact with water.
- Never put the motor unit or the adapter in water or in a dishwasher, as this causes permanent damage to these parts.
- Keep the adapter and the silicone tubes away from heated surfaces to avoid overheating and deformation of these parts.

SAVE THESE INSTRUCTIONS

Before first use


Disassemble, clean and disinfect all parts that come into contact with the breast and breast milk before you use the breast pump for the first time and after each use.

Fully charge the device before you use it for the first time or after a long period of non-use.

Cleaning and disinfecting


Overview

Clean and disinfect the parts that come into contact with the breast and breast milk as described below:

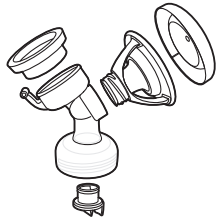	When	How
	Clean and disinfect before first use and after every use.*	Disassemble all parts and then clean them as described in chapter 'Cleaning' and disinfect them as described in chapter 'Disinfecting'.
		

* ⓘ Note: Disinfecting is especially important if your baby is less than 2 months old, was born prematurely, or has a weakened immune system due to illness or medical treatment. To ensure hygiene, perform the disinfecting step at least once a day.


The following parts do not come into contact with the breast and breast milk; clean them as described below:

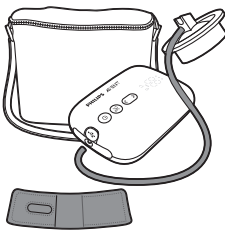	When	How
	Clean when needed.	Wipe with a clean, damp cloth with water and mild detergent. Clean the breast pump belt by hand wash with warm water and mild detergent or machine wash by selecting a cold, gentle cycle in a mesh laundry bag. Wipe the motor unit with a clean, soft cloth moistened with low-level disinfectant such as ethyl or isopropyl alcohol (70-90 %). Do not disinfect the tube.

⚠ **Caution:** Inspect the tubing. If you find condensation in the air tube, wash and dry it immediately or replace the tubing. If you see breast milk inside the tubing or tubing

connector, do not attempt to wash, clean or disinfect the tubing or tubing connector. Instead, contact the manufacturer.

Step 1: Disassembling

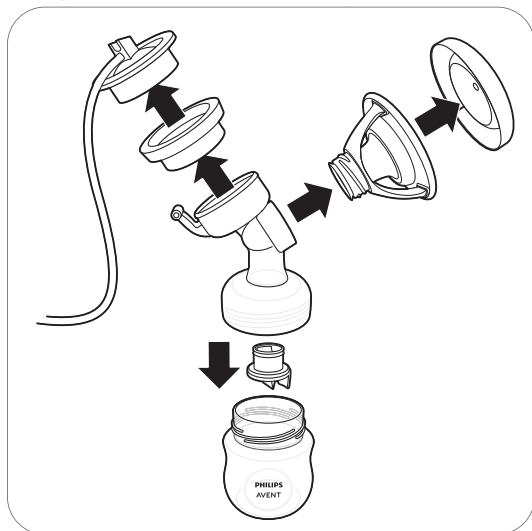Disassemble the breast pump kit, bottle and storage cup completely. Also remove the white valve from the pump body (Fig. B3).

Step 2: Cleaning

The parts that come into contact with the breast and breast milk can be cleaned manually or in the dishwasher.

⚠ Warning: Do not use antibacterial or abrasive cleaning agents to clean the breast pump parts as this may cause damage.

⚠ Caution: Never put the motor unit or the adapter in water or in a dishwasher, as this causes permanent damage to these parts.

⚠ Caution: Be careful when you remove and clean the white valve. If it gets damaged, your breast pump does not function properly. To remove the white valve, gently pull at the ribbed tab on the side of the valve. To clean the white valve, rub it gently between your fingers in warm water with some dishwashing liquid. Do not insert objects into the white valve, as this may cause damage.

Step 2A: Manual cleaning

Supplies needed:

- Mild dishwashing liquid
- Drinking-quality water
- Soft, clean brush used only for infant feeding items
- Clean tea towel or drying rack
- Clean bowl used only for infant feeding items


1. Rinse all parts under a running tap with lukewarm water. Do not place the parts in a sink to rinse.

2. Soak all parts for 5 minutes in warm water with some mild dishwashing liquid.

3. Clean all parts with a cleaning brush. Clean the feeding nipple with a brush as carefully as possible, to avoid damage.

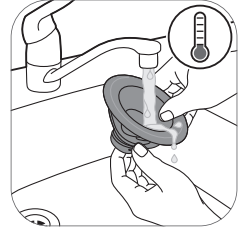4. Thoroughly rinse all parts under a running tap with cold clear water by holding them. Clean the bowl and the brush. Rinse them well and allow them to air-dry after each use. Wash them by hand or in a dishwasher at least every few days.

5. Leave all parts to air-dry on a clean, unused tea towel or drying rack in an area protected from dirt and dust. Do not use a tea towel to rub or to pat the items dry.

Step 2B: Cleaning in the dishwasher

Supplies needed:

- Mild dishwashing liquid or a dishwashing tablet
- Drinking-quality water

Note: Food colorings may discolor parts.

1. Place all parts on the top rack of the dishwasher.

2. Put detergent or dishwashing liquid, or a tablet in the machine. Run a standard program using hot water and a heated drying cycle (or sanitizing setting).

3. Remove the parts from dishwasher with clean hands. If items are not completely dry, place items on a clean, unused tea towel or paper towel to air-dry thoroughly before storing. Do not use a tea towel to rub or pat the items dry.

Step 3: Disinfecting


Supplies needed:

- A household pot
- Drinking-quality water

Caution: During disinfection with boiling water, prevent the bottle or other parts from touching the sides of the pot. This can cause irreversible product deformation or damage that Philips cannot be held liable for.

1. Fill a household pot with enough water to cover all parts and place all parts in the pot. Let the water boil for 5 minutes. Make sure the parts do not touch the sides of the pot.

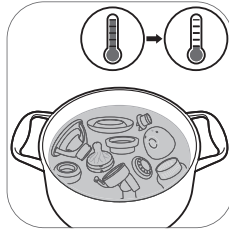2. Allow the water to cool down.

3. Gently remove the parts from the water. Place the parts on a clean surface/tea towel or drying rack to air-dry.

4. Store the dry items safely until needed. Ensure the clean pump parts, bottle brushes and wash bowls have air-dried thoroughly before storing. Items must be completely dry to help prevent germs and mold from growing. Store the dry items in a clean, protected area.


Using the breast pump

Flange size


The Philips Avent single/double electric breast pump Advanced has a soft, active flange. There is only one flange size. It gently stimulates your nipple to trigger milk flow. The flange is made of flexible silicone, fitting nipple sizes up to 30 mm/1.18 in.

Charging the breast pump


Fully charge the battery before you use the breast pump for the first time and when the battery is low. Fully charging the battery takes up to 4 hours.

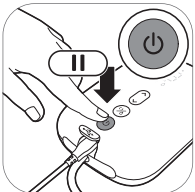
1. Switch off the device.

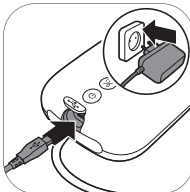
2. Insert the adapter into the power outlet and insert the small plug at the other end into the motor unit.

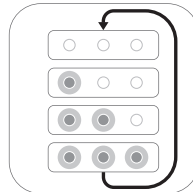3. The battery lights light up one by one and then all go out: the battery is charging.

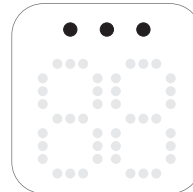
4. All 3 lights are steady on: the battery is fully charged.

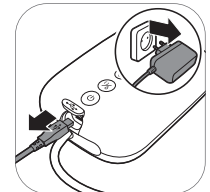5. Remove the adapter from the power outlet and pull the small plug out of the motor unit.


Battery indicator light description

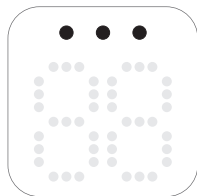
All lights are steady on. The battery is fully charged. Charged battery operates approximately 3 pumping sessions.*


Two lights are steady on. The battery is charged sufficiently for maximum two sessions.*

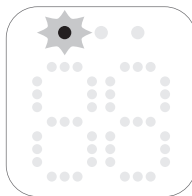First light flashes slowly. The battery is charged sufficiently for maximum one session.*

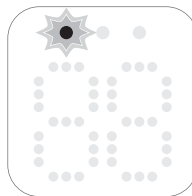First light flashes quickly. The battery is very low, the device will soon stop operating, charging is needed.

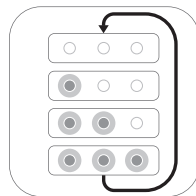Battery lights light up one by one and go out all together. The battery is charging.

* One session is the equivalent of 15 minutes (double electric breast pump) or 30 minutes (single electric breast pump).

Assembling the breast pump

⚠ Warning: Be careful, when you have disinfected the breast pump parts by boiling them, they can be very hot. To prevent burns, only start assembling the breast pump when disinfected parts have cooled down.

⚠ Warning: Wash your hands thoroughly with soap and water before you touch breast pump parts and breasts to prevent contamination. Avoid touching the inside of containers or lids.

⚠ Caution: Before each use, visually inspect the individual components for cracks, chips, tears, discoloration or deterioration. In the event that damage to the device is observed, discontinue use until the parts have been replaced.

⚠ Caution: Inspect the tubing. If you find condensation in the air tube, wash and dry it immediately or replace the tubing. If you see breast milk inside the tubing or tubing connector, do not attempt to wash, clean or disinfect the tubing or tubing connector. Instead, contact the manufacturer.

Make sure you have cleaned and disinfected the parts of the breast pump that come into contact with the breast and breast milk.

Check the breast pump kit parts for wear or damage before use and replace them if necessary.


ⓘ Note: Correct placement of the white valve, flange, silicone diaphragm and silicone tube is essential for the breast pump to form a proper vacuum.

1. Wash your hands thoroughly with soap and water.

2. Insert the white valve into the pump body from underneath. Push the white valve in the pump body as far as possible.

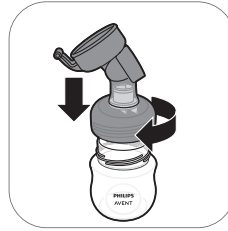
3. Screw the pump body clockwise onto the bottle until it is securely fixed.


4. Screw the flange onto the pump body.


5. Place the silicone diaphragm into the pump body.

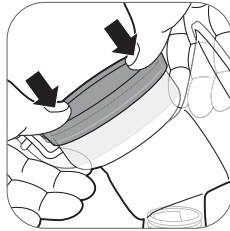
6. Make sure that the silicone diaphragm fits securely around the rim by pressing it down with your thumbs.

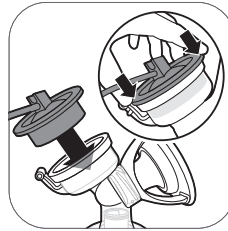
7. Attach the silicone tube and cap onto the diaphragm. Push down the cap until it is securely in place.

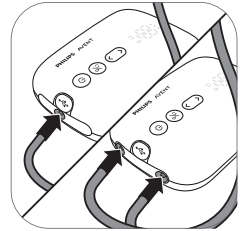8. Connect the tube(s) to the motor unit.

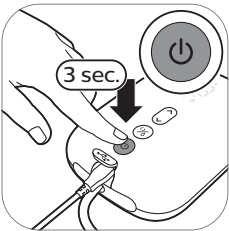
9A. When the breast pump is charged, you can use it cordless. Press and hold the on/off button until the device switches on.

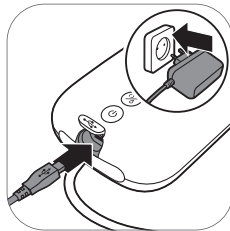9B. You can also use the breast pump connected to the power outlet. If the battery is low, insert the adapter into the power outlet and insert the small plug at the other end into the motor unit.

10. The breast pump is now ready for use.

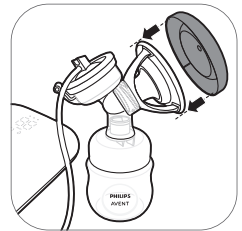① Note: You can place the cover over the flange to keep the breast pump clean when you take it with you and while you prepare for expressing.

Breast pump belt

The breast pump belt is intended to carry the motor unit of the electric breast pump while it is used to express and collect breast milk. The breast pump belt allows the electric breast pump to be used hands-free. The breast pump belt is intended to be worn over clothes.


Note: Pumping bra is not included in this product.

Using the breast pump belt


Charge the motor unit before using it in combination with the breast pump belt. Assemble and use the breast pump belt as follows:


1. Place the breast pump belt around your waist and close it with the hook and loop fastener.

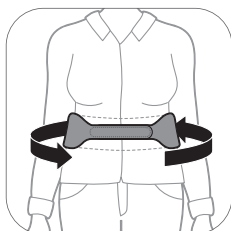
2. Turn and adjust the belt as necessary.

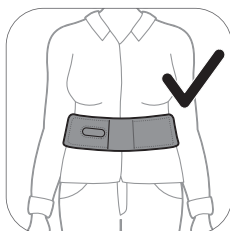
3A. Make sure the pockets of the breast pump belt are in the front.

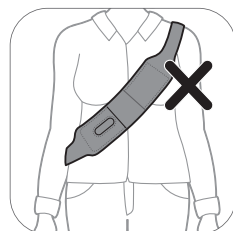
3B. Do not wear the breast pump belt across your body.

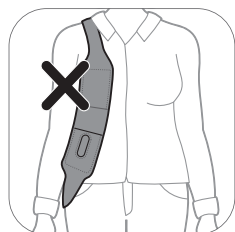
3C. Do not wear the breast pump belt across your shoulder.

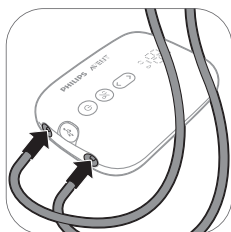
4. Connect the tube(s) to the motor unit.

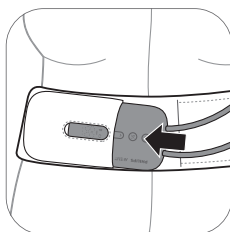5. Place the motor unit in the pocket of the breast pump belt with the buttons facing forwards.

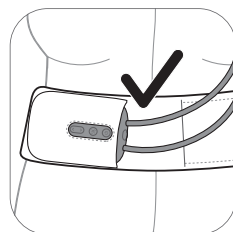
6. Make sure that the motor unit is placed in the pocket fully.

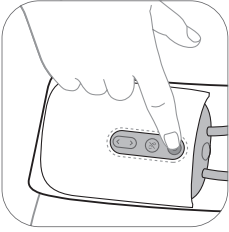7. You can now operate the motor unit from the belt while you are expressing.


⚠ Warning: If you experience skin irritation and/or skin redness when wearing the breast pump belt, stop using it immediately.

Motor unit part description

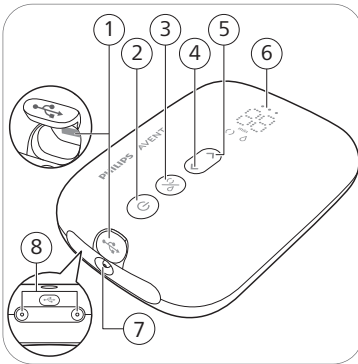- 1 USB power inlet
- 2 On/off button with pause/play function
- 3 Mode selection button
- 4 'Level down' button
- 5 'Level up' button
- 6 LED display
- 7 Tubing port for single breast pump
- 8 Tubing ports for double breast pump

LED display description

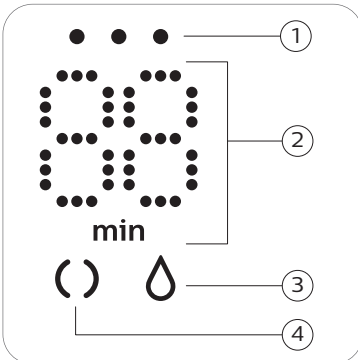- 1 Battery indicator lights
- 2 Timer (minutes) and vacuum level indication
- 3 Expression icon
- 4 Stimulation icon

Explanation of modes

The breast pump has two modes. Below you can find an explanation of these modes.

Mode icons	Mode	Explanation	No. of suction levels
	Stimulation mode 	Mode to stimulate the breast to start up the milk flow.	8 suction levels
	Expression mode 	Mode for efficient milk removal after milk has started flowing. After 90 sec. of stimulation, the device automatically switches to the expression mode.	16 suction levels
			The number on the display shows the running time in minutes.
			When you change the suction level, the display shows the selected suction level for a few seconds and then shows the running time again.

Pumping instructions

Warning: Always switch off the breast pump before you remove the pump body from your breast, in order to release the vacuum.

Warning: Do not continue pumping for more than five minutes at a time if you do not succeed in expressing any milk. Try to express at another time during the day.

Warning: If the process becomes very uncomfortable or painful, stop using the breast pump and consult your healthcare professional.

Caution: If discomfort is felt, then suction can be broken by inserting a finger between the breast and the flange.

Note: If you regularly express more than 125 ml/4 fl oz per session, you can purchase and use a 260 ml/9 fl oz Philips Avent Natural bottle to prevent overfilling and spillage.


Note: The breast pump can be used cordless if it is sufficiently charged, but you can also use it connected to the power outlet.

1. Wash your hands thoroughly with soap and water and make sure that your breast is clean.

2. Place the breast pump kit onto your breast. Make sure that you center your nipple.

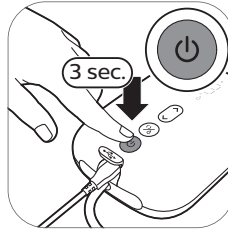
3. Press and hold the on/off button until the device switches on. The indicator lights go on to indicate that the device is in stimulation mode.

4. The breast pump starts in the stimulation mode and slowly increases its vacuum to the last used stimulation suction level (or level 5 when used for the first time).


5. The display shows the running time.

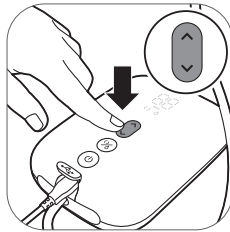
6. Use the 'level up' and 'level down' buttons to change the suction level to your preference.


7. After 90 seconds, the breast pump will switch smoothly to expression mode and slowly increase its vacuum to the last used expression suction level (or level 11 when used for the first time).

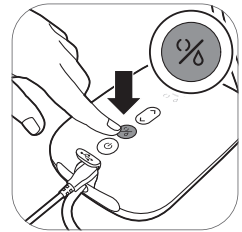
8. If your milk starts flowing earlier or if you want to switch from expression mode, you can switch modes with the mode selection button.

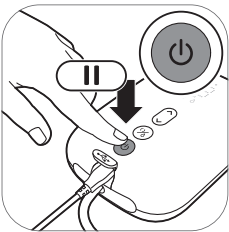9. Briefly press the on/off button if you want to pause. To continue, briefly press the on/off button again.

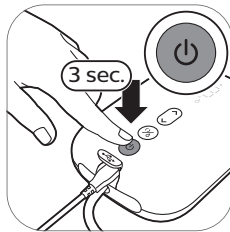10. Press and hold the on/off button to switch off.

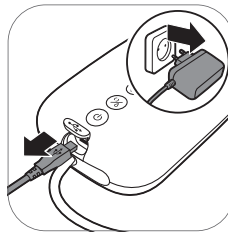11. Unplug the adapter from the power outlet and pull the small plug out of the motor unit.


- ① Tip: For optimal expression, choose the highest setting that still feels comfortable. This setting may differ across different sessions.
- ① Tip: You can use the double breast pump as a single breast pump by connecting only the tube of one pump kit to the motor unit.
- ① Note: The breast pump will switch off automatically after 30 minutes in order to prevent overheating.

Using the breast pump when you have symptoms of engorgement


Breast engorgement means your breasts are overfilled with milk. This is a painful and unpleasant condition. Frequent extraction of milk and completely emptying the breasts of milk can reduce engorgement. When you have symptoms of engorgement, make sure you use a lower vacuum level to minimize stretching of breast tissue.

- ① Note: When symptoms of engorgement are present, make sure you empty your breasts completely during each expression session.


After use

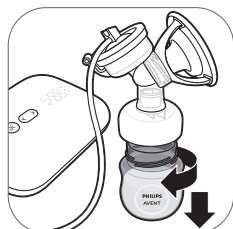1. Unscrew the bottle.

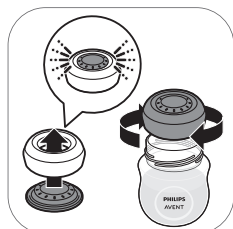2A. To store the breast milk: close the bottle with the sealing disc and screw ring.

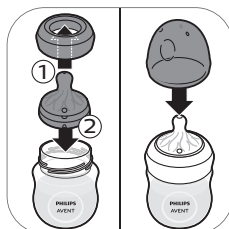
2B. To prepare the bottle for feeding your baby: place a nipple in a screw ring and screw the assembled screw ring onto the bottle. Seal the nipple with the bottle cap (see 'Assembling the bottle').

3. Remove the silicone tube and cap from the breast pump kit.

4. Disassemble, clean and disinfect all parts that have been in contact with the breast and breast milk (see 'Cleaning and disinfecting').

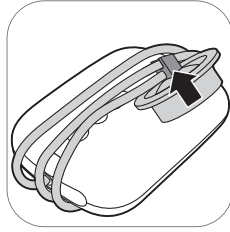5. Clean all parts that do not come into contact with the breast and breast milk (see 'Cleaning and disinfecting'). For easy storage, wrap the silicone tube around the motor unit and clip the cap onto the tube.

⚠ Warning: Refrigerate or freeze expressed breast milk immediately, label with date and time, or keep it at a room temperature of up to 77 °F (25 °C) for a maximum of 4 hours before you feed your baby.

ℹ Note: For more information about the breast pump and tips for expressing milk, visit www.philips.com/support.

Instructions for storing the breast milk

⚠ Warning: Never refreeze thawed breast milk to avoid a decrease in milk quality.

⚠ Warning: Never add fresh breast milk to frozen breast milk to avoid a decrease in milk quality and unintended thawing of the frozen milk.

- Write the date on the storage container. Include your child's name if you are giving the milk to a child care provider.
- Gently swirl the container to mix the cream part of the breast milk that may rise to the top back into the rest of the milk. Do not shake the milk. This can make some of the milk's valuable parts break down.
- Refrigerate or chill milk right after it is pumped, if possible. You can put it in the refrigerator, place it in a cooler or insulated cooler pack, or freeze it in small (2 to 4 ounce) batches for later feedings. Pumped milk at room temperature is fine without being refrigerated for up to 4 hours after pumping.

Below you find guidelines for storing breast milk:

Type of breast milk	Storage location and temperatures		
	Countertop 77 °F (25 °C) or colder (room temperature)	Refrigerator 40 °F (4 °C)	Freezer 0 °F (-18 °C) or colder
Freshly expressed or pumped	Up to 4 hours	Up to 4 days	Within 6 months is best Up to 12 months is acceptable
Thawed, previously frozen	1–2 hours	Up to 1 day (24 hours)	NEVER refreeze human milk after it has been thawed
Leftover from a feeding (baby did not finish the bottle)	Use within 2 hours after the baby is finished feeding		

Source: https://www.cdc.gov/breastfeeding/recommendations/handling_breastmilk.htm

Feeding breast milk with the bottle

Choosing the right nipple for your baby

Finding the right nipple is important:

If your newborn baby is consistently not taking enough milk throughout the feeding sessions or has complications in getting milk, switch to a nipple with a higher flow rate.

If persistent feeding issues occur, consult a healthcare professional.

Philips Avent Natural Response nipples are available with different flow rates. The breast pump comes with nipple number 2 (Slow). Use a higher flow rate if your baby falls asleep during feeding, gets frustrated, or plays with the nipple instead of drinking. Use a lower flow rate if your baby is gulping the milk down or if milk drips out of its mouth during feeding.

Philips Avent Natural feeding nipples are clearly numbered on the side, to indicate flow rate. For more information, see www.philips.com/avent.

For your child's safety and health

WARNING!

- Continuous and prolonged sucking of fluids will cause tooth decay.
- Always check food temperature before feeding.
- Throw away at the first signs of damage or weakness.
- Keep components not in use out of the reach of children.
- Never attach to cords, ribbons, laces or loose parts of clothing. The child can be strangled.
- Never use feeding nipples as a pacifier, to prevent choking hazard.
- Always use this product with adult supervision.
- Keep the bottle cap away from children to avoid suffocation.
- Do not allow children to play with small parts or walk/run while using bottles.
- Inspect all parts before use and pull the feeding nipple in all directions. Throw away any part at the first signs of damage or weakness.
- Always discard any breast milk that is left over at the end of a feed.

Caution: To avoid damage to the bottle


- Do not place in a heated oven, plastic can melt.
- Plastic material properties may be affected by disinfection and high temperatures. This can affect the fit of the bottle cap.
- Do not leave a feeding nipple in direct sunlight or heat, or leave in disinfectant for longer than recommended because this may damage the product.

Before using the bottle

- Inspect the bottle and feeding nipple before each use and pull the feeding nipple in all directions to prevent a choking hazard. Throw away at the first signs of damage or weakness.
- Clean and disinfect the bottle before first use and after every use (see 'Cleaning and disinfecting').


Assembling the bottle

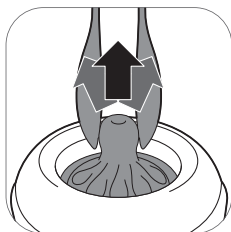1. The nipple is easier to assemble if you move it back-and-forth while pulling it upwards instead of pulling it up in a straight line.

2. Make sure you pull the nipple through until its lower part is aligned with the screw ring.

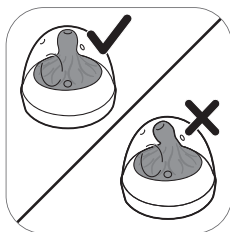3. Make sure you place the cap vertically onto the bottle so that the nipple sits upright.

4. Screw the screw ring with nipple and cap on the bottle.

5. To remove the cap, place your hand over the cap and your thumb in the dent of the cap.

Instructions for safe thawing and heating of breast milk

- If you use frozen breast milk, let it defrost completely before you heat it.
- Always thaw the oldest breast milk first. Remember first in, first out. Over time, the quality of breast milk can decrease.
- There are several ways to thaw your breast milk:
 - In the refrigerator overnight.
 - Set in a container of warm or lukewarm water.
 - Under lukewarm running water.
- Never thaw or heat breast milk in a microwave. Microwaving can destroy nutrients in breast milk and create hot spots, which can burn a baby's mouth.
- Use breast milk within 24 hours of thawing in the refrigerator (this means from the time it is no longer frozen or completely thawed, not from the time when you took it out of the freezer).
- Once breast milk is brought to room temperature or warmed after storing in the refrigerator or freezer, it should be used within 2 hours.
- Never refreeze breast milk once it has been thawed.

- Heat the bottle with defrosted or refrigerated breast milk in a bowl of hot water or in a bottle warmer. Remove the screw ring and sealing disc from the bottle.
- Always stir heated breast milk to ensure even heat distribution and test the temperature before serving.

Storing the bottles

- Store all parts in a dry, clean and covered container.
- Do not leave a feeding nipple in direct sunlight or heat, or leave in disinfectant for longer than recommended because this may damage the product.
- For hygiene reasons, we recommend replacing the feeding nipples after 3 months.

Compatibility

The Philips Avent single/double electric breast pump Advanced is compatible with all baby bottles in the Philips Avent Natural range and Philips Avent storage cups with adapter. We do not recommend using the Philips Avent Natural glass bottles with the breast pump.

Ordering accessories and spare parts

To buy accessories or spare parts, visit **www.philips.com/support** or go to your Philips dealer. You can also contact the Philips Consumer Care Center in your country.

When you order spare feeding nipples, make sure that you use a nipple with the correct flow rate for your baby (see 'Choosing the right nipple for your baby') and do not mix Philips Avent Anti-colic bottle parts and nipples with parts of the Philips Avent Natural bottles. They do not fit and could cause leakage or other issues.

Disposal

- This product contains a rechargeable lithium-ion battery which must be disposed of properly.

- Contact your local town or city officials for battery disposal information. In the US, you can also call 1-800-822-8837 or visit **www.call2recycle.org** for battery drop-off locations. In Canada, you can also call 1-888-224-9764 or visit **www.call2recycle.ca** for battery drop-off locations.
- For assistance, visit our website **www.philips.com/support** or call 1-800-542-8368 toll-free.
- Your product is designed and manufactured with high quality materials and components, which can be recycled and reused. For recycling information, please contact your local waste management facilities or visit: **www.recycle.philips.com**.

Full Two-Year Warranty

Philips North America LLC (USA) and Philips Electronics Ltd. (Canada) warrant this product against defects in materials or workmanship for a period of two years from the date of purchase, and agree to repair or replace any defective product without charge.

Important: This warranty does not cover damage resulting from accident, misuse or abuse, lack of reasonable care, the affixing of any attachment not provided with the product or loss of parts or subjecting the product to any but the specified voltage or batteries.

No responsibility is assumed for any special, incidental or consequential damages. In order to obtain warranty service, simply call toll-free 1-800-54-AVENT (1-800-542-8368).

Your appliance has been designed and developed with the greatest possible care and has an expected service life of 500 hours.

Support

For all product support, power consumption details and the EU Declaration of Conformity, please visit www.philips.com/support.

Troubleshooting

This chapter summarizes the most common problems you could encounter with the breast pump. If you are unable to solve the problem with the information below, visit www.philips.com/support for a list of frequently asked questions or contact the Consumer Care Center in your country.

Problem	Solution
I experience pain when I use the breast pump.	- Choose a suction level that is comfortable for you (see 'Pumping instructions'). - Check if the silicone diaphragm is not damaged (does not have small cracks, holes etc.). - Make sure that your nipple is placed in the center of the flange. - If the pain persists, stop using the breast pump and consult your healthcare professional.
The breast pump does not switch on while it is connected to power outlet.	- Only use the adapter supplied. - Make sure that the adapter has been inserted correctly into the power outlet and that the small plug has been inserted correctly into the motor unit. - Connect the adapter to another power outlet. Then press and hold the on/off button to check if the breast pump can be switched on. - If the problem persists, contact the Consumer Care Center via www.philips.com/support.
The breast pump does not switch on (indicator lights do not go on).	- Make sure that the battery is fully charged (see 'Charging the breast pump') or connect the breast pump to the power outlet.

Problem	Solution
The breast pump is not charging.	- Only use the adapter supplied. - Make sure that the adapter has been inserted correctly into the power outlet and that the small plug has been inserted correctly into the motor unit. - The breast pump is too cold or too hot for charging. Let the breast pump acclimatize to room temperature for 30 minutes. - If the problem persists, contact the Consumer Care Center via www.philips.com/support.
Charging the breast pump takes a very long time.	- Only use the adapter supplied. - Make sure that the adapter has been inserted correctly into the power outlet and that the small plug has been inserted correctly into the motor unit. - The breast pump is too cold or too hot. It is recommended to charge the breast pump at room temperature. - If the problem persists, contact the Consumer Care Center via www.philips.com/support.
I do not feel any suction/the suction level is too low.	- If possible, try a higher suction level. - Check if the silicone diaphragm is assembled correctly (see 'Assembling the breast pump'). - Check if the white valve is not damaged and is assembled correctly (see 'Assembling the breast pump'). - Check if the other parts of breast pump are intact and assembled correctly (see 'Assembling the breast pump'). - Make sure that the breast pump is positioned correctly on the breast to allow a seal between breast pump and breast. - If you still do not feel any suction, please contact the Consumer Care Center via www.philips.com/support.
The suction of the breast pump is too strong.	- If possible, try a lower suction level. When you express for the first time, start with the default suction level (level 5 in stimulation mode, level 11 in expression mode) and increase/decrease the level when needed. During different expression sessions, different levels may feel most comfortable for you. - Make sure that you only use the breast pump parts provided with this breast pump. - Check if the silicone diaphragm is not damaged (does not have small cracks, holes etc.). - If the problem persists, contact the Consumer Care Center via www.philips.com/support.
I express little/no milk when I use my breast pump.	- If you do not feel enough suction, please check Troubleshooting item 'I do not feel any suction/the suction level is too low'. - If you are using the Philips Avent breast pump for the first time, you may need some practice before you are able to express any milk. More expressing tips can be found on the Philips website www.philips.com/mybreastpump.

Problem	Solution
When I switch on the breast pump after storage, it does not start working immediately.	- Make sure that the battery is fully charged (see 'Charging the breast pump') or connect the breast pump to the power outlet. - Perhaps you stored the breast pump at a temperature below the minimum storage temperature or above the maximum storage temperature. Let the breast pump acclimatize in a room of 68 °F (20 °C) for 30 minutes to adjust its temperature to the operating range (41 °F to 104 °F/5 °C to 40 °C).
The display shows 'Er' and the device switches off automatically.	- If you are operating the breast pump on battery power, let it cool down for 30 minutes and try to switch it on again. - If you are operating the breast pump with the adapter connected to the power outlet, check if you are using the adapter that was supplied with the device. - If the problem persists, contact the Consumer Care Center via www.philips.com/support.

Undesirable side effects

When using the breast pump, the below undesirable side effects may occur. If you experience any of these symptoms, contact a healthcare professional or breastfeeding specialist.

Undesirable side effect	Description
Pain sensation	Pain felt on the breast or nipple, or felt as a response to the suction being above a comfortable level.
Sore nipples	Persistent pain in the nipples at the beginning of the pumping session, or that lasts throughout the entire pumping session, or pain between sessions.
Erythema (redness)	Redness of the skin of the breast.
Bruise, thrombus	A reddish-purple discoloration that does not fade when pressed. When a bruise fades, it may become green and brown.
Injured tissue on nipple (nipple trauma)	- Fissures or cracked nipples. - Skin tissue peeling off the nipple. Normally occurs in combination with cracked nipples and/or blisters. - Nipple tear. - Blister. Looks like small bubbles on the surface of the skin. - Bleeding. Cracked or torn nipples can lead to bleeding of the affected area.

Storage conditions

Keep the breast pump out of direct sunlight as prolonged exposure may cause discoloration. Store the breast pump and its accessories in a safe, clean and dry place.

If the breast pump has been stored in a hot or cold environment, place it in an environment with a temperature of 68 °F (20 °C) for 30 minutes to let it reach a temperature within the usage conditions (41 °F to 104 °F/5 °C to 40 °C) before you use it.

Fully charge the battery before storing your breast pump for an extended period of time.

Electromagnetic fields (EMF)

This Philips device complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Technical information

	Single electric breast pump	Double electric breast pump
Vacuum level	Stimulation: -60 to -200 mbar (-45 to -150 mmHg) Expression: -60 to -360 mbar (-45 to -270 mmHg) Adjustable in steps of 20 mbar Intermittent vacuum	
Cycle speed	Stimulation: 105-120 cycles/min, depending on vacuum level Expression: 53-85 cycles/min, depending on vacuum level	
Motor unit rated input	5 V direct current / 1.1 A	5 V direct current / 1.8 A
Mains adapter rated input	100-240 V  50/60 Hz 400 mA	
Mains adapter rated output	5 V direct current / 1.1 A	5 V direct current / 1.8 A
Adapter type number	S009AHz050yyyy The letters "yyyy" represent the output current from 0100 (1000 mA) to 0180 (1800 mA), increments in steps of 100 mA. The letter "z" denotes type of plug, it can be: U, V, B, S, C, A, D, K, T, I, E, F.	
Battery	Rechargeable battery: 3.6 V 2600 mAh (non-replaceable)	
Power port	Micro-USB	
Operating conditions	Temperature: 41 °F to 104 °F (5 °C to 40 °C) Relative humidity: 15 to 90 % (non-condensing)	
Storage and transportation conditions	Temperature: -4 °F to 113 °F (-20 °C to 45 °C) Relative humidity: up to 90 % (non-condensing)	
Operating pressure	700-1060 hPa (<3000 m/9842 ft. altitude)	
Net weight	Approximately 275 g	Approximately 355 g

	Single electric breast pump	Double electric breast pump
External dimensions	Motor unit: 145 mm x 95 mm x 45 mm (length x width x height)	
Device classification	Applicable classifications: Degree of protection against electric shock; - Class II ME equipment - Type BF applied part The breast pump and the adapter are specified as medical electrical system according to IEC 60601-1.	
Materials	ABS, silicone (motor unit) Polyester, nylon, polypropylene, rubber (breast pump belt) Polypropylene, PBT, silicone (other parts)	
Service life	500 hours	
Mode of operation	Continuous operation	Non-continuous operation: 55 min on, 180 min off at 72 °F (22 °C) (duty cycle 23 %) 20 min on, 90 min off at 104 °F (40 °C) (duty cycle 18 %)
Automatic power-off	After 30 minutes (after last user interaction)	
Ingress protection	IP22 (motor unit), IP20 (adapter)	

This breast pump contains a battery that is non-replaceable.

Electromagnetic compatibility - compliance information

This device has been tested and found to comply with IEC 60601-1-2:2014 Ed 4.0 for electromagnetic compatibility (EMC) as Class B according to CISPR 11:2009. Tests levels are listed in the accompanying tables. Only use this Philips Avent single/double electric breast pump Advanced in combination with the adapter supplied.

Emissions	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The device does not use RF energy for its intended function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The device is suitable for use in locations in residential environments and in establishments directly connected to a low voltage power supply network which supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Complies Class A Category	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Complies	

Guidance and manufacturer's declaration – electromagnetic immunity

The Philips Avent single/double electric breast pump Advanced is intended for use in the electromagnetic environment specified below. The customer or the user of the breast pump should assure that it is used in such an environment. The essential performance functions are the accuracy of the vacuum pressure and the backflow protection to withstand milk leakage into the motor unit. The Philips Avent single/double electric breast pump Advanced was tested for immunity to electromagnetic disturbances and passed using the criteria of no change in modes, vacuum accuracy ($\pm 20\%$), backflow protection and operation of the breast pump.

Immunity test	Immunity test level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 15 kV air	Floors can be a source of static electricity. When the air is very dry ($< 20\%$ RH) the user is advised to prefer to use the device in a room with concrete or wooden floor (instead of using it in a room with a synthetic floor).
Radiated RF EM fields IEC 61000-4-3	10 V/m 80 MHz – 2.7 GHz 80 % AM at 1 kHz	
Proximity fields from RF wireless communications equipment IEC 61000-4-3	385 (27), 450 (28), 710 (9), 745 (9), 780 (9), 810 (28), 870 (28), 930 (28), 1720 (28), 1845 (28), 1970 (28), 2450 (28), 5240 (9), 5500 (9), 5785 (9) MHz (V/m)	The device is immune against any RF wireless communication device such as smart phone, Wi-Fi or Bluetooth device at a distance of 30 cm (11.8 in) or more. Avoid a closer distance as this may result in improper operation.
RATED power frequency magnetic fields IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	
Electrical fast transient/burst IEC 61000-4-4	± 2 kV 100 kHz repetition frequency	Mains power quality expected of a typical commercial or residential environment.
Surges IEC 61000-4-5	± 1 kV Line-to-line ± 2 kV Line-to-ground	Mains power quality expected of a typical commercial or residential environment.
Conducted disturbances induced by RF fields IEC 61000-4-6	3 Vrms, before modulation is applied 0.15 MHz – 80 MHz 6 Vrms, before modulation is applied, in ISM and amateur radio bands between 0.15 MHz and 80 MHz 80 % AM at 1 kHz	

Immunity test	Immunity test level	Electromagnetic environment - guidance
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles single phase at 0° 0 % UT; 250/300 cycle	Mains power quality expected of a typical commercial or residential environment. At power interruptions the breast pump will automatically switch-over to battery mode without interruption of the function.

Note: UT is the AC mains voltage prior to application of the test level. If in the unlikely event the Philips Avent single/double electric breast pump Advanced is locked inoperative in a function by disturbances from other devices beyond the levels stated in the table above, the device will shut-off completely and restart. In that case make sure the distance to the other device causing this disturbance is enlarged.

FCC compliance information

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and the receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Explanation of symbols

The warning signs and symbols are essential to ensure that you use this device safely and correctly and to protect you and others from injury. Below you find the meaning of the warning signs and symbols on the label and in the user manual.

Symbol	Explanation
	Indicates to follow the instructions for use. ¹
	Indicates important information such as warnings and cautions. ²

Symbol	Explanation
	Indicates usage tips, additional information or a note. ³
	Indicates the Global Trade Identification Number and the Lot Code of the breast pump. ⁴
	Indicates the manufacturer. ⁵
	Indicates the manufacturing date. ⁶
	Indicates that the part of the appliance that comes into physical contact with the user (also known as the applied part) is of type BF (Body Floating) according to IEC 60601-1. The applied part is the breast pump kit. ⁷
	Indicates 'Class II Equipment'. The adapter is double-insulated (Class II). ⁸
	Indicates manufacturer's catalog number of the breast pump. ⁹
	Indicates 'direct current'. ¹⁰
	Indicates 'alternating current'. ¹¹
IP22	Indicates the protection against ingress of solid foreign objects greater than 12.5 mm and against harmful effects due to the ingress of dripping water when tilted at 15 degrees. ¹²
IP20	Indicates the protection against ingress of solid foreign objects greater than 12.5 mm and no protection against ingress of water. ¹³
	Indicates the batch number of the medical device. ¹⁴
	Indicates that the product is certified by UL. ¹⁵
	Indicates USB. ¹⁶
	Indicates the relative humidity upper limits to which the device can be safely exposed: up to 90 %. ¹⁷
	Indicates the storage and transportation temperature limits to which the medical device can be safely exposed: -4 °F to 113 °F/-20 °C to 45 °C. ¹⁸

Symbol	Explanation
	Indicates the Forest Stewardship Council (FSC). ¹⁹
	Indicates that the breast pump must be kept out of direct sunlight. ²⁰
	Indicates that the breast pump must be kept dry. ²¹
	Indicates the operating pressure range: 700 to 1060 hPa. ²²
	Indicates that this device contains a rechargeable battery which must be disposed of properly. See chapter 'Disposal' for more information. ²³
	Indicates that the adapter is certified by Intertek. ²⁴
	Identifies the UDI carrier, including the AIDC and human readable information. ²⁵
	Indicates that the material used for the packaging is corrugated cardboard/fiberboard. ²⁶
	Indicates that the material used for the packaging is non-corrugated cardboard/fiberboard. ²⁷
	Indicates that the material used for the plastic bag is low density polyethylene. ²⁸
	Indicates indoor use only. ²⁹
	Indicates that the part is certified with six level energy efficiency rating. ³⁰
	Indicates that the product should not be used after the period accompanying the symbol. The symbol indicates the expiry date of x years from the manufacturing date (FAB). ³¹

1 IEC 60601-1, Medical electrical equipment – Part 1: General requirements for basic safety and essential performance, Table D.2 Symbol 10 Refer to instruction manual/booklet

2 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.4.4 Caution

3 This symbol is not contained in a standard that is recognized by the FDA

4 ISO/IEC 16022, Information technology - Automatic identification and data capture technologies - Data Matrix barcode symbology specification

- 5 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.1.1 Manufacturer
- 6 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.1.3 Date of manufacture
- 7 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.1 Symbol 20 Type BF applied parts
- 8 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.1 Symbol 9 Class II equipment
- 9 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.1.6 Catalogue number
- 10 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.1 Symbol 1 Direct current
- 11 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.1 Symbol 4 Alternating current
- 12 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.3 Symbol 2 IP Code IEC 60529:1989+A1:1999+A2:2013, Degrees of protection provided by enclosures (IP Code)
- 13 IEC 60601-1, Medical electrical equipment – Part 1: General Requirements for basic safety and essential performance, Table D.3 Symbol 2 IP Code IEC 60529:1989+A1:1999+A2:2013, Degrees of protection provided by enclosures (IP Code)
- 14 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.1.5 Batch Code
- 15 The 'UL Certified' mark indicates that the product is certified by UL to meet the applicable safety standards for the USA and Canada
- 16 This symbol is not contained in a standard that is recognized by the FDA
- 17 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.3.8 Humidity limitation
- 18 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.3.7 Temperature limit
- 19 The FSC trademarks enable consumers to choose products that support forest conservation, offer social benefits, and enable the market to provide an incentive for better forest management
- 20 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.3.2 Keep away from sunlight
- 21 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.3.4 Keep dry
- 22 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.3.9 Atmospheric pressure limitation
- 23 This symbol is not contained in a standard that is recognized by the FDA
- 24 The 'ETL Verified' mark indicates that the adapter is certified by Intertek to meet the applicable safety standards for the USA and Canada
- 25 ISO 15223-1, Medical devices – Symbols to be used with medical device labels, labelling and information to be supplied, Part 1: General requirements, Clause 5.7.10 Unique device identifier
- 26 This symbol is not contained in a standard that is recognized by the FDA
- 27 This symbol is not contained in a standard that is recognized by the FDA
- 28 This symbol is not contained in a standard that is recognized by the FDA
- 29 This symbol is not contained in a standard that is recognized by the FDA
- 30 This symbol is not contained in a standard that is recognized by the FDA
- 31 This symbol is not contained in a standard that is recognized by the FDA

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Trademarks are the property of Koninklijke Philips N.V. / Les marques de commerce sont la propriété de Koninklijke Philips N.V.
All rights reserved. / Tous droits réservés.
Philips Consumer Lifestyle B.V.
a division of / une division de Philips North America LLC
Philips Electronics Ltd. / Philips Électronique Ltée