Visit us

www.philips.com

There may be certain differences between the content of this manual and your product. Please refer to the actual product. We reserve the right to update the product at any time without prior notice.

© 2023 AquaShield

All rights reserved

Philips and the Philips Shield Emblem are registered trademarks of Koninklijke Philips N.V. and are used under license.

This product has been manufactured by and is sold under the responsibility of Guangdong Deerma Technology Company Limited and Guangdong Deerma Technology Company Limited is the warrantor in relation to this product.

Rev A MAR 23

PHILIPS

Leg Massager

Knee Massager

PPM5521

User Manual

User Manual

1 Readme

Symbol Description

People Suitable to Use

Precautions for Use

Precautions for Maintenance

2 Product Features

3 Components

4 How to Use

Preparations Before Use

Operation

Charging

Indicator Description

5 Maintenance

Cleaning

Storage

Disposal

6 Troubleshooting

7 Technical Parameters

8 Service Guide for Product Warranty

Warranty and Service

Warranty Card

9 Packing List

1 Readme

Please carefully read the operation manual before use and keep it properly.
The manual is available at www.philips.com

Symbol Description

To avoid personal injury or property damage, be sure to observe the contents below.

- The following symbols are used to illustrate what you should observe.

 Indicates prohibition

 Indicates the actions that need to be subject to instructions

People Suitable to Use

 The product should be used with caution by the following groups/ under the following circumstances

- People who have a pacemaker or other related medical electrical products implanted in their bodies
- People who have metal, silicone or plastic materials injected into their bodies as a result of bone fractures or plastic surgery
- People who have ever had skin trauma
- People whose body temperature is above 38 degrees Celsius

- People who are easy to suffer from numbness of body
- People who have difficulty with respirators
- People suffering from hypertension
- People having acute symptoms such as visceral diseases (gastritis, hepatitis, gastroenteritis)
- People who have disturbances of skin sensation
- People suffering from body paresthesia
- People suffering from lumbar disc herniation, spondylitis
- People suffering from deformed arthritis, arthritis, and gout
- People suffering from spinal fractures and swelling of bones

 The product is not allowed to be used by the following groups/ under the following circumstances

- People who are unable to express their will or who cannot operate the machine themselves to use the product
- This product is not suitable for people who have physical, sensory or mental capacity decline and lack experience and knowledge (including children)
- Please stop using immediately when you feel lower back pain, dizziness, vomiting, palpitation and other abnormal body reactions during use.
- This product is set for 15-minute timing at a time. Do not use it for no more than 30 minutes each time. After each use, stop the device for 15 minutes to ensure the extension of its service life.
- Do not use the product under the following circumstances:
 - When you are severely tired
 - Within one hour before and after meals
 - When you are drunk

Precautions for Use

- ❗ The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- ❗ This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- ❗ Children should be supervised to ensure that they do not play with the appliance.
- ❗ If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- ❗ This appliance contains batteries that are only replaceable by skilled persons.
- ❗ Must be only supplied at safety extra low voltage corresponding to the marking on the appliance.
- ❗ Check whether the product is intact before use. If any damage is found, regardless of its extent, stop using the product immediately.
- ❗ Remove the dust on the charger on a regular basis. The accumulated dust on the plug is likely to cause poor insulation due to humidity, and eventually lead to fire.
- ❗ If you cannot use it or feel abnormal, stop using it immediately
- ❗ In case of power cut, please unplug the USB charging cable.
- ❗ When the product is not in use, the USB charging cable must be unplugged.
- ⊘ Do not disassemble, modify or repair the product by yourself, nor unzip the product by yourself.
- ⊘ Do not use the device when your body is wet, otherwise it is easy to get an electric shock.
- ⊘ Do not pull the product with force, nor stab it with a knife or a sharp object.


Precautions for Maintenance

- ❗ Be sure to unplug the charging cable during maintenance.
- ❗ Please properly dispose of the packaging bag of the product and keep it away from the infants.
- ⊘ The product has built-in lithium batteries. Do not disassemble or modify the product by yourself, immerse the product (controller) in water or toss it into the fire, otherwise the product is likely to be damaged.
- ⊘ Do not damage, modify, bend or pull the power cord with force. Do not place heavy objects on the power cord or squeeze the power cord.
- ⊘ Do not use volatile oil, thinner, alcohol, strong alkaline detergent, bleach or other cleaning agents for maintenance.
- ⊘ Keep the product out of reach of children and infants. Do not let them touch or operate the product, otherwise it is likely to lead to burns, electric shocks or personal injuries.
- ❗ Make the device cool down before storing.

2 Product Features

Through the airbag and warm compress functions, and equipped with multiple optional massage modes and levels, this product is suitable for different people for knee relaxation.

3 Components

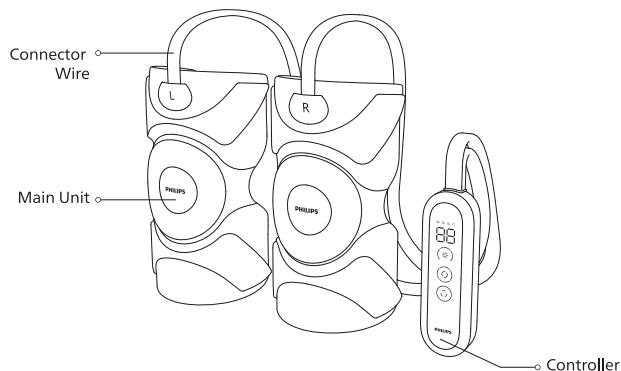
4 How to Use

Preparations Before Use

1. Check whether the device and its accessories are complete.
2. Wear the massager on the knees and adjust it to a comfortable state.

! As each person has different sensitivity to temperature, it is recommended to use the massager with clothes on for fear of scalding if choosing 45°C and above.

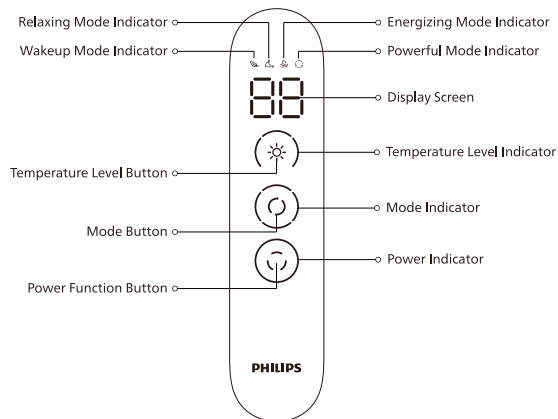Operation

Power Function Button

1. Press and hold the power function button for 2 seconds to turn on/off the massager. After power-on, the massager will enter the wakeup mode by default, enabling the warm compress at 40°C with a voice prompt "Massage On".
2. Double click the power function button to lock/unlock the screen with the voice prompts "Lock screen" and "Turn on screen". After the screen is locked, the status indicator is in the state of the orange breathing lamp, the screen displays the current warm compress temperature (it displays 00 when the warm compress is off), and the mode lamp turns on to indicate the current mode. Press and hold the power function button to turn off the massager.
3. The default timing of the product is 15 minutes. It is set to shut down automatically after 15 minutes of operation and enter the standby state with the voice prompt: "Turn off massage".

Mode Button

Press the mode button to switch the massage mode (Wakeup Mode - Relaxing Mode - Energizing Mode - Powerful Mode - Turn off air pressure) circularly. During adjustment, the mode level indicator turns blue progressively, the corresponding mode indicator lights on, and there is a voice prompt of the corresponding mode for each adjustment.

Temperature Level Button

Press the temperature level button to switch the temperature (40°C, 45°C, 50°C, heating off) circularly, and the display screen in the middle will display the current selected temperature. When adjusting, the corresponding temperature indicator is illuminated progressively (orange), and there is a voice prompt (40°C, 45°C, 50°C) of the corresponding temperature for each adjustment.

Low Battery Reminder

When the battery remaining capacity is less than 10%, the indicator lamp ring flashes orange for 5 seconds to prompt, and then the massager shuts down. At the same time, there is a voice prompt "The battery power is low, please charge it in time".

Charging State

When charging, the indicator lamp ring flashes blue like breathing until the massager is fully charged, then the blue lamp remains on.

Charging

1. Insert the charging plug into the charging port on the side of the battery holder;
2. Plug the charger plug into a 100-240V~, 50/60Hz power socket;
3. When charging, the indicator lamp ring flashes blue slowly until the massager is fully charged, then the blue lamp remains on;
4. It takes about 5 hours to fully charge the product.

Indicator Description

Indicator Type	State Indication	Voice Prompt
Power indicator	It turns orange by default after power-on. If the battery is low, it flashes quickly for 5 seconds to prompt and then shuts down. While charging, it flashes blue slowly, and remains blue after fully charged.	Massage On / Turn off massage / The battery power is low, please charge it in time (No voice prompt while charging)
Temperature level indicator	Display in orange for three levels, the first level light on by default after power-on. Press the button again, the second level light on.	40°C/45°C/50°C/ Turn off heating
Wakeup mode indicator	Enter this mode by default after power-on, the first green icon light on. Press the mode button to switch to the next mode.	Corresponding voice prompts of Wakeup Mode/ Relaxing Mode/ Energizing Mode/ Powerful Mode/ Turn off air pressure
Relaxing mode indicator	Press the mode button to enter relaxing mode, and the second icon turns blue.	
Energizing mode indicator	Press the mode button to enter energizing mode, and the third icon turns orange	

Indicator Type	State Indication	Voice Prompt
Powerful mode indicator	Press the mode button to enter powerful mode, and the fourth icon turns red	Corresponding voice prompts of Wakeup Mode/ Relaxing Mode/ Energizing Mode/ Powerful Mode/ Turn off air pressure
Dual 8 tube display	The set temperature is displayed in the heating state, and 00 is displayed after the heating is off	/

5 Maintenance

Cleaning

- 1.Remove the dust from the cloth cover of the product with a vacuum cleaner rather than by patting.
- 2.Wipe off the dirt on the surface with a clean rag.
- 3.For the heavy dirt, dampen a rag with the neutral detergent and then gently wipe it. Prevent the water from entering the product. Let it air dry.

! During maintenance, be sure to unplug the power cord, but do not unplug the cord with wet hands.

Storage

- Keep out of reach of children.
- Do not put it in places with direct sunlight, high temperature and humidity.
- Keep it in a dry and ventilated place.

Disposal

The lithium-ion battery used in this product is a valuable resource that can be recycled. Take out the lithium-ion battery when it is used up and discarded, and then hand it over to a professional institution that can recycle the rechargeable batteries for disposal.

1. Do not replace the battery by yourself during normal use.
2. The battery must be removed from the appliance before the it is discarded.
3. The appliance must be powered off before the battery is removed.
4. The battery must be disposed of in a safe manner.

Steps:

1. Pry open the cover with a tool and then unscrew with a screwdriver;
2. Take out the PCB screws on the lamp board support and main control board, and remove all the corresponding terminals;
3. Remove the screws on the cover with a screwdriver, and take out the cover;
4. Remove the battery from the lower case for safe disposal.

6 Maintenance

This chapter lists the most common problems to the product. If the problems you encounter are not included in the following tables, please call our service hotline or contact the authorized service center in your region.

Fault Description	Possible Cause	Solution
Startup Failure	Low battery	Connect the charger and use the product after charging
	The power function button is not on	Press and hold the power function button to power on
	Short circuit of the circuit board	Contact the service center
Stop working during use	The battery is low	Connect the charger and use the product after charging
	The 15-minute timed massage ends, and it shuts off.	Press and hold the power function button to power on
There is a sound coming from the product	The sound of air in□ation and bleeding of the internal air pump, which is a normal working state	No treatment required

Fault Description	Possible Cause	Solution
Product smokes	Short circuit of internal wires	Stop using it immediately and contact the service center

If the product has been left unused for a long time, be sure to confirm that its functions are normal before using it again		
Regular check	Are there any of the following faults?	In the event of the above faults, stop using the device immediately and cut off the power supply to prevent accidents. Please send it to our maintenance center.
	• Whether the power cord or the plug is hot • Whether the power cord is deeply scratched or deformed, or has a burnt smell • Whether there is abnormal sound or vibration during work • Whether the product can work properly after switch on • Whether there are other abnormalities or faults	

Note:
If the fault cannot be eliminated or the product still cannot operate satisfactorily, contact the after-sale service personnel, and do not disassemble or attempt to repair it by yourself.

7 Technical Parameters

Product Name	Philips Knee Massager
Product Model	PPM5521
Operating Voltage	3.7V ---
Operating Power	15W Max
Charging Voltage	5V ---
Charging Current	2A
Battery Capacity	5000mAh
Interface Type	Type-C
Charging Time	≤5hours
Net Weight	About 0.78kg
Product Dimension	1020×543×40mm

8 Service Guide for Product Warranty

Warranty and Service

If you need information or have any questions, please visit www.philips.com/-support or contact consumer service center in your country/region. If there is no consumer service center in your country/region, please contact your local distributor.

Within two years from the date of purchase, you will receive free warranty service for any damage caused by manufacturing process problems, or any damage to components under normal use conditions confirmed by our maintenance service personnel.

Warranty service excludes frequently replaced consumables, auxiliary equipment, shipping cost and door-to-door service. Please show proof of purchase to service personnel when servicing.

Malaysia
Consumer Care: +60-128868274
Email: philips.message@realx.tmy.net

Singapore
Consumer Care: +65-89420233
Email: philips.message@realx.tsg.net

Philippines
Consumer Care: +63-9604111314
Email: Philips.massager.ph@orohome.com

Thailand
Consumer Care: +66-20051888
Email: Philips_massager@orohome.co.th

Indonesia
Consumer Care: +021-39715012 / +021-39715013
Email: Philips_Massager.ID@orohomeindonesia.com

Vietnam
Consumer Care: +84-909800829
Email: Philips.massager.VN@orohome.vn

Warranty Card

Product Model		User Name	
Sales Agency		Phone Number	
Date of Purchase		Serial Number	
Special Technical Service Center		User Address	
Items of Maintenance	Date	Service Personnel	

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