

Spectral CT 7500

Table of Contents	
Drawing Nr	Drawing Title
01	Cover Sheet
02	Typica Site Layout
03	Typical Sections
04	3D View

Document Purpose

This document was prepared by Philips experts to demonstrate an example of correct equipment positioning of a specific Medical Equipment Solution offered by Philips within the space of a theoretical set of rooms. Information provided in this document was prepared for the convenience of facility designers and must be treated only as inspiration for technical design of medical facilities. Philips must be contacted to deliver support in preparing site-specific design that will facilitate the offered Medical Equipment Solution.

Responsibility

Philips does not take responsibility for any changes done to drawings by other parties, any damages resulting from use of information that is presented in this document, any damages caused by not involving Philips representatives in facility design coordination activities related to correct functional planning, room design and designing the layout of equipment that was configured by Philips for the Customer. Philips cannot accept responsibility for any damage caused by introducing facility design changes to rooms that are part of the Philips medical equipment solution that aren't discussed and aligned with Philips. Philips cannot accept responsibility for the partial use of Philips final issue drawings, however caused.

Scaling


All dimensions are indicated either in the Imperial (feet and fractional inches), Metric (millimeters) or both measurement systems. Do not scale from printed pdf files. Philips accepts no responsibility or liability for defective work due to scaling of these drawings.

Design Support

We encourage Customers and their Designers to contact Philips for medical facility design support related to Philips Medical Solutions. Philips will organize delivery of design and documentation support, that is related to the Philips offer and adapted to our Customer's site.

Legend

- Building Element
- Warning Light
- Elements installed on the Ceiling
- Elements installed on the Floor, Wall or unattached

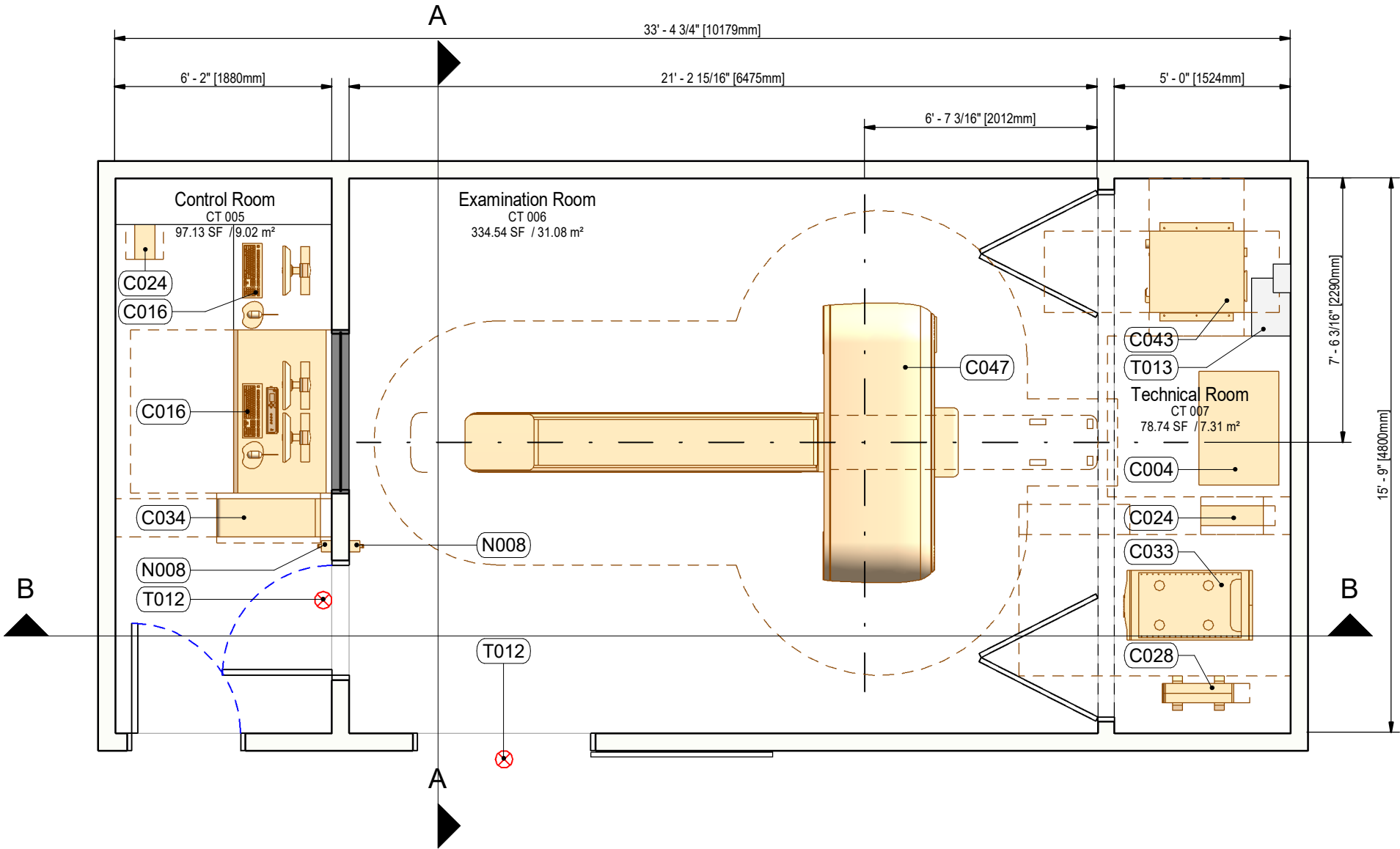

Code	Description
------	-------------

Control Room


C016	CT Operator Console
N008	Remote EPO 1 Switch
C024	Host Tower Cabinet
C016	CT Operator Console
C034	iMR Host Cabinet 12 Inch

Examination Room

N008	Remote EPO 1 Switch
T012	Warning Light
T012	Warning Light
C047	Spectral CT 7500

Technical Room

T013	Earth Reference and Dynamic Main Switch
C043	PDU Cabinet X4626
C028	Console UPS for Flex 3kVA
C004	Air Compressor Atlas-Copco SF1FF
C024	Host Tower Cabinet
C033	Hypersight Spectral Rack

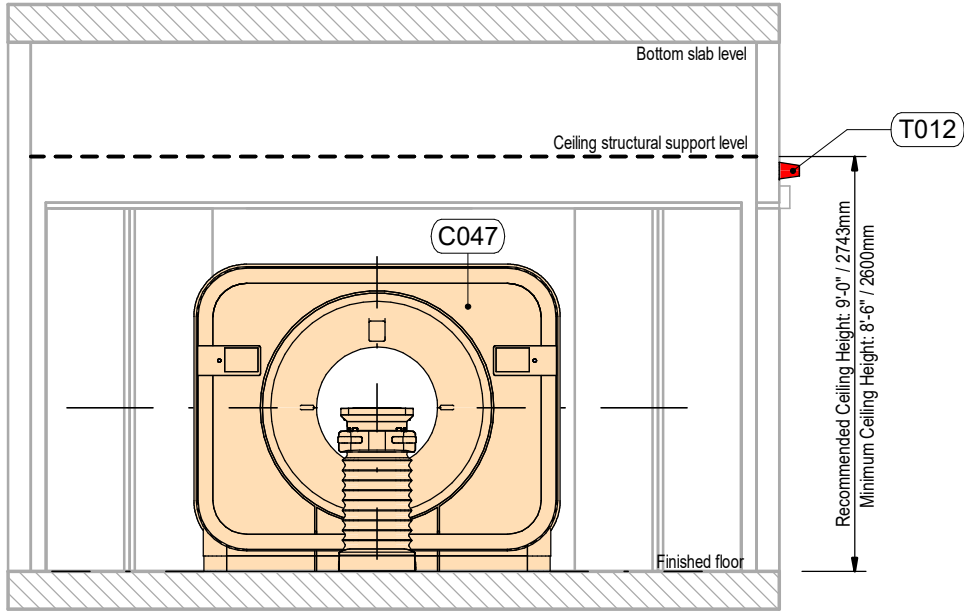


1: Section A
 □ 1 : 50

Code	Description
------	-------------


Control Room	
C016	CT Operator Console
N008	Remote EPO 1 Switch
C024	Host Tower Cabinet
C016	CT Operator Console
C034	iMR Host Cabinet 12 Inch

Examination Room	
N008	Remote EPO 1 Switch
T012	Warning Light
T012	Warning Light
C047	Spectral CT 7500

Technical Room	
T013	Earth Reference and Dynamic Main Switch
C043	PDU Cabinet X4626
C028	Console UPS for Flex 3kVA
C004	Air Compressor Atlas-Copco SF1FF
C024	Host Tower Cabinet
C033	Hypersight Spectral Rack

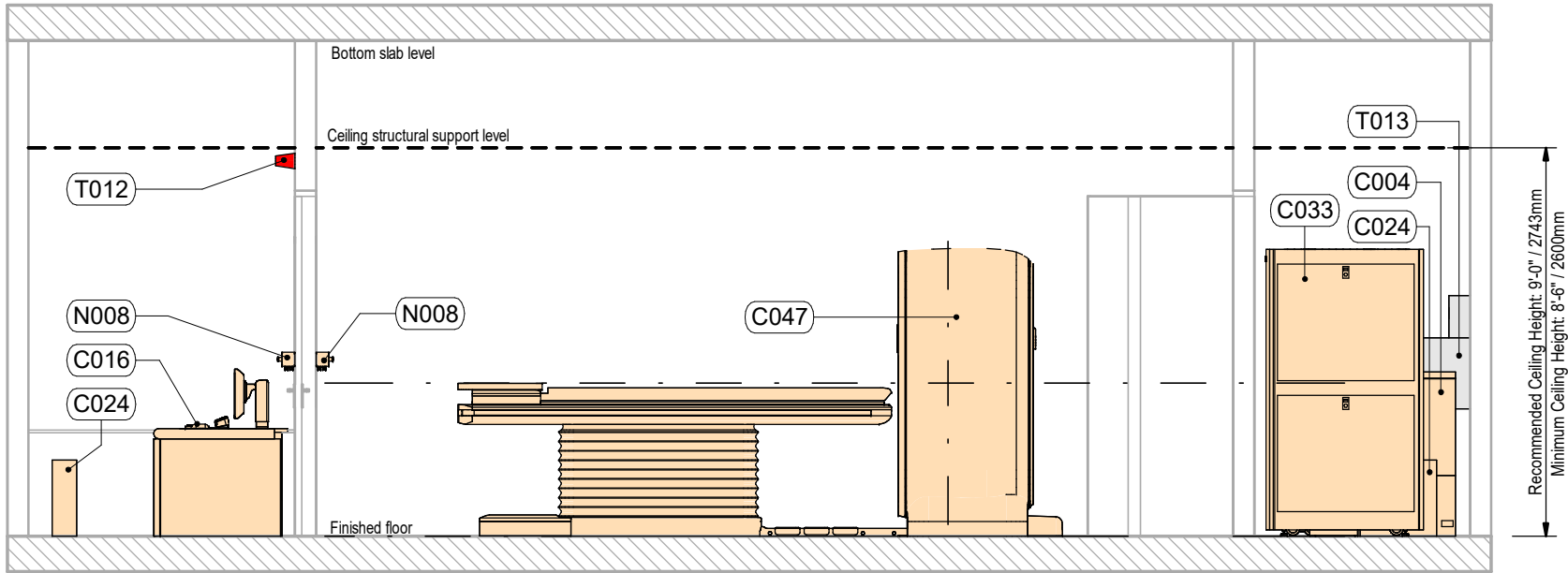

2: Section B
 □ 1 : 50

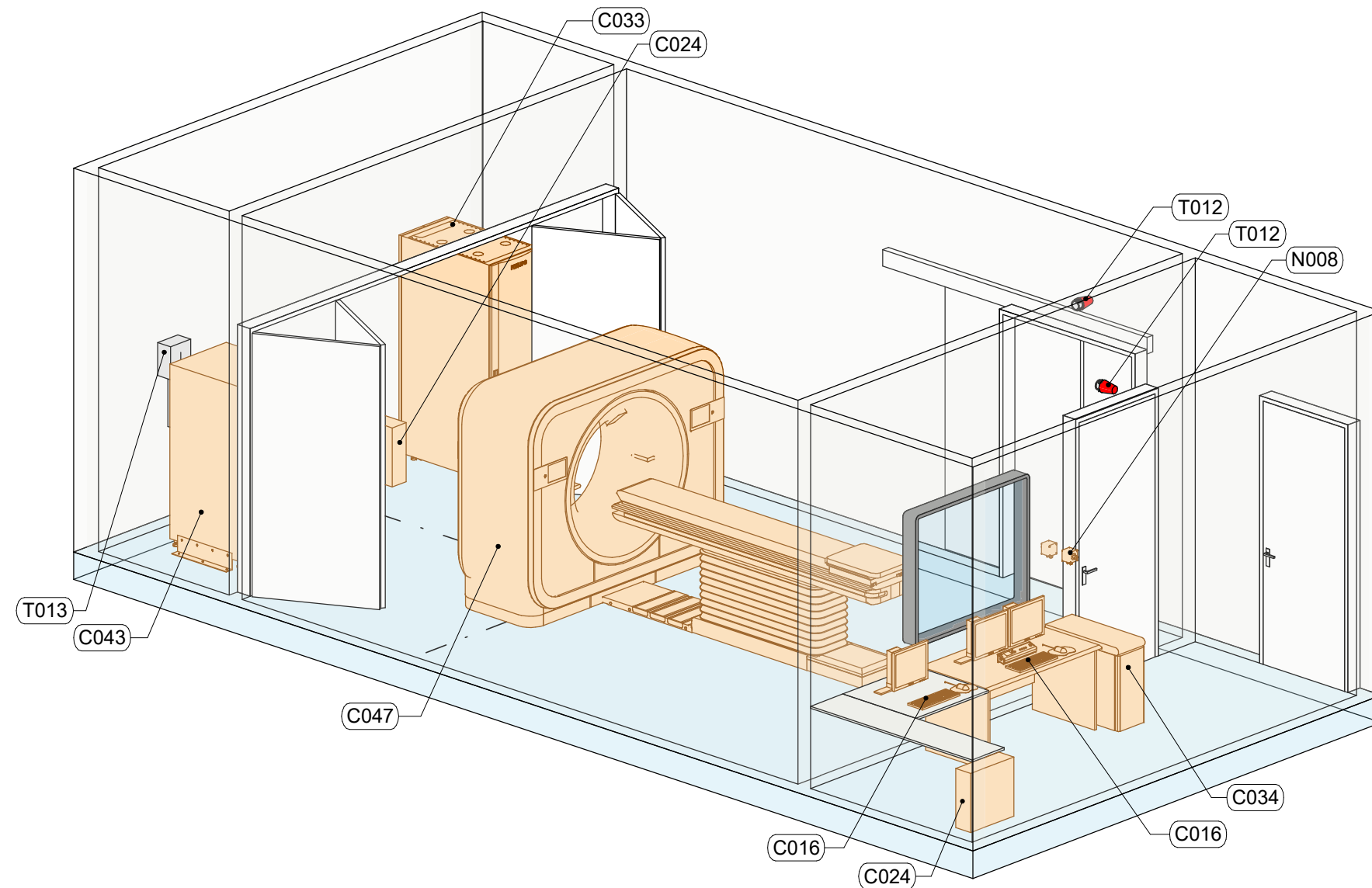

Code	Description
------	-------------

Control Room

C016	CT Operator Console
N008	Remote EPO 1 Switch
C024	Host Tower Cabinet
C016	CT Operator Console
C034	iMR Host Cabinet 12 Inch

Examination Room

N008	Remote EPO 1 Switch
T012	Warning Light
T012	Warning Light
C047	Spectral CT 7500

Technical Room

T013	Earth Reference and Dynamic Main Switch
C043	PDU Cabinet X4626
C028	Console UPS for Flex 3kVA
C004	Air Compressor Atlas-Copco SF1FF
C024	Host Tower Cabinet
C033	Hypersight Spectral Rack