
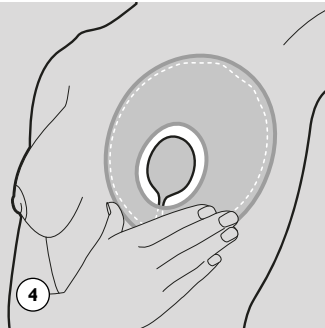
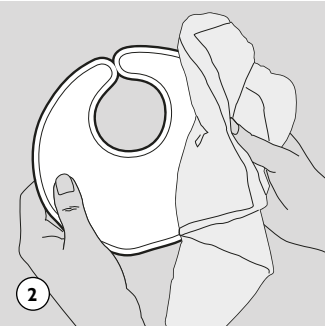
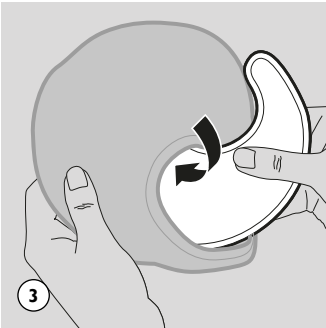
AVENT

A

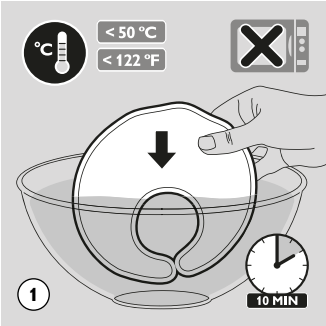GB	0844 33 80 489 national call rate	CH	0800 002 050 kostenfrei / 0848 000 292
IRL	0818 21 01 41 national call rate	KO	02-709-1200
NO	0800-023 0076 (mobiles vary)	NL	0800-023 0076
PL	22 397 15 06	PT	70 77 80 226 tarifa para chamadas nacionais
RU	8 800-200-08-80; 8 800-200-08-83 (para telefonos con varna)	SV	085 79 29 067 (ollika mobiltelefoner)
TR	444 7 445	SC	400 880 0008
TC	+86 (2) 8245 7399 / 0800-231-099	AR	(+71) 800AVENT
CA	1.800.54.AVENT	MX	1.800.54.AVENT
HU	06 1 700 81 52	IT	024 52 79 074 (variabile del cellulant)
BE	0800 80 190	FR	01 57 32 40 51 tarif appel local (tarif pour telephone portables variables)
ES	90 28 81 082 tarifa telefonica nacional	DE	0800 000 7520 kostenfrei / 0180 501 06 71 (mobiles vary)
EL	2111983029	LU	0800 26 550
AT	0800 180 016 kostenfrei / 820 401 132	AE	0800 180 016 kostenfrei / 820 401 132 (mobiles vary)
DA	35 44 41 30 pris ved indtæringsskald	SE	286 854 441-3
BG	+359 2 955 5115	GR	0818 21 01 41 national call rate
AE	0800 000 7520 kostenfrei / 0180 501 06 71	GB	0844 33 80 489 national call rate

GB
Note: Wash the sleeve before first use.
Before you apply the thermopad to the breast, make sure its temperature is acceptable by testing it on the back of your hand. Always use the thermopad together with the sleeve provided.

A Warm use: to stimulate milk flow before breastfeeding
• Immerse the thermopad in a bowl of hot water and let it warm to the optimum temperature (Fig. 1).
• Dry the thermopad with a cloth or towel (Fig. 2).
• Insert the thermopad into the sleeve supplied (Fig. 3).
• Apply the thermopad to the breast (Fig. 4).

B Cold use: to soothe sore or engorged breasts
• Place the thermopad in the refrigerator or freezer and let it cool to the optimum temperature for 24 hours (Fig. 5).
• Alternatively, immerse the thermopad in a bowl of ice water (Fig. 6).
• Insert the thermopad into the sleeve supplied (Fig. 7).
• Apply the thermopad to the breast (Fig. 8).

Cleaning and storage
Before you use the thermopad for the first time and after each use, gently clean the pad in warm water with some washing-up liquid. Do not use a brush, scourer or abrasive cleaning agent. • You can clean the sleeve by hand or in the washing machine at a low temperature. • Store the thermopad and sleeve in a cool, dry place out of direct sunlight and out of the reach of children.

Safety tips and troubleshooting
Do not heat the thermopad in a microwave. • Do not use the thermopad without the sleeve provided. Never put the frozen thermopad directly onto your skin. • Do not put the thermopad while you are feeding your baby or while expressing milk with a breast pump. • If the thermopad becomes damaged or appears to be leaking, stop using it and throw it away. • Do not swallow the contents of the thermopad. • After freezing, the gel inside the thermopad has a snow-like appearance. This is normal. The gel in the pad returns to its normal condition in a few hours. • If you have any concerns about using the thermopad, please consult your doctor or midwife.

Safety tips and troubleshooting
Do not heat the thermopad in a microwave. • Do not use the thermopad without the sleeve provided. Never put the frozen thermopad directly onto your skin. • Do not put the thermopad while you are feeding your baby or while expressing milk with a breast pump. • If the thermopad becomes damaged or appears to be leaking, stop using it and throw it away. • Do not swallow the contents of the thermopad. • After freezing, the gel inside the thermopad has a snow-like appearance. This is normal. The gel in the pad returns to its normal condition in a few hours. • If you have any concerns about using the thermopad, please consult your doctor or midwife.

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US
Note: Wash the sleeve before first use.
Before you apply the thermopad to the breast, make sure its temperature is acceptable by testing it on the back of your hand. Always use the thermopad together with the sleeve provided.

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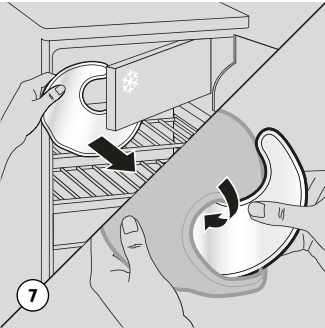
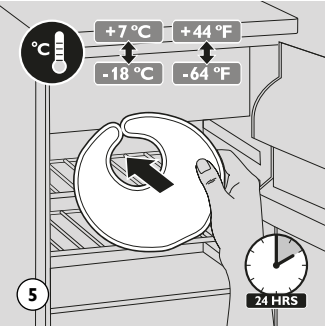
Safety tips and troubleshooting
Do not heat the thermopad in a microwave. • Do not use the thermopad without the sleeve provided. Never put the frozen thermopad directly onto your skin. • Do not put the thermopad on damaged skin. • Do not use the thermopad while you are feeding your baby or while expressing milk with a breast pump. • If the thermopad becomes damaged or appears to be leaking, stop using it and throw it away. • Do not swallow the contents of the thermopad. • After freezing, the gel inside the thermopad has a snow-like appearance. This is normal. The gel in the pad returns to its normal condition in a few hours. • If you have any concerns about using the thermopad, please consult your doctor or midwife.

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