



Breathe easier

Patient guide to oxygen therapy

PHILIPS

RESPIRONICS

Table of contents

The facts about COPD	3
Understanding oxygen therapy	4
Choices in oxygen therapy equipment	5
Safely living with oxygen	8
Traveling with oxygen	9
To learn more	10
Trip list for traveling with oxygen	11



The facts about COPD

You have been prescribed oxygen by your physician because you have been diagnosed with chronic obstructive pulmonary disease (COPD), which is a combination of diseases that make breathing difficult. Emphysema and chronic bronchitis are the two main forms of the disease that inhibit your ability to breathe.

When you are unable to exhale, there is not enough room for oxygen-rich air to enter the next time you take a breath. Causes of COPD can include smoking, chronic exposure to irritants, indoor and outdoor pollution, chronic infection, genetic predisposition, or other factors.


Symptoms include:

- shortness of breath
- chronic coughing
- increased mucus production

COPD symptoms may be subtle at first. Many people may think they are just getting older. As a result, they might make small changes in their lives to compensate for the onset of symptoms.




Emphysema

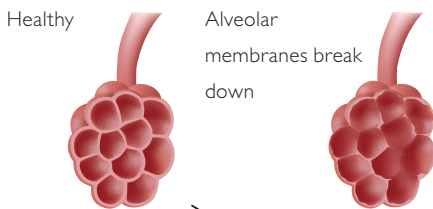

This illustration demonstrates the changes that can occur to the airways of the lungs due to COPD.

Chronic bronchitis

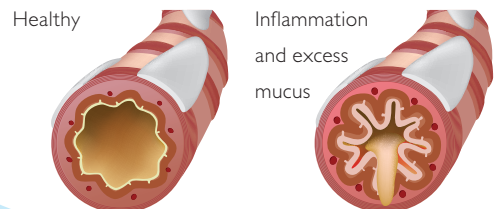

Understanding oxygen therapy

COPD is a progressive disease that can be treated with oxygen to supplement the good air that is entering the lungs. In addition to the other medications prescribed by the physician, the extra oxygen helps protect the body, helps you to function better, and allows you to stay more active.

Low levels of oxygen in the blood may be determined by a pulse oximetry or an arterial blood gas test. When oxygen blood levels are low it is called hypoxemia and can:

- constrict airways within the lungs, straining the heart to pump harder and eventually becoming larger and weaker
- increase the creation of red blood cells that carry oxygen through the blood stream which may cause blood clots, headaches, and high blood pressure
- harm the brain with memory and speech problems
- reduce physical activity

Oxygen therapy may help by reducing the effects of low oxygen levels by improving:

- energy and breathing
- sleep and mood
- mental alertness
- body's ability to carry out normal functions

Or by decreasing:

- morning headaches
- sleepiness
- irritability
- breathlessness



How many hours a day will I need to use the oxygen?

In most cases, oxygen is prescribed for continuous use. If your oxygen level is found to be low, using less than 15 hours a day has not been shown to provide a benefit, and does not protect your heart, brain and other organs of the body. If you are instructed to use continuous oxygen and choose to stop using it temporarily, it is best to do so only while resting quietly and not while sleeping, walking, or exerting yourself.*

Sometimes, oxygen is required only part of the day such as during exercise when you use more energy and therefore need more oxygen. The immediate benefits may be relief of breathlessness and an improvement in your ability to walk or perform daily activities.

Will I need oxygen when I sleep?

During sleep, your breathing slows. People who have low oxygen levels while awake are usually also lacking oxygen during sleep. In some cases, people who do not require oxygen while awake may require extra oxygen while sleeping. Your healthcare professional will determine if and how much oxygen you should receive at night.

*Source: the American Thoracic Society www.thoracic.org/clinical/copd-guidelines/for-patients/index.php

Choices in oxygen therapy equipment

You and your physician will establish the type of oxygen delivery system that is best for your lifestyle and your oxygen requirements. The supplemental oxygen liter flow prescribed can be determined by the oxygen saturation level in the blood found by the pulse oximetry test. Also taken into consideration is whether you will need a continuous flow of oxygen (gas is delivered throughout the breathing cycle) or what is called pulse dose (gas is only delivered while breathing in).

There are two types of oxygen delivery systems: stationary and portable. Stationary systems can provide a large amount of oxygen but are not mobile; portable systems allow you to be more mobile. A portable oxygen concentrator may be referred to as a POC. The following charts highlight some of the systems and their features.

“My POC has been a liberating machine. I cannot say enough good things about it.”

Vernon, POC user

- Equipment choice and modality can be generally based on the following:
- Prescription
 - Oxygen saturation level
 - Physical activity level
 - Amount of travel
 - Type of insurance



Types of oxygen therapy systems

Stationary	Features	Products
Concentrator	Uses a compressor that requires an electrical source	EverFlo
	Can weigh 8.5 to 50 pounds	EverFloQ
	Must be used in open, ventilated area	Millenium M10
Gas transfer	Generates oxygen and fill tanks	SimplyFlo
Liquid	Generates, stores liquid oxygen and refills a portable liquid oxygen device	UltraFill
Compressed gas	Generates, stores liquid oxygen and refills a portable liquid oxygen device	HomeLox
Compressed gas	Large steel or aluminum cylinders from devices that can generate and compress oxygen	E and H cylinders
Portable		
Concentrator	Lightweight, portable option that can run on batteries or AC and DC electricity	SimplyGo
Liquid oxygen	Easily filled from larger liquid storage or liquid generating devices	GoLox
Small compressed gas cylinders	Portable, small aluminum cylinders	UltraFill 2000 or 3000 psi cylinders

Choices in oxygen therapy equipment

Stationary oxygen systems

	Type	Users	Notable features
EverFlo 	5 LPM stationary concentrator	For users who require stationary oxygen of up to 5 LPM	<ul style="list-style-type: none"> • Easy to use and maintain • Compact • Lightweight (31 pounds) • Non-medical design • Transportable • Philips green product
EverFlo Q 	Ultra Quiet 5 LPM stationary concentrator	For users who require stationary oxygen of up to 5 LPM	<ul style="list-style-type: none"> • Ultra-quiet (40 dBA) sound level • Easy to use and maintain • Compact • Lightweight (31 pounds) • Non-medical design • Transportable
SimplyFlo 	2 LPM stationary concentrator	For users who require stationary oxygen up to 2 LPM	<ul style="list-style-type: none"> • Lightweight at 8.5 pounds • Smallest, lightest stationary concentrator in the world • Easily transportable
Millenium M10 	10 LPM stationary concentrator	Supports high-flow users	<ul style="list-style-type: none"> • Durable design • Supports a wider range of users with up to 10 LPM of oxygen output

Portable oxygen systems

	Type	Users	Notable features
<div>SimplyGo</div> <div></div>	Portable oxygen concentrator with continuous flow and pulse dose delivery	Meets the portable needs of nearly all oxygen users	<ul style="list-style-type: none">• Delivers up to 2 LPM continuous flow• Pulse settings 1 to 6• Weighs 10 pounds• Rugged• Requires no user maintenance
<div>UltraFill</div> <div></div>	Home oxygen system that fills gas cylinders in the home	Meets the needs of a wide range of oxygen users, including those who are highly active or require continuous flow oxygen	<ul style="list-style-type: none">• 3000 psi cylinders last >40% longer• One-handed connect and disconnect• Requires no user maintenance
<div>HomeLox</div> <div></div>	Home liquid oxygen system, that fits a small lightweight device or portable use	Provides users who would benefit from the long-lasting portability of traditional liquid oxygen with the ability to make and store liquid oxygen in the home	<div>HomeLox</div> <ul style="list-style-type: none">• Generates liquid oxygen in the home• Generates up to 3 liters per day• Hands-free filling and automatic shut-off• Requires no user maintenance
<div>GoLox</div> <div></div>			<div>GoLox</div> <ul style="list-style-type: none">• Weighs 3.8 pounds• Lasts up to 10 hours• Requires no user maintenance or batteries• Works with HomeLox and other home liquid oxygen systems

Safely living with oxygen

The potential for fire with oxygen is well-known as it make things burn faster. However, this is sometimes forgotten or overlooked. There are a few simple precautions that can be taken to create a safe environment in your home when using oxygen.

First and foremost is to make sure that you have working smoke detectors throughout your house. And, be sure that you have a working fire extinguisher as well. If you move around the house with your oxygen, you may need more than one fire extinguisher placed strategically throughout your living area.

Smoking can be very dangerous. No one, including yourself, should smoke in a room where you are using oxygen. Put a NO SMOKING sign in every room where oxygen is being used.

Be careful with your oxygen when you cook, keeping oxygen away from the stovetop and oven. Watch for splattering grease as it can catch fire.

Keep oxygen at least six feet away from the following:

- toys with electric motors
- electric baseboard or space heaters
- open flames, such as wood stoves or fireplaces
- electric blankets
- hairdryers, electric razors, and electric toothbrushes as they may spark

Other tips:

- Do not store oxygen in a trunk, box, or small closet.
- Storing oxygen under a bed is okay if air can move freely under the bed.
- Do not use oil, grease, or petroleum based products on or near the equipment.
- Avoid petroleum-based lotions and creams on your face and upper chest. Aloe vera and similar products can be used.
- Use caution with oxygen tubing so as not to trip or become entangled in furniture.
- Be familiar with the equipment and safety checks established by your medical equipment provider. Keep emergency service and other important numbers programmed into or posted near your phone.



Oxygen user spotlight—P.J.

P.J. is from Atlanta, Georgia, and uses the GoLox portable oxygen device for her supplemental oxygen needs. She is constantly on the go—working at her job at an advertising agency, selling real estate, entertaining friends, or volunteering her time.

Traveling with oxygen

Traveling with oxygen may seem intimidating. With a little planning, however, it can be easily—and safely—done. If you use oxygen on a regular basis, you will most likely need to use portable oxygen whenever you travel, whether it is across the street, a ride around town, or a trip across the country.

Traveling by car

If you are traveling by car and your vehicle is large enough to transport your oxygen concentrator, you can bring it with you and use wherever you have access to electrical outlets. You can use portable oxygen, which has become quite light and easy to use, while in the car and wherever you don't have access to electricity.

Air travel

When traveling by airplane, it is best to investigate the airlines and their policies on oxygen use during the flight. While booking, and again 24 to 48 hours in advance of your flight, notify the airline about your oxygen requirements as you may be required to sit in a certain seat. A portable oxygen concentrator (POC) can be used as long as it is an FAA-approved device such as the SimplyGo POC.

Bus or train travel

Bus and train lines vary in their rules and regulations. Most bus lines will allow you to carry on portable oxygen, but you'll need to check with them in advance of your trip. Some train lines may permit you to bring oxygen on board but require you use your own generator for power. Always call the bus or train line you're planning to use as far in advance of your trip so you can plan accordingly.

Cruise ships

If you're feeling adventurous enough to try a cruise, you might be surprised to find how easy most cruise lines make traveling with oxygen. Many cruise lines will provide oxygen while others will allow you to bring your portable oxygen concentrator. Again, be sure to call ahead.

“With my POC, I can refill my oxygen when I need it... without waiting.”

George, HomeLox/GoLox user



Oxygen user spotlight—Mark

Mark was diagnosed with COPD in 2003 and uses supplemental oxygen. He is an avid bicyclist in his hometown of Cheyenne, Wyoming and has ridden nearly 9,000 miles while relying on supplemental oxygen.

To learn more

For today's supplemental home oxygen user there are more choices than ever before. We're working with your care team to help them provide you with the right equipment to help make your life better. To learn more about the oxygen therapy equipment featured here, talk to your healthcare professional.



Allies in Better Sleep and Breathing

Philips Respironics is working in harmony with the care team to produce solutions that can lead to healthier lives.

Supplemental oxygen equipment is intended for prescription use by patients requiring high concentrations of oxygen. It may be considered for continuous use in home, institutional, and travel/mobile environments. It is not intended to be life-sustaining or life-supporting. And is not intended for pediatric use. You should talk to your doctor to determine if it is right for you.

Trip list for traveling with oxygen

- ☐ Bring current oxygen prescription that includes liter flow and duration
- ☐ Confirm arrangements with oxygen provider 24 hours in advance
- ☐ Pack power cords and batteries for all electrical equipment
- ☐ Carry extra nasal cannulas and tubing
- ☐ Pack medications in a carry-on bag
- ☐ Bring complete medication record
- ☐ Locate and bring insurance cards
- ☐ Secure summary of medical history from your physician

Medicine Name	Dose	Time	Purpose	Special Directions

Other Information

Known allergies	
Known problems with medication	
Your doctor	Phone
Your oxygen provider	Phone
Your pharmacy	Phone
In case of emergency	Phone

Respironics, GoLox, EverFlo, HomeLox, SimplyFlo, SimplyGo, and UltraFill are trademarks of Respironics, Inc. and its affiliates. All rights reserved.

Please visit www.philips.com/respironics



© 2013 Koninklijke Philips N.V.
All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.
Geyer WMB 07/28/13 MCI 4105753 PN 1112007