

Philips Avent comfort manual breast pump

Philips Avent is dedicated to helping mothers breastfeed for longer, because it is proven that breast milk is vitally important for healthy infant development

Research based product development. We have collaborated with leading scientists and healthcare professionals for over 25 years and conducted research with more than 15,000 mothers worldwide.

What makes this breast pump unique? Recent studies have proven a clear correlation between comfort and milk flow during expressing. Relaxed mothers feel more comfortable and express significantly more milk more easily. Independant research also confirms that infant suckling is a synergistic combination of breast compression and intra-oral vacuum, originating from the peristaltic action of the infant's tongue. The Philips Avent Comfort breast pump is based on exactly these insights.



Unique, natural expressing position

The breast pump has a unique breast shield design. Mothers can sit in a comfortable position when expressing, with no need to lean forward yet milk still flows directly into the container. This improved position helps mothers feel more comfortable and, as a result, express more milk, naturally.

94% of mothers (strongly) agreed that being comfortable when expressing is important for good milk flow.*

Clinically proven to be effective

The breast pump mimics the infant's natural peristaltic suckling behaviour by combining breast compression with gentle suction. It features our clinically proven massage cushion with soft petals that gently compress the breast to stimulate milk flow. This innovative cushion has a distinctive, velvety texture that feels warm on the breast, further helping mothers to relax whilst expressing. The removable cushion also prevents milk pooling, aiding hygiene and eliminating the risk of lost milk. 90% of mothers were (very) satisfied with the effectiveness of the breast pump.*



Effortless expressing

The breast pump is compact with an ergonomically shaped handle for effortless finger-tip control. This allows mothers to keep a tight seal on the breast, making pumping considerably easier. 97% of mothers (strongly) agreed that the new breast pump is easy to use.*

*Independent home placement test among 73 mums in the UK, August 2011



Product components

Comfort breast pump	Specifications	Benefits
Pumping mode	Manual	Full control of vacuum level and pumping speed
Key features	 Unique compact breast shield design with integrated funnel Clinically proven massage cushion with distinctive, velvety texture Compact design 	 Unique, natural expressing position for easier milk flow Gentle stimulation of natural let down and milk flow Effortless expressing, storage and transport
Vacuum levels	As determined by manual compression of handle. Range: 0-250mmHg	Full control of vacuum level and speed
What's in the pack	1x breast pump expression unit (including 1x new Philips Avent Natural bottle) 1x newborn flow teat 1x hygiene cover for cushion 1x sealing disc 1x instruction manual for Comfort breast pump 1x instruction manual for Natural bottle 1x day breastpad sample pack 1x night breastpad sample pack	All mothers need to start expressing, feeding and storing breast milk
Materials	Polypropylene: pump body, hygiene cover for funnel, handle, diaphragm stem, Natural bottle. Silicone: cushion, diaphragm, valve, newborn flow teat TPE: sealing disc. Entirely BPA free	Safe and hygienic to use. All parts can be sterilised
Dimensions	Total height: 17.1cm (when fully assembled) Total width: 13.5cm (when fully assembled)	Compact size for easy expressing. Reduced storage space
Weight	142g (when fully assembled)	Lightweight for easy transport
Cleaning	All parts can be immersed in water and sterilised	Easy to keep the pump hygienic
Assembly	Easy visual matching of parts. Few small pieces	Intuitive assembly
Instructions of use	Printed instruction manual with easy-reference troubleshooting guide	All a mother needs to know about starting to express in one leaflet
Warranty	2 years	Philips Avent consumer care service available for any support
Accessories	Large soft massage cushion available for purchase separately. Funnel size: Ø2.5cm (Ø1.95cm funnel size on regular cushion that comes standard with the Philips Avent Comfort breast pump)	Comfortable expressing for larger than average nipples
Compatibility	Compatibility across the Philips Avent range (including e.g. Classic bottles and milk storage containers)	Flexibility in combining different feeding and storage options
Available spare parts	Pump body, massage cushion, massage cushion large, diaphragm and stem, valve, handle	Can be ordered online or via the Philips Avent consumer care service

Other products to help mothers breastfeed for longer



2 in 1 thermopads: can be used warm, to stimulate milk flow, or cold to soothe sore and engorged breasts.



Day and night breast pads: ultra dry, with silky soft feel. Designed to be extra discreet


keeping mothers comfortable between

feeds. Available for day and night.



Moisturising nipple cream: softens dry or sensitive nipples and helps prepare the delicate breast area for feeding.



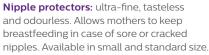



Milk storage containers: for easy and

safe storage of breast milk in the fridge or freezer. Compatible with all Philips Avent breast pumps.

© 2016 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. (Royal Philips) or their respective owners.



Philips Avent, a business of Philips Electronics UK Ltd. Philips Centre, Guildford Business Park, Guildford, Surrey, GU2 8XH GB: 0844 33 80 489 5 pence/minute (mobile rates may vary) IRL: 0818 21 01 41 national call rate (mobile rates may vary)