



**PHILIPS**

Prodiva 1.5T CS

dS MSK M

MR Coil

## Enhance your throughput

### Faster MSK exams with Breeze Workflow


Breeze Workflow for MSK applications is powered by the dS MSK M (dStream MSK Medium) coil. It enables efficient workflow and high throughput through a lightweight coil design and ease of positioning for MSK and other applications. MSK exams with the dS MSK M coil can be performed while the posterior coils (head/neck base and spine) are still present on the table, connected and actively decoupled. Feet first exams for most MSK studies are supported with the dS MSK M coil.

# dS MSK M

Coverage	23 cm
Maximum number of elements	Up to 8 elements (standalone) Up to 23 elements (combined with neuro/vascular/spine coil)
Main applications	Knee, foot, ankle, long bones, wrist, shoulder, pediatric, neuro, vascular
Data sampling	dStream digital broadband
Coil connection	dS Interface
SENSE parallel imaging	Up to 8
Weight	0.98 kg




Simplify exams with SmartSelect: Automatically detects and selects the coil and coil elements which maximize the SNR matching the area to be scanned.

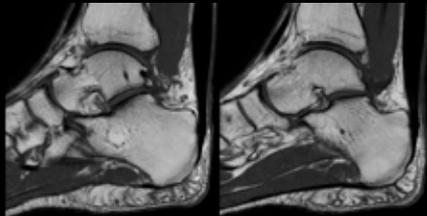

**Anatomy: Ankle**  
**Sequence: Sag T1**  
**Scan Time: 3.18 m**  
**Thickness: 2.5mm**



**Anatomy: Knee**  
**Sequence: PDFS**  
**Scan Time: 3.20m**  
**Thickness: 3mm**

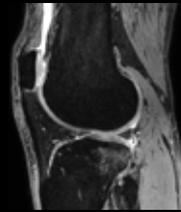

**Anatomy: Knee**  
**Sequence: FFE**  
**Scan Time: 4.20m**  
**Thickness: 1mm**

This material is not for distribution in the USA.

© 2017 Koninklijke Philips N.V. All rights reserved.  
Specifications are subject to change without notice.  
Trademarks are the property of Koninklijke Philips N.V.  
or their respective owners.

4522 991 29701 \* SEP 2017



**How to reach us**  
[www.philips.com/healthcare](http://www.philips.com/healthcare)  
[healthcare@philips.com](mailto:healthcare@philips.com)

**Product information**  
[www.philips.com/mrcoils](http://www.philips.com/mrcoils)