

## Enhance your throughput

## Faster body exams with Breeze Workflow

The dS (dStream) TorsoCardiac Coil 1.5T is an integral part of the Breeze Workflow when used in combination with the dS NeuroVascularSpine (NVS) 1.5T coil. This thin, flat, lightweight anterior coil easily fits around each patient's body shape for fast, efficient, and comfortable scanning. Thanks to short cables, mini connectors and compact electronics, the anterior coil is ultralight. Breeze connect reduces the number of positioning steps and enables you to accelerate your daily exam workflow.


## dS TorsoCardiac

42 cm (landscape mode) 55 cm (portrait mode)
Up to 20 elements (8 anterior + 12 posterior elements)
Thorax, abdomen, pelvis, lower peripherals and cardiac (combined with optional posterior NVS Spine coil). Body oncology
dStream digital broadband
dS Interface
Up to 16
~1.30 kg




Together with the dS NVS coil, the dS TorsoCardiac coil delivers:


- Fast patient set up with Breeze Workflow thanks to ease of coil handling
- Easy exams with SmartSelect. It automatically detects and selects the right coil and coil elements to create the highest SNR for the selected area.
  - Simplifies patient positioning and coil placement
  - No need for manual coil or element selection
- Fast imaging with proven dS SENSE parallel imaging offering up to 16 times faster imaging in core applications

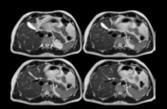


Anatomy: Abodomen Sequence: Multivane XD Scan Time: 2.48m Thickness: 3mm

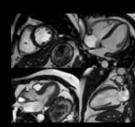

Anatomy: Cardiac Sequence: Cine bTFE

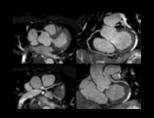

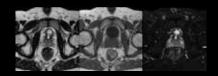

Anatomy: Prostate Sequence: T2/T1/SPAIR Scan Time: 3.17m/1.56m/3.21m Thickness: 3mm

This material is not for distribution in the USA.

© 2017 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.



How to reach us www.philips.com/healthcare healthcare@philips.com

**Product information** www.philips.com/mrcoils

4522 991 30071 \* SEP 2017