

dStream WholeBody

MR coils



Integrated solution for whole body MR imaging

Key benefits

- Full anatomic coverage up to 200 cm
- Conforming and strapless design to promote patient comfort
- Digitization in the coils for increased SNR
- dS SENSE to speed up acquisition in all imaging directions


An integrated coil solution for whole body and peripheral vascular related imaging. It includes two dStream Torso coils. Combined with the dStream Posterior, HeadNeck and dStream Base it enables 200 cm coverage with a maximum of 108 channels. The flexible, lightweight easyto-position dStream Torso coil is designed to conform both in right-left and foot-head directions to support virtually all patients. This allows large coverage and comfortable strap-free operation.

dStream WholeBody coil

Systems	1.5T and 3.0T systems with dStream
Coverage	200 cm
Maximum number of channels	108
Main applications	Whole Body, Abdomen, Pelvis, Cardiac, Peripheral vascular
Coil solution type	Integrated
Data sampling	dStream digital broadband
Coil connection	Triple FlexConnect
SENSE parallel imaging	dS SENSE enhanced parallel imaging performance



T2w TSE Resolution: 1.2 x 1.5 x 6.0 mm Scan time: 0:24 min /station Achieva dStream 1.5T

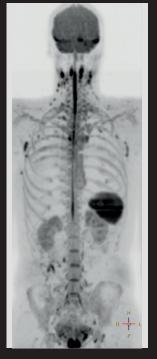

DWIBS Resolution: 5.0 x 5.0 x 6.0 mm Scan time: 1:58 min /station Achieva dStream 1.5T



T1w TSE Resolution: 1.6 x 2.1 x 6.0 mm Scan time: 1:32 min /station Ingenia 1.5T



STIR Resolution: 1.6 x 2.6 x 6.0 mm Scan time: 1:20 min /station Ingenia 1.5T

© 2018 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.



How to reach us healthcare@philips.com www.philips.com/healthcare

More information www.philips.com/mrcoils

4522 991 34621 * APR 2018