



Philips OptiChamber
Diamond
Valved Holding Chamber
and Spacer

Compact design 140ml volume
For asthma inhaler medications
For children and adults
Optional mask accessories



HH1305

Helping you stay in control

Designed to improve medication delivery

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.(1,2)

Effective treatment

- Get medicine to the lungs**

Comfort*

- One size mouthpiece works for all
- Optional mask contours to face to make it easier to fit

Designed for performance

- Special features make it easy to use and maintain
- Low resistance valves let you breathe easily
- Anti-static valved holding chamber gives you more time

For commonly prescribed medications

- For commonly prescribed medications

PHILIPS

Valved Holding Chamber and Spacer

Compact design 140ml volume For asthma inhaler medications, For children and adults, Optional mask accessories

Highlights


Get medicine to the lungs



OptiChamber Diamond is designed to improve aerosol delivery to the lungs and can help you maximize the effectiveness of your inhaled therapies.

Breathe easier

Low resistance valves

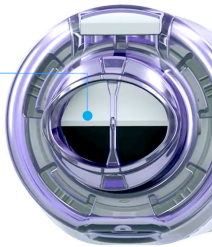

Low resistance valves open freely, even under low flow rates from children, allowing you or your child to breathe easily through the chamber.*

Time to breathe



Anti-static chamber

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The

anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.*1

Easy use, easy maintenance



High-flow whistle

The adapter keeps your inhaler securely in place, while the built-in whistle lets you know if you are breathing in too fast. Both the mouthpiece and adapter are easily removed for cleaning.


For adults and children



Designed with a stepped mouthpiece, the OptiChamber Diamond can be used by both children and adults.

HH1305/00

Comfort

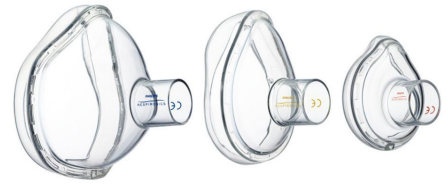

The detachable LiteTouch facemask has a soft contoured cushion that fits easily and rests gently on the face to provide a more comfortable treatment. The SoftTouch seal helps reduce leakage.**

Compatible



OptiChamber Diamond is designed for use with all your commonly prescribed inhaler medications.

Valved Holding Chamber and Spacer

Compact design 140ml volume For asthma inhaler medications, For children and adults, Optional mask accessories

Specifications

Maintenance

- Cleaning: Warm water and liquid soap
- Life span data: Replace after one year

Product details

- Length: 14.2 cm (5.6")

- Mouthpiece: Interfaces w 22 mm connectors
- Volume: 140 ml

Material

- Chamber: Acrylonitrile Butadiene Styrene
- LiteTouch mask: Polyester, Silicone
- Valves: Silicone

HH1305/00



Issue date 2018-06-07

Version: 9.1.1

12 NC: 0000 010 79820
UPC: 3 83735 79820 7

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* Slator L., von Hollen D., Sandell D., Hatley R.H.M. In vitro comparison of the Effect of Inhalation Delay and flow rate on the emitted dose from three valved holding chambers. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 37-43.

* (1) Berlinski A., von Hollen D., Hatley R.H.M., Hardaker L.E.A., Nikander K. Drug delivery in asthmatic children following coordinated and uncoordinated inhalation maneuvers: a randomized crossover trial. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2017, 30(3): 182-189.

* ** Hatley R.H.M., von Hollen D., Sandell D., Slator L. In-vitro Characterization of the OptiChamber Diamond valved holding chamber. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 24-36.

* *** Tong, K. et al. An instrumented Valved Holding Chamber with facemask to measure application forces and flow in young asthmatic children. *Journal of Aerosol Med and Pulm Drug Del.* 2014; 27 (Suppl 1): S55-62.