



Philips OptiChamber
Diamond
Valved Holding Chamber &
Mask

Large LiteTouch Mask included
Compact design 140ml volume
For common asthma medications



HH1308

Helping you stay in control

Designed to improve medication delivery

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.(1,2)

Effective treatment

- Get medicine to the lungs(3)

Comfort

- Optional mask contours to face to make it easier to fit
- One size works for all

Designed for performance

- Special features make it easy to use and maintain
- Low resistance valves let you breathe easily
- Anti-static valved holding chamber gives you more time

For commonly prescribed medications

- For commonly prescribed medications

PHILIPS

Highlights


Get medicine to the lungs



OptiChamber Diamond is designed to improve aerosol delivery to the lungs and can help you maximize the effectiveness of your inhaled therapies.

Breathe easier

Low resistance valves

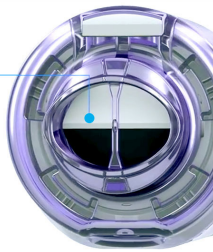

Low resistance valves open freely, even under low flow rates from children, allowing you or your child to breathe easily through the chamber.(1)

Compatible



OptiChamber Diamond is designed for use with all your commonly prescribed inhaler medications.

Easy use, easy maintenance



High-flow whistle

The adapter keeps your inhaler securely in place, while the built-in whistle lets you know if you are breathing in too fast. Both the mouthpiece and adapter are easily removed for cleaning.

Time to breathe




Anti-static chamber

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The

anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.(1,2)

Comfort

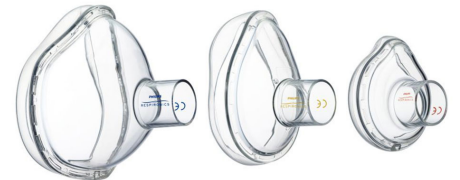

The detachable LiteTouch facemask has a soft contoured cushion that fits easily and rests gently on the face to provide a more comfortable treatment. The SoftTouch seal helps reduce leakage.(4)

For adults and children



Designed with a stepped mouthpiece, the OptiChamber Diamond can be used by both children and adults.

Specifications

Product details

- For use with: Adults
- Length: 14.2 cm (5.6")
- Mouthpiece: Interfaces w 22 mm connectors
- Volume: 140 ml

Package contents

- Include: Large LiteTouch Mask, OptiChamber Diamond

Maintenance

- Cleaning: Warm water and liquid soap
- Life span data: Replace after one year

Material

- Chamber: Acrylonitrile Butadiene Styrene
- LiteTouch mask: Polyester, Silicone
- Valves: Silicone



Issue date 2018-03-14

Version: 11.1.1

12 NC: 0000 010 79828
UPC: 3 83735 79828 3

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* (1) Slatore L., von Hollen D., Sandell D., Hatley R.H.M. In vitro comparison of the Effect of Inhalation Delay and flow rate on the emitted dose from three valved holding chambers. Journal of Aerosol Medicine and Pulmonary Drug Delivery, 2014, 27(S1): 37-43.

* (2) Berlinski A., von Hollen D., Hatley R.H.M., Hardaker L.E.A., Nikander K. Drug delivery in asthmatic children following coordinated and uncoordinated inhalation maneuvers: a randomized crossover trial. Journal of Aerosol Medicine and Pulmonary Drug Delivery, 2017, 30(3): 182-189.

* (3) Hatley R.H.M., von Hollen D., Sandell D., Slatore L. In-vitro Characterization of the OptiChamber Diamond valved holding chamber. Journal of Aerosol Medicine and Pulmonary Drug Delivery, 2014, 27(S1): 24-36.

* (4) Tong, K. et al. An instrumented Valved Holding Chamber with facemask to measure application forces and flow in young asthmatic children. Journal of Aerosol Med and Pulm Drug Del. 2014; 27 (Suppl 1): S55-62.