



Valved Holding Chamber with Mask

- Small LiteTouch Mask included
- For children aged 0–18 m
- For common asthma medications
- Compact design 140 ml volume

HH1306/00

Helping you stay in control

Designed to improve medication delivery

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended for longer, giving you more time to inhale.(1,2)

Benefits

Comfort

- Optional mask contours to face to make it easier to fit
- One-size-fits-all mouthpiece

Effective treatment

- Get medicine to the lungs**

Designed for performance

- Low-resistance valves let you breathe easily
- Special features make it easy to use and maintain
- Anti-static valved holding chamber gives you more time

For commonly prescribed medications

- For commonly prescribed medications

PHILIPS

OptiChamber Diamond

Features

Get medicine to the lungs



OptiChamber Diamond is designed to improve aerosol delivery to the lungs and can help you maximise the effectiveness of your inhaled therapies.

Breathe easier



Low resistance valves open freely, even under low flow rates from children, allowing you or your child to breathe easily through the chamber.*

Comfort



The detachable LiteTouch facemask has a soft contoured cushion that fits easily and rests gently on the face to provide a more comfortable treatment. The SoftTouch seal helps reduce leakage.***

Compatible



OptiChamber Diamond is designed for use with all of your commonly prescribed inhaler medications.

Specifications

Please note that this is a pre-sales leaflet. The contents of this leaflet reflect the best of our knowledge per date and country mentioned above. The contents of this leaflet are subject to change without notice. Philips does not accept any liability as to the contents of this leaflet.


Product details

For use with	Children aged 1–5 yr
Volume	140 ml
Length	14.2 cm (5.6")
Mouthpiece	Interfaces w 22 mm connectors

Package contents

Include	Small LiteTouch Mask OptiChamber Diamond
---------	---

Easy use, easy maintenance

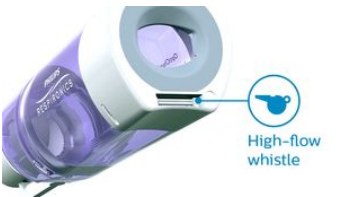

The adapter keeps your inhaler securely in place, while the built-in whistle lets you know if you are breathing in too fast. Both the mouthpiece and adapter are easily removed for cleaning.

For adults and children



Designed with a stepped mouthpiece, the OptiChamber Diamond can be used by both children and adults.

Time to breathe



OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.*1

Material

Chamber	Acrylonitrile Butadiene Styrene
LiteTouch mask	Polyester, Silicone
Valves	Silicone

Maintenance

Cleaning	Warm water and liquid soap
Life span data	Replace after one year

* **Slator L., von Hollen D., Sandell D., Hatley R.H.M. In vitro comparison of the Effect of Inhalation Delay and flow rate on the emitted dose from three valved holding chambers. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 37-43.

* (1)Berlinski A., von Hollen D., Hatley R.H.M., Hardaker L.E.A., Nikander K. Drug delivery in asthmatic children following coordinated and uncoordinated inhalation maneuvers: a randomized crossover trial. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2017, 30(3): 182-189.

* **Hatley R.H.M., von Hollen D., Sandell D., Slator L. In-vitro Characterization of the OptiChamber Diamond valved holding chamber. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 24-36.

* ***Tong, K. et al. An instrumented Valved Holding Chamber with facemask to measure application forces and flow in young asthmatic children. *Journal of Aerosol Med and Pulm Drug Del.* 2014; 27 (Suppl 1): S55-62.



data subject to change
2018, July 31

Version: 12.1.1
EAN: 383735798221

© 2018 Koninklijke Philips N.V.
All rights reserved.

Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

www.philips.com