

Philips DreamWear
Nasal mask

Fit pack -frame and 4 cushions
INTL

HH1111



DreamWear®

Nasal mask with mask cushion underneath the nose

The innovative design* of DreamWear offers you the freedom and comfort to sleep in the way you would like.

Freedom of movement

- Innovative mask design for more freedom in sleep therapy
- More comfort with DreamWear
- For a more natural sleep
- Soft, flexible mask frame

No red marks on the bridge of the nose

- Cushion under the nose
- Slim headgear design

Unobstructed view

- Hose connection on the head
- Read without disabilities or watch TV

PHILIPS

Nasal mask
Fit pack -frame and 4 cushions INTL

Highlights

Cushion under the nose



The mask cushion lies gently under the nose, which causes no pressure points on the bridge of the nose and prevents unpleasant skin irritations on the nostrils.

A more natural sleep



DreamWear combines the advantages of minimal contact masks and nasal masks, provides high comfort and can contribute to successful sleep therapy.

Hose connection on the head




The exhaled air is not discharged directly at the nose, but at the top of the head. This hose positioning allows more freedom of movement during sleep. Users say they no longer need to sleep in a specific position with DreamWear, but have more freedom.*

Read or watch TV



Thanks to the open design, DreamWear offers a very large field of view, so you can easily wear glasses and read, watch TV or use a tablet PC in bed before going to sleep.

Soft, flexible mask frame

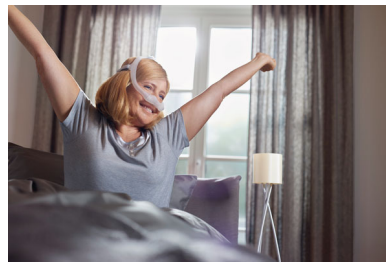

The flexible material, through which the air flows, gently rests on the skin, prevents pressure points and provides a comfortable fit.

Innovative mask design

The innovative design* of DreamWear prevents pressure sores and skin irritation on the nostrils and on the nose bridge, and allows the user, in addition to a clear field of view, more freedom of movement through the hose connection on the head.

Slim headgear design



DreamWear unites comfort and aesthetics - the headgear is unobtrusive and sleek in color, making DreamWear as discrete and unobtrusive as possible.

More comfort with DreamWear

Users said DreamWear was more comfortable to wear than their previously prescribed mask, as it is barely noticeable.

HH1111/01

Specifications

Freedom of movement

- Clear field of vision: Individual sleep position

Innovative design

- Comfort and performance: Reduced skin irritation

Service

- Warranty: 90-days

Intent of use

- To be used on patients >30kg: with CPAP prescription

Specifications

- Operating pressure: 4 -20 cmH2O
- Sound pressure level: 18 dBA
- N° of parts: Five
- Headgear material: Urethane foam/nylon/Spandex
- Frame material: Polycarbonate
- Cleaning instructions: mask: Hand wash with mild detergent, Rinse thoroughly, air dry
- Cleaning: headgear, tubing: Hand wash with mild detergent, Rinse thoroughly, hang to dry
- Frame/Elbow/Swivel material: Polycarbonate
- Fabric sleeves material: Polyester/Spandex



Issue date 2018-09-22

Version: 6.3.1

12 NC: 0000 011 16700

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

www.philips.com

* Source: Philips User Survey 2015

* The costs of DreamWear are private and are usually not reimbursed by the statutory health insurance.

* The scope of supply includes a standard mask frame (size M) and four cushions for all sizes (S, M, MW, L). The sleep therapy device is not included.