



Philips SideStream
Reusable Nebulizer

Treatment in 6 minutes*

Can be used with mask
Suitable for adults & children
Replace after 6 months

HH1320/00

Treatment in 6-8 minutes*

Proven technology makes treatment fast and easy

The SideStream Reusable Nebulizer can be used for 6 months before replacing. With its unique design, SideStream offers faster drug delivery and shorter treatment times, and can be used with most commonly prescribed respiratory medications..

Fast treatment times

- Designed to draw in air

Proven technology

- Treatment delivered in 6-8 minutes*

Consistent performance

- Rely on consistent treatment

PHILIPS

Reusable Nebulizer

Treatment in 6 minutes* Can be used with mask, Suitable for adults & children, Replace after 6 months

Specifications

Maintenance

- Replace: After 6 months

Product details

- Aerosol flow to patient: 16 L/min
- Driving gas: 6 L/min
- Venturi air entrainment: 10 L/min

HH1320/00

Highlights

Effective delivery over time




Staying compliant with your breathing treatment is easier when taking your medication is consistent and reliable. SideStream nebulizers keep the driving air and the liquid drug separate for ongoing performance.

Fast, effective treatment*



The nebulizer top, combined with an InnoSpire compressor, SideStream delivers faster aerosol drug delivery and shorter treatment times, so you can get back to your daily life.

Fast, effective treatment time

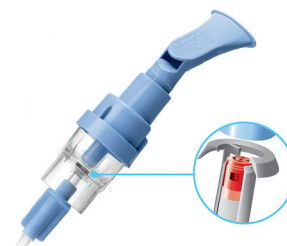

The SideStream nebulizer's venturi system boosts the airflow supplied by the compressor, delivering medication more efficiently and making your overall treatment faster.



Issue date 2018-07-26

Version: 1.0.1

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

** Using a 2.5ml Salbutamol dose