

Philips OptiChamber
Diamond
Valved Holding Chamber
and Spacer

Compact design 140ml volume

For asthma inhaler medications
For children and adults
Optional mask accessories

HH1329



Helping you stay in control

Designed to improve medication delivery

OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.(1,2)

Comfort

- Optional mask contours to face to make it easier to fit
- One size mouthpiece works for all

Effective treatment

- Get medicine to the lungs**

Designed for performance

- Special features make it easy to use and maintain
- Low resistance valves let you breathe easily
- Anti-static valved holding chamber gives you more time

For commonly prescribed medications

- For commonly prescribed medications

PHILIPS

Valved Holding Chamber and Spacer

Compact design 140ml volume For asthma inhaler medications, For children and adults, Optional mask accessories

Highlights

Breathe easier



Low resistance valves open freely, even under low flow rates from children, allowing you or your child to breathe easily through the chamber.*

Comfort




The detachable LiteTouch facemask has a soft contoured cushion that fits easily and rests gently on the face to provide a more comfortable treatment. The SoftTouch seal helps reduce leakage.**

Compatible



OptiChamber Diamond is designed for use with all your commonly prescribed inhaler medications.

Easy use, easy maintenance

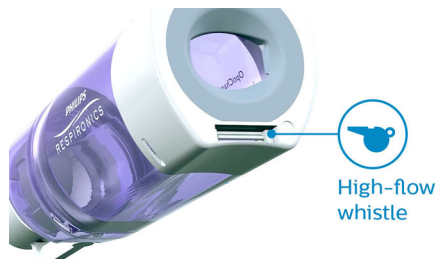

The adapter keeps your inhaler securely in place, while the built-in whistle lets you know if you are breathing in too fast. Both the mouthpiece and adapter are easily removed for cleaning.

For adults and children



Designed with a stepped mouthpiece, the OptiChamber Diamond can be used by both children and adults.

HH1329/00

Get medicine to the lungs



OptiChamber Diamond is designed to improve aerosol delivery to the lungs and can help you maximize the effectiveness of your inhaled therapies.

Time to breathe



OptiChamber Diamond is designed to allow medication to be delivered to the lungs. The anti-static material inside the tube allows the aerosol to be suspended longer, giving you more time to inhale.*1

Valved Holding Chamber and Spacer

Compact design 140ml volume For asthma inhaler medications, For children and adults, Optional mask accessories

Specifications

Maintenance

- Cleaning: Warm water and liquid soap
- Life span data: Replace after one year

Product details

- Length: 14.2 cm (5.6")

- Mouthpiece: Interfaces w 22 mm connectors
- Volume: 140 ml

Material

- Chamber: Acrylonitrile Butadiene Styrene
- LiteTouch mask: Polyester, Silicone
- Valves: Silicone

HH1329/00



Issue date 2018-08-23

Version: 3.0.1

12 NC: 8670 001 46071
EAN: 38 37300 90597 4

© 2018 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* (1)Berlinski A., von Hollen D., Hatley R.H.M., Hardaker L.E.A., Nikander K. Drug delivery in asthmatic children following coordinated and uncoordinated inhalation maneuvers: a randomized crossover trial. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2017, 30(3): 182-189.

* *Slator L., von Hollen D., Sandell D., Hatley R.H.M. In vitro comparison of the Effect of Inhalation Delay and flow rate on the emitted dose from three valved holding chambers. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 37-43.

* **Hatley R.H.M., von Hollen D., Sandell D., Slator L. In-vitro Characterization of the OptiChamber Diamond valved holding chamber. *Journal of Aerosol Medicine and Pulmonary Drug Delivery*, 2014, 27(S1): 24-36.

* ***Tong, K. et al. An instrumented Valved Holding Chamber with facemask to measure application forces and flow in young asthmatic children. *Journal of Aerosol Med and Pulm Drug Del.* 2014; 27 (Suppl 1): S55-62.