

PHILIPS

IPL - Hair removal device

Lumea Prestige

For use on body, face, bikini
15 minutes to treat lower legs
>250,000 lamp flashes
Cordless design



SC2009/00



Feel continuously hair free*

Philips Lumea Prestige is the most effective cordless IPL solution and works wonders to prevent the reappearance of body and facial hair. When applied regularly, the gentle pulses of light keep skin silky-smooth every day.

Long-lasting hair removal with IPL technology

- Enjoy smooth skin every day
- Effortlessly effective
- Adapted for safe and effective use in your own home
- Suitable for a wide variety of hair and skin types

Gentle and effective on body, face and bikini

- Safe and effective even on sensitive areas
- Large treatment window for fast application
- Facial attachment for safe facial treatment
- Attachment for specialised treatment of bikini area hair
- Gentle treatment

Maximum freedom, minimum hassle

- Cordless convenience
- No replacement parts, no hidden costs
- Five adjustable light energy settings



asimpleswitch.com


Highlights

IPL technology



Philips Lumea uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in professional beauty salons. Philips has adapted this technology for safe and effective use in the comfort of your own home. Philips worked closely with leading dermatologists to develop our breakthrough hair removal system. For more than 10 years we carried out extensive consumer research with more than 2000 volunteers.

Enjoy smooth skin every day

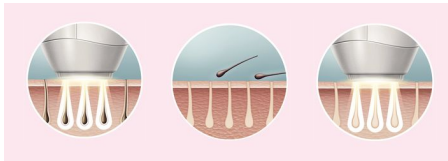

Philips Lumea applies gentle pulses of light to the hair root. As a consequence, the amount of hair your body grows gradually decreases. Repeating the treatment leaves your skin beautifully hair free and touchably smooth.

Effortlessly effective

Clinical results after treatment: **example 1** (leg)



Clinical results after treatment: **example 2** (bikini line)



Our clinical studies have shown significant hair reduction in as little as four fortnightly treatments resulting in smooth hair free skin. To maintain these results, simply repeat the treatment when needed. Time between

treatments may vary based on your individual hair regrowth. - To enlarge this picture, please click on the image in the image gallery on the top of this page

Suitable hair and skin types




Philips Lumea works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Philips Lumea is NOT effective on red, light-blonde or white/grey hair. Philips Lumea is also NOT suitable for very dark skin.

Application areas



Use with confidence to remove hair on the face (upper lip, chin and sideburns) and body parts including the legs, underarms, bikini area, belly and arms.

Fast treatment on legs

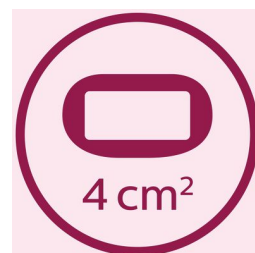

Big treatment window for fast application on large areas like legs.

Extra safety for the face



The facial attachment has its own specially designed integrated light filter for safe and convenient use on the upper lip, chin and sideburns. It can also be used on other sensitive areas of the body, e.g. underarms and other hard-to-reach areas.

Bikini area attachment



Bikini area attachment for effective treatment of bikini hair. Hair in this area tends to be typically stronger and thicker than leg hair.



asimpleswitch.com

Philips Green Logo

Philips Green Products can reduce costs, energy consumption and CO2 emissions. How? They offer a significant environmental improvement in one or more of the Philips Green Focal Areas — Energy efficiency, Packaging, Hazardous substances, Weight, Recycling and disposal and Lifetime reliability.

Specifications

Technical specifications

High performance lamp: Delivers > 250 000 flashes

Application areas

Body areas: Arms, Legs, Belly, Bikini, Underarms

Face areas: Chin, Sideburns, Upper lip

Attachments

Body attachment (4 cm²): For use below the neckline

Bikini area attachment (2 cm²): Special treatment bikini hair

Facial attachment (2 cm²): For safe application on face

Items included

Charger: Multi voltage charger

Instructions for use: User manual

Storage: Luxurious pouch

Accessory: Cleaning cloth

Power

Charging time: 100 min

Running time: Minimum 370 flashes

Voltage: 100-240 V

Safety and adjustable settings

Integrated UV filter: Protects skin from UV light

5 light energy settings: Adjustable to your skin type

Integrated safety system: Prevent unintentional flashing

Service

Warranty: 2 years global warranty

Application time

Lower leg: 8 min

Underarm: 1 min

Bikini line: 1 min

Face areas: 1 min

Application mode

Slide and Flash: For convenient application

Corded / cordless use: Cordless

Stamp and Flash: For treatment on curvy areas

