

PHILIPS

IPL hair removal
system

Lumea Precision
Plus

For use on body, face, bikini

15 minutes to treat lower legs

Lifetime >140.000 light pulses

Cordless design



SC2008/11



The ideal way to prevent hair regrowth for smooth skin at home

Philips Lumea IPL is the most effective cordless IPL solution and works wonders to prevent the reappearance of hair on face and body. Gentle pulses of light, applied regularly, keep skin silky-smooth every day.

Long lasting smoothness

- Enjoy smooth skin every day
- Effortlessly effective

Hair removal with IPL technology

- Adapted for safe use in the comfort of your own home
- Suitable for a wide variety of hair and skin types

Gentle and effective on body, face and bikini

- Safe and effective even on sensitive areas
- Bigger body attachment for fast treatment
- Facial attachment for safe facial treatment
- Attachment for specialised treatment of bikini area hair
- Five adjustable light energy settings

Maximum freedom, minimum hassle

- Cordless convenience
- No replacement parts, no hidden costs



asimpleswitch.com

Highlights

Long lasting smoothness



Philips Lumea applies gentle pulses of light to the hair root. As a consequence the hair sheds naturally and hair regrowth is inhibited. Repeating this treatment regularly will keep your skin smooth every day.

Effortlessly effective



Our clinical studies have shown significant hair reduction in as little as four bi-weekly treatments resulting in smooth hair free skin. To maintain these results, simply repeat the treatment when needed. Time between treatments may vary based on your individual hair regrowth.

IPL technology



Philips Lumea uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in professional beauty salons. Philips has adapted this technology for safe and effective use in the comfort of your own home. Philips worked closely with leading dermatologists to develop our breakthrough hair removal system. For more than 10 years we carried out extensive consumer research with more than 2000 women.

Suitable hair and skin types

Suitable hair and skin tones

Skin Tone	White	Beige	Light brown	Medium brown	Dark brown	Brownish black and darker
Black	⊗	⊗	⊗	⊗	⊗	⊗
Dark brown	✓	✓	✓	✓	✓	⊗
Brown	✓	✓	✓	✓	✓	⊗
Dark blond	✓	✓	✓	✓	✓	⊗
Light blond	⊗	⊗	⊗	⊗	⊗	⊗
Red	⊗	⊗	⊗	⊗	⊗	⊗
White/grey	⊗	⊗	⊗	⊗	⊗	⊗

● This product is suitable for you to use ⊗ This product is neither suitable nor safe for your skin and hair type


Philips Lumea works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Philips Lumea is NOT effective on red, light-blond or white/grey hair. Philips Lumea is also NOT suitable for very dark skin.

Application areas



Use with confidence to remove hair on face (upper lip, chin and sideburns) and body parts including legs, underarm, bikini area, belly and arms.

Fast treatment on legs

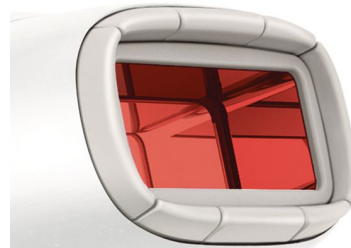

Bigger body attachment for fast treatment on large areas like legs.

Extra safety for the face



The facial attachment has its own specially designed integrated light filter for safe and convenient use on the upper lip, chin and

sideburns. It can also be used on other sensitive areas on the body, eg underarms and other hard-to-reach areas.

Bikini area attachment



Bikini area attachment for effective treatment of bikini hair. Hair in this area tends to be typically stronger and thicker than leg hair.

Five light energy settings



Philips Lumea has five adjustable light energy settings to ensure gentle but effective treatment. Used correctly, Philips Lumea's IPL treatment is safe and gentle to use even on sensitive skin and sensitive body areas.



asimpleswitch.com

Philips Green Logo

Philips Green Products can reduce costs, energy consumption and CO2 emissions. How? They offer a significant environmental improvement in one or more of the Philips Green Focal Areas – Energy efficiency, Packaging, Hazardous substances, Weight, Recycling and disposal and Lifetime reliability.

Specifications

Technical specifications

High performance lamp: Delivers > 140,000 flashes

Safety and adjustable settings

5 light energy settings: Adjustable to your skin type

Integrated UV filter: Protects skin from UV light

Integrated safety system: Prevent unintentional flashing

Application areas

Body areas: Legs, Arms, Belly, Bikini, Armpits

Face areas: Chin, Sideburns, Upper lip

Application mode

Slide and Flash: For convenient application

Step and Flash: For treatment on curvy areas

Corded / cordless use: Cordless

Attachments

Body attachment (4cm2): For use below the neckline

Bikini area attachment (2cm2): Special treatment bikini hair

Facial attachment (2cm2): For safe application on face

Application time

Lower leg: 8 min

Armpit: 1 min

Bikini line: 1 min

Face areas: 1 min

Items included

Storage: Pouch

Instructions for use: Quick Start Guide, User manual

Charger: Multi voltage charger

Power

Voltage: 100-240 V

Charging time: 100 min

Running time: Minimum 370 flashes

Service

Warranty: 2 years global warranty

