

# PHILIPS

IPL - Hair removal device

Lumea Advanced

For body and face

With skin tone sensor

15 minutes to treat lower legs

+ Satin Compact pen trimmer



BRI923/60



## Be hair free for at least 3 months\*

### At least 85% hair reduction in 3 treatments\*\*

Philips Lumea Advanced works wonders to prevent the reappearance of hair on the body. Gentle pulses of light, applied regularly, keep skin silky-smooth every day.

#### **Long lasting hair removal with IPL technology**

- Expert IPL technology at-home, developed with dermatologists
- Suitable for a wide variety of hair and skin types
- FREE Lumea App for a personal coach at your fingertips

#### **Safe, fast and effective on body and face**


- Proven safe and effective treatment
- Integrated skin tone sensor

#### **Specialized treatment of different body areas**

- Safe and effective even on sensitive areas
- Large treatment window for fast application
- Precision attachment for extra safety on the face
- Attachment for specialised treatment of bikini area hair
- Special Beauty Edition including Satin Compact pen trimmer


# Highlights

## Long lasting smoothness

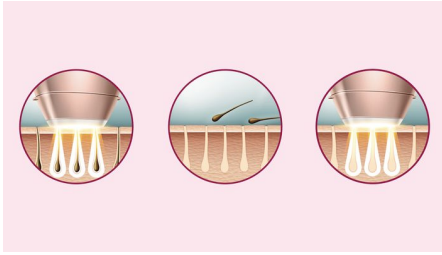

IPL stands for Intense Pulsed Light. Philips Lumea applies gentle pulses of light to the hair root, putting the follicle into a resting phase. As a consequence, the amount of hair your body grows gradually decreases. Repeating the treatment leaves your skin beautifully hair free and touchably smooth. Treatment to prevent hair regrowth is safe and gentle, even on sensitive areas. Philips Lumea is clinically tested and developed with dermatologists for easy and effective treatment, from the comfort of your own home.

## Proven effectiveness

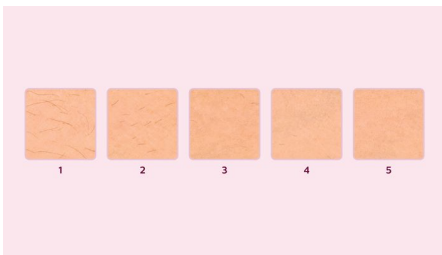

Objective studies show at least 85% hair reduction in as little as 3 treatments\*\*. The first 4 treatments must be carried out every 2 weeks, after which you should already be able to see the desired results. To maintain your results, simply touch up every 4 weeks. After just 8 touch up treatments you can be hair-free for at least 3 months\*. To enlarge this picture, please click on the image in the image gallery at the top of this page.

## Suitable hair and skin types



Philips Lumea Prestige works effectively, easily and safely on a wide variety of hair and skin types. It works on naturally dark blonde, brown and black colored hairs and on skin tones from very white to dark brown. As with other IPL based treatments, Philips Lumea cannot be used to treat white / grey, light blonde or red hair and is not suitable for very dark skin. This is due to the high contrast required between the pigment in hair color and pigment in the skin tone.

## Skin Tone Sensor



The integrated skin tone sensor measures the treated skin complexion at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for this product it will automatically stop emitting pulses.

## Application areas



Use with confidence to remove hair on face (upper lip, chin and sideburns) and body parts including legs, underarm, bikini area, stomach and arms.

## Fast treatment on legs



Big treatment window for fast application on large areas like legs.

## Extra safety for the face



With the additional integrated light filter, it can be used for safe treatment of facial hair on the upper lip, chin, cheeks and other sensitive areas.

## Bikini area attachment



Bikini area attachment for effective treatment of bikini hair. Hair in this area tends to be typically stronger and thicker than leg hair.

# Specifications

## Application areas

**Body areas:** Arms, Legs, Belly, Armpits, Bikini

**Face areas:** Chin, Sideburns, Upper lip

## Attachments

**Body attachment (4cm<sup>2</sup>):** For use below the neckline

**Face attachment (2cm<sup>2</sup>):** For use below the eyes

**Bikini area attachment (2cm<sup>2</sup>):** Special treatment bikini hair

## Application mode

**Slide and Flash:** For easy motion application

**Corded / cordless use:** Corded

**Stamp and Flash:** For treatment on curvy areas

## Application time

**Bikini:** 4min

**Armpit:** 2,5 min

**Face areas:** 2 min

**Lower leg:** 15 min

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**Skin Tone sensor:** Detects your skin tone

## Power

**Voltage:** 100-240 V

**Corded use**

## Items included

**Instructions for use:** User manual

**Storage:** Pouch

**Adapter:** 24V / 1500mA

**Accessory:** Cleaning cloth

## Pen trimmer

**Accessories:** Eyebrow comb, Cleaning brush

**Ease of use:** Compact handle

**Power:** AAA battery

## Service

**Warranty:** 2 years global warranty

