

PHILIPS

CPAP with
Humidifier

DreamStation Go

Thinnest travel humidifier¹

Uses tap or bottled water³

Works with any mask you like²

DreamMapper compatibility

HH1540/00



At-home CPAP performance and comfort on the go

Keep your comfortable, all-night therapy going, everywhere you go. The DreamStation Go system includes our powerful, lightweight travel CPAP and innovative heated humidifier to deliver the same therapy enjoyed by over 5 million users.

Heated Humidifier

- Our smallest CPAP with the thinnest travel humidifier⁵
- Use any mask you like
- Freedom to use tap water³
- The same features and connectivity you enjoy at home
- Stay in control with DreamMapper, even on the go
- Manages water usage evenly for comfort throughout the night
- Smart humidification. All night long.
- Micro-flexible standard 12mm tubing

Highlights

Pack performance and comfort



At half the size of our trusted System one and DreamStation platforms, the compact DreamStation Go is our smallest, lightest PAP⁵. The DreamStation Go system combines the humidifier and trusted device together in one compact, convenient package that fits easily into your suitcase. The system complies with FAA requirements⁶ so now where your travel takes you, you can take your therapy on the go.

Works with any mask



When you've found a mask you're comfortable with, the last thing you want to do is be forced to switch your mask while travelling. The DreamStation Go heated humidifier is designed to work with any mask that works for you, so you can enjoy on-the-go heated humidification from the comfort of your own mask.²

Freedom to use tap water³




With the freedom to use tap, bottled or distilled water³, the DreamStation Go heated humidifier reduces the inconvenience of taking along heavy, bulky containers or searching for distilled water when you arrive. That means you can choose to travel light and still enjoy the same comfort of heated humidification that you experience at home.

Connectivity on the go



Bring all the benefits of your home CPAP along, every time you travel. The powerful DreamStation Go features an intuitive color touchscreen for quick, simple setup and easy navigation, and includes Bluetooth connectivity for tracking and monitoring of your progress across both home and travel CPAPs.

DreamMapper capable

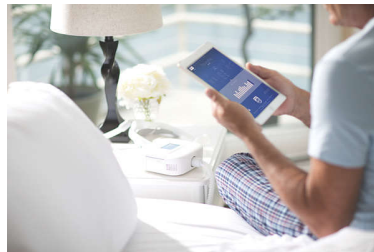

Sleep apnea doesn't stay behind when you're away from home. That's why the DreamStation Go is compatible with our

DreamMapper self-management tool. DreamMapper combines your home and travel data in one place so you can see how well your night's sleep went, how you're progressing, and track your improvements towards your goals — every day, anywhere you are.

Water saving technology



The DreamStation Go heated humidifier's water-saving technology adapts to a wide range of ambient conditions to deliver humidification evenly throughout the full night's sleep.⁴ You can enjoy the same all-night comfort you experience at home, no matter where your travels take you.

Monitors and adjusts therapy

DreamStation Go features water-saving technology designed to give you a full night of heated humidification⁴ by actively monitoring ambient conditions and therapy settings and automatically adjusting moisture delivery when humidifier target time is enabled. Also, in the event you underfill or forget to fill the humidifier, our heater plate shutdown feature will take over to prevent dry, uncomfortable air from being delivered to you for the remainder of the night.

Micro-flexible tubing

The DreamStation Go is the only heated humidifier with micro-flexible standard 12mm tubing. The 12mm tubing is Philips' smallest and lightest tubing ever offered⁷, providing enhanced packability, flexibility and freedom of movement while you sleep. Fast and easy click connections on each end of the tubing. Features a slim connection to our DreamWear mask, and includes a 22mm ISO adapter for other masks.

Specifications

Product Details

Size: 25.3 x 15.1 x 5.9 cm

Weight: 1150.4 g

Warranty: 2 years

Humidifier Details

Size: 14.7 x 15.1 x 5.9 cm

Weight: 344.7 g

Humidifier Capacity: 240 ml

Tubing compatibility: Micro-flexible 12mm, 15mm and 22mm standard tubing

Water compatibility: Distilled, tap, bottled³

Humidifier usage while on Battery power: No

Technical information

Pressure range: 4 to 20 cm H₂O

Pressure mode(s): CPAP, APAP

Flex modes: C-Flex, C-Flex+, A-Flex

Device setup: Touchscreen color display

Data storage capacity: MicroSD card (optional), On-board: > 1 year

Altitude compensation: Automatic

Electrical requirements: 100-240 VAC, 50/60 Hz

Battery technology: Lithium Ions

Output power: 50W



- * 1. When compared against Transcend, HDM Z1, AirMini (does not offer heated humidifier), Devilbiss
- * 2. DreamStation Go is equipped with a universal adapter designed to fit all standard 22 mm OD connections.
- * 3. In very hard water areas, it is recommended that you use distilled or bottled water.
- * 4. Adjusts humidification output when water-saver target time is enabled. Refer to the User Manual for more information.
- * 5. 2016 Internal measurement of DreamStation. Includes standard power supply.
- * 6. RTCA/DO-160G section 21, category M
- * 7. When compared to Philips 22mm Performance tubing