

# PHILIPS

IPL – Hair removal  
device

Lumea Essential

For use on body and face

22 minutes to treat lower legs

>200,000 lamp flashes

Compact corded design

SC1983/00

 Product Picture

## Easy, long term hair reduction

Philips Lumea Essential works wonders to prevent the reappearance of hair on face and body. Gentle pulses of light, applied regularly, keep skin silky-smooth every day.

### **Long lasting hair removal with IPL technology**

- Enjoy smooth skin every day
- Effortlessly effective
- Adapted for safe use in the comfort of your own home
- Suitable for different hair and skin types

### **Safe and effective on face and body**

- Safe and effective even on sensitive areas
- Skin Tone Sensor
- Gentle treatment

### **Maximum freedom, minimum hassle**

- No replacement parts, no hidden costs
- The smallest and lightest device within the Lumea range
- Slide and Flash mode for a more convenient application
- Five adjustable light energy settings


# Highlights

## IPL technology




Philips Lumea uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in professional beauty salons. Philips has adapted this technology for effective use in the safety and comfort of your own home. Philips worked closely with skin experts to develop this technology. For more than 10 years Philips carried out extensive consumer research on IPL technology with more than 2000 women.

## Enjoy smooth skin every day

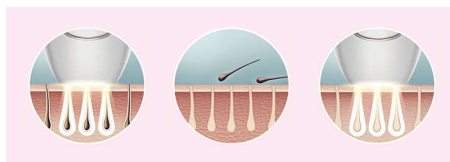

Philips Lumea applies gentle pulses of light to the hair root. As a consequence, the amount of hair your body grows gradually decreases. Repeating the treatment leaves your skin beautifully hair free and touchably smooth.

## Effortlessly effective

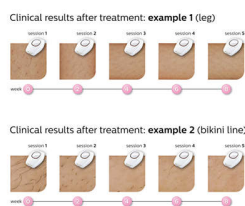

Our clinical studies have shown significant hair reduction in as little as four fortnightly treatments resulting in smooth hair free skin. To maintain these results, simply repeat the treatment when needed. Time between treatments may vary based on your individual

hair regrowth. - To enlarge this picture, please click on the image in the image gallery on the top of this page

## Suitable hair and skin types



Philips Lumea works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Philips Lumea is NOT effective on red, light-blond or white/grey hair. This product is also NOT suitable for dark skin.

## Application areas

Use with confidence to remove hair on face (upper lip, chin ) and body parts including legs, arm pits, bikini area, belly and arms.

## Skin Tone Sensor



For extra safety the integrated skin tone sensor measures the treated skin complexion at the beginning of each session and occasionally

during the session. If it detects a skin tone that is too dark for this product it will automatically stop emitting pulses.

## Gentle treatment



Gentle to your skin.

## Smallest and lightest Lumea


The Philips Lumea Essential is the smallest and lightest in the Lumea range. It is easy to hold and easy to reach all the areas you want to treat.

## No replacement parts needed



Philips Lumea is a complete solution that comes ready to work right out of the box. More importantly, it does not require any replacement lamps or gels.

## Slide and Flash mode

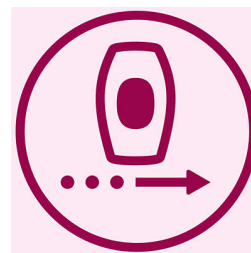

Main application mode enables sliding motion for a more convenient and continuous application.

# Specifications

**Application areas**

**Body areas:** Arms, Legs, Belly, Bikini, Armpits

**Face areas:** Chin, Sideburns, Upper lip

**Safety and adjustable settings**

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**Skin Tone sensor:** Detects your skin tone

**Application mode**

**Slide and Flash:** For easy motion application

**Corded/cordless use:** Corded

**Stamp and Flash:** For treatment on curvy areas

**Power**

**Voltage:** 100-240 V

**Technical specifications**

**High performance lamp:** Delivers <gt/> 200,000 flashes

**Service**

**Warranty:** 2 years global warranty

**Items included**

**Instructions for use:** User manual

**Adapter:** 24V / 1500mA

**Application time**

**Lower leg:** 11 min

**Armpit:** 1,5 min

**Bikini line:** 1,5 min

**Face areas:** 2,5 min

